

Method Validation Study for the Determination of Residues of XDE-208 and its Major
Metabolites in Water using Offline Solid-Phase Extraction and Liquid Chromatography with
Tandem Mass Spectrometry Detection

INTRODUCTION

Scope

This method is applicable for the quantitative determination of residues of XDE-208 (sulfoxaflor), and its major metabolites in water. The method was validated over the concentration range of 0.05-50.0 µg/L with a validated limit of quantitation of 0.05 µg/L. A procedure without solid phase extraction was validated also to allow the efficient analysis of samples which do not need to reach the levels mentioned above. The procedure without SPE was validated over the concentration ranges of 0.25-50.0 µg/L with a validated limit of quantitation of 0.25 µg/L. Common and chemical names, molecular formulas, and the nominal masses for the analyte and related compounds are given in Table 1.

This study was conducted to fulfill data requirements outlined in the U. S. EPA Residue Chemistry Test Guidelines, OPPTS 850.1700 (1), the European Commission Guidance Document on Residue Analytical Methods, SANCO/825/00 rev.7 (2) and SANCO/3029/99 (3), and PMRA Residue Chemistry Guidelines as Regulatory Directive Dir98-02 (4).

Method Principle

Residues of XDE-208 (sulfoxaflor), and its major metabolites are prepared for analysis by adding 1.0 N hydrochloric acid to the sample aliquot. If quantitating at levels below 0.25 µg/L is unneeded or if expected residue levels in samples are above 0.25 µg/L, a stable isotope internal standard mixture is added and the sample is filtered. The sample is analyzed by liquid chromatography with positive-ion electrospray ionization tandem mass spectrometry (LC/MS/MS). If quantitating at levels below 0.25 µg/L is needed or if residue levels in samples are unknown or expected below 0.25 µg/L, the solution is purified using a reverse-phase polymeric solid-phase extraction (SPE) column. Samples are eluted from the SPE column with an acetonitrile/water (80:20) solution with 0.1 % formic acid. Next, a stable isotope internal standard mixture is added and the eluate is diluted with a water/acetonitrile (95:5) with 0.01 % formic acid. The sample is analyzed by liquid chromatography with positive-ion electrospray ionization tandem mass spectrometry (LC/MS/MS).

Safety Precautions

Each analyst must be acquainted with the potential hazards of the equipment, reagents, products, solvents, and procedures used in this method before commencing laboratory work. SOURCES

OF INFORMATION INCLUDE: OPERATION MANUALS, MATERIAL SAFETY DATA SHEETS, LITERATURE, AND OTHER RELATED DATA. Safety information should be obtained from the supplier. Disposal of waste materials, reagents, reactants, and solvents must be in compliance applicable governmental requirements.

Acetonitrile and methanol are flammable and should be used in well-ventilated areas away from ignition sources. Formic acid and hydrochloric acid are corrosive and can cause severe burns. It is imperative that proper eye and personal protection equipment be worn when handling these reagents

Test Substance/Analytical Standard and Internal Standard

Test Substance/ Analytical Standard	TSN Number	Percent Purity	Certification Date	Reference
XDE-208 (sulfoxaflor) ^a	105878	99.7	10-Feb-2009	ML-AL-2008- 003623REV
X11519540 ^a	106498	98	05-Oct-2009	FA&PC 09-228627
X11579457 ^a	030941-0001	97	18-Jun-2008	FA&PC 08-184890
X11719474 ^a	030626-0003	99.5	14-May-2008	FA&PC 08-181812
X11843864 ^b	030721-0002	100	24-Feb-2010	FA&PC 09-236732
X11944782 ^b	031118-0001	97	16-Dec-2009	FA&PC 09-203293

^a Test substance/analytical standard

^b Internal standard

The above standards may be obtained from the Test Substance Coordinator, Dow AgroSciences LLC, 9330 Zionsville Road, Building 304, Indianapolis, IN 46268-1054.

Equipment, Glassware, and Materials

Equipment, glassware, materials, reagents, and chemicals considered to be equivalent to those specified may be substituted with the understanding that their performance must be confirmed by appropriate tests/proper scientific judgment. Common laboratory glassware and supplies are assumed to be readily available. Unless specified otherwise, class A volumetric glassware is used to prepare analytical standards, fortification solutions, and calibration standards. Suggested equipment and materials are listed in the following sections.

Laboratory Equipment

Balance, analytical, Model AE100, [Mettler-Toledo, Inc.](#)

Balance, pan, Model BB2440, Mettler-Toledo, Inc.

Pipet, positive-displacement, 20-50 μL capacity, catalog number M50, [Gilson Inc.](#)

Pipet, positive-displacement, 100-1000 μL capacity, catalog number M1000, Gilson Inc.

Vacuum manifold, Model spe-12G, [Mallinckrodt Baker, Inc.](#)

Vortex mixer, Model G-560, [Scientific Industries, Inc.](#)

Chromatographic System

Column, analytical, Zorbax SB-C8, 4.6 mm x 75 mm, 3.5- μm particle size, catalog number 866953-906, [Agilent Technologies.](#)

Liquid chromatograph, Symbiosis Pharma, [Spark Holland Inc.](#)

Mass spectrometer, Model API QTRAP 5500, [MDS/Sciex.](#)

Mass spectrometer data system, Model Analyst 1.5.1, MDS/Sciex.

Glassware and Materials

Column, Oasis HLB SPE, 60-mg sorbent, 3-mL reservoir, catalog number WAT094226, [Waters.](#)

Pipet tip, positive-displacement, 50- μL capacity, catalog number CP50, Gilson Inc.

Pipet tip, positive-displacement, 1000- μL capacity, catalog number CP1000, Gilson Inc.

Vial, autosampler, 2-mL, catalog number C4000-1W, [National Scientific Company.](#)

Vial, syringeless filter, 0.45- μm PTFE. catalog number 09-923-28, [Fisher Scientific.](#)

Vial, 7-mL, with PTFE-lined screw cap, catalog number 03-340-60B, Fisher Scientific.

Vial, 11-mL, with PTFE-lined screw cap, catalog number 03-340-60C, Fisher Scientific.

Vial cap, for autosampler vial, catalog number C4000-54B, National Scientific Company.

Reagents

Acetonitrile, HPLC grade, catalog number 2856-10, Mallinckrodt Baker.

Formic acid, 96%, ACS grade, catalog number 251364-500G, [Sigma-Aldrich.](#)

Hydrochloric acid, 0.01 N, certified concentration, catalog number SA62-1, Fisher Scientific.

Hydrochloric acid, 1.0 N, certified concentration, catalog number SA48-1, Fisher Scientific.

Methanol, HPLC grade, catalog number 3041-10, Mallinckrodt Baker.

Water, OmniSolv grade, catalog number WX-0004-1, [EMD Chemicals](#).

Prepared Solutions

acetonitrile containing 0.01% formic acid (mobile phase B)

Pipet 0.20 mL of formic acid into a 2000-mL graduated mixing cylinder containing approximately 1900 mL of acetonitrile. Dilute to volume with acetonitrile. Cap the cylinder and invert it multiple times to mix well prior to use.

acetonitrile/water (80:20) (v/v) containing 0.1% formic acid

Measure 800 mL of acetonitrile using a 1-L graduated cylinder. Add 1.0 mL of formic acid and approximately 150 mL of water into the cylinder and allow to equilibrate to room temperature. Dilute to volume with water.

water containing 0.01% formic acid (mobile phase A)

Pipet 0.20 mL of formic acid into a 2000-mL graduated mixing cylinder containing approximately 1900 mL of water. Dilute to volume with water. Cap the cylinder and invert it multiple times to mix well prior to use.

water/ acetonitrile (95/5) (v/v) containing 0.01% formic acid

Measure 50 mL of acetonitrile containing 0.01% formic acid (prepared above as mobile phase B) into a 1000-mL graduated cylinder. Add 900 mL of water containing 0.01% formic acid (prepared above as mobile phase A) and allow to equilibrate to room temperature. Dilute to volume with water containing 0.01% formic acid.

EXPERIMENTAL

Instrumental Conditions

Typical HPLC Operating Conditions

Instrumentation:	Spark Holland Symbiosis Pharma
Column:	
Primary	Agilent Zorbax SB-C8 4.6 x 75 mm, 3.5- μ m
Alternate	Phenomenex Synergi Hydro-RP 80A 4.6 x 75 mm, 4- μ m
Column Temperature:	ambient (approximately 22 °C)
Injection Volume:	10 μ L
Injection Wash	
Wash Port 1	700 μ L of acetonitrile/water/formic acid (80:20:0.1)
Wash Port 2	700 μ L of methanol
Wash Port 3	700 μ L of water
Run Time:	approximately 14 minutes
Mobile Phase:	A – water containing 0.01% formic acid B – acetonitrile containing 0.01% formic acid
Mobile Phase Split:	approximately 200 μ L/min split to source

Gradient:	Time (min)	Flow Rate (mL/min)	Solvent A (percent)	Solvent B (percent)
	0:01	1.00	100	0
	3:00	1.00	100	0
	8:00	1.00	0	100
	10:00	1.00	0	100
	11:00	1.00	100	0
	14:00	1.00	100	0

Flow Diverter	
Flow to Waste	0.0 min \rightarrow 5.0 min
Flow to Source	5.0 min \rightarrow 8.6 min
Flow to Waste	8.6 min \rightarrow end of run

Typical Mass Spectrometry Operating Conditions

Instrumentation: Applied Biosystems QTRAP 5500 MS System
 Applied Biosystems Analyst 1.5.1 data system

Ionization Mode: electrospray
 Polarity: positive
 Scan Type: MRM
 Resolution: Q1 – unit, Q3 – unit
 Curtain Gas (CUR): 40 psi
 Collision Gas (CAD): medium
 Ion Source Gas 1 (GS1): 40 psi
 Ion Source Gas 2 (GS2): 60 psi
 Temperature (TEM): 475 °C
 Entrance Potential (EP): 10 volts
 IonSpray Voltage (IS): 5000 volts

Acquisition Time Delay: 0.00 minutes
 Period Duration: 8.50 minutes
 Dwell Time: 50 ms

Analytes:	Precursor Ion, Q1	Product Ion, Q3	Declustering Potential, v	Collision Energy, v	Cell Exit Potential, v
X11422208					
quantitation	278.1	174.1	60	13	22
confirmation	278.1	154.1	60	38	18
X11719474					
quantitation	296.1	174.1	51	14	20
confirmation	296.1	105.1	51	42	19
X11519540					
quantitation	254.1	175.1	86	26	22
confirmation	254.1	154.1	86	55	20
X11579457					
quantitation	253.1	174.1	52	12	22
confirmation	253.1	154.1	52	35	20
X11843864 (M+3 ISTD) (XDE-208 stable isotope)					
quantitation	281.1	177.1	63	13	24
X11944782 (M+4 ISTD) (X11719474 stable isotope)					
quantitation	300.1	178.1	61	16	22

Representative spectra, calibration curves, and chromatograms are shown in Figures 1-21. Typical chromatograms for the determination and confirmation of XDE-208 (sulfoxaflor), and its major metabolites in water are illustrated in Figures 11-21.

Preparation of Standard Solutions

Preparation of XDE-208 (Sulfoxaflor) and Metabolite Stock Solutions

Weigh 0.1000 g of X11422208 (XDE-208) analytical standard and quantitatively transfer to a 100-mL volumetric flask with acetonitrile. Dilute to volume to obtain a 1000- $\mu\text{g}/\text{mL}$ stock solution.

Weigh 0.1000 g of X11519540 analytical standard and quantitatively transfer to a 100-mL volumetric flask with acetonitrile. Dilute to volume to obtain a 1000- $\mu\text{g}/\text{mL}$ stock solution.

Weigh 0.1000 g of X11579457 analytical standard and quantitatively transfer to a 100-mL volumetric flask with acetonitrile. Dilute to volume to obtain a 1000- $\mu\text{g}/\text{mL}$ stock solution.

Weigh 0.1000 g of X11719474 analytical standard and quantitatively transfer to a 100-mL volumetric flask with acetonitrile. Dilute to volume to obtain a 1000- $\mu\text{g}/\text{mL}$ stock solution.

Pipet 5.0 mL of the 1000- $\mu\text{g}/\text{mL}$ solutions above into a single 200-mL volumetric flask and dilute to volume with acetonitrile to obtain a mixed solution containing 25.0 $\mu\text{g}/\text{mL}$ of each compound.

Pipet 20.0 mL of the 25.0- $\mu\text{g}/\text{mL}$ solution above into a 200-mL volumetric flask and dilute to volume with acetonitrile to obtain a mixed solution containing 2.50 $\mu\text{g}/\text{mL}$ of each compound.

Pipet 20.0 mL of the 2.50- $\mu\text{g}/\text{mL}$ solution above into a 200-mL volumetric flask and dilute to volume with acetonitrile to obtain a mixed solution containing 0.250 $\mu\text{g}/\text{mL}$ of each compound.

Pipet 20.0 mL of the 0.250- $\mu\text{g}/\text{mL}$ solution above into a 200-mL volumetric flask and dilute to volume with acetonitrile to obtain a mixed solution containing 0.025 $\mu\text{g}/\text{mL}$ of each compound.

Pipet 20.0 mL of the 0.025- $\mu\text{g}/\text{mL}$ solution above into a 200-mL volumetric flask and dilute to volume with acetonitrile to obtain a mixed solution containing 0.0025 $\mu\text{g}/\text{mL}$ of each compound.

Prepare solutions for spiking samples by diluting the above stock solutions with acetonitrile as follows:

Concentration of Stock Soln. µg/mL	Aliquot of Stock Soln. mL	Final Soln. Volume mL	Spiking Soln. Final Conc. µg/mL	Equivalent Sample Conc. ^a µg/L
0.0025	60.0	100	0.0015	0.015
0.0025	--	--	0.0025	0.025
0.025	20.0	100	0.0050	0.050
0.025	50.0	100	0.0125	0.125
0.025	--	--	0.025	0.250
0.250	20.0	100	0.050	0.500
0.250	50.0	100	0.125	1.25
0.250	--	--	0.250	2.50
2.50	20.0	100	0.500	5.00
2.50	50.0	100	1.25	12.50
2.50	--	--	2.50	25.0
25.0	20.0	100	5.00	50.0
25.0	50.0	100	12.50	125.0
25.0	--	--	25.0	250.0

^a The equivalent sample concentration is based on fortifying a 10.0-mL water sample with 100 µL of spiking solution.

Preparation of the Mixed Stable-Isotope Internal Standard Solution

Weigh 0.0050 g of the X11843864 (XDE-208 M+3) stable-isotope internal standard and quantitatively transfer to a 50-mL volumetric flask with acetonitrile. Dilute to volume with acetonitrile to obtain a 100-µg/mL stock solution.

Weigh 0.0050 g of the X11944782 (X11719474 M+4) stable-isotope internal standard and quantitatively transfer to a 50-mL volumetric flask with acetonitrile. Dilute to volume with acetonitrile to obtain a 100-µg/mL stock solution.

Pipet 1.00 mL each of the 100-µg/mL solutions above into a single 100-mL volumetric flask and dilute to volume with acetonitrile to obtain a mixed solution containing 1.00 µg/mL of each compound.

Preparation of Calibration Standards for Samples

Prepare calibration standards by dispensing 250 µL of the 1.00 µg/mL mixed stable-isotope internal standard solution and 1000 µL of the 0.0015-1.25 µg/mL spiking solutions into a series of 20-mL volumetric flasks and dilute to volume with a water/acetonitrile solution (95:5) containing 0.01% formic acid.

The concentrations of the calibration standards are as follows:

Concentration of Spiking Soln. µg/mL	Aliquot of Spkg. Soln. mL	Final Soln. Volume mL	Calib Soln. Final Conc. ng/mL	Equivalent Sample Conc. µg/L
0.0015	1.00	20.0	0.075	0.015
0.0025	1.00	20.0	0.125	0.025
0.0050	1.00	20.0	0.250	0.050
0.0125	1.00	20.0	0.625	0.125
0.0250	1.00	20.0	1.25	0.250
0.0500	1.00	20.0	2.50	0.500
0.125	1.00	20.0	6.25	1.25
0.250	1.00	20.0	12.5	2.50
0.500	1.00	20.0	25.0	5.00
1.25	1.00	20.0	62.5	12.50

^a Conversion from equivalent sample concentration to final concentration is outlined in the example calculations (Figures 22-23).

Preparation of XDE-208 (X11422208) Standards to Determine Isotopic Crossover

Using a 1000-µL syringe or positive-displacement pipet, dispense 1000 µL of the 0.250 µg/mL X11422208 mixed spiking solution into a 20-mL volumetric flask and dilute to volume with a water/acetonitrile solution (95:5) containing 0.01% formic acid. The resulting solution contains 12.5 ng/mL of a mixed X11422208 solution.

Preparation of X11843864 Internal Standards to Determine Isotopic Crossover

Using a 250-µL syringe or positive-displacement pipet, dispense 250 µL of the 1.00 µg/mL XDE-208 mixed internal standard solution into a 20-mL volumetric flask and dilute to volume with a water/acetonitrile solution (95:5) containing 0.01% formic acid. The resulting solution contains 12.5 ng/mL of a mixed X11422208 internal standard solution.

Sample Origin, Numbering, Preparation and Storage

Untreated control samples of the water were obtained from the Dow AgroSciences LLC Sample Management Group. All samples were tracked in the Dow AgroSciences LLC Regulatory Labs Information Management System (RLIMS) database. Unique sample numbers were assigned to the samples to track them during receipt, storage, and analysis. The chemical and physical properties of the water samples used in the validation study are listed below. Complete source documentation is included in the study file.

Sample Group Number	Water Type	pH	Hardness (mg equiv. CaCO ₃ /L)	Total Suspended Solids (ppm)	Alkalinity (mg CaCO ₃ /L)	Total Organic Carbon (ppm)	Dissolved Organic Carbon (ppm)
189	(Tap) Drinking Water	7.9	351	10	252	9.4	9.1
190	(Well) Ground Water	7.9	441	6	340	1.6	1.3
191	(Pond) Surface Water	8.3	255	6	137	11.6	7.3
192	(Pond) Surface Water	7.9	262	10	188	10.6	9.4
193	(Well) Ground Water	8.0	324	24	298	7.6	5.3
194	(Tap) Drinking Water	8.2	249	4	202	2.5	2.3

No sample preparation was required for the water samples prior to analysis. Samples were stored refrigerated at approximately 4 °C after their time of sampling and during the course of the method validation study, except when they were removed for taking aliquots for sample analysis. For the purposes of this study, tap water was

Analysis Procedure

If quantitation at levels below 0.25 µg/L is unnecessary or residue levels in samples are expected above 0.25 µg/L, the procedure without SPE purification may be used for decreased analysis time. If quantitation at levels below 0.25 µg/L is needed or residue levels in samples are unknown or expected below 0.25 µg/L, the procedure with SPE purification must be followed.

Analysis Procedure Without SPE Purification

1. Pipet 4.0 mL or weigh 4.0-gram portions of the control water sample into a series of 7-mL glass vials.
2. Add 40 µL of 1.0 N hydrochloric acid to the sample vial.
3. For preparing fortified samples, dispense 40 µL of the appropriate spiking solutions into the sample vial to encompass the necessary concentration range.
4. Using a positive-displacement pipet, dispense 50 µL of the 1.00 µg/mL mixed internal standard solution into the sample vial.
5. Cap the sample vial with a PTFE-lined cap and then mix the sample for 1-2 seconds.

6. If necessary, filter the samples with a Whatman syringeless 0.45- μ m filter device.
7. Analyze the samples and calibration standards by LC/MS/MS with positive-ion electrospray tandem mass spectrometry.

Analysis Procedure With SPE Purification

1. Pipet 10.0 mL or weigh 10.0-gram portions of the control water sample into a series of 11-mL glass vials.
2. Add 100 μ L of 1.0 N hydrochloric acid to each sample vial.
3. For preparing fortified samples, dispense 100- μ L aliquots of the appropriate spiking solutions into the sample vial to encompass the necessary concentration range.
4. Cap the sample vial with a PTFE-lined cap and then mix the sample for 1-2 seconds.
5. Purify the sample using the following SPE procedure:
 - a. Place an OASIS[®] HLB SPE column (3-mL, 60 mg) on the vacuum manifold.
 - b. Condition the SPE column with 1 mL of acetonitrile followed by 1 mL of 0.01 N hydrochloric acid. Dry the SPE column under full vacuum (\approx -10 inches Hg) for 5 seconds between solvents.
 - c. Transfer the sample from Step 9.2.4 to the SPE column. Draw the sample through the column at a flow rate of approximately 1 mL/min, discarding the eluate.
 - d. Rinse the sample vial with 1 mL of 0.01 N hydrochloric acid and transfer to the SPE column. Draw the solvent through the column at a flow rate of approximately 1 mL/min, discarding the eluate. Dry the column under full vacuum (\approx -10 inches Hg) for approximately 60 seconds.
 - e. Elute the XDE-208 and metabolites from the SPE column with two 500- μ L aliquots of an acetonitrile/water solution (80:20) containing 0.10% formic acid, collecting the eluate in a 7-mL vial.
6. Using a positive-displacement pipet, dispense 25 μ L of the 1.00 μ g/mL mixed internal standard solution into the sample vial.
7. Adjust the volume in the sample vial to 2.0 mL with a water/acetonitrile solution (95:5) containing 0.01% formic acid by comparing the height of the solution in the sample vial to identical vials containing 2.0 mL of the water/acetonitrile solution (95:5) containing 0.01% formic acid.
8. Transfer the sample to a 2-mL autosampler vial or a 96-well plate.
9. Analyze the samples and calibration standards by LC/MS/MS with positive-ion electrospray tandem mass spectrometry.

Calculations

For XDE-208 and X11719474

Inject the series of calibration standards described in the calibration standard preparation section using the conditions listed in the instrument section and determine the peak areas for the analyte and internal standard.

For each sample and standard calculate the quantitation ratio (analyte peak area/internal standard peak area).

Prepare a standard curve using linear regression analysis with 1/x weighting by plotting the equivalent analyte concentration on the abscissa (x-axis) and the respective quantitation ratio on the ordinate (y-axis) as shown in Figures 7-8. Power regression or quadratic curve fit may also be used if appropriate.

Determine the concentration ($\mu\text{g/L}$) and/or recovery (%) from the sample as described in the example calculation outlined in Figure 22.

For X11519540 and X11579457

Inject the series of calibration standards described in the calibration standard preparation section using the conditions listed in the instrument section and determine the peak areas for the analyte.

Prepare a standard curve using linear regression analysis with 1/x weighting by plotting the equivalent analyte concentration on the abscissa (x-axis) and the respective peak area on the ordinate (y-axis) as shown in Figures 9-10. Power regression or quadratic curve fit may also be used if appropriate.

Determine the concentration ($\mu\text{g/L}$) and/or recovery (%) from the sample as described in the example calculation outlined in Figure 23.

Confirmation of Residue Identity

The method is specific for the determination of XDE-208 (sulfoxaflor), and its major metabolites X11519540, X11579457, and X11719474 by virtue of the chromatographic separation and selective detection system used. To demonstrate confirmation, a total of two MS/MS ion

transitions are monitored for each analyte. The following example is given for the analyte XDE-208.

Inject the series of calibration standards described in the calibration standard preparation section using the conditions listed in the instrument section and determine the peak areas for the analyte as indicated below.

XDE-208 (X1142208)	Q1/Q3 m/z 278/174 (quantitation)
	Q1/Q3 m/z 278/154 (confirmation)

For each standard, calculate the respective confirmation ratios.

$$\text{Confirmation Ratio} = \frac{\text{peak area of confirmation ion transition}}{\text{peak area of quantitation ion transition}}$$
$$\text{Confirmation Ratio} = \frac{\text{XDE - 208 peak area at } m/z \text{ 278/154}}{\text{XDE - 208 peak area at } m/z \text{ 278/174}}$$

For example, using the data for XDE-208 (X1142208) confirmation from the 0.125-ng/mL standard, from set 091185 S04, found in Figure 12:

$$\text{Confirmation Ratio} = \frac{4208}{7369} = 0.5710$$

Confirmation of the presence of the analyte is indicated when the retention time of the samples matches that of the standards and the confirmation ratio is in the range of $\pm 20\%$ of the average found for the standards. The confirmation MS/MS ion transition area ratios for each recovery sample at or above the LOQ were all within the range of $\pm 20\%$ of the average found for the standards within the sets with only one exception (20.6%). Example chromatograms in Figures 11-21 include both quantitation and confirmation transitions.

Determination of Isotopic Crossover

In this assay, analytes and internal standards are quantitated using MS/MS transitions characteristic of each compound. When using stable-isotope labeled internal standards, there is a possibility that isotopic contributions will occur between the transitions used for quantitation of the unlabeled and labeled compounds. This isotopic overlap between the analyte and the internal standard is determined empirically by analyzing standard solutions of each compound and should be addressed for accurate determination of concentrations (5, 6).

To determine the isotopic crossover for each analyte and respective stable isotope internal standard, inject the crossover standards described in the Preparation of Standard Solutions section, and determine the peak areas for the analyte and internal standard.

For example, to determine the contribution of the unlabeled XDE-208 to X11843864 (XDE-208-[M+3] stable isotope) using the XDE-208 standard crossover data from Figure 11 for the quantitation ion:

$$\text{Crossover Factor (analyte} \rightarrow \text{ISTD)} = \frac{\text{peak area of internal standard transition}}{\text{peak area of analyte transition}}$$

$$\text{Crossover Factor (analyte} \rightarrow \text{ISTD)} = \frac{\text{peak area at } m/z \text{ 281/177}}{\text{peak area at } m/z \text{ 278/174}}$$

$$\text{Crossover Factor (analyte} \rightarrow \text{ISTD)} = \frac{0}{656387} = 0$$

In a similar manner, to determine the contribution of X11843864 to the unlabeled XDE-208 using the internal standard crossover data from Figure 11 for the quantitation ion:

$$\text{Crossover Factor (ISTD} \rightarrow \text{analyte)} = \frac{\text{peak area of analyte transition}}{\text{peak area of internal standard transition}}$$

$$\text{Crossover Factor (ISTD} \rightarrow \text{analyte)} = \frac{\text{peak area at } m/z \text{ 278/174}}{\text{peak area at } m/z \text{ 281/177}}$$

$$\text{Crossover Factor (ISTD} \rightarrow \text{analyte)} = \frac{0}{656132} = 0$$

During the validation study, no significant mass spectral isotopic crossover was observed for quantitation or confirmation transitions for the two analytes with stable isotope internal standards and therefore no correction was performed. If isotopic crossover is encountered it should be assessed and the respective quantitation ratios corrected for accurate determination of concentrations.

Statistical Treatment of Data

Statistical treatment of data included the calculation of regression equations, coefficients of determination (r^2) for describing the linearity of calibration curves, and means, standard deviations, and relative standard deviations of the results for the fortified recovery samples.

Standardization of SPE Elution Profile

There is a possibility that variation in the Oasis HLB SPE columns may influence the elution profile of the analytes. If it is necessary to obtain an elution profile for each lot of SPE columns used to ensure optimum recovery and clean-up efficiency, the following procedure can be used:

1. To a 11-mL vial containing 10.0 mL of 0.01 N hydrochloric acid, add 100 μ L of the 0.025- μ g/mL mixed spiking solution and vortex mix.
2. Place an OASIS[®] HLB SPE column (3-mL, 60 mg) on the vacuum manifold.

3. Condition the SPE column with 1 mL of acetonitrile followed by 1 mL of 0.01 N hydrochloric acid. Dry the SPE column under full vacuum (\approx -10 inches Hg) for 5 seconds between solvents.
4. Transfer the sample from Step 1 to the SPE column. Draw the sample through the column at a flow rate of approximately 1 mL/min, discarding the eluate.
5. Rinse the sample vial with 1 mL of 0.01 N hydrochloric acid and transfer to the SPE column. Draw the solvent through the column at a flow rate of approximately 1 mL/min, discarding the eluate. Dry the column under full vacuum (\approx -10 inches Hg) for approximately 60 seconds.
6. Elute the XDE-208 and metabolites from the SPE column with four separate 500- μ L aliquots of an acetonitrile/water solution (80:20) containing 0.10% formic acid, collecting the each fraction in a separate 7-mL vial.
7. For each fraction collected, dispense 25 μ L of the 1.00 μ g/mL mixed internal standard solution into the sample vial.
8. Adjust the volume in the sample vial to 2.00 mL with a water/acetonitrile solution (95:5) containing 0.01% formic acid by comparing the height of the solution in the sample vial to identical vials containing 2.00 mL of the water/acetonitrile solution (95:5) containing 0.01% formic acid.
9. Transfer the sample to a 2-mL autosampler vial or a 96-well plate.
10. Analyze the samples and calibration standards by LC/MS/MS with positive-ion electrospray tandem mass spectrometry.
11. Calculate the percent recovery as described in the Calculations section.

A typical elution profile is illustrated in Figure 24. If the elution profile differs from that shown, adjust the volume of the acetonitrile/water solution (80:20) containing 0.10% formic acid solution to be collected in Step 5 of the Analysis Procedure with SPE Purification section.

Supplemental Notes

1. The instrumental conditions may be modified to obtain optimal chromatographic separation and sensitivity.
2. Aliquot volumes may differ if different stock solution concentrations are available. Final volume of prepared solutions can vary accordingly with use demand and/or reagent availability.
3. Each instrument should be tested for carryover by injecting a reagent blank after the highest calibration standard solution. After evaluation of the result, the analyst may decide to avoid

injecting high calibration standard solutions or treated samples before untreated samples.

4. The type of regression model can be chosen to give the best fit (coefficient of determination) for the data.

5. If the sample extracts contain analyte concentrations that exceed the linear range of the standard calibration curve (approximately 62.5 ng/mL, equivalent to 62.50 µg/L in the original sample), dilute those samples with an appropriate amount of water/acetonitrile solution (95:5) containing 0.01% formic acid to obtain responses at least 20% less than the highest calibration standard.

CONCLUSION

This method is applicable for the quantitative determination of residues of XDE-208 (sulfoxaflor), and its major metabolites X11519540, X11579457, and X11719474 in water. The method was validated over the concentration ranges of 0.05-50.0 µg/L with a validated limit of quantitation of 0.05 µg/L. A procedure without solid phase extraction was validated also to allow the efficient analysis of samples which do not need to reach the levels mentioned above. The procedure without SPE was validated over the concentration ranges of 0.25-50.0 µg/L with a validated limit of quantitation of 0.25 µg/L.

ARCHIVING

The protocol, raw data, and the original version of the final report are all filed in the Dow AgroSciences LLC archives at 9330 Zionsville Road in Indianapolis, Indiana 46268-1054.

REFERENCES

1. *Residue Chemistry Test Guidelines, OPPTS 850.7100, Data Reporting for Environmental Chemistry Methods*; EPA-712-C-96-348; U.S. Environmental Protection Agency. U.S. Government Printing Office: Washington, DC, 1996.
2. *Guidance Document on Residue Analytical Methods-Sanco.825/00 rev. 7*. Directorate General Health and Consumer Protection, March 17, 2004.
3. *Guidance for Generating and Reporting Methods of Analysis in Support of Pre-registration Data Requirements for Annex II (part A, Section 4) and Annex III (part A, Section 5) of Directive 91/414-Sanco.3029/99 rev. 4*. Directorate General Health and Consumer Protection, July 11, 2000.
4. *Residue Chemistry Guidelines, Regulatory Directive Dir98-02, Section 3, Residue Analytical Method*, Submission Management and Information Division, Pest Management Regulatory Agency, Health Canada: Ottawa, Ontario, June 1, 1998.
5. Jenden, D. J.; Roch, M.; Booth, R. A. *Anal. Biochem.* **1973**, 55 438-448.
6. Barbalas, M. P.; Garland, W. A. *J. Pharm. Sci.* **1991**, 80(10), 922-927.
7. Keith, L. H.; Crummett, W.; Deegan, J., Jr.; Libby, R. A.; Taylor, J. K.; Wentler, G. *Anal. Chem.* 1983, 55, 2210-2218.

The information herein is presented in good faith, but no warranty, express or implied, is given nor is freedom from any patent owned by Dow AgroSciences LLC or by others to be inferred. In the hands of qualified personnel, the procedures are expected to yield results of sufficient accuracy for their intended purposes, but recipients are cautioned to confirm the reliability of their techniques, equipment, and standards by appropriate tests. Anyone wishing to reproduce or publish the material in whole or in part should request written permission from Dow AgroSciences LLC.

Table 1. Identity and Structures of XDE-208 and Related Compounds

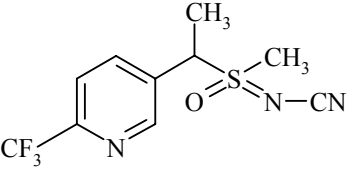
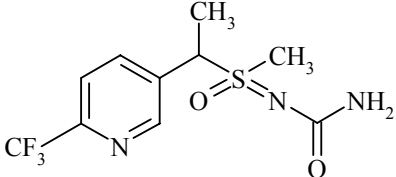
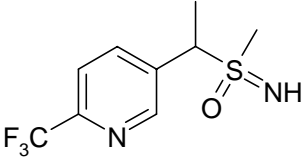
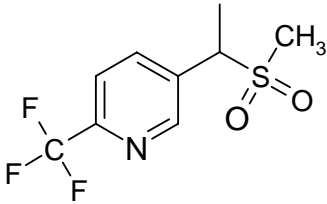
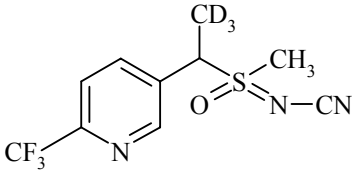
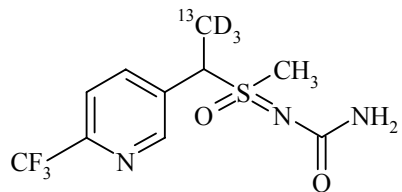
Common Name of Compound	Structural Formula and Chemical Name
<p>XDE-208 (sulfoxaflor)</p> <p>Molecular Formula: C₁₀H₁₀F₃N₃OS</p> <p>Formula Weight: 277.27</p> <p>Nominal Mass: 277</p> <p>CAS Number 946578-00-3</p>	 <p>[1-(6-trifluoromethylpyridin-3-yl)ethyl](methyl)-oxido-λ⁴-sulfanylidene cyanamide</p>
<p>X11719474</p> <p>Molecular Formula: C₁₀H₁₂F₃N₃O₂S</p> <p>Formula Weight: 295.29</p> <p>Nominal Mass: 295</p> <p>CAS Number not available</p>	 <p>N-(methyl)oxido{1-[6-(trifluoromethyl)pyridine-3-yl]ethyl}-λ⁴-sulfanylidene) urea</p>
<p>X11579457</p> <p>Molecular Formula: C₉H₁₁F₃N₂OS</p> <p>Formula Weight: 252.26</p> <p>Nominal Mass: 252</p> <p>CAS Number not available</p>	 <p>[5-[1-(S-methylsulfonimidoyl)ethyl]-2-(trifluoromethyl)pyridine</p>
<p>X11519540</p> <p>Molecular Formula: C₉H₁₀F₃NO₂S</p> <p>Formula Weight: 253.24</p> <p>Nominal Mass: 253</p> <p>CAS Number not available</p>	 <p>5-(1-methylsulfonyl)ethyl)-2-(trifluoromethyl)pyridine</p>

Table 1 (Cont.) Identity and Structures of XDE-208 and Related Compounds

Common Name of Compound	Structural Formula and Chemical Name
<p>X11843864 (XDE-208 M+3) stable-isotope</p> <p>Molecular Formula: C₁₀D₃H₇F₃N₃OS Formula Weight: 280.29 Nominal Mass: 280</p> <p>CAS Number not available</p>	 <p>[1-(6-trifluoromethylpyridin-3-yl)-2,2,2-²H₃-ethyl](methyl)-oxido-λ⁴-sulfanylidene cyanamide</p>
<p>X11944782 (X11719474 M+4) stable-isotope</p> <p>Molecular Formula: ¹³CC₉D₃H₉F₃N₃O₂S Formula Weight: 299.30 Nominal Mass: 299</p> <p>CAS Number not available</p>	 <p>N-(methyl)oxido{2-¹³C-2,2,2-²H₃-1-[6-(trifluoromethyl)pyridine-3-yl]ethyl}-λ⁴-sulfanylidene urea</p>