

ATTACHMENT A

LABORATORY STANDARD OPERATING PROCEDURES

ATTACHMENT A-1
COLUMBIA ANALYTICAL SERVICES, INC.
STANDARD OPERATING PROCEDURES

EXT-3550BPines, Rev. 0	Ultrasonic Extraction of Solids for 8270C SIM Analysis for Indiana Pines
8270CSIMPines, Rev. 0	Analysis of PAHs by Gas Chromatography/Mass Spectroscopy Using Selective Ion Monitoring (GC/MS/SIM) for Indiana Pines Site
HRMS-8290, Rev. 5.2	Analysis of Polychlorinated Dibenzo-p-dioxins and Polychlorinated Dibenzofurans by High-Resolution Gas Chromatography/High-Resolution Mass Spectroscopy (HRGC/HRMS)
MET-3050Pines, Rev. 0	Metals Digestion, Soils, Sediments, and Sludge for ICP Analysis for the Pines Indiana Site
MET-6010BPines, Rev. 1	Determination of Metals and Trace Elements by Inductively Coupled Plasma Atomic Emission Spectroscopy (ICP) for Indiana Pines Site
MET-7471APines, Rev. 0	Determination of Mercury in Solid or Semisolid Waste by Cold Vapor Atomic Absorption Spectroscopy
MET-3050B, Rev. 3	Metals Digestion, Soils, Sediments, and Sludge for ICP and GFAA Analysis
MET-GFAA, Rev. 3	Determination of Trace Metals by Graphite Furnace Atomic Absorption Spectroscopy (GFAAS)
MET-ICSPines, Rev. 0	Total Sulfur for Ion Chromatography for Indiana Pines Site
GEN-300Pines, Rev. 1	Determination of Anions Using Ion Chromatography for Indiana Pines Site
ADM-MDL, Rev. 5	The Determination of Method Detection Limits
SMO-GEN, Rev. 2	Sample Receiving

STANDARD OPERATING PROCEDURE
for
**ULTRASONIC EXTRACTION OF SOLIDS
FOR 8270C SIM ANALYSIS
FOR INDIANA PINES**

SOP No.: EXT-3550BPines

Revision: 0

December 3, 2004

Approved by: *Michael D. Calk*
Supervisor
Vicky St. Aub
QA Coordinator
Michael K. Perry
Laboratory Manager

12/6/04
Date
12/3/04
Date
12/9/04
Date

© Columbia Analytical Services, Inc., 2004
One Mustard St., Suite 250
Rochester, NY 14609

Annual review of this SOP has been performed
and the SOP still reflects current practice.

Initials: _____ Date: _____
Initials: _____ Date: _____
Initials: _____ Date: _____

NON-CONTROLLED COPY
Will Not Be Updated

Table of Contents	<u>Page</u>
1. Scope and Applicability.....	3
2. Summary of Method.....	3
3. Definitions.....	3
4. Health and Safety Warnings.....	4
5. Cautions.....	4
6. Interferences.....	4
7. Personnel Qualifications.....	4
8. Equipment and Supplies.....	5
9. Procedure.....	6
9.1. Calibration and Standardization.....	6
9.2. Sample Collection.....	6
9.3. Sample Handling and Preservation.....	7
9.4. Sample Preparation	7
9.5. Troubleshooting.....	8
9.6. Data Acquisition, Calculations, and Data Reduction Requirements.....	9
10. Data and Records Management.....	9
11. Quality Control and Quality Assurance.....	9
12. References.....	9
Attachments	
Extraction Benchsheet.....	10
Flow Chart of Extraction.....	11

1. SCOPE AND APPLICABILITY

- 1.1. This SOP uses EPA SW-846 Method 3550B for extracting PAHs from solids such as soils, sludges, and wastes. The ultrasonic process ensures intimate contact of the sample matrix with the extraction solvent.
- 1.2. This SOP was modified specifically for PAH analysis by 8270C in SIM mode for the Indiana Pines Site project.

2. METHOD SUMMARY

A 30 g sample is mixed with anhydrous sodium sulfate to form a free flowing powder. This is solvent extracted three times using ultrasonic extraction. The extract is then concentrated and a portion of the concentrate is removed for analysis.

3. DEFINITIONS

- 3.1. **Extraction batch** – a group of no more than 20 field samples extracted on the same day with the same reagents under the same conditions.
- 3.2. **Laboratory Control Sample (LCS)** - An aliquot of sodium sulfate to which a known quantity of the method analyte is added in the laboratory. The LCS is analyzed exactly like a sample, and its purpose is to determine whether the methodology is in control, and whether the laboratory is capable of making accurate and precise measurements.
- 3.3. **Matrix** - the predominant material, component, or substrate (e.g., soil, sludge, etc.) of which the sample to be analyzed is composed.
- 3.4. **Matrix Spike (MS/MSD)** - In the matrix spike analysis, a predetermined quantity of a standard solution of the analytes of interest is added to a sample matrix prior to sample extraction. The purpose of the matrix spike is to evaluate the effects of the sample matrix on the methods used for the analyses. Percent recovery is calculated for the analytes detected to measure accuracy. The RPD between the MS and MSD is calculated to measure precision.
- 3.5. **Method Blank (MB)** – an artificial sample known to be free of the analytes of interest. Used to measure contamination introduced during extraction and analysis.
- 3.6. **Sample** - a portion of material to be analyzed that is contained in single or multiple containers and identified by a unique sample number.
- 3.7. **Surrogates (Surrogate Standards)** - an organic compound which is similar to the target analyte(s) in chemical composition and behavior in the analytical process. Surrogate compounds are added to every blank, sample, matrix spike, matrix spike duplicate, LCS,

matrix spike blank, and standard. These are used to evaluate analytical efficiency by measuring recovery. Surrogates are not expected to be detected in environmental media.

4. HEALTH AND SAFETY WARNINGS

The toxicity or carcinogenicity of each reagent used in this method has not been precisely determined; however, each chemical should be treated as a potential health hazard. Exposure to these reagents should be reduced to the lowest possible level. The laboratory is responsible for maintaining a current awareness file of OSHA regulations regarding the safe handling of the chemicals specified in this method. Concentrated sulfuric acid and the 50% sodium hydroxide solution are moderately toxic and extremely irritating to skin and mucous membranes. Use these reagents in a fume hood whenever possible. If eye or skin contact occurs, flush with large volumes of water. Always wear safety glasses or a shield for eye protection, protective clothing, and observe proper mixing when working with these reagents.

5. CAUTIONS

- *An open container of sodium sulfate may become contaminated during storage in the laboratory.*
- *Avoid concentrating the sample to less than 1 ml. When the volume of solvent is reduced below 1 ml, analytes may be lost*

6. INTERFERENCES

- 6.1. Solvents, reagents, glassware, and other sample processing hardware may yield artifacts and/or interferences to sample analysis. All these materials must be demonstrated to be free from interferences under the conditions of the analysis by analyzing method blanks.
- 6.2. Interferences coextracted from the samples will vary considerably from source to source.
- 6.3. Phthalate esters contaminate many types of products commonly found in the laboratory. Plastics, in particular, must be avoided because phthalates are commonly used as plasticizers and are easily extracted from plastic materials.

7. PERSONNEL QUALIFICATIONS

At a minimum, personnel must have attained at least a 4-year degree (or 2-yr degree plus one year experience) in a science-related field and have successfully completed an Initial Demonstration of Capability and the Training Plan Form (attached). Training and Demonstration of Capability are in accordance with NELAC 2002 standard.

8. EQUIPMENT AND SUPPLIES

- 8.1. Ultrasonic Disrupter - A horn type device equipped with a titanium tip, or a device that will give equivalent performance. The horn should be tuned prior to sample extraction. The disrupter must have a minimum power wattage of 300 watts, with pulsing capability. Follow the manufacturers instructions for preparing the disrupter for extraction of samples with low concentration. Use a 1/2" horn.
- 8.2. Beakers - 400 ml.
- 8.3. Vacuum filtration apparatus.
 - 8.3.1. Buchner funnel.
 - 8.3.2. Filter paper - Whatman No. 41 or equivalent.
 - 8.3.3. Side- Arm Flask - 1000 mL Pyrex.
- 8.4. Erlenmeyer flasks - 250 mL and 500 mL.
- 8.5. Kuderna-Danish (K-D) apparatus.
 - 8.5.1. Concentrator tube - 10 ml, graduated (Kontes K-570050-1025 or equivalent).
 - 8.5.2. Evaporation flask - 250 mL and 500 mL (Kontes K-570001-500 or equivalent). Attach to concentrator tube with clamps.
 - 8.5.3. Snyder column - Three ball macro (Kontes K-503000-0121 or equivalent).
 - 8.5.4. Plastic Clips - 19/24 Joint (KT675300-0019 or equivalent).
- 8.6. Glass boiling beads - 3 mm glass beads.
- 8.7. Water bath - Heated, with concentric ring cover, capable of temperature control ($\pm 5^{\circ}\text{C}$). The bath should be used in a hood.
- 8.8. Vials and Caps - 20 mL scintillation vials with screw caps and aluminum foil liner, 2 mL capacity glass with crimp tops for GC auto sampler.
- 8.9. Pipets - glass volumetric 1 mL or 2 mL.
- 8.10. 5 inch Pyrex funnel with a small pad of Pyrex glass wool.
- 8.11. Volumetric flasks - 10 mL and 1 and 2 mL

- 8.12. Balance - Top loading, capable of accurately weighing to the nearest 0.01g.
- 8.13. Tongue Depressors - 6 inch standard.
- 8.14. Oven - drying.
- 8.15. LabConco RapidVap Evaporation System
- 600mL Sample Tubes – End point Stem 1.5mL
 - 8 place rack for 600ml Sample tubes
- 8.16. Standards, LCS, MS/MSD, and surrogate spiking solutions – concentrations and recipes are discussed in the analytical SOP.
- 8.17. Sodium sulfate - granular - anhydrous reagent grade, heated at 120°C for 16 hours, and stored in a glass bottle. Store at room temperature. Expires upon manufacturer's indications or 3 years from receipt if no indication is provided.
- CAUTION: An open container of sodium sulfate may become contaminated during storage in the laboratory.*
- 8.18. Extraction solvent - .Methylene chloride - pesticide quality or equivalent. Store at room temperature. Expires upon manufacturer's indications or 3 years from receipt if no indication is provided. Boiling point 39°C.
- 8.19. Exchange solvent - .Hexane, - Pesticide quality or equivalent. Store at room temperature. Expires upon manufacturer's indications or 3 years from receipt if no indication is provided. Boiling point 68.7°C.

9. PROCEDURE

- 9.1. **Calibration and Standardization** – Tune the horn according to manufacturer's instructions prior to sample extraction.
- 9.2. **Sample Collection**

Purchased, precleaned, certified sample containers should be glass or Teflon, and have screw-caps with Teflon lined septa. In situations where Teflon is not available, solvent-rinsed aluminum foil may be used as a liner. However, acidic or basic samples may react with the aluminum foil, causing eventual contamination of the sample. Plastic containers or lids may NOT be used for the storage of samples due to the possibility of sample contamination from the phthalate esters and other hydrocarbons within the plastic. Sample containers should be filled with care so as to prevent any portion of the collected sample coming in contact with the sampler's gloves, thus causing contamination. Samples shall be stored at 0-6°C and shipped to the laboratory within 48 hrs.

9.3. Sample Handling and Preservation

Sample collection, preservation, and custody management is in accordance with NELAC 2002 Standard.

Soil samples must be iced or refrigerated at 0-6°C from the time of collection until extraction. Extract samples within 14 days from collection. Store sample extracts in extract coolers at 0-6°C and analyze within 40 days of extraction.

9.4. Sample Preparation

- 9.4.1. Sediment/soil samples - Decant and discard any water layer on a sediment sample. Mix sample thoroughly, especially composited samples. Discard any foreign objects such as sticks, leaves, and rocks.
- 9.4.2. Cut, shred, or otherwise break down gummy, fibrous or clay-like materials to allow mixing and maximum exposure of the sample surfaces for extraction. The professional judgment of the analyst is required for handling of these difficult matrices.
- 9.4.3. The following step should be performed rapidly to avoid loss of the more volatile extractables: Weigh approximately 30 g of sample into a 400 ml beaker. Record the weigh to the nearest 0.1 g. Add 60 g of sodium sulfate and mix thoroughly. Nonporous or wet samples (gummy or clay type) that do not have a freeflowing sandy texture may be mixed with more sodium sulfate if needed. After addition of sodium sulfate, the sample should be free flowing. Add 1 mL of 1 ppm BN surrogate in methanol to all samples, spikes, LCS, and blanks. Add 1 mL of 1 ppm PAH spike in methanol to the LCS, MS, and MSD. Immediately add 100 mL of Methylene chloride.
- 9.4.4. Place the bottom surface of the tip of the disrupter horn about 1/2 in. below the surface of the solvent, but above the sediment layer.
- 9.4.5. Extract ultrasonically for 3 minutes, with output control knob set at 10 (full power) and with mode switch on Pulse (pulsing energy rather than continuous energy) and percent-duty cycle knob set at 50%. Do not use microtip probe.
- 9.4.6. Decant and filter extracts through a Buchner funnel, using Whatman No. 41 filter paper and vacuum filtration, into a sidearm flask.
- 9.4.7. Repeat the extraction two more times with two additional 100 ml portions of solvent. On the final ultrasonic extraction, pour the entire sample into the Buchner funnel and rinse with extraction solvent. Keep the solvent extract and discard the sample on the filter.

9.4.8. Dry and Concentrate Extract with LabConco RapidVap Nitrogen Evaporation System

9.4.8.1. Add extract to a 600 mL Sample tube

9.4.8.2. Turn on power switch

9.4.8.3. Set parameters – Input the % Speed, Temperature, Time and Number of Samples to be concentrated. (See users manual for Time, Temperature and Speed Setting Guidelines per volume per solvent.)

9.4.8.4. Turn on Nitrogen gas.

9.4.8.5. Load sample tubes into RapidVap chamber. Close lid and swing both latched over the lid and tighten knobs.

9.4.8.6. Press Run to begin concentration. Be sure to check the samples at the beginning of the program to ensure samples are not rotating so quickly that they are splashing out of the sample tube. If this occurs, press the Stop button and adjust the speed. To resume press Run.

9.4.8.7. Allow samples to evaporate to a final volume of 1 mLs. When the time in the program is complete the instrument will beep and stop. Check sample volume. If additional time is needed, reset program and press Run.

CAUTION: Avoid concentrating the sample to less than 1 ml. When the volume of solvent is reduced below 1 ml, analytes may be lost.

9.4.8.8. The final volume is transferred to a vial with a Teflon lined screw-cap or crimp top, and label appropriately. Store extract in extract cooler at temperature of 0-6°C.

9.4.9. At this point samples are ready to be analyzed for the target analytes using the SOP 8270CSIMPines.

9.5. **Troubleshooting**

9.5.1. Ultrasonic probes - Pitting of the probe tip will occur with use. When the bottom begins to cup, the probe needs to be sent out to be machined. Eventually, the probe will need to be replaced.

9.5.2. Bath / RapidVap temperature should be monitored. It is not to exceed 15-20°C of the boiling point of the solvent being evaporated.

9.6. Data Acquisition, Calculations, and Data Reduction Requirements

9.6.1. Be sure all documentation is complete and legible.

9.6.2. All benchsheets shall be maintained in the Extraction log binder for reference.
On a routine basis, the benchsheets are bound and archived.

10. DATA AND RECORDS MANAGEMENT

10.1. Responsibilities - It is the responsibility of the analyst to perform the analysis according to this SOP and to complete all documentation required for data review. Analysis and interpretation of the results are performed by personnel in the laboratory who have demonstrated the ability to generate acceptable results utilizing this SOP. Final review and sign-off of the data is performed by the department supervisor or designee.

10.2. Data Review - Data must be reviewed by the analyst and a peer (supervisor or qualified analyst) using a Data Quality Checklist before the results are validated and reported to the client. This checklist is found in SOP 8270CSIMPines.


11. QUALITY CONTROL AND QUALITY ASSURANCE

11.1. Extract a MB, LCS, MS and MSD (or LCSD if there is insufficient volume) for every batch of 20 or fewer samples. These QC samples are subjected to exactly the same analytical procedures as those used on actual samples. QC criteria and corrective action for these and the surrogates are specified in SOP 8270CSIMPines.


12. REFERENCES

“Test Methods For Evaluating Solid Waste Physical/Chemical Methods,” USEPA SW-846, December 1996.

[illegible]

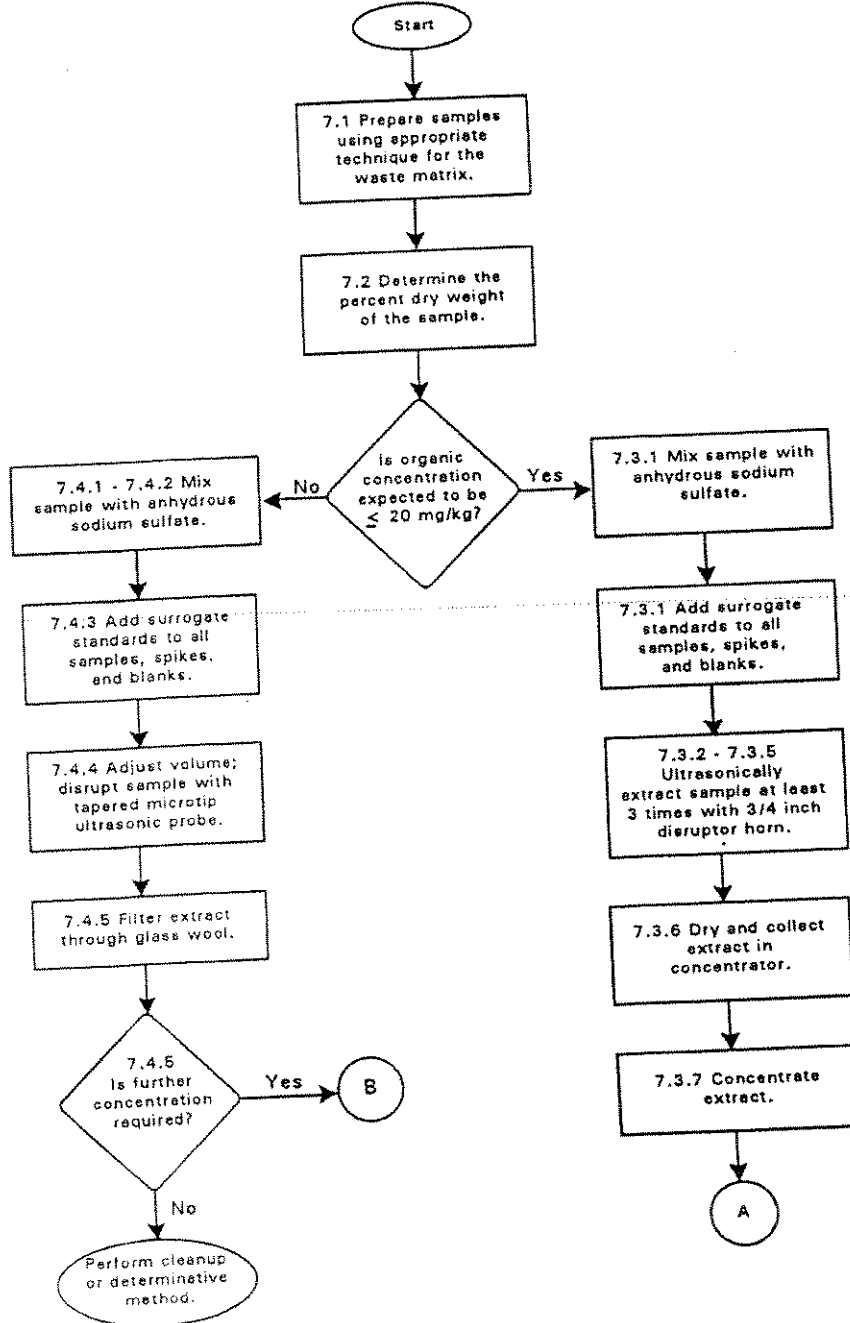METHOD 3550B
ULTRASONIC EXTRACTION

METHOD 3550B
continued

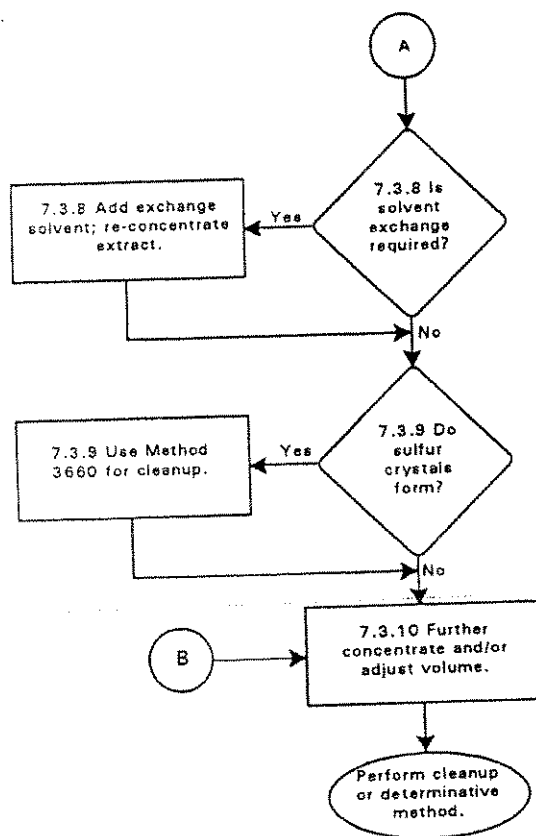

SOP Code: 8270CSIMPines
Revision: 0
Date: 12/21/04
Page 1 of 25

STANDARD OPERATING PROCEDURE

for

ANALYSIS OF PAHs BY GAS CHROMATOGRAPHY/MASS SPECTROMETRY USING SELECTIVE ION MONITORING (GC/MS/SIM) FOR INDIANA PINES SITE

SOP Code: 8270CSIMPines

Revision:0

December 21, 2004


Approved by:


Department Manager

12/21/04
Date


Quality Assurance

12/21/04
Date

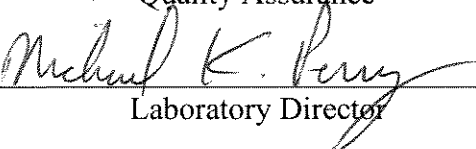
Laboratory Director

12/21/04
Date

© Columbia Analytical Services, Inc., 2004
1 Mustard Street, Suite 250
Rochester, NY 14609

Annual review of this SOP has been performed
and the SOP still reflects current practice.

Initials: _____ Date: _____
Initials: _____ Date: _____
Initials: _____ Date: _____

DOCUMENT CONTROL

NUMBER: _____

Initials: _____ Date: _____

Table of Contents	<u>Page</u>
1. Scope and Applicability.....	3
2. Summary of Method.....	3
3. Definitions.....	3
4. Health and Safety Warnings.....	5
5. Interferences.....	6
6. Personnel Qualifications.....	6
7. Equipment and Supplies.....	7
8. Procedure.....	10
8.1. Calibration and Standardization.....	10
8.2. Sample Collection.....	13
8.3. Sample Handling and Preservation.....	13
8.4. Sample Preparation	14
8.5. Sample Analysis.....	14
8.6. Troubleshooting.....	14
8.7. Data Acquisition, Calculations, and Data Reduction Requirements.....	16
8.8. Computer Hardware and Software.....	17
9. Data and Records Management.....	17
10. Quality Control and Quality Assurance.....	17
11. References.....	20

Attachments

Table 1: Target Compound List with MRLs.....	21
Table 2: Decafluorotriphenylphosphine Ion Abundance Criteria.....	22
Table 3: Characteristic Ions of Target Analytes.....	23
Data Review Checklist.....	24
Training Plan Form.....	25

1 SCOPE AND APPLICABILITY

This SOP uses EPA SW-846 Method 8270C (in the SIM mode) for the determination of low concentration levels of poly aromatic hydrocarbons(PAHs) in aqueous, soil, sludge, sediment, and various types of waste samples. Table 1 lists the compounds that are routinely determined by this procedure along with their method reporting limits (MRLs) in water and soil matrixes.

This procedure gives gas chromatographic/mass spectrometric (GC/MS) conditions for the detection of parts per billion (ppb) levels of PAHs.

This SOP was specifically modified for the Indiana Pines Site Project.

2 SUMMARY OF METHOD

Samples are extracted with solvent by EPA methodologies. A two μ L sample aliquot is injected into the gas chromatograph (GC) by splitless injection. The PAHs are separated by a fused silica capillary column, and the compounds are detected by a mass selective detector (MSD) using the SIM mode. The retention and the ratio of two characteristic ions of each analyte are used for identification. The response of either the primary ion or the secondary ion is used for quantitation.

3 DEFINITIONS

- 3.1 **Analysis Window** – Samples are analyzed in a set referred to as a “window”. The window begins with the injection of the DFTPP tune verification standard. Standards, required QC samples, and samples may be run for 12 hours in this window. A new window must be opened to continue analysis.
- 3.2 **Retention Time Window:** the time period established within which a target analyte is qualitatively determined to be present in the sample.
- 3.3 **Initial Calibration** - analysis of analytical standards for a series of different specified concentrations; used to define the linearity and dynamic range of the response of the detector to the target compounds.
- 3.4 **Matrix** - the predominant material, component, or substrate (e.g., surface water, drinking water, etc.) of which the sample to be analyzed is composed.

- 3.5 **QA/QC Samples:** Samples added to a sample preparation batch, or an analytical batch to provide quality assurance checks on the analysis.
- 3.5.1 **Laboratory Control Sample** - a matrix spiked sample with compounds representative of the target analytes. This is used to document laboratory performance.
- 3.5.2 **Matrix Spike/MSD** - an aliquot of sample spiked with a known concentration of target analyte(s). The spiking occurs prior to sample preparation and analysis. A matrix spike is used to document the bias of a method in a given sample matrix.
- 3.5.3 **Method Blank** - an analyte-free matrix to which all reagents are added in the sample volumes or proportions as used in sample processing. The method blank should be carried through the complete sample preparation and analytical procedure. The method blank is used to document contamination resulting from the preparation and analytical process.
- 3.6 **Percent Drift or Difference (%D)** - Used to compare two values, the percent difference indicates both the direction and the magnitude of the comparison, i.e., the percent difference may be either negative, positive, or zero. (In contrast, see relative percent difference).
- 3.7 **% Relative Standard Deviation (%RSD):** statistical measure of variation. Used in this method to measure the relative variation of initial calibration standards. Calculated by dividing the standard deviation of the individual calibration factors by the average calibration factor and multiplying by 100 to express as a percentage.
- 3.8 **Relative Percent Difference (RPD)** – The absolute value of the difference of two values divided by the average of the same two values. Used to compare the precision of the analysis. The result is always a positive number.
- 3.9 **Sample** - a portion of material to be analyzed that is contained in single or multiple containers and identified by a unique sample number.
- 3.10 **Surrogates** (Surrogate Standards) – an organic compound which is similar to the target analyte(s) in chemical composition and behavior in the analytical process. For semivolatiles, surrogate compounds are added to every blank, sample, matrix spike, matrix spike duplicate, LCS, matrix spike blank, and standard. These are used to evaluate analytical efficiency by measuring recovery. Surrogates are not expected to be detected in environmental media.

- 3.11 Internal Standards** – Internal standards are organic compounds which are similar to the analytes of interest but which are not found in the samples. The chosen internal standards are used to calibrate the instrument's response.
- 3.12 Organic Free Reagent Water** - ASTM Type II Deionized Water.
- 3.13 Batch** - Samples processed together as a unit, not to exceed 20 investigative samples. See ADM-BATCH for further discussion.
- 3.14 Independent (or Initial) Calibration Verification (ICV also known as Reference Check)** - A standard from a different source as the calibration standards used to verify the calibration curve.
- 3.15 Continuing Calibration Verification (CCV)** - A standard from the same source as the calibration standards used to verify the curve with each daily run and throughout the run at specified intervals.
- 3.16 Method Detection Limit (MDL):** a statistically derived value representing the lowest level of target analyte that may be measured by the instrument with 99% confidence that the value is greater than zero
- 3.17 Method Reporting Limit (MRL):** The minimum amount of a target analyte that can be measured and reported quantitatively. The MRL is equivalent to Practical Quantitation Level (PQL) and Estimated Quantitation Level (EQL). Typically, the MRL is calculated as five times the MDL (although this is a rule of thumb and not intended to be a strict policy of establishing the MRL for a compound).
- 3.18 Neat Stock Standard** - A purchased, single component assayed reference material having a stated purity used to prepare working calibration standards.

4 HEALTH AND SAFETY WARNINGS

- 4.1** Methylene chloride have been tentatively classified as a known or suspected human or mammalian carcinogen. The toxicity or carcinogenicity of the remaining chemicals used in this method has not been precisely defined. Each compound, mixture of compounds, and surrogates, as well as the samples, should be treated as a potential health hazard. Exposure to these chemicals should be reduced to the lowest level possible through the use of gloves (to minimize absorption through the skin) and hoods (to minimize inhalation). Material safety data sheets (MSDS) are available for all these reagents and solvents.

- 4.2 Samples may contain high concentrations of polynuclear aromatics, amines, phenols, and pesticides/PCBs; therefore, exposure to samples via contact must be minimized. Wear lab coat, gloves, and glasses when handling samples and reagents.
- 4.3 All applicable safety and compliance guidelines set forth by CAS, and by federal, state, and local regulations must be followed during performance of this procedure. All work must be stopped in the event of known or potential compromise to the health or safety of any CAS employee, and must be reported immediately to a laboratory supervisor.

5 INTERFERENCES

- 5.1 Sources of interference in this method can be grouped into three broad categories.
 - 5.1.1 Supply: Contaminated solvents, reagents, or GC carrier gas.
 - 5.1.2 Hardware: Contaminated glassware, analytical equipment; e.g., syringe, injection port, column surfaces, and/or detector surfaces.
 - 5.1.3 Matrix: Compounds extracted from the sample to which the detector will respond.
- 5.2 Interference from phthalate esters can best be minimized by avoiding contact with any plastic materials and checking all solvents and reagents for phthalate contamination. Exhaustive cleanup of reagents and glassware may be required to eliminate background phthalate contamination.

6 PERSONNEL QUALIFICATIONS

At a minimum, personnel must have attained at least a 4-year degree (or 2-yr degree plus one year experience) in a science-related field and have successfully completed an Initial Demonstration of Capability and the Training Plan Form (attached). Training and Demonstration of Capability are in accordance with NELAC 2002 standard. This method is restricted to use by or under the supervision of analysts experienced in the use of gas chromatograph/mass spectrometers and skilled in the interpretation of mass spectra.

7 EQUIPMENT AND SUPPLIES

7.1 Gas Chromatograph/Mass Spectrometer System (GC/MS)

An analytical system complete with a temperature programmable gas chromatograph suitable for splitless injection and all required accessories, including autosampler, analytical column, and carry gas. The capillary column should be directly coupled to the source of the mass spectrometer.

Column - RTx-5Sil MS with Integra guard (or equivalent)
30 m x 0.25 mm ID
0.5 μ m film thickness
silicone-coated fused capillary

7.2 Mass Spectrometer (MS)

A MS capable of scanning from 35 to 500 amu every second or less, using 70 volts (nominal) electron energy in the electron impact ionization mode. The mass spectrometer must be capable of producing a mass spectrum for Decafluorotriphenylphosphine (DFTPP) which meets all of the criteria in Table 2 when 50 ng of DFTPP is injected onto the GC/MS system.

7.3 GC/MS Interface

Any GC-to-MS interface that gives acceptable calibration points at 50 ng per injection for each compound of interest and achieves acceptable tuning performance criteria may be used.

7.4 In-line Gas Purifier (optional):

Agilent Technologies part no. 5182-9705 (or equivalent) to remove water, oxygen, and hydrocarbons.

7.5 Microsyringes - 10, 50, 100, 250, 500, and 1000- μ L

7.6 Balance

7.6.1 Top-loading Scale: Capable of weighing to 0.1 g.

7.6.2 Analytical Balance: Capable of weighing to 0.1 mg.

7.7 Spatula - Stainless steel.

7.8 Helium: High purity grade (99.99%) used as carrier gas for the gas chromatograph.

7.9 Analyte-free Reagent Water: Water for which no target analytes are observed at or about the MRL, or the MDL, depending upon the project.

7.10 Solvents: acetone, Methylene chloride, methanol, and other appropriate standards that is pesticide grade or better. Purchased commercially. Store at room temperature. Expires upon manufacturer's indications or 3 years from receipt when no indication is given.

7.11 Stock Standards for Initial Calibration and Continuing Calibration Verification Standards Store all standards at -10-20°C. All stock standards expire upon the manufacturer's indications or one year if no indication is given.

7.11.1 SemiVolatile Internal Standards: 4000 µg/mL in Methylene Chloride, NSI.

7.11.2 8270SIM Stock Standards

7.11.2.1 Ready Stock (CLP surrogates): 200ppm in Methylene Chloride, NSI.

7.11.2.2 1-Methylnaphthalene: 2000 µg/mL in Methylene Chloride, ChemService.

7.11.3 Tuning Mix: 2500 µg/mL in Methylene Chloride, Restek.

7.12 Intermediate Calibration Standards

Intermediate Calibration Standards have a 6 months expiration date and must be stored between -10 and -20 °C.

7.12.1 Semi-volatile Organic Compounds (50 µg/mL in Methylene Chloride): Take 2.50 mL of 200ppm ReadyStock (CLP surrogates) and 0.25 mL of 2000µg/mL 1-Methylnaphthalene standard and dilute to 10 mL with Methylene Chloride.

7.12.2 Semi-volatile Organic Compounds (10 µg/mL in Methylene Chloride): Take 2.0 mL of the 50 ug/mL in Methylene Chloride solution and dilute to 10 mL with Methylene Chloride.

7.13 Working Calibration Standard

The following calibration standards are prepared from the 10 µg/mL intermediate standard to a final volume of 1 mL with Methylene Chloride. These standards expire in 6 months and must be stored between -10 and -20 °C.

Final Concentration (µg/mL)	Aliquot of 10 µg/mL (µL)	Final Volume (mL)
0.1	10	1
0.2	20	1
0.5	50	1
1.0	100	1
2.0	200	1
5.0	500	1

7.14 Decafluorotriphenylphosphine (DFTPP) Tune Standard

Decafluorotriphenylphosphine Intermediate Standard (25 ppm in Methylene Chloride) - Dilute a 100 µL aliquot of 2500 ppm solution to a 10-mL volumetric flask with Methylene Chloride. This dilution expires in 6 months and must be stored between -10 and -20 °C.

7.15 Pentafluorotributylamine (PFTBA) - Calibration gas

7.16 Laboratory Control Sample (LCS) and Matrix Spiking Standards –

7.16.1 8270 LCS mix 1 (100 ppm), Supelco - Store at -10-20°C. Standard expires upon the manufacturer's indications or one year if no indication is given.

7.16.2 LCS Working standard solution (1.0 µg/mL) – Take 0.5 mL of 100 ppm 8270LCS mix 1 and dilute to 50 mLs with Methanol. This dilution expires in 6 months and is stored between -10 and -20 °C.

7.17 ICV standards –

7.17.1 8270/625 Mid-stock (200ppm) in Methylene Chloride, Supelco. Store at -10-20°C. Standard expires upon the manufacturer's indications or 6 months if no indication is given.

7.17.2 Intermediate Solution for ICV (10 µg/mL) Take 0.05 mL of 200ppm 8270/625 Midstock and dilute to 1.0 mL with Methylene Chloride. Expires in 6 months. Store between -10 and -20 °C.

7.17.3 Working Standard Solution for ICV (1.0 µg/mL) Take 0.10 mL of the intermediate ICV solution and dilute to 1.0 mL with Methylene Chloride. Expires in 6 months. Store between -10 and -20 °C.

8 PROCEDURE

8.1 Calibration and Standardization

8.1.1 MSD Tuning: Tune the MSD to meet the criteria in Table 2 for a 50 ng injection of DFTPP. Acquire the mass spectrum of DFTPP as follows: Average three scans (the peak apex scan and the scans immediately preceding and following the apex). Alternatively, the apex scan, or the scan preceding the apex, or the scan following the apex can be used. When background subtraction is required, subtract using a single scan no more than 20 scans prior to the elution of DFTPP. Do not use any part of the DFTPP peak to background subtract. Use identical MSD instrument conditions for all subsequent standards, samples, spiked samples, and QC samples associated with a DFTPP analysis. Use the analysis time for the DFTPP to define the start of the 12-hour analysis window.

8.1.2 Initial Calibration (ICAL)

Prior to conducting any sample analyses, prepare a multi-point calibration. See Section 7 for preparation of the Working Calibration Standards. Select a primary and a secondary ion from each analyte for identification. Tabulate the area response of the primary quantitation ion (see Table 3) versus the concentration for each internal standard, analyte, and surrogate. Calculate the response factor (RF) for each analyte and surrogate relative to the associated internal standard using the following formula.

$$RF = (A_X)(C_{IS}) / (A_{IS})(C_X)$$

Where

- A_X = Peak area of the analyte or surrogate's characteristic ion;
 A_{IS} = Peak area of the associated internal standard's characteristic ion;
 C_X = Concentration of the analyte or surrogate in the calibration sample; and
 C_{IS} = Concentration of the associated internal standard.

8.1.2.1 Calculate the average response factor, RF_{AVE} , for each analyte and surrogate from RFs of each of the calibration levels.

$$RF_{AVE} = \sum RF_i / n$$

8.1.2.2 Calculate the standard deviation (SD) and the percent relative standard deviation (% RSD) for each analyte.

$$\% \text{ RSD} = (SD)(100) / RF_{AVE}$$

8.1.2.3 If the % RSD of any analyte or surrogate is $\leq 15\%$, assume linearity over the calibration range. Use the RF_{AVE} for the analyte or surrogate to quantitate sample analytes. Alternatively, calibrate the analytes by linear regression or quadratic regression. If the coefficient correlation of linear or quadratic regression is < 0.990 , then corrective action must be taken to eliminate the problem prior to re-attempting calibration. If the calibration criteria are not met, check standards with bad injection and re-analyzed standard. If bad injection is not evident, perform maintenance and calibration. Demonstrate that the calibration is in-control before proceeding with the analysis.

8.1.2.4 Supervisory review and approval of the ICAL is required.

8.1.2.5 Create an ICAL File. Place the following ICAL documents in the file.

- ICAL Checklist filled out and approved
- Sequence report
- DFTPP Tune analysis Report
- Blank analysis Quantitation Report
- Calibration Status Report - Initial
- Response Factor Report
- Data Analysis Parameters Report
- Copy of the calibration curve for any compound that uses a curve (instead on the average RF)
- Quantitation Report for each calibration standard (including manual integration documentation)
- ICV Quantitation Report and Evaluate Continuing Calibration Report
- Calibration Status Report - Final

8.1.3 Initial Calibration Verification: Following the initial calibration, analyze an ICV standard to verify that there are no systematic errors in any of the analyte calibrations. The ICV standard is at a concentration near the midpoint of the calibration curve, but it is made from a source other than the source used for the preparation of the calibration standards. The ICV standard must contain all the analytes that are in the calibration standards. The concentration of each analyte is calculated using the procedure for quantitation determined during the ICAL. The percent recovery for each analyte must be within 70 - 130%. If an analyte's percent recovery is outside the limits, take corrective action to determine why. Some compounds may exceed this criteria and the initial calibration may still be valid. Use professional judgment when evaluating reactive compounds or those exhibiting poor chromatographic behavior.

8.1.4 Continuing Calibration Verification:

8.1.4.1 Check the MSD's tune by injecting 50 ng of DFTPP as described above at the start of a 12-hour analysis window. If the criteria found in Table 2 are met, then continue the check the initial calibration curve by analyzing the CCV. If the first run of the DFTPP fails, retry. If the second run also fails, inspect the system for potential maintenance needs. Take corrective action before attempting to retune.

8.1.4.2 After the tuning criteria have been verified, verify the initial calibration by analyzing a midrange continuing calibration verification (CCV) standard. The 1 ppm level standard is recommended. Compare the results to the ICAL.

- The analytes using RF_{AVE} for quantitation must have a percent difference $\leq 20\%$.
- The analytes using an equation for quantitation must have a percent recovery of 80 to 120%.

8.1.4.3 If the tune criteria and the continuing calibration criteria are met, check the retention times of all compounds, surrogates, and internal standards against the initial calibration. If the retention time for any internal standard changes by more than 30 seconds from the retention time from the mid-point standard of the most recent initial calibration, inspect the system for malfunctions and make corrections as required. If the area for any of the internal standards changes by a factor of 2 (-50% to +100%) from the response of the mid-point standard of the most recent initial calibration, make corrections to the system.

8.1.4.4 If the CCV standard analysis meets the continuing calibration verification criteria, begin sample analysis. If the CCV criteria are not met, evaluate if the problem is due to internal standard response, to surrogate response, or to the response of one of the target analytes. If the problem is due to one or more of the internal standards and/or one or more of the surrogates, re-prepare the internal standard/surrogate solution and re-analyze. If the problem persists, perform maintenance. If the results for one or more of the target analytes do not meet criteria for continuing calibration verification, perform maintenance. Demonstration of in-control is required before proceeding with the analysis.

8.2 Sample Collection

- 8.2.1** Containers used to collect samples are to be purchased precleaned and certified. The sample containers should be 1-liter amber glass with Teflon lined screw-top cap for water and 16 oz. glass jar, or metal sleeve for soil sample.
- 8.2.2** Sample containers should be filled with care so as to prevent any portion of the sample coming in contact with the sampler's gloves, thus potentially contaminating the sample. Samples should not be collected or stored in the presence of exhaust fumes. If the sample comes in contact with the sampler (e.g., if an automatic sampler is used), run reagent water through the sampler and use the rinsate as a field blank.

8.3 Sample Handling and Preservation

- 8.3.1** Water and soil samples are iced or refrigerated at 0-6 °C in the dark from time of collection until extraction.
- 8.3.2** Extract water samples within 7 days of collection and analyze the extracts within 40 days of preparation of the extract.
- 8.3.3** Extract Soil and sludge samples within 14 days of collection and analyze the extract within 40 days of preparation of the extract.

8.4 Sample Preparation

Method Selection: Prepare samples by one of the following methods prior to GC/MS analysis.

<u>Matrix</u>	<u>Methods</u>
Water	3510
Soil/sediment	3550

8.5 Sample Analysis

Identify analytes by retention time and ion ratios, using two ions per analyte. Quantitate that compound based on the integrated abundance from the primary or the secondary characteristic ion (see Table 3). Use the internal standard nearest the retention time of that given analyte. Quantify all compounds based on the initial, multi-point calibration using either RF_{AVE} or the calibration equation that was determined during the ICAL.

8.6 Troubleshooting

8.6.1 Carrier Gas Purifier - If in-line purifiers or scrubbers are in place, these purifiers should be changed as recommended by the supplier.

8.6.2 Gas Chromatograph

8.6.2.1 Chromatographic performance can often be improved by clipping off a small portion of the front of the capillary column. The cut needs to be straight and clean (uniform, without fragmentation) by using the proper column cutting tool.

8.6.2.2 Over time, the capillary column will exhibit poorer overall performance as peak resolution deteriorates due to analysis of contaminated samples. The length of time for this to occur will depend on the samples analyzed. When a noticeable decrease in column performance is evident and other maintenance options do not result in improvement, the replace the column. This is especially evident when difficulties are experienced in conjunction with calibration.

8.6.3 Mass Selective Detector (MSD)

- 8.6.3.1 Tune the MSD as needed to achieve acceptable and consistent performance.
 - 8.6.3.2 Check the foreline pump oil level weekly. Add pump fluid until the oil level in the window is near, but not above, the upper line.
 - 8.6.3.3 Replace the foreline pump oil every six months. This may be done by the chemist or by an authorized service engineer.
 - 8.6.3.4 Check the level of PFTBA in the calibration vial every six months. Refill the vial if necessary.
 - 8.6.3.5 Check the diffusion pump fluid level annually. Replace the fluid if the level is low or dark or cloudy. This may be done by the chemist or by an instrument service person.
 - 8.6.3.6 Clean the ion source as needed, depending upon the performance of the MSD. This may be done by the chemist or by an authorized service engineer.
 - 8.6.3.7 Other maintenance is performed as recommended by Hewlett-Packard.
- 8.6.4 Maintenance log – Document all Preventive maintenance, as well as instrument repair, in the appropriate instrument maintenance log. Most routine maintenance and troubleshooting are performed by CAS staff. Other maintenance or repairs may, or may not require factory service, depending upon the nature of the task. Any maintenance performed by outside services must also be documented – either through notes in the log or through documents provided by the service. The log entries will include the date maintenance was performed, symptoms of the problem, serial numbers of major equipment upgrades or replacements. The datafile name of the first acceptable run after maintenance is to be documented in the maintenance log.

8.7 Data Acquisition, Calculations, and Data Reduction Requirements

8.7.1 Introduction: The GC/MS data processing software uses the Hewlett-Packard RTE Integrator to generate the raw data used to calculate each analyte's RF_{AVE} values, the sample amounts, and the spike values. The software does three passes through each data file. The first two identify and integrate each internal standard and surrogate. The third pass uses the time-drift information from the first two passes to search for all the calibrated analytes in the proper retention times and with the proper characteristic quantitation ions.

8.7.2 Target Analyte Concentrations in Water Samples: The concentration of an analyte, C_X , in a water sample is calculated using the RF_{AVE} as follows.

$$C_X \text{ in } \mu\text{g/L} = (A_X)(C_{IS})(DF) / (A_{IS})(RF_{AVE})$$

Where

A_X = Peak area of the analyte's (or surrogate's) characteristic ion;
 A_{IS} = Peak area of the associated internal standard's characteristic ion;
 C_{IS} = Concentration of the associated internal standard;
 RF_{AVE} = Average response factor for analyte from the ICAL; and
 DF = Dilution factor, if the sample was diluted prior to analysis. If no dilution was made, $DF = 1$.

8.7.3 Target Analyte Concentration in Soil Samples The concentration of an analyte, C_X , in a soil sample is calculated using the RF_{AVE} as follows.

$$C_X \text{ in } \mu\text{g/Kg} = (A_X)(C_{IS})(DF) / (A_{IS})(RF_{AVE})(W_S)(D)$$

Where

A_X = Peak area of the analyte's (or surrogate's) characteristic ion;
 A_{IS} = Peak area of the associated internal standard's characteristic ion;
 C_{IS} = Concentration of the associated internal standard;
 RF_{AVE} = Average response factor for analyte from the ICAL;
 W_S = Weight of sample extracted in gram
 D = Percent dry weight of sample $\div 100$ (if the result is to be reported on a dry weight basis; and
 DF = Dilution factor

8.8 Computer Hardware and Software

IBM-compatible PC with HP Chemstation software including EnviroQuant with Extracted Ion Current Profile (EICP), or equivalent.

9 DATA AND RECORDS MANAGEMENT

9.1 Responsibilities - It is the responsibility of the analyst to perform the analysis according to the instructions in this SOP and to complete all documentation required for data review. Analysis and interpretation of the results are only to be performed by personnel in the laboratory who have demonstrated the ability to generate acceptable results utilizing this SOP. This demonstration is in accordance with the training program of the laboratory. Final review and sign-off of the data is performed by the department supervisor/manager or designee.

9.2 Data Review – Data will be reviewed by the instrument analyst and a qualified peer using a Data Review Checklist (attached) and validated by a supervisor.

10 QUALITY ASSURANCE/QUALITY CONTROL REQUIREMENTS

10.1 Tune and Calibration: The acceptance criteria for MSD tuning verification, initial calibration, and continuing calibration verification are detailed in the procedure. If tune does not meet criteria, perform target-tune and then re-analyze DFTPP. If tune check with DFTPP is still out-of-control, perform maintenance. Demonstrate that the system is in-control before proceeding with the analysis of CCV standard.

10.2 Internal Standards: The area counts of the internal standards of the initial daily CCV must be within ± 2 times (i.e., -50 % to +100 %) of the responses of the internal standards of the mid-point ICAL standard. Compare the internal standard responses for all analyses in the 12-hour window with the initial daily CCV. Internal standard area counts of subsequent samples and standards must be within ± 2 times (i.e., -50% to +100%) of the responses of the internal standards in the CCV analyzed at the start of the 12-hour window. Internal standards must have RT ± 0.5 min from the ICAL for the CCV and ± 0.5 from the CCV for the samples. CCV internal standard area counts must be -50% to 100% of the initial calibration mid-point standard area counts, otherwise, a new curve is required. Sample internal standard area counts must be -50% to 100% of the CCV area counts. Re-analyze any internal standard outlier unless a matrix interference can be clearly demonstrated by the first analysis. Flag as estimated any associated results from each outlying Internal Standard.

10.3 Surrogate Recoveries

Surrogate compound	Water Recovery Limits (%)	Soil Recovery Limits (%)
2-Fluorobiphenyl	27-114	23-120
Nitrobenzene -d5	22-124	18-125
Terphenyl- d14	23-139	19-145

If a surrogate fails acceptance, the sample must be evaluated for matrix interferences and “historical results”. Reanalyze the sample to confirm the interference. If confirmed, re-extract the sample unless confirmed by MS/MSD or there is insufficient sample volume. If needed contact client and flag the data in the report. If surrogates are diluted more than 10 times, report as “D”, diluted below calibration. For package reports, include initial and confirmation analysis results.

- 10.4 Method Blank The method blank should not have target analytes detected at or above the method reporting limit. If there is MB contamination, determine whether the contamination is from the instrument or due to contamination from the extraction. If the contamination is due to extraction, re-extract the batch with clean glassware or flag the data appropriately.

10.5 Laboratory Control Sample Recoveries

Compound	Water Recovery Limits (%)	Soil Recovery Limits (%)
Acenaphthene	49-116	42-112
Acenaphthylene	45-122	44-114
Anthracene	54-120	49-113
Benzo(a)anthracene	61-116	47-116
Benzo(a)pyrene	60-118	41-122
Benzo(b)fluoranthene	60-116	48-117
Benzo(g,h,i)perylene	50-125	34-126
Benzo(k)fluoranthene	54-120	41-123
Chrysene	60-117	45-117
Dibenzo(a,h)anthracene	31-139	29-129
Fluoranthene	56-121	36-122
Fluorene	48-122	40-113
Indeno(1,2,3-cd)pyrene	50-125	40-122
Naphthalene	39-109	44-101
Phenanthrene	47-128	51-110
Pyrene	60-113	35-128

If the LCS fails acceptance limits for any target compounds, the analyst must evaluate the system and calibration. If no problems are found, then the batch QC must be evaluated to determine what corrective action must be taken. This may involve the Project manager or Department Supervisor. Corrective action will depend on specific project, client, or state agency. The batch shall be re-extracted or the data may be flagged in the final report.

10.6 Matrix Spike Recoveries

Compound	Water Recovery Limits (%)	Soil Recovery Limits (%)
Acenaphthene	49-116	42-112
Acenaphthylene	45-122	44-114
Anthracene	54-120	49-113
Benzo(a)anthracene	61-116	47-116
Benzo(a)pyrene	60-118	41-122
Benzo(b)fluoranthene	60-116	48-117
Benzo(g,h,i)perylene	50-125	34-126
Benzo(k)fluoranthene	54-120	41-123
Chrysene	60-117	45-117
Dibenzo(a,h)anthracene	31-139	29-129
Fluoranthene	56-121	36-122
Fluorene	48-122	40-113
Indeno(1,2,3-cd)pyrene	50-125	40-122
Naphthalene	39-109	44-101
Phenanthrene	47-128	51-110
Pyrene	60-113	35-128

If the matrix spike fails acceptance, the sample must be evaluated for matrix interferences requirements. Evaluate the recovery of the duplicate MS and/or batch LCS. If the LCS is acceptable, continue with the analysis and assume matrix interferences.

11 REFERENCES

- 11.1 *Semi-volatile Organic Compounds by Gas Chromatography/Mass Spectrometry (GC/MS): Capillary Column Technique*, Method 8270C, Revision 3, December 1996 in *Test Methods for Evaluating Solid Waste, Physical/Chemical Methods*, U. S. EPA, SW-846, Final Update III
- 11.2 *Determinative Chromatographic Separations*, Method 8000B, Revision 2, December 1996 in *Test Methods for Evaluating Solid Waste, Physical/Chemical Methods*, U. S. EPA, SW-846, Final Update III

Table 1
Target Compound List

Compound	MRL for Water (µg/L)	MRL for Soil (µg/Kg)
Acenaphthene	0.2	6.6
Acenaphthylene	0.2	6.6
Anthracene	0.2	6.6
Benzo (a) Anthracene	0.1	3.3
Benzo (a) Pyrene	0.2	6.6
Benzo (b) Fluoranthene	0.2	6.6
Benzo (g,h,i) Perylene	0.2	6.6
Benzo (k) Fluoranthene	0.2	6.6
Chryene	0.2	6.6
Dibenzo (a,h) Anthracene	0.2	6.6
Fluoranthene	0.2	6.6
Fluorene	0.2	6.6
Indeno (1,2,3-cd) Pyrene	0.2	6.6
Naphthalene	0.2	6.6
Phenanthrene	0.2	6.6
Pyrene	0.2	6.6

Table 2
Decafluorotriphenylphosphine Ion Abundance Criteria

Ion Mass (m/z)	Relative Abundance Criteria
51	30 to 60% of m/z 198
68	<2 % of m/z 69
70	<2 % of m/z 69
127	40 to 60 % of m/z 198
197	<1 % of mass 198
198	Base peak, 100 % relative abundance
199	5 to 9% of m/z 198
275	10 to 30 % of m/z 198
365	>1 % of mass 198
441	Present but less than m/z 443
442	>40 % of m/z 198
443	17 to 23 % of m/z 442

Table 3
Characteristic Ions of Target Analytes

Compound	Primary Ion	Secondary Ion(s)
<i>Internal Standards</i>		
Naphthalene-d8	136	68, 108
Acenaphthene-d10	164	162, 160
Phenanthrene-d10	188	94, 80
Chrysene-d12	240	120, 236
Perylene-d12	264	260, 265
<i>Surrogates</i>		
Nitrobenzene-d5	80	82, 54
2-Fluorobiphenyl	172	171
Terphenyl-d14	244	122, 212
<i>Analytes</i>		
Naphthalene	128	127
Acenaphthylene	152	151, 153
Acenaphthene	154	153, 152
Fluorene	166	165, 167
Phenanthrene	178	179, 176
Anthracene	178	176, 179
Fluoranthene	202	101, 203
Pyrene	202	200, 203
Benz(a)anthracene	228	229, 226
Chrysene	228	226, 229
Benzo(b)fluoranthene	252	253, 125
Benzo(k)fluoranthene	252	253, 125
Benzo(a)pyrene	252	253, 125
Indeno(1,2,3-cd)pyrene	276	138, 278
Dibenz(a,h)anthracene	278	276, 138
Benzo(g,h,i)perylene	276	138, 277



GC/MS SEMI-VOLATILES DATA QUALITY CHECKLIST

Method: _____

Instrument: _____
Analysis Date: _____

Yes	No	NA			Yes	No	NA
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	1.	Extraction Benchsheet complete?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	2.	Holding Times met method requirements?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	3.	Instrument Run Log complete?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	4.	DFTPP Tune met method requirements?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	5.	CCVs acceptable (win, RT, %D)?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	6.	Quantitation reports present for all samples, blk, lcs, ms, msd?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	7.	Has the analyst initialed and dated each quantitation report?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	8.	Are spectral match details for all sample hits present? (packages only)	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	9.	Graphically verifies spectral matches for all sample hits? (routine report)	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	10.	Has current ICAL been used to quantitate all sample results?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	11.	Are all analyses within 12-hr window?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	12.	Are all Internal standard retention times within 30 sec. of opening CCV?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	13.	Are all IS responses within -50% to +100% of opening CCV?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	14.	Surrogate recoveries are within QC limits for samples and QC?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	15.	Method Blank results < PQL?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	16.	LCS recoveries within QC limits?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	17.	MS/MSD recoveries within QC limits?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>		• RPDs between MS/MSD within QC limits?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	18.	All sample concentrations within Linear Range?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	19.	Dilution factors verified and calculated correctly?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	20.	All peak integrations acceptable?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	21.	Are all manual integrations flagged, initialed and dated?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	22.	Are internal COCs included in package (if applicable)?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

Analyst: _____
Date: _____

Peer Review: _____
Date: _____

COMMENTS:

****Comments must be provided for any items noted above as "No"**

Training Plan for Analysis of SVOCs by GC/MS

Procedure: _____

SOP: _____ Revision: _____ Date: _____

Trainee: _____

1. Read SOP Trainer: _____ Trainee: _____ Date: _____

2. Read appropriate EPA Method Trainer: _____ Trainee: _____ Date: _____

3. Demonstrated understanding of the scientific basis of the analysis
Gas chromatography
Mass spectrometry Trainer: _____ Trainee: _____ Date: _____

4. Demonstrated familiarity with related SOPs

ADM-BATCHSEQ

ADM-SIGFIG

ADM-DREV

ADM-DATAENTRY

ADM-NCAR

ADM-TRANDOC

ADM-INT

ADM-MDL

Trainer: _____ Trainee: _____ Date: _____

5. Observe performance of SOP
- sample preparation (soil, water, other matrices) and sample loading
 - analytical sequence setup
 - DFTPP tuning evaluation
 - initial calibration and continuing calibration verification
 - sample analysis
 - EnviroQuant use
 - data reduction and reporting

Trainer: _____ Trainee: _____ Date: _____

6. I have read, understood and agree to perform the most recent version of the SOP:

Signature: _____ Date: _____

7. Perform SOP with supervision
- including all items listed in 5.

Trainer: _____ Trainee: _____ Date: _____

8. Independent performance of the SOP
- including all items listed in 5.
 - IDC (4 mid-range standards performed before client samples are analyzed)
 - attach IDC certificate, raw data, and summary spreadsheet.

Trainer: _____ Trainee: _____ Date: _____

SOP Code: HRMS-8290
Date: 12/07/04
Authors: XL, JF, AV, RD, KV


STANDARD OPERATING PROCEDURE


ANALYSIS OF POLYCHLORINATED DIBENZO-*p*-DIOXINS AND POLYCHLORINATED DIBENZOFURANS BY HIGH-RESOLUTION GAS CHROMATOGRAPHY/HIGH-RESOLUTION MASS SPECTROMETRY (HRGC/HRMS)

SOP Code: HRMS-8290

Revision: 5.2klv

12/07/04

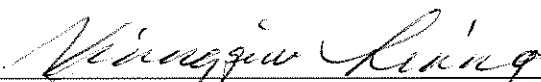Approved by: 
Xiangqiu Liang, Laboratory Director

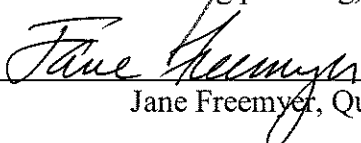
Jane Freemyer, Quality Assurance Manager

12/8/04
Date

12/08/04
Date

© Columbia Analytical Services, Inc., 2003
10655 Richmond Ave., Suite 130A
Houston, TX 77042

TABLE OF CONTENTS

	Page
1.0 Application.....	3
2.0 Method Summary.....	4
3.0 Definitions.....	6
4.0 Health and Safety Warnings.....	7
5.0 Cautions.....	7
6.0 Interferences.....	8
7.0 Personnel Qualifications	8
8.0 Equipment and Supplies	10
9.0 Instrument Calibration and Standarization.....	16
10.0 Sample Collection.....	22
11.0 Sample Handling and Preservation.....	22
12.0 Sample Preparation and Analysis	22
13.0 Data Analysis and Calculations	33
14.0 Data and Records Management	38
15.0 Quality Assurance Requirements.....	40
16.0 References.....	48

Standard Operating Procedure

ANALYSIS OF POLYCHLORINATED DIBENZO-*p*-DIOXINS AND POLYCHLORINATED DIBENZOFURANS BY HIGH-RESOLUTION GAS CHROMATOGRAPHY/HIGH-RESOLUTION MASS SPECTROMETRY (HRGC/HRMS)

1.-1 APPLICATION

- 1.1 This method provides procedures for the detection and quantitative measurement of polychlorinated dibenzo-*p*-dioxins (PCDDs), and polychlorinated dibenzofurans (PCDFs) in a variety of environmental matrices and at part-per-trillion to part-per-quadrillion concentrations. The following compounds can be determined by this method:

<u>Analyte</u>	<u>CAS Registry No</u>
2,3,7,8-Tetrachlorodibenzo- <i>p</i> -dioxin (TCDD)	1746-01-6
1,2,3,7,8-Pentachlorodibenzo- <i>p</i> -dioxin (PeCDD)	40321-76-4
1,2,3,4,7,8-Hexachlorodibenzo- <i>p</i> -dioxin (HxCDD)	39227-28-6
1,2,3,6,7,8-Hexachlorodibenzo- <i>p</i> -dioxin (HxCDD)	57653-85-7
1,2,3,7,8,9-Hexachlorodibenzo- <i>p</i> -dioxin (HxCDD)	19408-74-3
1,2,3,4,6,7,8-Heptachlorodibenzo- <i>p</i> -dioxin (HpCDD)	35822-46-9
1,2,3,4,5,6,7,8-Octachlorodibenzo- <i>p</i> -dioxin (OCDD)	3268-87-9
2,3,7,8-Tetrachlorodibenzofuran (TCDF)	51207-31-9
1,2,3,7,8-Pentachlorodibenzofuran (PeCDF)	57117-41-6
2,3,4,7,8-Pentachlorodibenzofuran (PeCDF)	57117-31-4
1,2,3,4,7,8-Hexachlorodibenzofuran (HxCDF)	70648-26-9
1,2,3,6,7,8-Hexachlorodibenzofuran (HxCDF)	57117-44-9
1,2,3,7,8,9-Hexachlorodibenzofuran (HxCDF)	72918-21-9
2,3,4,6,7,8-Hexachlorodibenzofuran (HxCDF)	60851-34-5
1,2,3,4,6,7,8-Heptachlorodibenzofuran (HpCDF)	67562-39-4
1,2,3,4,7,8,9-Heptachlorodibenzofuran (HpCDF)	55673-89-7
1,2,3,4,5,6,7,8-Octachlorodibenzofuran (OCDF)	39001-02-0
Total Tetrachlorodibenzo- <i>p</i> -dioxin (TCDD)	41902-57-5
Total Pentachlorodibenzo- <i>p</i> -dioxin (PeCDD)	36088-22-9
Total Hexachlorodibenzo- <i>p</i> -dioxin (HxCDD)	34465-46-8
Total Heptachlorodibenzo- <i>p</i> -dioxin (HpCDD)	37871-00-4
Total Tetrachlorodibenzofuran (TCDF)	55722-27-5
Total Pentachlorodibenzofuran (PeCDF)	30402-15-4
Total Hexachlorodibenzofuran (HxCDF)	55684-94-1
Total Heptachlorodibenzofuran (HpCDF)	38998-75-3

- 1.2 The analytical method calls for the use of high-resolution gas chromatography and high-resolution mass spectrometry (HRGC/HRMS) on purified sample extracts. Table 1 lists the various sample types covered by this analytical protocol, the 2,3,7,8-TCDD-based method calibration limits (MCLs), and other pertinent information. Samples containing concentrations of specific congeneric analytes (PCDDs and PCDFs – see Table 4) considered within the scope of this method that are greater than ten times the upper MCLs must be analyzed by a protocol designed for such concentration levels, e.g., Method 8280. An optional method for reporting the analytical results using a 2,3,7,8 TCDD toxicity equivalency factor (TEF) is described.

The sensitivity of this method is dependent upon the level of interferences within a given matrix (see Table 1.) The calibration range of the method for a 1L water sample is 10 to 2000ppq for TCDD/TCDF and PeCDD/PeCDF, and 1.0 to 200ppt for a 10g soil, sediment, fly ash, or tissue sample for the same analytes. Analysis of a one tenth aliquot of the sample permits measurement of concentrations up to 10 times the upper MCL. The actual limits of detection and quantitation will differ from the lower MCL, depending on the complexity of the matrix.

- 1.3 This method is designed for use by analysts who are experienced with residue analysis and skilled in HRGC/HRMS.
- 1.4 Because of the extreme toxicity of many of these compounds, the analyst must take the necessary precautions to prevent exposure to materials known or believed to contain PCDDs or PCDFs. It is the responsibility of the laboratory personnel to ensure that safe handling procedures are employed.

2.0 METHOD SUMMARY

- 2.1 This procedure uses matrix specific extraction, analyte specific cleanup, and HRGC/HRMS analysis techniques.
- 2.2 If interferences are encountered, the method provides selected cleanup procedures to aid the analyst in their elimination. A simplified analysis flow chart is presented at the end of this method.
- 2.3 A specified amount (see Table 1) of soil, sediment, fly ash, water, sludge (including paper pulp), still bottom, fuel oil, chemical reactor residue, fish tissue, or human adipose tissue is spiked with a solution containing specified amounts of each of the nine isotopically labeled PCDDs/PCDFs. The sample is then extracted according to a matrix specific extraction procedure. Aqueous samples that are judged to contain 1 percent or more solids, and solid samples that show an aqueous phase, are filtered. The solid phase (including the filter) and the aqueous phase are extracted separately, and the extracts combined before extract cleanup.

The extraction procedures are:

- 2.3.1 Toluene: Soxhlet extractions for soil, sediment, fly ash, and paper pulp samples;
 - 2.3.2 Methylene chloride: Jar/Separatory Funnel extraction for water samples;
 - 2.3.3 Toluene: Dean-Stark extraction for fuel oil and chemical samples
 - 2.3.4 Toluene extraction for still bottom samples and sludge samples
 - 2.3.5 Hexane/methylene chloride: Soxhlet extraction for fish tissue samples
 - 2.3.6 Methylene chloride extraction for human adipose tissue samples.
 - 2.3.7 As an option, all solid samples (wet or dry) may be extracted with toluene using a Soxhlet/Dean Stark extraction system. The decision for the selection of an extraction procedure for chemical reactor residue samples is based on the appearance (consistency, viscosity) of the samples.
- 2.4 The extracts are submitted to a sulfuric acid washing treatment and dried. Following a solvent exchange step, the extracts can be cleaned up by column chromatography on silica gel, alumina or activated carbon.
- 2.5 The preparation of the final extract for HRGC/HRMS analysis is accomplished by adding 10 to 50 μ L (depending on the matrix) of a nonane solution containing 50pg/ μ L of the recovery standards $^{13}\text{C}_{12}$ -1,2,3,4-TCDD and $^{13}\text{C}_{12}$ -1,2,3,7,8,9-HxCDD. The former is used to determine the percent recoveries of tetra- and pentachlorinated PCDD/PCDF congeners, while the latter is used to determine the percent recoveries of the hexa-, hepta- and octachlorinated PCDD/PCDF congeners.
- 2.6 One μ L of the concentrated extract are injected into an HRGC/HRMS system capable of performing selected ion monitoring at resolving powers of at least 10,000 (10 percent valley definition).
- 2.7 The identification of OCDD and nine of the fifteen 2,3,7,8-substituted congeners, for which a $^{13}\text{C}_{12}$ -labeled standard is available in the sample fortification and recovery standard solutions, is based on their elution at their exact retention time (within 0.005 retention time units measured in the routine calibration) and the simultaneous detection of the two most abundant ions in the molecular ion region. The remaining six 2,3,7,8-substituted congeners (i.e., 2,3,4,7,8-PeCDF; 1,2,3,4,7,8-HxCDD; 1,2,3,6,7,8-HxCDF; 1,2,3,7,8,9-HxCDF; 2,3,4,6,7,8-HxCDF, and 1,2,3,4,7,8,9-HpCDF), for which no carbon-labeled internal standards are available in the sample fortification solution, and all other PCDD/PCDF congeners are identified when their relative retention times fall within their respective PCDD/PCDF retention time windows, as established from the routine calibration data, and the simultaneous detection of the two most abundant ions in the molecular ion region. The identification of OCDF is based on its retention time relative to $^{13}\text{C}_{12}$ -OCDD and the simultaneous detection of the two most abundant ions in the molecular ion region. Identification is also based on a comparison of the ratios of the integrated ion abundance of the molecular ion species to their theoretical abundance ratios.

2.8 Quantitation of the individual congeners, total PCDDs and total PCDFs is achieved in conjunction with the establishment of a multipoint (five points) calibration curve for each homologue, during which each calibration solution is analyzed once.

3.0 DEFINITIONS

3.1. Abbreviations

PCDD	=	Polychlorinated dibenzo- <i>p</i> -dioxin
PCDF	=	Polychlorinated dibenzofuran
TCDD	=	Tetrachlorodibenzo- <i>p</i> -dioxin
PeCDD	=	Pentachlorodibenzo- <i>p</i> -dioxin
HxCDD	=	Hexachlorodibenzo- <i>p</i> -dioxin
HpCDD	=	Heptachlorodibenzo- <i>p</i> -dioxin
OCDD	=	Octachlorodibenzo- <i>p</i> -dioxin
TCDF	=	Tetrachlorodibenzofuran
PeCDF	=	Pentachlorodibenzofuran
HxCDF	=	Hexachlorodibenzofuran
HpCDF	=	Heptachlorodibenzofuran
OCDF	=	Octachlorodibenzofuran
HxCDPE	=	Hexachlorodiphenyl ether
HpCDPE	=	Heptachlorodiphenyl ether
OCPE	=	Octachlorodiphenyl ether
NCDPE	=	Nonachlorodiphenyl ether
DCDPE	=	Decachlorodiphenyl ether
PFK	=	Perfluorokerosene
CS	=	Cleanup standard
IS	=	Internal standard
RS	=	Recovery standard
HRGC	=	High-resolution gas chromatography
HRMS	=	High-resolution mass spectrometry
TEF	=	Toxicity equivalence factor
TEQ	=	Toxicity equivalent

3.2 Isotope-Labeled Standards

- 3.2.1 Recovery standards are added to the sample extract immediately before analysis by HRGC/HRMS. Recoveries of the cleanup and internal standards are determined by comparing the peak areas of the cleanup and internal standards with the peak areas of the recovery standards.
- 3.2.2 The cleanup standard is added to the sample extract prior to initiation of cleanup procedures. Loss of cleanup standard reflects losses occurring during cleanup. Generally, greater losses of cleanup standard occur when exhaustive cleanup techniques are indicated.
- 3.2.3 Internal standards are added to the sample before extraction. Losses of internal standards reflect losses occurring during both extraction and cleanup. It is difficult to recover the internal standard (and the analytes) from some sample matrices. Severe losses of internal standards may result in impaired detection limits.

3.3 Toxicity Equivalence Factors

- 3.3.1 Not all of the 17 regulated dioxin/furan isomers have the same degree of toxicity. Toxicity Equivalence Factors (TEF), relative to the most toxic dioxin/furan isomer, 2,3,7,8-TCDD, have been established. The concentration of each 2,3,7,8-congener is multiplied by the appropriate toxicity equivalence factor (TEF). The individual results of these calculations are summed to determine the 2,3,7,8-TCDD toxicity equivalent (TEQ). Table 8

4.0 HEALTH and SAFETY WARNINGS

- 4.1. This method is to be used only by analysts experienced with residue analysis and skilled in mass spectral analytical techniques.
- 4.2 Safety training for working with PCDDs and PCDFs is required prior handling these chemicals.
- 4.3. Because of the extreme toxicity of these compounds, the analyst must take necessary precautions to prevent exposure to self, or to others, of materials known or believed to contain PCDDs or PCDFs.

5.0 CAUTIONS

- 5.1 Low-level contamination is always a possibility for HRGC/HRMS analysis, due to the chemical properties of dioxins/furans.

- 5.1.1. When HRGC/HRMS glassware becomes contaminated, it is discarded.
- 5.1.2 Disposable glassware and material is used whenever possible during extraction and cleanup.

6.0 INTERFERENCES

- 6.1 Solvents, reagents, glassware and other sample processing hardware may yield discrete artifacts or elevated baselines that may cause misinterpretation of the chromatographic data. All of these materials must be demonstrated to be free from interferences under the conditions of analysis by performing laboratory method blanks. Analysts should avoid using PVC gloves.
- 6.2 The use of high purity reagents and solvents helps minimize interference problems. Purification of solvents by distillation in all-glass systems may be necessary.
- 6.3 Interferences coextracted from the sample will vary considerably from matrix to matrix. PCDDs and PCDFs are often associated with other interfering chlorinated substances such as polychlorinated biphenyls (PCBs), polychlorinated diphenyl ethers (PCDPEs), polychlorinated naphthalenes, and polychlorinated alkyldibenzofurans. These analytes may be found at concentrations several orders of magnitude higher than the analytes of interest. Retention times of target analytes must be verified using reference standards. These values must correspond to the established retention time windows. While cleanup techniques are provided as part of this method, unique samples may require additional cleanup steps to achieve lower detection limits.
- 6.4 A high-resolution capillary column (60m DB-5, J&W Scientific, or equivalent) is used in this method. However, no single column is known to resolve all isomers. The 60m DB-5 GC column is capable of 2,3,7,8-TCDD isomer specificity. In order to determine the concentration of the 2,3,7,8-TCDF (if detected on the DB-5 column), the sample extract must be reanalyzed on a column capable of 2,3,7,8-TCDF isomer specificity (e.g., DB-225, SP-2330, SP-2331, or equivalent).

7.0 PERSONNEL QUALIFICATIONS

All analysts must demonstrate proficiency in the method by completing the following Training Plan:

Training Plan for Analysis of Method 8290 by HRMS

SOP File: HRMS-8290 (Rev. 4/01-09-01)

Trainee: _____

- | | | | |
|---|----------------|----------------|-------------|
| 1. Read and study SOP | Trainer: _____ | Trainee: _____ | Date: _____ |
| 2. Read Methods 8000B and 8290 | Trainer: _____ | Trainee: _____ | Date: _____ |
| 3. Demonstrated scientific understanding of the analysis | Trainer: _____ | Trainee: _____ | Date: _____ |
| Sample preparation | | | |
| HR-Gas chromatography | | | |
| HR-Mass spectrometry | | | |
| 4. Demonstrated familiarity with related SOPs | Trainer: _____ | Trainee: _____ | Date: _____ |
| SOP for Analytical Batches and Analytical Sequences | | | |
| SOP for Making Entries into Logbooks and onto Benchsheets | | | |
| SOP for Manual Integration of Chromatographic Peaks | | | |
| SOP for Significant Figures | | | |
| SOP for Nonconformity and Corrective Action Documentation | | | |
| SOP for Determination of Method Detection Limits | | | |
| 5. Observe performance of SOP | Trainer: _____ | Trainee: _____ | Date: _____ |
| - sample preparation (soil, water, other matrices) and sample loading | | | |
| - analytical sequence setup | | | |
| - initial calibration and continuing calibration verification | | | |
| - sample analysis | | | |
| - software introduction | | | |
| - data reduction and reporting | | | |
| 6. Perform SOP with supervision | Trainer: _____ | Trainee: _____ | Date: _____ |
| - sample preparation (soil, water, other matrices) and sample loading | | | |
| - analytical sequence setup | | | |
| - initial calibration and continuing calibration verification | | | |
| - sample analysis | | | |
| - software use | | | |
| - data reduction and reporting | | | |
| 7. Independent performance of the SOP | Trainer: _____ | Trainee: _____ | Date: _____ |
| - sample preparation (soil, water, other matrices) and sample loading | | | |
| - analytical sequence setup | | | |
| - initial calibration and continuing calibration verification | | | |
| - sample analysis | | | |
| - software proficiency | | | |
| - data reduction and reporting | | | |
| - initial demonstration of competency | | | |
| - IPR study | | | |

- single blind PE sample

8. Instrument operation and maintenance
- autosampler
 - gas chromatograph and capillary column installation
 - mass spectrometer
 - data system

Trainer: _____ Trainee: _____ Date: _____

8.0 EQUIPMENT AND SUPPLIES

- 8.1 High-Resolution Gas Chromatograph/High-Resolution Mass Spectrometer/Data System (HRGC/HRMS/DS) - The GC must be equipped for temperature programming, and all required accessories must be available, such as syringes, gases, and capillary columns.

8.1.1 GC Injection Port - The GC injection port must be designed for capillary columns. The use of splitless injection techniques is recommended. On column 1 µL injections can be used on the 60m DB-5 column. The use of a moving needle injection port is also acceptable. When using the method described in this protocol, a 2 µL injection volume is used consistently (i.e., the injection volumes for all extracts, blanks, calibration solutions and the performance check samples are 2 µL). One-µL injections are allowed; however, laboratories must remain consistent throughout the analyses by using the same injection volume at all times.

8.1.2 Gas Chromatograph/Mass Spectrometer (GC/MS) Interface – The GC/MS interface components should withstand 350°C. The interface must be designed so that the separation of 2,3,7,8-TCDD from the other TCDD isomers achieved in the gas chromatographic column is not appreciably degraded. Cold spots or active surfaces (adsorption sites) in the GC/MS interface can cause peak tailing and peak broadening. It is recommended that the GC column be fitted directly into the mass spectrometer ion source without being exposed to the ionizing electron beam. Graphite ferrules should be avoided in the injection port because they may adsorb PCDD and PCDF compounds. Vespel or equivalent, ferrules are recommended.

8.1.3 Mass Spectrometer - The static resolving power of the instrument must be maintained at a minimum of 10,000 (10 percent valley). Figure 4

8.1.4 Data System - A dedicated data system is employed to control the rapid multiple-ion monitoring process and to acquire the data. Quantitation data (peak areas or peak heights) and SIM traces (displays of intensities of each ion signal being monitored including the lock-mass ion as a function of time) must be acquired

during the analyses and stored. Quantitations may be reported based upon computer generated peak areas or upon measured peak heights (chart recording). The data system must be capable of acquiring data at a minimum of 10 ions in a single scan. It is also recommended to have a data system capable of switching to different sets of ions (descriptors) at specified times during an HRGC/HRMS acquisition. The data system should be able to provide hard copies of individual ion chromatograms for selected gas chromatographic time intervals. It should also be able to acquire mass spectral peak profiles and provide hard copies of peak profiles to demonstrate the required resolving power. The data system should permit the measurement of noise on the base line.

8.2 GC Columns

8.2.1 In order to have an isomer specific determination for 2,3,7,8-TCDD and to allow the detection of OCDD/OCDF within a reasonable time interval in one HRGC/HRMS analysis, use of the 60m DB-5 fused silica capillary column is recommended. Minimum acceptance criteria must be demonstrated and documented. At the beginning of each 12 hour period (after mass resolution and GC resolution are demonstrated) during which sample extracts or concentration calibration solutions will be analyzed, column operating conditions must be attained for the required separation on the column to be used for samples. Isomer specificity for all 2,3,7,8-substituted PCDDs/PCDFs cannot be achieved on the 60m DB-5 GC column alone. In order to determine the proper concentrations of the individual 2,3,7,8-substituted congeners, the sample extract must be reanalyzed on another GC column that resolves the isomers.

8.2.2 30m DB-225 fused silica capillary column, (J&W Scientific) or equivalent.

8.3 Miscellaneous Equipment and Materials - The following list of items does not necessarily constitute an exhaustive compendium of the equipment needed for this analytical method.

8.3.1 Nitrogen evaporation apparatus with variable flow rate.

8.3.2 Balances capable of accurately weighing to 0.001g.

8.3.3 Centrifuge.

8.3.4 Water bath, capable of being temperature controlled within $\pm 2^{\circ}\text{C}$.

8.3.5 Stainless steel or glass container large enough to hold contents of one pint sample container.

- 8.3.6 250mL polypropylene beaker.
- 8.3.7 Drying oven.
- 8.3.8 20mL scintillation vials.
- 8.3.9 Laboratory hoods.
- 8.3.10 Pipets, disposable, Pasteur, 150mm long x 5mm ID.
- 8.3.11 Pipets, disposable, serological, 10mL, for the preparation of the carbon columns.
- 8.3.12 2mL screw-top vials.
- 8.3.13 Electric meat grinder with a 3 to 5mm hole size inner plate.
- 8.3.14 Carbon, 120-440 mesh, activated just before use.
- 8.3.15 Teflon boiling chips (or equivalent).

NOTE: Teflon boiling chips may float in methylene chloride, may not work in the presence of any water phase, and may be penetrated by nonpolar organic compounds.
- 8.3.16 6mL and 60mL polypropylene reservoirs.
- 8.3.17 Frits that fit in the polypropylene reservoirs.
- 8.3.18 Glass fiber filters, 0.50 μ m, Whatman GFF, or equivalent.
- 8.3.19 Dean-Stark trap, 5 or 10mL.
- 8.3.20 All glass Soxhlet apparatus, 500mL flask.
- 8.3.21 Soxhlet/Dean Stark extractor (optional), all glass, 500mL flask.
- 8.3.22 Glass funnels, sized to hold 170mL of liquid.
- 8.3.23 Desiccator.
- 8.3.24 2-liter Separatory funnels/Screw top jars
- 8.3.25 Rotary evaporator with a temperature controlled water bath.

- 8.3.26 Glass wool.
- 8.3.27 Extraction jars, glass, 250mL, with teflon lined screw cap.
- 8.3.28 Volumetric flasks, Class A - 5mL to 100mL.
- 8.3.29 Auto sampler vials with 150uL inserts.

NOTE: Reuse of glassware should be minimized to avoid the risk of contamination. All glassware that is reused must be scrupulously cleaned as soon as possible after use, according to the following procedure: Rinse glassware with water. Soak glassware overnight in a mixture of equal parts of water, Alconox and RBS-35. Rinse with water and allow drying. Rinse with high purity acetone, toluene, dichloromethane and hexane. Store in a clean environment

- 8.4 All standards should be logged in and numbered when delivered. Details of all dilutions of standards should also be entered in the PCDD/PCDF Standards Logbook. *Note:* store at room temperature in the dark—do not refrigerate.
- 8.5 Organic-free reagent water - All references to water in this method refer to organic-free reagent water
- 8.6 Column Chromatography Reagents
- 8.7 Carbon 120-440 mesh (Supelco Envi-Carb or equivalent); weigh 0.5g into 20mL scintillation vials and activate overnight at a minimum of 100 C.
- 8.8 Silica gel, high purity grade, type 60, 70-230 mesh. Activate overnight at a minimum temperature of 120⁰C. Store in a glass bottle sealed with Teflon lined screw cap.
- 8.9 Silica gel impregnated with sodium hydroxide. Add one part (by weight) of 1M NaOH solution to two parts (by weight) silica gel (activated) in a screw cap bottle and mix with a glass rod until free of lumps. Store in desiccator in a glass bottle sealed with Teflon lined screw cap.
- 8.10 Silica gel impregnated with 40 percent (by weight) sulfuric acid. Add two parts (by weight) concentrated sulfuric acid to three parts (by weight) silica gel (activated), mix with a glass rod until free of lumps. Store in desiccator in a glass bottle sealed with Teflon lined screw cap.
- 8.11 Quartz sand.

- 8.12 Extraction thimbles, 43mm x 123mm.
- 8.13 Sulfuric acid, concentrated, ACS grade, specific gravity 241.84.
- 8.14 Sodium chloride, analytical reagent, 5 percent (w/v) in organic-free reagent water.
- 8.15 Desiccating agent
- 8.16 Sodium Sulfate (powder, anhydrous), Na_2SO_4 .
- 8.17 Solvents
 - 8.17.1 Methylene Chloride. High purity distilled in glass or highest available purity.
 - 8.17.2 Hexane. High purity distilled in glass or highest available purity.
 - 8.17.3 Methanol. High purity distilled in glass or highest available purity.
 - 8.17.4 Nonane. High purity distilled in glass or highest available purity.
 - 8.17.5 Toluene. High purity distilled in glass or highest available purity.
 - 8.17.6 Tridecane, High purity distilled in glass or highest available purity.
 - 8.17.7 Acetone. High purity distilled in glass or highest available purity.
- 8.18 High-Resolution Concentration Calibration Solutions. Five Nonane solutions containing unlabeled (17) and carbon-labeled (11) PCDDs and PCDFs at known concentrations are used to calibrate the instrument. The concentration ranges are homologue dependent, with the lowest values for the tetrachlorinated dioxin and furan (1.0 pg/ μL) and the highest values for the octachlorinated congeners (1000pg/ μL). Table 4
 - 8.18.1 Calibration solutions may be obtained from the Environmental Monitoring Systems Laboratory, U.S. EPA Cincinnati, Ohio. However, additional secondary standards must be obtained from commercial sources, and solutions should be prepared in the analyst's laboratory. It is the responsibility of the laboratory to ascertain that the calibration solutions received (or prepared) are indeed at the appropriate concentrations before they are used to analyze samples.
 - 8.18.2 Store the concentration calibration solutions in 1mL minivials at room temperature in the dark.
 - 7.18.3 Calibration solutions, Set contains one 0.2ml ampule each of HRCC1, HRCC2, HRCC3, HRCC4 and HRCC5 solutions. Table 4
- 8.19 The Window Defining/GC Column Performance Mix Solution - This solution contains the first and last eluting isomers for each homologous series from tetra- through heptachlorinated congeners. The solution also contains a series of other TCDD isomers for the purpose of documenting the chromatographic resolution. The $^{13}\text{C}_{12}$ -2,3,7,8-TCDD and TCDF are also present. Table 6

- 7.19.1 The three nonspecific TCDD isomers (ng/ml) in this solution used to verify column performance. Table 6
- 8.20 Internal standard solution- Individual labeled analytes are combined to make an internal compound stock. This stock solution is diluted with Nonane to make internal standard spiking solution. Table 9
- 8.22 Recovery standard solution- Individual labeled analytes are combined to make a recovery compound stock. This stock solution is diluted with Nonane to make recovery standard spiking solution. Table 9
- 8.23 Cleanup standard solution- Individual labeled analyte is used to make a stock. This stock solution is diluted with Nonane to make cleanup standard spiking solution
- 8.24 Matrix spiking solution (natives) - Individual natives analytes are combined to make a matrix reference compound stock. This stock solution is diluted with Nonane to make matrix reference spiking solution. Table 9

9.0 INSTRUMENT CALIBRATION AND STANDARDIZATION

- 9.1 The total cycle time for data acquisition must be < 1 second. The total cycle time includes the sum of all the dwell times and voltage reset times.
- 9.2 Acquire SIM data for all the ions listed in the five descriptors. Table 5
- 9.3 Initial Calibration - Initial calibration is required before any samples analyzed for PCDDs and PCDFs. Initial calibration is also required if any routine calibration does not meet the required criteria. All five high-resolution concentration calibration solutions must be used for the initial calibration.
- 9.4 Tune the instrument with PFK.
- 9.5 Inject 1 μ L of the GC column performance check solution and acquire SIM mass spectral data. The total cycle time must be <1 second. The laboratory must not perform any further analysis until it is demonstrated and documented that the column performance check criterion was met.
- 9.6 By using the same GC and MS conditions that produced acceptable results with the column performance check solution, analyze a 1 μ L portion of each of the five concentration calibration solutions once with the following mass spectrometer operating parameters.

- 9.7 The ratio of integrated ion current (homologous series quantitation ions) must be within the indicated control limits (set for each homologous series) for all unlabeled calibration standards.
- 9.8 The ratio of integrated ion current for the ions belonging to the carbon-labeled internal and recovery standards must be within the acceptable control limits.
- 9.9 For each selected ion current profile (SICP) and for each GC signal corresponding to the elution of a target analyte and of its labeled standards, the signal-to-noise ratio (S/N) must be better than or equal to 2.5. Measurement of S/N is required for any GC peak that has an apparent S/N of less than 5:1. The result of the calculation must appear on the SICP above the GC peak in question. Figure 3
- 9.10 Calculate the 17 relative response factors (RF) for unlabeled target analytes [RF(n); n = 1 to 17] relative to their appropriate internal standards and the nine RFs for the labeled $^{13}\text{C}_{12}$ internal standards [RF(m); m = 18 to 26] relative to the recovery standards according to the following formula:

$$RF_n = \frac{(A_n^1 + A_n^2) \times Q_{is}}{(A_{is}^1 + A_{is}^2) \times Q_n} \qquad RF_{is} = \frac{(A_{is}^1 + A_{is}^2) \times Q_{re}}{(A_{rs}^1 + A_{rs}^2) \times Q_{is}}$$

Where:

A_n^1 and A_n^2 = sum of the integrated ion abundances of the quantitation ions for unlabeled PCDDs/PCDFs.

A_{is}^1 and A_{is}^2 = sum of the integrated ion abundances of the quantitation ions for the labeled internal standard PCDDs/PCDFs.

A_{rs}^1 and A_{rs}^2 = sum of the integrated ion abundances of the quantitation ions the labeled recovery standards.

Q_{is} = quantity of the internal standard injected is (pg)

Q_{rs} = quantity of the recovery standard injected (pg).

Q_n = quantity of the unlabeled PCDD/PCDF analyte injected (pg).

The RF_n and RF_{is} are dimensionless quantities; the units used to express Q_{is} , Q_{rs} and Q_n must be the same.

- 9.11 The mean response factor (RF) and the percent relative standard deviation (%RSD) are determined by the following equations

$$\overline{RF}_n = \frac{\sum_{j=1}^5 RF_{n(j)}}{5}; \quad \%RSD = \frac{SDCR}{MVCS}$$

Where:

n represents a particular PCDD/PCDF (2,3,7,8-substituted) congener (n = 1 to 17), and j is the injection number (or calibration solution number; j = 1 to 5).

SDCR= Standard Deviation of the Calibration response

MVCS= Mean Value of the Calibration response

- 9.12 The relative response factors to be used for the determination of the concentration of total isomers in a homologous series are calculated as follows:

- 9.12.1 For congeners that belong to a homologous series with only one isomer (e.g., OCDD and OCDF) or only one 2,3,7,8-substituted isomer (TCDD, PeCDD, HpCDD, and TCDF), the mean RF used will be the same as the mean RF.

NOTE: The calibration solutions do not contain $^{13}\text{C}_{12}$ -OCDF as an internal standard. This is because a minimum resolving power of 12,000 is required to resolve the [M+6] ion of $^{13}\text{C}_{12}$ -OCDF from the [M+2] ion of OCDD (and [M+4] from $^{13}\text{C}_{12}$ -OCDF with [M] of OCDD). Therefore, the RF for OCDF is calculated relative to $^{13}\text{C}_{12}$ -OCDD.

- 9.12.2 For congeners that belong to a homologous series containing more than one 2,3,7,8-substituted isomer, the mean RF used for those homologous series will be the mean of the RFs calculated for all individual 2,3,7,8-substituted congeners using the equation below:

$$\overline{RF}_k = \frac{1}{t} \sum_{n=1}^t RF_n$$

where:

k = 27 to 30; with 27 = PeCDF; 28 = HxCDF;
29 = HxCDD; and 30 = HpCDF (Table 11)

t = total number of 2,3,7,8-substituted isomers present in the calibration solutions for each homologous series (e.g., two for PeCDF, four for HxCDF, three for HxCDD, two for HpCDF).

NOTE: Presumably, the HRGC/HRMS response factors of different isomers within a homologous series are different. However, this analytical protocol will make the assumption that the HRGC/HRMS responses of all isomers in a homologous series that do not have the 2,3,7,8-substitution pattern are the same as the responses of one or more of the 2,3,7,8-substituted isomer(s) in that homologous series.

- 9.12.3 Relative response factors [RF_m] to be used for the determination of the percent recoveries for the nine internal standards are calculated as follows:

$$RF_m = \frac{A_{is}^m \times Q_{rs}}{Q_{is}^m \times A_{rs}}$$

$$\overline{RF} = 1/5 \sum_{j=1}^5 RF_m(j)$$

$m = 18$ to 26 (congener type) and $j = 1$ to 5 (injection number),

A_{is}^m = sum of the integrated ion abundances of the quantitation ions for a given internal standard ($m = 18$ to 26),

A_{rs} = sum of the integrated ion abundances of the quantitation ions for the appropriate recovery standard, Q_{rs} , Q_{is}^m = quantities of, respectively, the recovery standard (rs) and a particular internal standard ($is = m$) injected, (pg),

RF_m = relative response factor of a particular internal standard (m) relative to an appropriate recovery standard, as determined from one injection, and RF = calculated mean relative response factor of a particular internal standard (m) relative to an appropriate recovery standard, as determined from the five initial calibration injections (j).

9.12.4 Criteria for Acceptable Calibration - The criteria listed below for acceptable calibration must be met before sample analyses are performed.

- 9.12.4.1 The percent relative standard deviations for the mean response factors [RF_n and RF_m] from the 17 unlabeled standards must not exceed ± 20 percent, and those for the nine labeled reference compounds must not exceed ± 30 percent.
- 9.12.4.2 The S/N for the GC signals present in every SICP (including the ones for the labeled standards) must be ≥ 10 .
- 9.12.4.3 The ion abundance ratios must be within the specified control limits.

NOTE: If the criterion for acceptable calibration is met, the analyte specific RF can then be considered independent of the analyte quantity for the calibration concentration range. The mean RFs will be used for all calculations until the routine calibration criteria are no longer met. At such time, new mean RFs will be calculated from a new set of injections of the calibration solutions.

9.12.5 Routine Calibration (Continuing Calibration Check) – Routine calibrations must be performed at the beginning of a 12-hour period after successful mass resolution and GC resolution performance checks. A routine calibration is also required at the end of a 12-hour shift. Inject 1 μ L of the concentration calibration solution HRCC-3 standard. By using the same HRGC/HRMS conditions, determine and document an acceptable calibration.

9.12.6 Criteria for Acceptable Routine Calibration - The following criteria must be met before further analysis is performed.

- 9.12.6.1 The measured RFs [RF_n for the unlabeled standards] obtained during the routine calibration runs must be within ± 20 percent of the mean values established during the initial calibration.
- 9.12.6.2 The measured RFs [RF_m for the labeled standards] obtained during the routine calibration runs must be within ± 30

percent of the mean values established during the initial calibration.

- 9.12.6.3 The ion abundance ratios must be within the allowed control limits.
- 9.12.6.4 If either one these criteria, is not satisfied, repeat one more time. If these criteria are still not satisfied, the entire routine calibration process must be reviewed. It is realized that it may not always be possible to achieve all RF criteria. For example, it has occurred that the RF criteria for $^{13}\text{C}_{12}$ -HpCDD and $^{13}\text{C}_{12}$ -OCDD were not met, however, the RF values for the corresponding unlabeled compounds were routinely within the criteria established in the method. In these cases, 24 of the 26 RF parameters have met the QC criteria, and the data quality for the unlabeled HpCDD and OCDD values were not compromised as a result of the calibration event. In these situations, the analyst must assess the effect on overall data quality as required for the data quality objectives and decide on appropriate action. Corrective action would be in order, for example, if the compounds for which the RF criteria were not met included both the unlabeled and the corresponding internal standard compounds.

NOTE: An initial calibration must be carried out whenever the HRCC-3, the sample fortification, or the recovery standard solution is replaced by a new solution from a different lot.

10.0 SAMPLE COLLECTION

- 10.1. Samples are typically delivered by a express shipping service or by the client.

11.0 SAMPLE HANDLING AND PRESERVATION

- 11.1 Samples are logged into the Sample Receiving Logbook, and labeled. A database entry is also made, which triggers the printing of a set of tracking forms.
- 11.2 Storage and Holding Times - All samples, except fish and adipose tissue samples, must be stored at 4⁰C in the dark, extracted within 30 days and completely analyzed within 45 days of extraction. Fish and adipose tissue samples must be stored at -10⁰C \pm 2⁰C in the dark, extracted within 30 days and completely analyzed within 45 days of collection. Whenever samples are analyzed after the holding time expiration date, the results should be considered to be minimum concentrations and should be identified as such.

NOTE: The holding times listed in this method are recommendations. PCDDs and PCDFs are very stable in a variety of matrices, and holding times may be as high as a year for certain matrices. Sample extracts, however, should always be analyzed within 45 days of extraction.

- 11.3. Fish and tissue samples must be shipped on dry ice and arrive at 4 \pm 2 degrees C

12.0 SAMPLE PREPARATION AND ANALYSIS

- 12.1. Before extraction, a visual inspection of a sample is done. If it is not homogeneous, the entire contents of the sample container is transferred to a sheet of aluminum foil and then mixed thoroughly with a polystyrene spoon prior to removing an aliquot for analysis.
- 12.2 If a soil, sediment, or paper pulp sample contains more than 25 percent water, the two phases are separated, the aqueous phase is discarded, and the remaining solid phase is analyzed.
- 12.2.1 Transfer an estimated 50g aliquot of the sample to a centrifuge tube. Centrifuge for 15 minutes at 1000rpm.
- 12.2.2. Thoroughly mix the solid with a polystyrene spoon, and weigh out aliquots for analysis and dry weight determination.
- 12.2.3. Return the remaining sample to the original sample container.

- 12.3. If a water sample has a solids content greater than or equal to 1 percent, then the solid phase is separated from the aqueous phase, and both phases are extracted separately.

12.3.1. Filter the sample through a 0.5 μ m filter.

12.3.2. Extract the filter and solids as for a soil sample. Extract the aqueous filtrate as for a water sample.

12.3.3 Combine the two extracts and proceed with the cleanup and analysis.

12.4 Addition of Internal Standard

12.4.1 Use a portion of 10g of the sample to be analyzed. Transfer the sample portion to a tared extraction thimble and determine its weight.

12.4.2 Except for adipose tissue, add 100 μ L of internal standard solution to the sample(s).

12.4.3 For water samples, mix the internal standard solution with 10.0mL acetone in a pre-labeled scintillation vial and add to the corresponding sample container.

12.5 Extraction of Fish and Tissue

12.5.1 To a tared thimble, add approximately 30g anhydrous sodium sulfate to a 20g portion of a homogeneous fish samples and mix thoroughly with a pasteur pipet. After breaking up any lumps, place the spiked 200 μ L of internal standard into the fish/Sodium Sulfate mixture in the Soxhlet apparatus. Add approximately 270mL Hexane/Methylene Chloride (1:1) to the Soxhlet apparatus and reflux for 16 hours. The solvent must cycle completely through the system five times per hour.

NOTE: As an option, a rotary evaporator may be used in place of the KD apparatus for the concentration of the extracts.

12.5.2 Evaporate the extract to 5-10 mL. Do not allow the extract to go to dryness.

12.5.3 To a tarred 20 mL scintillation vial, transfer exactly one half of the extract. Weigh the scint vial containing the extract, record the weight and allow to dry in a gravity oven overnight before weighing a second time.

12.5.4 Add 500 μ L tridecane to the 500mL flask.. Concentrate the extract on a rotary evaporator to an apparent volume of 0.5mL.

- 12.5.5 Add 20mL hexane to the 500mL flask. Again concentrate to 0.5mL, using either a rotary evaporator or by air-drying in the hood.
- 12.5.6 Add 30mL Hexane and decant the contents of the 500mL flask into a 250mL screw top flask. Rinse the 500mL flask with an additional 30mL portion of hexane in the sonicator and add the rinse to the jar.
- 12.5.7 Spike with 100 μ L cleanup standard and proceed with sulfuric acid clean up.

12.6 Fuel Oil/Chemical

- 12.6.1 Extract sample by mixing 1g of sample with 60mL hexane in a 250mL jar. Spike with 100 μ L internal spiking standard.
- 12.6.2 The sample extract volume should be in 60 mL hexane. Partition the Hexane extract against approximately 10mL of concentrated sulfuric acid. Shake for 30 seconds. Allow a minimum of 1 hour for separation. Remove and discard the sulfuric acid layer (bottom). Repeat the acid washing until little color is visible in the acid layer (perform a maximum of four acid washings).
- 12.6.3 Partition the extract against 20mL of 5 percent (w/v) sodium chloride/HPLC water. Shake for 30 seconds. Allow a minimum 20 minutes for separation. Remove and discard the aqueous layer (bottom).
- 12.6.4 Add 500 μ L tridecane to each sample jar and allow to evaporate until near dryness, using a rotary evaporator or by air-drying under the hood. The tridecane will prevent the sample from going to complete dryness. The sample is now ready for column cleanup

12.7 Fly Ash

NOTE: Because of the tendency of fly ash to "fly", all handling steps should be performed in a hood in order to minimize contamination.

- 12.7.1 Weigh about 10g fly ash to two decimal places in a tared 250 mL extraction jar. Add 100 μ L internal standard solution to the sample(s). Add 150mL of 1M HCl to the fly ash sample. Seal the jar with the Teflon lined screw cap and shake for 3 hours at room temperature.

12.7.1 Rinse a 0.5µm filter with organic-free reagent water. Assemble a Buchner funnel into a 1L flask and place the filter in the funnel. Pour the sample on the filter and filter the fly ash cake with approximately 500mL organic-free reagent water. A disposable filter apparatus can be used.

12.7.2 Place the sample and the filter paper into an extraction thimble. Add 10g anhydrous powdered sodium sulfate.

12.7.3 Extract in a Soxhlet extraction apparatus charged with approximately 270mL toluene for 16 hours using a five cycle/hour schedule.

NOTE: As an option, a Soxhlet/Dean Stark extractor system may be used, with toluene as the solvent. No sodium sulfate is added when using this option.

12.7.4 Cool and add 0.5mL tridecane to the 500mL round bottom flask. Concentrate the extract to near dryness on a rotary evaporator at 40-50°C. Add 20mL Hexane and evaporate to near dryness, or alternatively, air-dry under the hood.

12.7.5 Transfer the concentrate to a 250 screw top jar using 20mL Hexane. Rinse the flask in a sonicator with two 20mL portions of hexane and add the rinses to the jar.

12.7.6 The sample is now ready for sulfuric acid clean up.

12.8 Aqueous samples

12.8.1 Allow the sample to come to ambient temperature. Spike 100µL internal standard solution into a 20 mL scintillation vial. Add 10 mL acetone and pour into the appropriate sample jar.

12.8.2 When the sample is judged to contain 1 percent or more solids, the sample must be filtered through a 0.5µm filter. If the suspended solids content is too great to filter through the filter, centrifuge the sample, decant, and then filter the aqueous phase.

12.8.3 Combine the solids from the centrifuge bottle(s) with the particulates on the filter and with the filter itself and proceed with the Soxhlet extraction.

12.8.4 Pour the aqueous filtrate into a 2L screw top jar/or 2L separatory funnel. Add 100mL Methylene Chloride to the sample bottle, seal and shake for 30 seconds to rinse the inner surface. Pour the solvent from the sample jar to the 2L jar/or 2L

separatory funnel. Close and shake the 2L jar/or 2L separatory funnel for 3 minutes. Periodic venting is necessary.

- 12.8.5 Allow the organic layer to separate from the water phase for a minimum of 10 minutes. If the emulsion interface between layers is more than one third the volume of the solvent layer, the analyst must employ mechanical techniques to complete the phase separation (e.g., centrifuge the extract).
- 12.8.6 Transfer the Methylene Chloride 250mL glass jar.
- 12.8.7 Repeat the extraction once more with fresh 50mL portions of methylene chloride. Transfer to a 250 mL screw top jar.
- 12.8.8 After extracting the water sample, pour into a 1L graduated cylinder and record the volume.
- 12.8.9 If the extract is colorless or has a light tint, add 500 μ L tridecane and evaporate to near 500 μ L and proceed with the silica gel/carbon column clean up. If the extract has a medium to dark tint, evaporate to near 500 μ L (using a rotary evaporator or by air-drying under the hood), and proceed to sulfuric acid clean up. Spike 100 μ L of cleanup standard before proceeding with cleanup.

12.9 Soil/Sediment/Paper Pulp

- 12.9.1 Add approximately 20g anhydrous powdered sodium to an extraction thimble. Tare the thimble. Transfer 10g of sample to the thimble. Spike with internal standard solution.

NOTE: As an option, a Soxhlet/Dean Stark extractor system may be used, with toluene as the solvent. See Tables 15, 16, 17.

- 12.9.2 Add approximately 270mL toluene to the Soxhlet apparatus and reflux for 16 hours. The solvent must cycle completely through the system five times per hour.
- 12.9.3 Cool the extraction apparatus. Add 500 μ L tridecane to the flask and concentrate to near dryness on a rotary evaporator at 45- 50°C. Add 20mL Hexane to the flask and evaporate to near dryness, with either a rotary evaporator or by air-drying under the hood.
- 12.9.4 Transfer the extract to a 250mL screw top jar with 30mL hexane. Rinse

the 500mL flat bottom in a sonicator with an additional 30mL rinse of Hexane and transfer to the screw top jar. Final volume should be approximately 60mL hexane.

12.9.5 If the extract is colorless or has a light tint, add 500 μ L tridecane and evaporate to near 500 μ L and proceed with the silica gel/carbon column clean up. If the extract has a medium to dark tint, proceed to sulfuric acid clean up.

12.9.6 Transfer approximately 1g of sample into a tarred 20 mL scint vial. Weigh and Record the weight of the vial containing the sample. Place in a gravity oven and allow drying overnight. Re-weigh the dried sample and record the weight.

12.10 Clean up Procedures

All samples must be spiked with 100 μ L cleanup standard solution before any clean up procedure is performed.

12.10.1 Sulfuric Acid Clean up

12.10.1.1 The sample extract volume should be in 60 mL hexane. Partition the Hexane extract against approximately 10mL of concentrated sulfuric acid. Shake for 30 seconds. Allow a minimum of 1 hour for separation. Remove and discard the sulfuric acid layer (bottom). Repeat the acid washing until little color is visible in the acid layer (perform a maximum of four acid washings).

12.10.1.2 Partition the extract against 20mL of 5 percent (w/v) sodium chloride/HPLC water. Shake for 30 seconds. Allow a minimum 20 minutes for separation. Remove and discard the aqueous layer (bottom).

12.10.1.3 Add 500 μ L tridecane to each sample jar and allow to evaporate until near dryness, by air-drying under the hood. The tridecane will prevent the sample from going to complete dryness.

12.10.1.4 The sample is now ready for column cleanup.

12.10.2 Silica/Carbon Column Cleanup

12.10.2.1 Pack a 6mL column as follows: Insert a frit at the bottom of the column. Remove the pre-weighed, activated carbon from the oven, and add approximately 2ml toluene. Cool down is not necessary.

Transfer the slurry to the 6mL column and allow the carbon to settle. Rinse the 20mL vial with addition rinses of toluene until most of the carbon has been transferred. Rinse the column sides of residual carbon. Allow the toluene to completely drain. Insert a frit on top and push down to the top of the carbon. Elute with 5mL dichloromethane then 5mL hexane. Discard the eluate. Check the column for channeling. If channeling is observed, discard the column. Do not tap a wetted column.

- 12.10.2.2 Pack a 60ml reservoir, with silica gel as follows: Insert a frit and push it to the bottom. Place 1 teaspoon activated neutral silica gel in the column and tap the column gently to settle the silica gel. Add 1 teaspoon sodium hydroxide-impregnated silica gel, 2 teaspoon sulfuric acid-impregnated silica gel, and 1 teaspoon sodium sulfate. Tap the column gently after each addition. Elute with 30mL hexane and stop the flow just before exposure of the top layer of sodium sulfate to air. Discard the eluate. Check the column for channeling. If channeling is observed, discard the column. Do not tap the wetted column.
- 12.10.2.3 Place the silica gel column on top of the vacuum manifold and the carbon column of the bottom opposite the silica gel column. The column is now ready for loading the sample extract.
- 12.10.2.4 Transfer the 2mL extract to the silica gel column and allow eluting until the extract level is at the top of the sodium sulfate. A slow vacuum can be used to facilitate the extract into the column.
- 12.10.2.5 Rinse the flask containing the extract with 5mL hexane and load into the Silica Gel column. Allow eluting until the extract level is at the top of the sodium sulfate. A slow vacuum can be used to facilitate the extract into the column.
- 12.10.2.6 Slowly add 30mL of Hexane to the Silica Gel column. Add another 30mL, then 20mL more of Hexane when space permits for a total of 80mL.

NOTE: At this point, neither pressure nor vacuum is necessary. Do not adjust the drip rate.

- 12.10.2.7 Remove the used silica gel column when the drip from the carbon column stops. Reverse the carbon column by placing on top of the

manifold. Elute the column with 5mL Dichloromethane, then 5mL of 1:1 Ethyl Acetate/Hexane.

12.10.2.8 Place a labeled 20mL scintillation vial under the carbon column and elute with a total of approximately 20mL Toluene. Be sure to collect the Toluene. Allow the extract to evaporate to dryness under the fume hood.

12.10.2.9 Add approximately 2mL hexane to the scintillation vial. Vortex and sonicate, then transfer the extract to a clear 2mL screw top vial. Allow evaporation to dryness under the fume hood.

NOTE: If possible, transfer the extract to the labeled autosampler vial when the volume reaches approximately 150µL.

12.10.2.10 Add 160uL Hexane to the 2mL vial. Rinse the sides of the 2mL vial with only the 160uL that was added. Transfer to an appropriately labeled autosampler vial. Allow sample to evaporate to dryness, under the hood.

12.10.2.11 The sample is ready to be spiked with 20uL recovery standard solution.

12.11 Chromatographic/Mass Spectrometric Conditions and Data Acquisition Parameters

12.11.1 Gas Chromatograph
Column coating: DB-5
Film thickness: 0.25µm
Column dimension: 60 m x 0.25mm
Injector temperature: 300⁰C
Splitless valve time: 1 min
Interface temperature: 300⁰C
Temperature program:

STAGE	INITIAL TEMP, C	INITIAL HOLD TIME, MIN	TEMPERATURE RAMP, C/MIN	FINAL TEMPERATURE C	FINAL HOLD TIME, MIN
1	150	5	35	215	5
2			1.5	230	6
3			7	315	5

12.11.2 Mass Spectrometer

12.11.2.1 The mass spectrometer must be operated in a selected ion monitoring (SIM) mode with a total cycle time (including the voltage-reset time) of one second or less. At a minimum, the ions listed in Table 5 for each of the five SIM descriptors must be monitored. Note that with the exception of the last descriptor (OCDD/OCDF), all descriptors contain 10 ions. The selection (Table 6) of the molecular ions M and M+2 for $^{13}\text{C}_{12}\text{-HxCDF}$ and $^{13}\text{C}_{12}\text{-HpCDF}$ rather than $^{13}\text{M}+2$ and $^{13}\text{M}+4$ (for consistency) was made to eliminate, even under high-resolution mass spectrometric conditions, interferences occurring in these two ion channels for samples containing high levels of native HxCDDs and HpCDDs. It is important to maintain the same set of ions for both calibration and sample extract analyses. The selection of the lock-mass ion is left to the performing laboratory.

NOTE: At the option of the analyst, the tetra- and pentachlorinated dioxins and furans can be combined into a single descriptor.

12.11.2.2 The recommended mass spectrometer tuning conditions are based on the groups of monitored ions shown in Table 5. By using a PFK molecular leak, tune the instrument to meet the minimum required resolving power of 10,000 (10% valley) at m/z 304.9824 (PFK) or any other reference signal close to m/z 303.9016 (from TCDF). Fig.4 By using peak matching conditions and the aforementioned PFK reference peak, verify that the exact mass of m/z 380.9760 (PFK) is within 5 ppm of the required value. Note that the selection of the Low Mass and High Mass ions must be such that they provide the largest voltage jump performed in any of the five mass descriptors.

12.12 Analysis

12.12.1 Remove the sample or blank extract from storage. With a stream of dry, purified nitrogen, reduce the extract volume to dryness.

12.12.12 Inject a 1µL aliquot of the extract into the GC, operated under the conditions that have been established to produce acceptable results with the performance check solution.

12.12.3 Acquire SIM data. Use the same acquisition and mass spectrometer operating conditions previously used to determine the relative response factors.

NOTE: The acquisition period must at least encompass the PCDD/PCDF overall retention time window previously determined. Selected ion current profiles (SICP) for the lock-mass ions (one per mass descriptor) must also be recorded and included in the data package. These SICPs must be true representations of the evolution of the lock-mass ions amplitudes during the HRGC/HRMS run. The analyst may be required to monitor a PFK ion, not as a lock-mass, but as a regular ion, in order to meet this requirement. It is recommended to examine the lock-mass ion SICP for obvious basic sensitivity and stability changes of the instrument during the GC/MS run that could affect the measurements. Report any discrepancies in the case narrative.

12.12.4 Identification Criteria - For a gas chromatographic peak to be identified as a PCDD or PCDF, it must meet all of the following criteria:

12.12.4.1 Retention Times

12.12.4.1.1 For 2,3,7,8-substituted congeners, which have an isotopically labeled internal or recovery standard present in the sample extract (this represents a total of 10 congeners including OCDD), the retention time (RRT; at maximum peak height) of the sample components must be within -1 to +3 seconds of the isotopically labeled standard.

- 12.12.4.1.2 For 2,3,7,8-substituted compounds that do not have an isotopically labeled internal standard present in the sample extract, the retention time must fall within 0.005 retention time units of the relative retention times measured in the routine calibration. Identification of OCDF is based on its retention time relative to $^{13}\text{C}_{12}$ -OCDD as determined from the daily routine calibration results.
- 12.12.4.1.3 For non-2,3,7,8-substituted compounds (tetra through octa; totaling 119 congeners), the retention time must be within the corresponding homologous retention time windows established by analyzing the column performance check solution.
- 12.12.4.1.4 The ion current responses for both ions used for quantitative purposes (e.g., for TCDDs: m/z 319.8965 and 321.8936) must reach maximum simultaneously (± 2 seconds).
- 12.12.4.1.5 The ion current responses for both ions used for the labeled standards (e.g., for $^{13}\text{C}_{12}$ -TCDD: m/z 331.9368 and m/z 333.9339) must reach maximum simultaneously (± 2 seconds).

NOTE: The analyst is required to verify the presence of 1,2,8,9-TCDD and 1,3,4,6,8-PeCDF in the SICPs of the daily performance checks. Should either one compound be missing, the analyst is required to take corrective action as it may indicate a potential problem with the ability to detect all the PCDDs/PCDFs.

12.12.4.2 Abundance Ratios

- 12.12.4.2.1 The integrated ion currents for the two ions used for quantitation purposes must have a ratio between the lower and upper limits established for the homologous series to which the peak is assigned. Table 7

12.12.4.3 Signal-to-Noise Ratio

12.12.4.3.1 All ion current intensities must be ≥ 2.5 times noise level for positive identification of a PCDD/PCDF compound or a group of coeluting isomers.

12.12.4.4 Polychlorinated Diphenyl Ether Interferences

12.12.4.4.1 In addition to the above criteria, the identification of a G peak as a PCDF can only be made if no signal having a S/N ≥ 2.5 is detected at the same retention time (± 2 seconds) in the corresponding polychlorinated diphenyl ether channel.

13.0 DATA ANALYSIS AND CALCULATIONS

13.1 For gas chromatographic peaks that have met the criteria, calculate the concentration of the PCDD or PCDF compounds using the formula:

$$C_x = \frac{A_x \times Q_{is}}{A_{is} \times W \times RF_n}$$

Where:

C_x = concentration of unlabeled PCDD/PCDF congeners (or group of coeluting isomers within an homologous series) in pg/g,

A_x = sum of the integrated ion abundances of the quantitation ions for unlabeled PCDDs/PCDFs,

A_{is} = sum of the integrated ion abundances of the quantitation ions for the labeled internal standards,

Q_{is} = quantity, in pg, of the internal standard added to the sample before extraction,

W = weight, in g, of the sample (solid or organic liquid), or volume in mL of an aqueous sample, and

RF_n = calculated mean relative response factor for the analyte
[RF_n with $n = 1$ to 17]

If the analyte is identified as one of the 2,3,7,8-substituted PCDDs or PCDFs, RF_n is the value calculated using the equation section

13.1.1

However, if it is a non-2,3,7,8-substituted congener, the $RF(k)$ value is calculated using the equation $RF(k) = 27$ to 30]. Table 11

- 13.2 Calculate the percent recovery of the nine internal standards measured in the sample extract, using the formula:

$$\% \text{ Recovery} = \frac{A_{is} \times Q_{rs}}{Q_{is} \times A_{rs} \times RF_n} \times 100$$

where:

A_{is} = sum of the integrated ion abundances of the quantitation ions for the labeled internal standard,

A_{rs} = sum of the integrated ion abundances of the quantitation ions for the labeled recovery standard; the selection of the recovery standard depends on the type of congeners,

Q_{is} = quantity, in pg, of the internal standard added to the sample before extraction,

Q_{rs} = quantity, in pg, of the recovery standard added to the cleaned-up sample residue before HRGC/HRMS analysis, and

RF_m = calculated mean relative response factor for the labeled internal standard relative to the appropriate recovery standard. This represents the mean [RF with $m = 18$ to 26]. (Table 10)

NOTE: For human adipose tissue, adjust the percent recoveries by adding 1 percent to the calculated value to compensate for the 1 percent of the extract diverted for the lipid determination.

- 13.3 If the concentration in the final extract of any of the fifteen 2,3,7,8-substituted PCDD/PCDF compounds exceeds the upper method calibration limits (MCL) listed in Table 1 (e.g., 200 pg/ μ L for TCDD in soil), the linear range of response versus concentration may have been exceeded, and a second analysis of the sample (using a one tenth aliquot) should be undertaken. The volumes of the internal and recovery standard solutions should remain the same as described for the sample preparation. For the other congeners (including OCDD), however, report the measured concentration and indicate that the value exceeds the MCL.
- 13.4 If a smaller sample size would not be representative of the entire sample, one of the following options is recommended:
- (1) Re-extract an additional aliquot of sufficient size to insure that it is representative of the entire sample. Spike it with a higher concentration of internal standard. Prior to GC/MS analysis, dilute the sample so that it has a concentration of internal standard equivalent to that present in the calibration standard. Then, analyze the diluted extract.

- (2) Re-extract an additional aliquot of sufficient size to insure that it is representative of the entire sample. Spike it with a higher concentration of internal standard. Immediately following extraction, transfer the sample to a volumetric flask and dilute to known volume. Remove an appropriate aliquot and proceed with cleanup and analysis.
- (3) Use the original analysis data to quantitate the internal standard recoveries. Respike the original extract (note that no additional cleanup is necessary) with 100 times the usual quantity of internal standards. Dilute the re-spiked extract by a factor of 100. Reanalyze the diluted sample using the internal standard recoveries calculated from the initial analysis to correct the results for losses during isolation and cleanup.

13.5 The total concentration for each homologous series of PCDD and PCDF is calculated by summing up the concentrations of all positively identified isomers of each homologous series. Therefore, the total should also include the 2,3,7,8-substituted congeners. The total number of GC signals included in the homologous total concentration value must be specified in the report. If an isomer is not detected, use zero (0) in this calculation.

13.6 Sample Specific Estimated Detection Limit - The sample specific estimated detection limit (EDL) is the concentration of a given analyte required to produce a signal with a peak height of at least 2.5 times the background signal level. An EDL is calculated for each 2,3,7,8-substituted congener that is not identified, regardless of whether or not other non-2,3,7,8-substituted isomers are present. Two methods of calculation can be used, as follows, depending on the type of response produced during the analysis of a particular sample.

13.6.1 Samples giving a response for both quantitation ions that are less than 2.5 times the background level.

13.6.2 Use the expression for EDL (specific 2,3,7,8-substituted PCDD/PCDF) below to calculate an EDL for each absent 2,3,7,8-substituted PCDD/PCDF (i.e., $S/N \leq 2.5$). The background level is determined by measuring the range of the noise (peak to peak) for the two quantitation ions of a particular 2,3,7,8-substituted isomer within an homologous series, in the region of the SICP trace corresponding to the elution of the internal standard (if the congener possesses an internal standard) or in the region of the SICP where the congener is expected to elute by comparison with the routine calibration data (those congeners that do not have a $^{13}\text{C}_{12}$ -labeled standard), multiplying that noise height by 2.5,

and relating the product to an estimated concentration that would produce that peak height.

Use the formula:

$$EDL = \frac{2.5 \times H_x \times Q_{is}}{H_{is} \times W \times RF_n}$$

where:

EDL = estimated detection limit for homologous 2,3,7,8-substituted PCDDs/PCDFs.

H_x = sum of the height of the noise level for each quantitation ion for the unlabeled PCDDs/PCDFs

H_{is} = sum of the height of the signal level for each quantitation ion for the labeled internal standard.

- 13.6.3 Samples characterized by a response above the background level with a S/N of at least 2.5 for both quantitation ions.
- 13.6.4 When the response of a signal having the same retention time as a 2,3,7,8-substituted congener has a S/N in excess of 2.5 and does not meet any of the other qualitative identification criteria, calculate the "Estimated Maximum Possible Concentration" (EMPC), except that A should represent the sum of the x area under the smaller peak and of the other peak area calculated using the theoretical chlorine isotope ratio.
- 13.6.5 The relative percent difference (RPD) of any duplicate sample results are calculated as follows:

$$RPD = \frac{|S_1 - S_2|}{\frac{(S_1 + S_2)}{2}} \times 100$$

S₁ and S₂ represent sample and duplicate sample results.

- 13.6.6 The 2,3,7,8-TCDD toxicity equivalents of PCDDs and PCDFs present in the sample are calculated, if requested by the data user, according to the method recommended by the Chlorinated Dioxins Workgroup (CDWG) of the EPA and the Center for Disease Control (CDC). This method assigns a 2,3,7,8-TCDD toxicity equivalency factor (TEF) to each of the fifteen 2,3,7,8-substituted PCDDs and PCDFs and to OCDD and OCDF.

(Table 8) The 2,3,7,8-TCDD equivalent of the PCDDs and PCDFs present in the sample is calculated by summing the TEF times their concentration for each of the compounds or groups of compounds. The exclusion of other homologous series such as mono-, di-, and tri-chlorinated dibenzodioxins and dibenzofurans does not mean that they are non-toxic. The above procedure for calculating the 2,3,7,8-TCDD toxicity equivalents is not claimed by the CDWG to be based on a thoroughly established scientific foundation. The procedure, rather, represents a "consensus recommendation on science policy". Since the procedure may be changed in the future, reporting requirements for PCDD and PCDF data would still include the reporting of the analyte concentrations of the PCDD/PCDF congeners.

13.7 Two GC Column TEF Determination

- 13.7.1 The concentration of 2,3,7,8-TCDD (see note below), is calculated from the analysis of the sample extract on the 60m DB-5 fused silica capillary column. The chromatographic separation between the 2,3,7,8-TCDD and its close eluters (1,2,3,7/1,2,3,8-TCDD and 1,2,3,9-TCDD) must be equal or less than 25 percent valley. Figure 5
- 13.7.2 The concentration of the 2,3,7,8-TCDF is obtained from the analysis of the sample extract on the 30m DB-225 fused silica capillary column. However, the GC/MS conditions must be altered so that:
- (1) only the first three descriptors (i.e., tetra-, penta-, and hexachlorinated congeners) are used;
 - (2) and the switching time between descriptor 2 (pentachlorinated congeners) and descriptor 3 (hexachlorinated congeners) takes place following the elution of $^{13}\text{C}_{12}$ - 1,2,3,7,8-PeCDD.
 - (3) The chromatographic separation between the 2,3,7,8-TCDF and its close eluters (2,3,4,7-TCDF and 1,2,3,9-TCDF) must be equal or less than 25 percent valley. Figure 6

NOTE: The confirmation and quantitation of 2,3,7,8-TCDD may be accomplished on the SP-2330 GC column instead of the DB-5 column.

- 13.7.3 For a gas chromatographic peak to be identified as a 2,3,7,8-substituted PCDD/PCDF congener, it must meet the ion abundance and signal-to-

noise ratio criteria, respectively. In addition, the retention time identification criterion applies here for congeners for which a carbon-labeled analogue is available in the sample extract. However, the relative retention time (RRT) of the 2,3,7,8-substituted congeners for which no carbon-labeled analogues are available must fall within 0.006 units of the carbon-labeled standard RRT. Experimentally, this is accomplished by using the attributions described and the results from the routine calibration run on the SP-2330.

14.0 DATA AND RECORDS MANAGEMENT

- client
- 14.1 CAS reports the analytical data produced in its laboratories to the client via the certified analytical report. This report typically includes a transmittal letter, a case narrative, project information, specific test results, quality control data, chain of custody information, and any other project-specific support documentation. The following procedures describe our data reduction, validation and reporting procedures.
- 14.2 All data is initially reviewed and processed by analysts using appropriate methods (e.g. chromatographic software, instrument printouts, hand calculation, etc.) A file of all raw data is printed, reviewed for completeness and quality criteria against an in-house checklist and signed off by the analyst. The operations manager reviews all reported data against the raw data; validating completeness and quality. The final report data package is then reviewed by the project manager for compliance with previously established project requirements. Typically, all data is reported in the units and MCLs listed in Appendix C.
- 14.3 Assessment of the analytical data includes a check on data consistency by looking for comparability of duplicate analyses, comparability of previous data from the same sampling location (if available), adherence to accuracy and precision control limits, and anomalous low or high parameter values. The results of this review will be discussed with either the departmental supervisor or lab director for resolution prior to final release of the package.
- 14.4 Once the data has been checked for accuracy and acceptability, the final report and raw data is forwarded to the lab director or quality assurance coordinator, who further reviews the data package for errors. When the entire data set has been found to be acceptable the lab director signs the report, the report is distributed and the raw data is filed for approximately one year, and then archived. All hard copy and electronic backups are archived in a secured file room for a period of at least 5 years from the date of the final report. It is not unusual to have various clients require 10-year retention of records, therefore, the archivist, project chemist, and possibly the client are consulted prior to destruction of the records.
- 14.5 The integrity of the data generated in the laboratory is primarily assessed by the analyst, supervisor and project chemist through the use of a variety of measures that may include

the reagent blanks, laboratory fortified blanks, duplicates, matrix spikes and QC samples. The numerical criteria for evaluation of these QC samples are listed in Appendix C; these various QC sample analyses are evaluated using the flow diagrams found in Figures 12-1 through 12-9. Other validation measures of the data include a check of the linearity of calibration curve, an accuracy check of the QC standards and a check of the system sensitivity. Data transcriptions and calculations are also reviewed. Specific calculations used for determining the concentration or value of the measured parameters from the raw data are given in each of the analytical methods or CAS SOPs.

case 14.6 When an analyst determines that the data has met the data quality objectives (and/or any client-specific data quality objectives) of the method and has qualified any anomalies in a clear, acceptable fashion, the data is validated by the supervisor. Prior to release of the report to the client, the project manager must also review the entire body of data for completeness and to ensure that any and all client-specified objectives were successfully achieved. If required, samples exceeding any established state/federal maximum contaminant level or reportable concentration level, must be reported to the client. A narrative may be written by the project manager to explain any unusual problems with a specific analysis or sample, client-specific objectives, exceedences, etc... The original raw data, along with a copy of the final report, is archived. CAS maintains control of analytical results by adhering to standard operating procedures and by observing sample custody requirements. All data are calculated and reported in units consistent with project specifications, to enable easy comparison of data from report to report. Typical qualifiers used to flag analytical results are listed in Appendix D.

14.7 A document control system ensures that all documents are accounted for when the project is complete. A service request number is assigned to each project for reporting and filing purposes. This number is associated with each order number (sample).

- 14.7.1 The archiving system includes all of the following items for each set of analyses performed:
- Chain-of-custody documentation
 - Benchsheets describing sample preparation
 - Sample analysis sequence
 - Analysis benchsheets and instrument printouts
 - Chromatograms and peak integration reports for all samples, standards, blanks, spikes and reruns
 - Log book ID number for the appropriate standards
 - Copies of report submitted to the client
 - Copies of Nonconformity and Corrective Action Report (NCAR) forms, if needed

- 14.7.2 Individual sets of analyses are indexed by analysis date and/or service request number. Since many analyses are performed with computer-based data systems, the final sample concentrations can be automatically calculated. If additional calculations are needed, they are written on the integration report or securely stapled to the chromatogram, if done on a separate sheet.
- 14.7.3 The archive room is an off-site file room in which files shall be maintained for a period of at least five years (from date of report issue). It is not unusual to have various clients require a 10-year retention of records, therefore, the archivist, project manager, and possibly the client are consulted prior to destruction of the records. The archive cabinet and/or off site storage area is kept locked and access keys are controlled. All documents must be signed out if needed outside of the archive room and returned in a timely manner. A designated archivist monitors filing, incoming, and outgoing data from the archive.

15.0 QUALITY CONTROL AND QUALITY ASSURANCE

- 15.1 Refer to Chapter One for specific quality control (QC) procedures. Quality control to validate sample extraction is covered in Method 3500. If extract cleanup was performed, follow the QC in Method 3600 and in the specific cleanup method.
- 15.2 System Performance Criteria - System performance criteria are presented below. The laboratory may use the recommended GC column. It must be documented that all applicable system performance criteria were met before analysis of any sample is performed. Figure 3 provides a typical 12-hour analysis sequence, whereby the response factors and mass spectrometer resolving power checks must be performed at the beginning and the end of each 12-hour period of operation. A GC column performance check is only required at the beginning of each 12-hour period during which samples are analyzed. An HRGC/HRMS method blank run is required between a calibration run and the first sample run. The same method blank extract may thus be analyzed more than once if the number of samples within a batch requires more than 12 hours of analyses.

15.2.1 GC Column Performance

- 15.2.1.1 Inject 1 μ L of the column performance check solution and acquire selected ion monitoring (SIM) data within a total cycle time of ≤ 1 second.
- 15.2.1.2 The chromatographic separation between 2,3,7,8-TCDD and the peaks representing any other unlabeled TCDD isomers must be resolved with a valley of ≤ 25 percent, where:

$$\text{Valley percent} = \frac{(X)}{(Y)} \times 100$$

x = measured from the 2,3,7,8-closest TCDD eluting isomer, and

y = the peak height of 2,3,7,8-TCDD

NOTE: It is the responsibility of the laboratory to verify the conditions suitable for the appropriate resolution of 2,3,7,8-TCDD from all other TCDD isomers. The GC column performance check solution also contains the known first and last PCDD/PCDF eluters under the conditions specified in this protocol. Their retention times are used to determine the eight homologue retention time windows that are used for qualitative and quantitative purposes. All peaks (that includes $^{13}\text{C}_{12}$ -2,3,7,8-TCDD) should be labeled and identified on the chromatograms. Furthermore, all first eluters of a homologous series should be labeled with the letter F, and all last eluters of a homologous series should be labeled with the letter L.

Any individual selected ion current profile (SICP) (for the tetras, this would be the SICP for m/z 322 and m/z 304) or the reconstructed homologue ion current (for the tetras, this would correspond to m/z 320 + m/z 322 + m/z 304 + m/z 306) constitutes an acceptable form of data presentation. An SICP for the labeled compounds (e.g., m/z 334 for labeled TCDD) is also required.

- 15.2.1.3 The retention times for the switching of SIM ions characteristic of one homologous series to the next higher homologous series must be indicated in the SICP. (see Figure 5) Accurate switching at the appropriate times is absolutely necessary for accurate monitoring of these compounds. Allowable tolerance on the daily verification with the GC performance check solution should be better than 10 seconds for the absolute retention times of all the components of the mixture. Particular caution should be exercised for the switching time between the last tetrachlorinated congener (i.e., 1,2,8,9-TCDD) and the first pentachlorinated congener (i.e., 1,3,4,6,8-PeCDF), as these two compounds elute within 15 seconds of each other on the 60m DB-5 column. A laboratory with a GC/MS system that is not capable of detecting both congeners (1,2,8,9-TCDD and 1,3,4,6,8-PeCDF) within one analysis must take corrective action. If the recommended column is not used,

then the first and last eluting isomer of each homologue must be determined experimentally on the column which is used, and the appropriate isomers must then be used for window definition and switching times.

15.2.2 Mass Spectrometer Performance

- 15.2.2.1 Chromatography time for PCDDs and PCDFs exceeds the long term mass stability of the mass spectrometer. Because the instrument is operated in the high-resolution mode, mass drifts of a few ppm (e.g., 5 ppm in mass) can have serious adverse effects on instrument performance. Therefore, a mass drift correction is mandatory. To that effect, it is recommended to select a lock-mass ion from the reference compound (PFK is recommended) used for tuning the mass spectrometer. The selection of the lock-mass ion is dependent on the masses of the ions monitored within each descriptor. However, an acceptable lock-mass ion at any mass between the lightest and heaviest ion in each descriptor can be used to monitor and correct mass drifts. The level of the reference compound (PFK) metered into the ion chamber during HRGC/HRMS analyses should be adjusted so that the amplitude of the most intense selected lock-mass ion signal (regardless of the descriptor number) does not exceed 10 percent of the full scale deflection for a given set of detector parameters. Under those conditions, sensitivity changes that might occur during the analysis can be more effectively monitored.

NOTE: Excessive PFK (or any other reference substance) may cause noise problems and contamination of the ion source resulting in an increase in downtime for source cleaning.

- 15.2.2.1 Documentation of the instrument resolving power must then be accomplished by recording the peak profile of the high-mass reference signal (m/z 380.9760) obtained during the above peak matching experiment by using the low-mass PFK ion at m/z 304.9824 as a reference. The minimum resolving power of 10,000 must be demonstrated on the high-mass ion while it is transmitted at a lower accelerating voltage than the low-mass reference ion, which is transmitted at full sensitivity. The format of the peak profile representation must allow manual determination of the resolution, i.e., the horizontal axis must be a calibrated mass scale (amu or ppm per division). The result of the peak width measurement (performed at 5 percent of the maximum, which

corresponds to the 10 percent valley definition) must appear on the hard copy and cannot exceed 100ppm at m/z 380.9760 (or 0.038amu at that particular mass).

15.3 Quality Control Samples

15.3.1 Performance Evaluation Samples - Included among the samples in all batches may be samples (blind or double blind) containing known amounts of unlabeled 2,3,7,8-substituted PCDDs/PCDFs or other PCDD/PCDF congeners.

15.3.2 Performance Check Solutions

- 15.3.2.1 At the beginning of each 12-hour period during which samples are to be analyzed, an aliquot of the 1) GC column performance check solution and 2) high-resolution concentration. calibration solution, HRCC-3, shall be analyzed to demonstrate adequate GC resolution and sensitivity, response factor reproducibility, and mass range calibration, and to establish the PCDD/PCDF retention time windows. A mass resolution check shall also be performed to demonstrate adequate mass resolution using an appropriate reference compound (PFK is recommended). If the required criteria are not met, remedial action must be taken before any samples are analyzed.
- 15.3.2.2 To validate positive sample data, the routine or continuing calibration and the mass resolution check must be performed also at the end of each 12-hour period during which samples are analyzed. Furthermore, an HRGC/HRMS method blank run must be recorded following a calibration run and the first sample run.
- 15.3.2.3 If the laboratory operates only during one period (shift) each day of 12 hours or less, the GC performance check solution must be analyzed only once (at the beginning of the period) to validate the data acquired during the period. However, the mass resolution and continuing calibration checks must be performed at the beginning as well as at the end of the period.
- 15.3.2.4 If the laboratory operates during consecutive 12-hour periods (shifts), analysis of the GC performance check solution must be performed at the beginning of each 12-hour period. The mass

resolution and continuing calibration checks from the previous period can be used for the beginning of the next period.

- 15.3.2.5 Results of at least one analysis of the GC column performance check solution and of two mass resolution and continuing calibration checks must be reported with the sample data collected during a 12-hour period.
 - 15.3.2.6 Deviations from criteria specified for the GC performance check or for the mass resolution check invalidate all positive sample data collected between analyses of the performance check solution, and the extracts from those positive samples shall be reanalyzed.
 - 15.3.2.7 If the continuing calibration check performed at the end of a 12-hour period fails by no more than 25 percent RPD for the 17 unlabeled compounds, and 35 percent RPD for the 9 labeled reference compounds, use the mean RFs from the two daily routine calibration runs to compute the analyte concentrations, instead of the RFs obtained from the initial calibration.
 - 15.3.2.8 A new initial calibration (new RFs) is required immediately (within two hours) following the analysis of the samples, whenever the RPD from the end-of-shift routine calibration exceeds 25 percent or 35 percent, respectively. Failure to perform a new initial calibration immediately following the analysis of the samples will automatically require reanalysis of all positive sample extracts analyzed before the failed end-of-shift continuing calibration check.
- 15.3.3 The GC column performance check mixture, high-resolution concentration calibration solutions, and the sample fortification solutions may be obtained from the EMSL-CIN. However, if not available from the EMSL-CIN, standards can be obtained from other sources, and solutions can be prepared in the laboratory. Concentrations of all solutions containing 2,3,7,8-substituted PCDDs/PCDFs, which are not obtained from the EMSL-CIN, must be verified by comparison with the EPA standard solutions that are available from the EMSL-CIN.
- 15.3.4 Field Blanks - Each batch of samples usually contains a field blank sample of uncontaminated soil, sediment or water that is to be fortified before analysis. In addition to this field blank, a batch of samples may include a rinsate, which is a portion of the solvent (usually trichloroethylene) that was used to rinse sampling equipment. The rinsate is analyzed to assure the sampling equipment has not contaminated any samples.

15.3.5 Fortified Field Blank

- 15.3.5.1 Weigh a 10g portion or use 1L (for aqueous samples) of the specified field blank sample and add 100 μ L of the solution containing the nine internal standards diluted.
- 15.3.5.2 Add 20 μ L of the recovery standard solution and analyze a 1 μ L aliquot of the concentrated extract.
 - 15.3.5.2.1 Calculate the concentration of 2,3,7,8-substituted PCDDs/PCDFs and the percent recovery of the internal standards.
 - 15.3.5.2.2 Extract and analyze a new simulated fortified field blank whenever new lots of solvents or reagents are used for sample extraction or for column chromatographic procedures.

15.3.6 Rinsate Sample

- 15.3.6.1 The rinsate sample must be fortified like a regular sample.
- 15.3.6.2 Take a 100mL (+ 0.5mL) portion of the sampling equipment rinse solvent (rinsate sample), filter, if necessary, and add 100 μ L of the solution containing the nine internal standards.
- 15.3.6.3 Using a KD apparatus, concentrate to approximately 5mL.

NOTE: As an option, a rotary evaporator may be used in place of the KD apparatus for the concentration of the rinsate.
- 15.3.6.4 Transfer the 5mL concentrate from the KD concentrator tube in 1mL portions to a 1mL minivial, reducing the volume in the minivial as necessary with a gentle stream of dry nitrogen.
- 15.3.6.5 Rinse the KD concentrator tube with two 0.5mL portions of hexane and transfer the rinses to the 1mL minivial. Blow down with dry nitrogen as necessary.
- 15.3.6.6 Just before analysis, add 20 μ L recovery standard solution and reduce the volume to its final volume. No column chromatography is required.

- 15.3.6.7 Analyze an aliquot following the same procedures used to analyze samples.
 - 15.3.6.8 Report percent recovery of the internal standard and the presence of any PCDD/PCDF compounds in µg/L of rinsate solvent.
- 15.3.7 Duplicate Analyses, if specified by the project instructions:
- 15.3.7.1 In each batch of samples, locate the sample specified for duplicate analysis, and analyze a second 10g soil or sediment sample portion or 1L water sample, or an appropriate amount of the type of matrix under consideration.
 - 15.3.7.2 The results of the laboratory duplicates (percent recovery and concentrations of 2,3,7,8-substituted PCDD/PCDF compounds) should agree within 25 percent relative difference (difference expressed as percentage of the mean). Report all results.
 - 15.3.7.3 Recommended actions to help locate problems:
 - 15.3.7.4 Verify satisfactory instrument performance.
 - 15.3.7.5 If possible, verify that no error was made while weighing the sample portions.
 - 15.3.7.6 Review the analytical procedures with the performing laboratory personnel.
- 15.3.8 Matrix Spike and Matrix Spike Duplicate
- 15.3.8.1 Locate the sample for the MS and MSD analyses (the sample may be labeled "double volume").
 - 15.3.8.2 Add an appropriate volume of the matrix spike fortification solution and of the sample fortification solution, adjusting the fortification level.
 - 15.3.8.3 Analyze the MS and MSD samples.

15.3.8.4 The results obtained from the MS and MSD samples (concentrations of 2,3,7,8-substituted PCDDs/PCDFs) should agree within 20 percent relative difference.

15.3.9 Percent Recovery of the Internal Standards

15.3.9.1 For each sample, method blank and rinsate, calculate the percent recovery. The percent recovery should be between 40 percent and 135 percent for all 2,3,7,8-substituted internal standards.

NOTE: A low or high percent recovery for a blank does not require discarding the analytical data but it may indicate a potential problem with future analytical data.

15.4 Identification Criteria

15.4.1 If either one of the identification criteria is not met for a homologous series, it is reported that the sample does not contain unlabeled 2,3,7,8-substituted PCDD/PCDF isomers for that homologous series at the calculated detection limit.

15.4.2 If the first initial identification criteria are met, but the criteria appearing are not met, that sample is presumed to contain interfering contaminants. This must be noted on the analytical report form, and the sample should be rerun or the extract reanalyzed.

15.5 Unused portions of samples and sample extracts should be preserved for 90 days after sample receipt to allow further analyses.

15.6 Reuse of glassware is to be minimized to avoid the risk of contamination.

16.0 REFERENCES

1. "Control of Interferences in the Analysis of Human Adipose Tissue for 2,3,7,8-Tetrachlorodibenzo-p-dioxin". D. G. Patterson, J.S. Holler, D.F. Grote, L.R. Alexander, C.R. Lapeza, R.C. O'Connor and J.A. Liddle. Environ. Toxicol. Chem. 5, 355-360 (1986).
2. "Method 8290: Analytical Procedures and Quality Assurance for Multimedia Analysis of Polychlorinated Dibenzo-p-Dioxins and Dibenzofurans by High-Resolution Gas Chromatography/High-Resolution Mass Spectrometry". Y. Tondeur and W.F. Beckert. U.S. Environmental Protection Agency, Environmental Monitoring Systems Laboratory, Las Vegas, NV.
3. "Carcinogens - Working with Carcinogens", Department of Health, Education, and Welfare, Public Health Service, Center for Disease Control. National Institute for Occupational Safety and Health. Publication No. 77-206, August 1977.
4. "OSHA Safety and Health Standards, General Industry", (29 CFR 1910), Occupational Safety and Health Administration, OSHA 2206 (revised January 1976).
5. "Safety in Academic Chemistry Laboratories", American Chemical Society Publication, Committee on Chemical Safety (3rd Edition, 1979.)
6. "Hybrid HRGC/MS/MS Method for the Characterization of Tetrachlorinated Dibenzodioxins in Environmental Samples." Y. Tondeur, W.J. Niederhut, S.R. Missler, and J.E. Campana, Mass Spectrom. 14, 449-456 (1987).
6. USEPA National Dioxin Study - Phase II, "Analytical Procedures and Quality Assurance Plan for the Determination of PCDD/PCDF in Fish", EPA-Duluth, October 26, 1987.

TABLE 1
MATRICES, SAMPLE SIZES, AND METHOD CALIBRATION LIMITS (PARTS PER TRILLIONS)

	Water	Soil Sediment Paper pulp^b	Fly ash	Fish tissue^c	Human adipose tissue	Sludges Fuel oil	Still- bottom
Lower MCL^a	0.01	1.0	1.0	1.0	1.0	5.0	10
Upper MCL^a	2	200	200	200	200	1000	2000
Weight (g)	1000	10	10	10	2	1	1
IS spiking levels (ppt)	1	100	100	100	100	500	1000
Final extr. vol. (μL)	10-50	10-50	50	10-50	10-50	50	50

^a For other congeners, multiply the values by 1 for TCDF/TCDD, by 2.5 for PeCDD/PeCDF, HxCDD/HxCDF/HpCDD/HpCDF, and by 5 for OCDD/OCDF.

^b Sample dewatered prior to analysis.

^c An additional 10g sample is used for determination of lipid content.

TABLE 2

THE FIFTEEN 2,3,7,8-SUBSTITUTE PCDD AND PCDF CONGENERS

PCDD	PCDF
2,3,7,8-TCDD(*)	2,3,7,8-TCDF(*)
1,2,3,7,8-PeCDD(*)	1,2,3,7,8-PeCDF(*)
1,2,3,6,7,8-HxCDD(*)	2,3,4,7,8-PeCDF
1,2,3,4,7,8-HxCDD	1,2,3,6,7,8-HxCDF
1,2,3,7,8,9-HxCDD(+)	1,2,3,7,8,9-HxCDF
1,2,3,4,6,7,8-HpCDD(*)	1,2,3,4,7,8-HxCDF(*)
	2,3,4,6,7,8-HxCDF
	1,2,3,4,6,7,8-HpCDF(*)
	1,2,3,4,7,8,9-HpCDF

* The ¹³C-labeled analogue is used as an internal standard

+ The ¹³C-labeled analogue is used as recovery standard

TABLE 3
ISOMERS OF CHLORINATED DIOXIN AND FURANS AS A FUNCTION OF THE NUMBER OF CHLORINE ATOMS

Number of Chlorine Atoms	Number of Dioxin Isomers	Number of 2378-Dioxin	Number of Furan Isomers	Number of 2378-furans
1	2	-	4	-
2	10	-	16	-
3	14	-	28	-
4	22	1	38	1
5	14	1	28	2
6	10	3	16	4
7	2	1	4	2
8	1	1	1	1
Total	75	7	135	10

TABLE 4
HIGH RESOLUTION CONCENTRATION CALIBRATION SOLUTIONS

Unlabeled Analytes	Concentrations (pg/uL)				
	1	2	3	4	5
2,3,7,8-TCDD	1	2.5	10	50	200
2,3,7,8-TCDF	1	2.5	10	50	200
1,2,3,7,8-PeCDD	2.5	6.25	25	125	500
1,2,3,7,8-PeCDF	2.5	6.25	25	125	500
2,3,4,7,8-PeCDF	2.5	6.25	25	125	500
1,2,3,4,7,8-HxCDD	2.5	6.25	25	125	500
1,2,3,6,7,8-HxCDD	2.5	6.25	25	125	500
1,2,3,7,8,9-HxCDD	2.5	6.25	25	125	500
1,2,3,4,7,8-HxCDF	2.5	6.25	25	125	500
1,2,3,6,7,8-HxCDF	2.5	6.25	25	125	500
1,2,3,7,8,9-HxCDF	2.5	6.25	25	125	500
2,3,4,6,7,8-HxCDF	2.5	6.25	25	125	500
1,2,3,4,6,7,8-HpCDD	2.5	6.25	25	125	500
1,2,3,4,6,7,8-HpCDF	2.5	6.25	25	125	500
1,2,3,4,7,8,9-HpCDF	2.5	6.25	25	125	500
OCDD	5	12.5	50	250	1000
OCDF	5	12.5	50	250	1000
Internal Standards					
13C12 – 2,3,7,8 – TCDD	50	50	50	50	50
13C12 – 2,3,7,8 – TCDF	50	50	50	50	50
13C12 – 1,2,3,7,8 – PeCDD	50	50	50	50	50
13C12 – 1,2,3,7,8 – PeCDF	50	50	50	50	50
13C12 – 1,2,3,6,7,8 – HxCDD	125	125	125	125	125
13C12 – 1,2,3,4,7,8 – HxCDF	125	125	125	125	125
13C12 – 1,2,3,4,6,7,8 – HpCDD	125	125	125	125	125
13C12 – 1,2,3,4,6,7,8 – HpCDF	125	125	125	125	125
13C12 – OCDD	250	250	250	250	250
Recovery Standards					
13C12 – 1,2,3,4 – TCDD	50	50	50	50	50
13C12 – 1,2,3,7,8,9 – HxCDD	125	125	125	125	125

TABLE 5
IONS MONITORED FOR HRGC/HRMS ANALYSIS OF PCDDs/PCDFs

DESCRIPTOR	EXACT M/Z	M/Z TYPE	ELEMENT COMPOSITION	COMPOUND
1	292.9825	LOCK	C_7F_{11}	PFK
	303.9016	M	$C_{12}H_4^{35}Cl_4O$	TCDF
	305.8987	M+2	$C_{12}H_4^{35}Cl^{37}O$	TCDF
	315.9419	M	$^{13}C_{12}H_4^{35}Cl_4O$	TCDF(S)
	317.9389	M+2	$^{13}C_{12}H_4^{35}Cl_3^{37}ClO$	TCDF(S)
	319.8965	M	$C_{12}H_4^{35}ClO_2$	TCDD
	321.8936	M+2	$C_{12}H_4^{35}Cl_3^{37}ClO_2$	TCDD
	327.8847	M	$C_{12}H_4^{37}Cl_4O_2$	TCDD(CS)
	330.9792	QC	C_7F_{13}	PFK
	331.9368	M	$^{13}C_{12}H_4^{35}Cl_4O_2$	TCDD(S)
	333.9339	M+2	$^{13}C_{12}H_4^{35}Cl^{37}ClO_2$	TCDD(S)
	375.8364	M+2	$C_{12}H_4^{35}Cl_5^{37}ClO$	HXCDPE
2	339.8597	M+2	$C_{12}H_3^{35}Cl_4^{37}ClO$	PECDF
	341.8567	M+4	$C_{12}H_3^{35}Cl_3^{37}Cl_2O$	PECDF
	351.9000	M+2	$^{13}C_{12}H_3^{35}Cl_4^{37}ClO$	PECDF(S)
	353.8970	M+4	$^{13}C_{12}H_3^{35}Cl_3^{37}Cl_2O$	PECDF(S)
	355.8546	M+2	$C_{12}H_3^{35}Cl_3^{37}ClO_2$	PECDD
	357.8516	M+4	$C_{12}H_3^{35}Cl_3^{37}Cl_2O_2$	PECDD
	367.8949	M+2	$^{13}C_{12}H_3^{35}Cl_4^{37}ClO_2$	PECDD(S)
	369.8919	M+4	$^{13}C_{12}H_3^{35}Cl_3^{37}Cl_2O_2$	PECDD(S)
	409.7974	M+2	$C_{12}H_3^{35}Cl_6^{37}ClO$	HPCDPE
3	373.8208	M+2	$C_{12}H_2^{35}Cl_5^{37}ClO$	HXCDF
	375.8178	M+4	$C_{12}H_2^{35}Cl_4^{37}Cl_2O$	HXCDF
	383.8639	M	$^{13}C_{12}H_2^{35}Cl_6O$	HXCDF(S)
	385.8610	M+2	$^{13}C_{12}H_2^{35}Cl_5^{37}ClO$	HXCDF(S)

TABLE 5 (con't)

DESCRIPTOR	EXACT M/Z	M/Z TYPE	ELEMENT COMPOSITION	COMPOUND
3	389.8157	M+2	$C_{12}H_2^{35}Cl_5^{37}ClO_2$	HXCDD
	391.8127	M+4	$C_{12}H_2^{35}Cl_4^{37}Cl_2O_2$	HXCDD
	392.9760	LOCK	C_9F_{15}	PFK
	401.8559	M+2	$^{13}C_{12}H_2^{35}Cl_5^{37}ClO_2$	HXCDD(S)
	403.8529	M+4	$^{13}C_{12}H_2^{35}Cl_4^{37}C_2IO$	HXCDD(S)
	430.9729	QC	C_9F_{17}	PFK
	445.7555	M+4	$C_{12}H_2^{35}Cl_5^{37}Cl_2O$	OCDPE
	407.7818	M+2	$C_{12}H^{35}Cl_6^{37}ClO$	HPCDF
4	409.7789	M+4	$C_{12}H^{35}Cl_5^{37}Cl_2O$	HPCDF
	417.8253	M	$^{13}C_{12}H^{35}Cl_7^{37}O$	HPCDF(S)
	419.8220	M+2	$^{13}C_{12}H^{35}Cl_6^{37}ClO$	HPCDF(S)
	423.7766	M+2	$C_{12}H^{35}Cl_6^{37}ClO_2$	HPCDD
	425.7737	M+4	$C_{12}H^{35}Cl_5^{37}Cl_2O_2$	HPCDD
	430.9729	LOCK	C_9F_{17}	PFK
	435.8169	M+2	$^{13}C_{12}H^{35}Cl_6^{37}ClO_2$	HPCDD(S)
	437.8140	M+4	$^{13}C_{12}H^{35}Cl_5^{37}Cl_2O_2$	HPCDD(S)
	479.7165	M+4	$C_{12}H^{35}Cl_7^{37}Cl_2O$	NCDPE
	441.728	M+2	$C_{12}H^{35}Cl_7^{37}ClO$	OCDF
	443.7399	M+4	$C_{12}H^{35}Cl_6^{37}Cl_2O$	OCDF
5	457.7377	M+2	$C_{12}H^{35}Cl_7^{37}ClO_2$	OCDD
	459.7348	M+4	$C_{12}H^{35}Cl_6^{37}Cl_2O_2$	OCDD
	469.7779	M+2	$^{13}C_{12}H^{35}Cl_7^{37}ClO_2$	OCDD(IS)
	471.7750	M+4	$^{13}C_{12}H^{35}Cl_6^{37}Cl_2O_2$	OCDD(IS)
	513.6775	M+4	$C_{12}H^{35}Cl_8^{37}Cl_2O_2$	DCDPE
	442.9728	QC	$C_{10}F_{17}$	PFK

Nuclidic masses used:

H= 1.007825

C= 12.000000

^{13}C = 13.003355

F= 18.9984

O= 15.994915

^{35}Cl = 34.968853

^{37}Cl = 36.965903

PFK= PERFLUOROKEROSENE

SS= INTERNAL STANDARD

CS=CLEANUP STANDARD (only one m/z)

TABLE 6
PCDD AND PCDF CONGENERS PRESENT IN THE GC
PERFORMANCE EVALUATION SOLUTION AND USED FOR DEFINING THE H GC HOMOLOGUE
RETENTION TIME WINDOWS ON A 60-M DB-5 COLUMN

# Chlorine Atoms	PCDD Positional Isomer		PCDF Positional Isomers	
	First Eluters	Last Eluters	First Eluters	Last Eluters
4 ^a	1,3,6,8	1,2,8,9	1,3,6,8	1,2,8,9
5	1,2,4,6,8/1,2,4,7,9	1,2,3,8,9	1,3,4,6,8	1,2,3,8,9
6	1,2,4,6,7,9/1,2,4,6,8,9	1,2,3,4,6,7	1,2,3,4,6,8	1,2,3,4,8,9
7	1,2,3,4,6,7,9	1,2,3,4,6,7,8	1,2,3,4,6,7,8	1,2,3,4,7,8,9
8	1,2,3,4,6,7,8,9		1,2,3,4,6,7,8,9	

^a In addition to these two TCDD isomers, the 1,2,3,4,-, 1,2,3,7-, 1,2,3,8-, 2,3,7,8-, ¹³C₁₂ – 2,3,7,8-, and 1,2,3,9-TCDD isomers must be present as a check of column resolution

TABLE 7
THEORICAL ION ABUNDANCE RATIOS AND THEIR CONTROL LIMITS FOR
PCDD AND PCDF

# Chlorine Atoms	Ion Type	Theoretical Abundance Ratio	Control Limits	
			Lower	Upper
4	M/M+2	0.77	0.65	0.89
5	M+2/M+4	1.55	1.32	1.78
6	M+2/M+4	1.24	1.05	1.43
6 ^a	M/M+2	0.51	0.43	0.59
7 ^b	M/M+2	0.44	0.37	0.51
7	M+2/M+4	1.04	0.88	1.20
8	M+2/M+4	0.89	0.76	1.02

^a Used only for ¹³C₁₂ – HxCDF (IS)

^b Used only for ¹³C₁₂ – HpCDF (IS)

TABLE 8
2,3,7,8-TCDD TOXICITY EQUIVALENCY FACTORS (TEFs) FOR THE POLYCHLORINATED
DIBENZODIOXINS AND DIBENZOFURANS

Analyte	TEF ^a
2,3,7,8-TCDD	1.00
1,2,3,7,8-PeCDD	0.50
1,2,3,6,7,8-HxCDD	0.10
1,2,3,7,8,9-HxCDD	0.10
1,2,3,4,7,8-HxCDD	0.10
1,2,3,4,6,7,8-HpCDD	0.01
1,2,3,4,6,7,8,9-OCDD	0.001
2,3,7,8-TCDF	0.1
1,2,3,7,8-PeCDF	0.05
2,3,4,7,8-PeCDF	0.5
1,2,3,6,7,8-HxCDF	0.1
1,2,3,7,8,9-HxCDF	0.1
1,2,3,4,7,8-HxCDF	0.1
2,3,4,6,7,8-HxCDF	0.1
1,2,3,4,6,7,8-HpCDF	0.01
1,2,3,4,7,8,9-HpCDF	0.01
1,2,3,4,6,7,8,9-OCDF	0.001

^a Taken from “Interim Procedure for Estimating Risk Associated with Exposures to Mixtures of Chlorinated Dibenzo-*p*-Dioxin and –Dibenzofurans 1989 Update”. (EPA/625/3-89/016, March 1989).

TABLE 9
CONCENTRATION OF STOCK AND SPIKING SOLUTION CONTAINING PCDDs/PCDFs

PCDDs/PCDFs	Internal Standard Stock Solution (ng/mL)	Internal Standard Spiking Solution (ng/mL)	MS Stock Fortification Solution (ng/mL)	MS Spiking Fortification Solution (ng/mL)
2,3,7,8-TCDD			100	2
2,3,7,8-TCDF			100	2
1,2,3,7,8-PeCDD			250	5
1,2,3,7,8-PeCDF			250	5
2,3,4,7,8-PeCDF			250	5
1,2,3,4,7,8-HxCDD			250	5
1,2,3,6,7,8-HxCDD			250	5
1,2,3,7,8,9-HxCDD			250	5
1,2,3,4,7,8-HxCDF			250	5
1,2,3,6,7,8-HxCDF			250	5
1,2,3,7,8,9-HxCDF			250	5
2,3,4,6,7,8-HxCDF			250	5
1,2,3,4,6,7,8-HpCDD			250	5
1,2,3,4,6,7,8-HpCDF			250	5
1,2,3,4,7,8,9-HpCDF			250	5
OCDD			500	10
OCDF			500	10
Internal Standard				
$^{13}\text{C}_{12}$ – 2,3,7,8-TCDD	100	10		
$^{13}\text{C}_{12}$ – 2,3,7,8-TCDF	100	10		
$^{13}\text{C}_{12}$ – 1,2,3,7,8-PeCDD	100	10		
$^{13}\text{C}_{12}$ – 1,2,3,7,8-PeCDF	100	10		
$^{13}\text{C}_{12}$ – 1,2,3,6,7,8-HxCDD	250	25		
$^{13}\text{C}_{12}$ – 1,2,3,4,7,8-HxCDF	250	25		
$^{13}\text{C}_{12}$ – 1,2,3,4,6,7,8-HpCDD	250	25		
$^{13}\text{C}_{12}$ – 1,2,3,4,6,7,8-HpCDF	250	25		
$^{13}\text{C}_{12}$ – OCDD	500	50		
Cleanup Standard				
$^{37}\text{Cl}_4$ – 2,3,7,8-TCDD		8.0		
Recovery Standard				
$^{13}\text{C}_{12}$ – 1,2,3,4-TCDD		50		
$^{13}\text{C}_{12}$ – 1,2,3,7,8,9-HxCDD		50		

TABLE 10
REFERENCE COMPOUNDS FOR QUANTITATION OF NATIVE AND LABELED
PCDDs AND PCDF

Compound Number	Type	Name	Reference Compound
1	Native	2,3,7,8-TCDD	18
2	Native	1,2,3,7,8-PeCDD	19
3	Native	1,2,3,4,7,8-HxCDD	20
4	Native	1,2,3,4,6,7,8-HxCDD	20
5	Native	1,2,3,7,8,9-HxCDD	20
6	Native	1,2,3,4,6,7,8-HpCDD	21
7	Native	OCDD	22
8	Native	2,3,7,8-TCDF	23
9	Native	1,2,3,7,8-PeCDF	24
10	Native	2,3,4,7,8-PeCDF	24
11	Native	1,2,3,4,7,8-HxCDF	25
12	Native	1,2,3,6,7,8-HxCDF	25
13	Native	1,2,3,7,8,9-HxCDF	25
14	Native	2,3,4,6,7,8-HxCDF	25
15	Native	1,2,3,4,6,7,8-HpCDF	26
16	Native	1,2,3,4,7,8,9-HpCDF	26
17	Native	OCDF	22
18	Internal Standard	¹³ C – 2,3,7,8-TCDD	27
19	Internal Standard	¹³ C – 1,2,3,7,8-PeCDD	27
20	Internal Standard	¹³ C – 1,2,3,6,7,8-HxCDD	28
21	Internal Standard	¹³ C – 1,2,3,6,7,8-HpCDD	28
22	Internal Standard	¹³ C – OCDD	28
23	Internal Standard	¹³ C – 2,3,7,8-TCDF	27
24	Internal Standard	¹³ C – 1,2,3,7,8-PeCDF	27
25	Internal Standard	¹³ C – 1,2,3,4,7,8-HxCDF	28
26	Internal Standard	¹³ C – 1,2,3,4,6,7,8-HpCDF	28
27	Recovery Standard	¹³ C – 1,2,3,4-TCDD	
28	Recovery Standard	¹³ C – 1,2,3,7,8,9-HxCDD	
29	Cleanup Standard	³⁷ Cl – 1,2,3,7,8-TCDD	27

TABLE 11
RELATIVE RESPONSE FACTOR [RF (NUMBER)] ATTRIBUTIONS

Number	Specific Congener Name
1	2,3,7,8-TCDD (and total TCDDs)
2	2,3,7,8-TCDF (and total TCDFs)
3	1,2,3,7,8-PeCDD (and total PeCDDs)
4	1,2,3,7,8-PeCDF
5	2,3,4,7,8-PeCDF
6	1,2,3,4,7,8-HxCDD
7	1,2,3,6,7,8-HxCDD
8	1,2,3,7,8,9-HxCDD
9	1,2,3,4,7,8-HxCDF
10	1,2,3,6,7,8-HxCDF
11	1,2,3,7,8,9-HxCDF
12	2,3,4,6,7,8-HxCDF
13	1,2,3,4,6,7,8-HpCDD (and totals HpCDD)
14	1,2,3,4,6,7,8-HpCDF
15	1,2,3,4,7,8,9-HpCDF
16	OCDD
17	OCDF
18	$^{13}C_{12}$ – 2,3,7,8-TCDD
19	$^{13}C_{12}$ – 2,3,7,8-TCDF
20	$^{13}C_{12}$ – 1,2,3,7,8-PeCDD
21	$^{13}C_{12}$ – 1,2,3,7,8-PeCDF
22	$^{13}C_{12}$ – 1,2,3,6,7,8-HxCDD
23	$^{13}C_{12}$ – 1,2,3,4,7,8-HxCDF
24	$^{13}C_{12}$ – 1,2,3,4,6,7,8-HpCDD
25	$^{13}C_{12}$ – 1,2,3,4,6,7,8-HpDf
26	$^{13}C_{12}$ – OCDD
27	Total PeCDFs
28	Total HxCDFs
29	Total HxCDDs
30	Total HpCDFs

FIGURE 1

GENERAL STRUCTURES OF DIBENZO-*p*-DIOXINS (TOP) AND DIBENZOFURANS (BOTTOM)

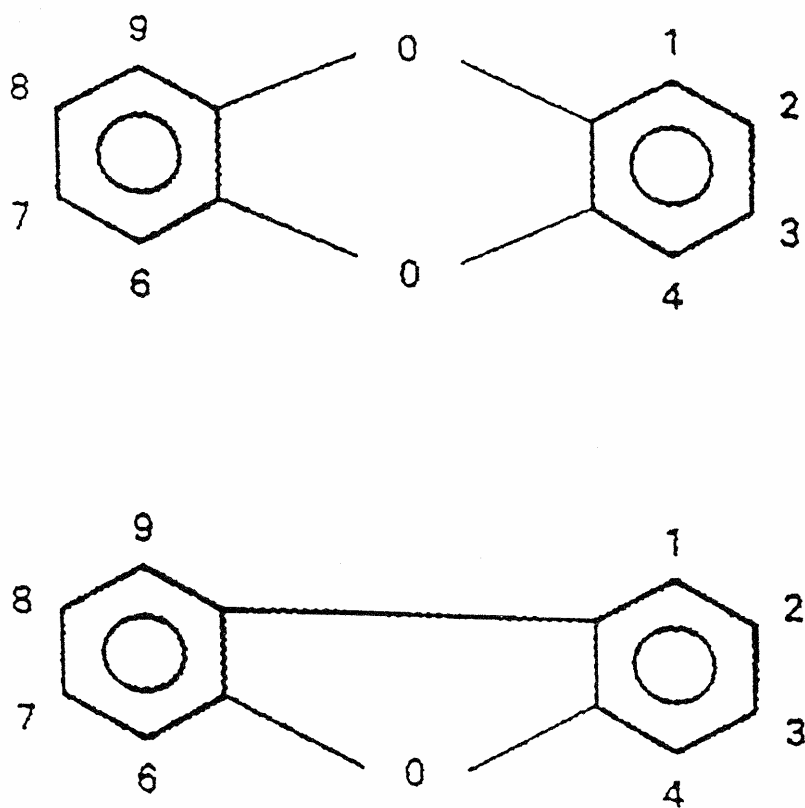


FIGURE 2
TYPICAL 12-HOUR ANALYSIS SEQUENCE OF EVENTS
ANALYTICAL PROCEDURE


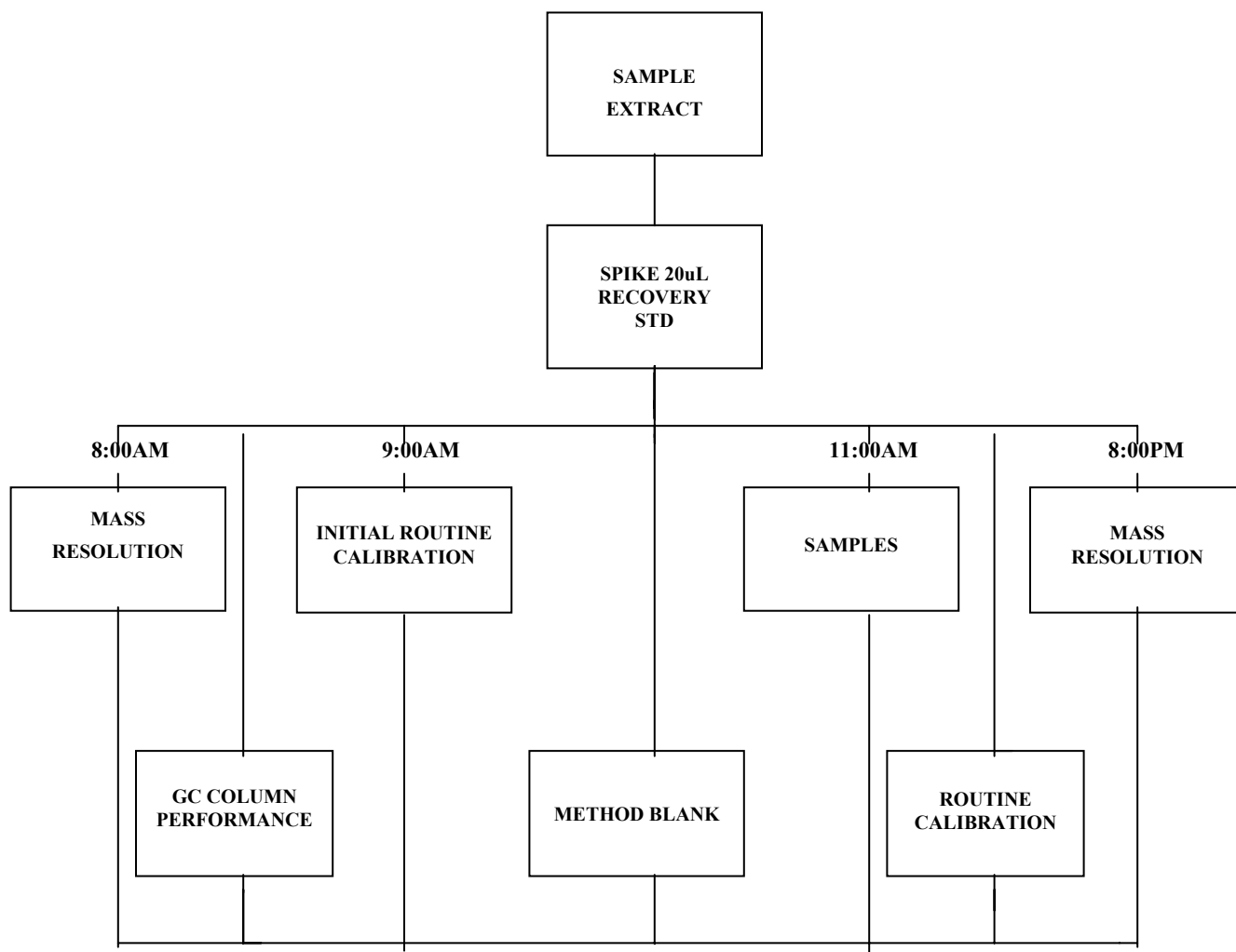


FIGURE 3
 MANUAL DETERMINATION OF S/N

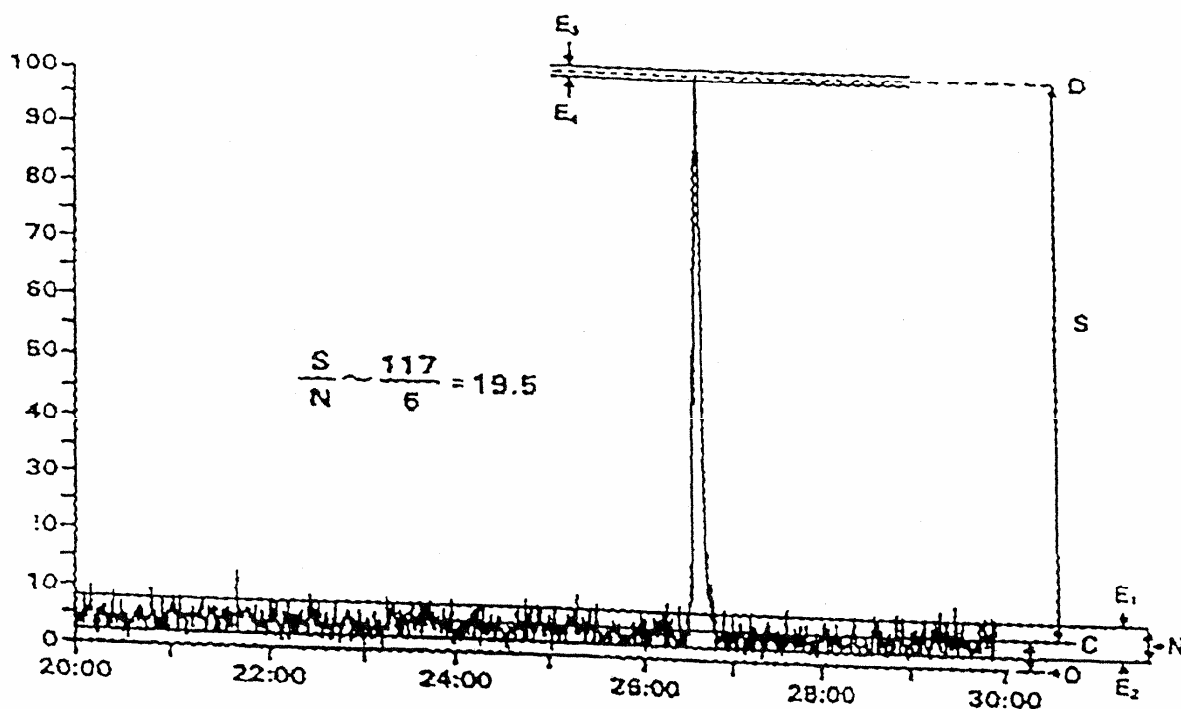

The peak height (S) is measured between the mean noise (lines C and D.) These mean signal values are obtained by tracing the line between the baseline average noise extremes, E1 and E2, and between the apex average noise extremes, E3 and E4, at the apex of the signal.

NOTE: It is imperative that the instrument interface amplifier electronic zero offset be set high enough so that negative going baseline noise is recorded.

FIGURE 4
MASS RESOLUTION CHECK

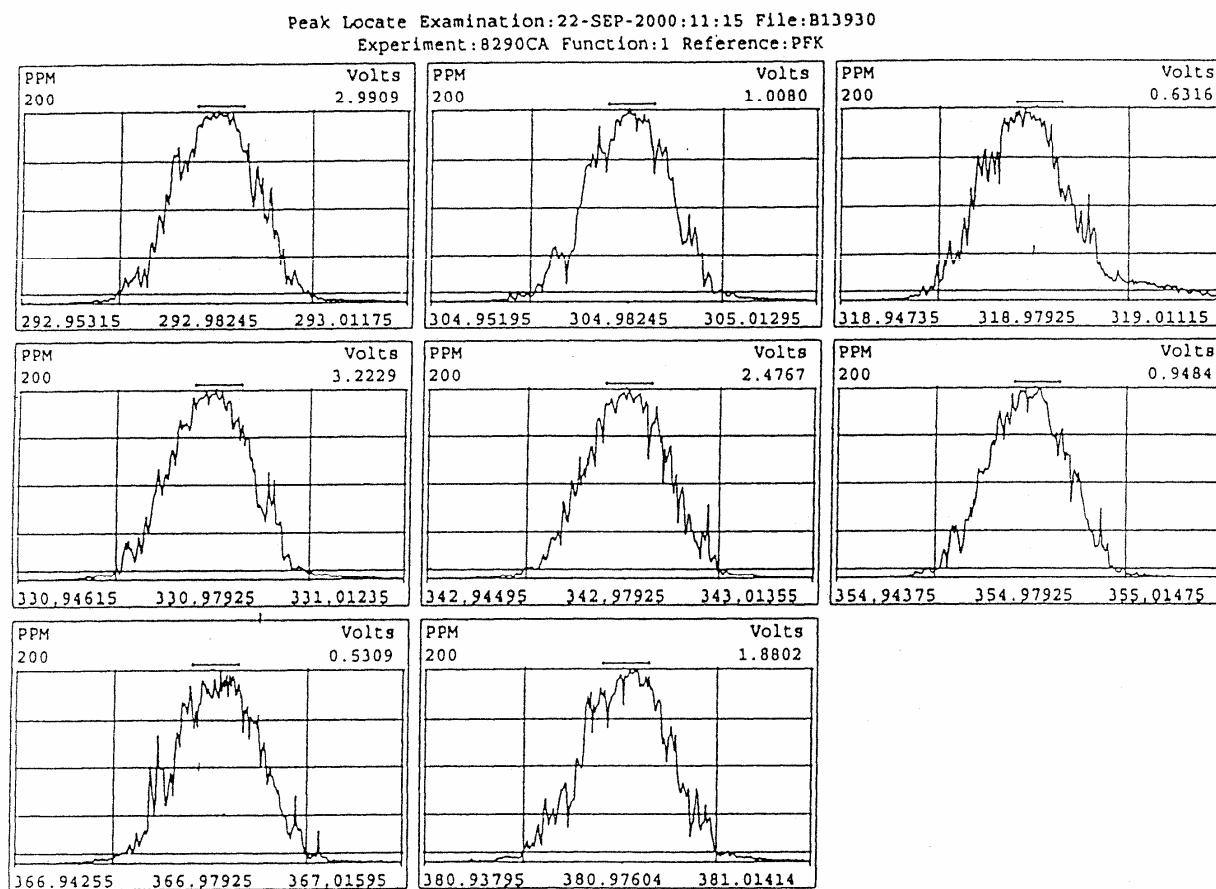


FIGURE 5
DB-5 WINDOW DEFINITION

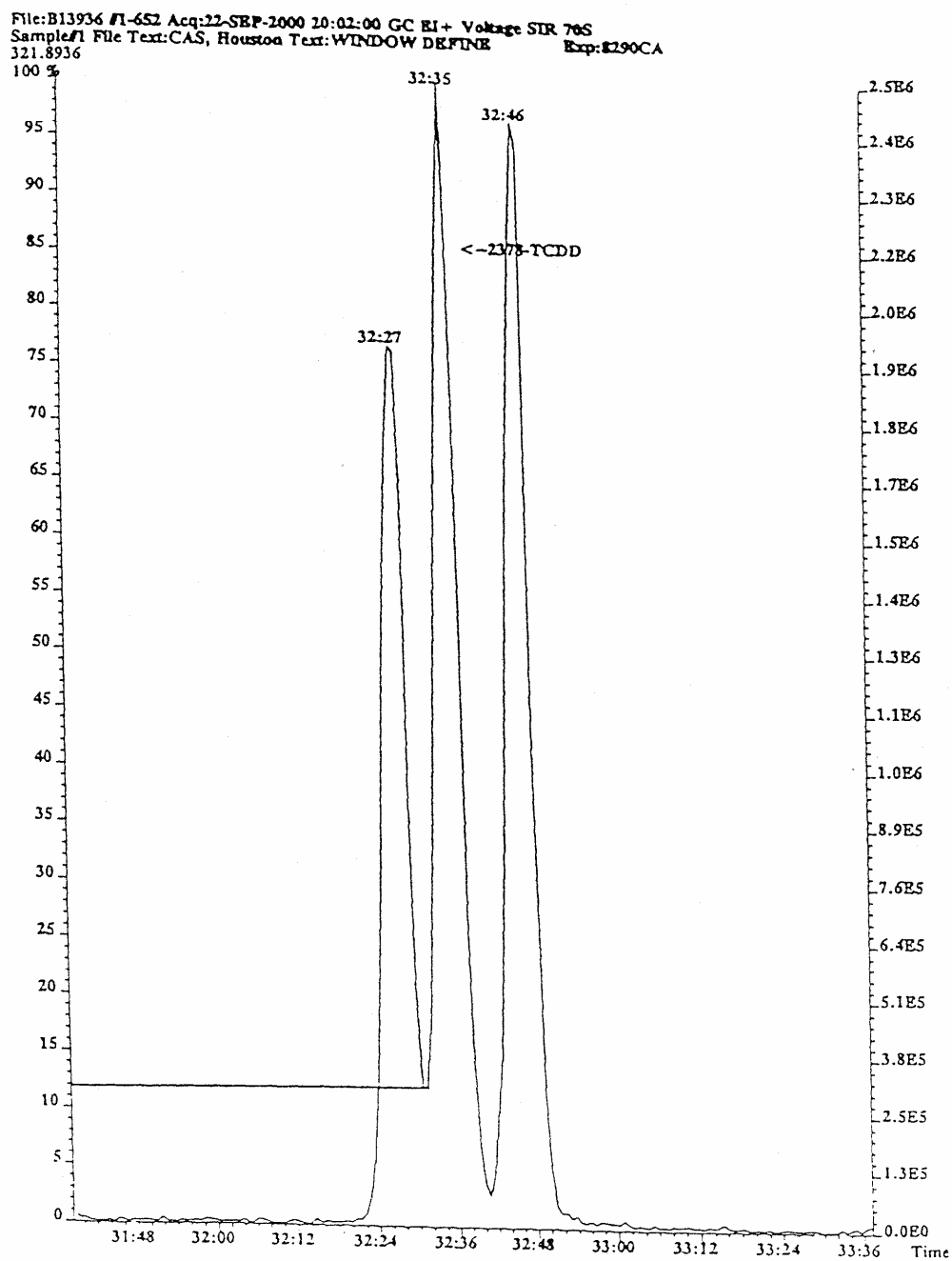


FIGURE 5 (cont.)

DB-5 WINDOW DEFINITION

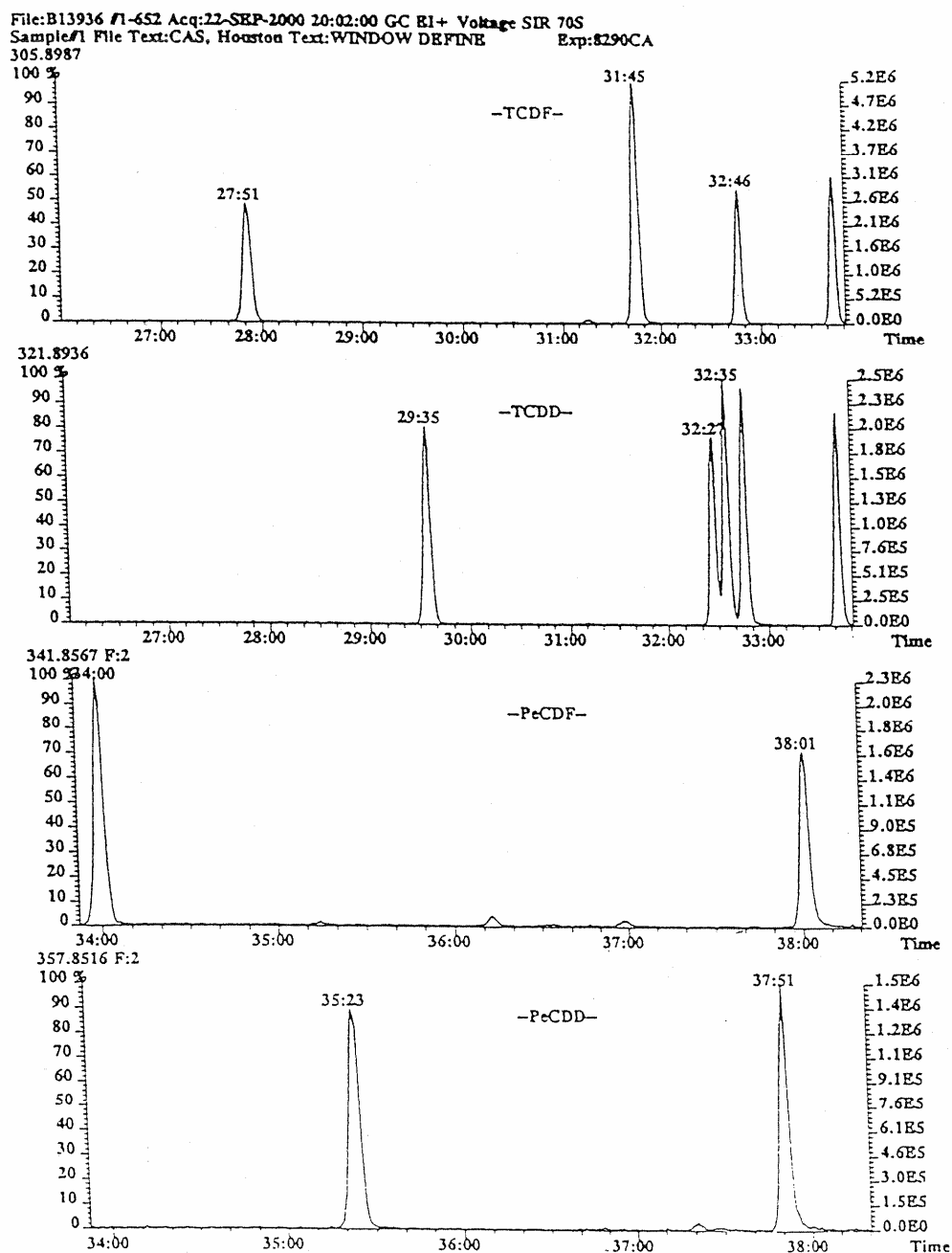


FIGURE 5 (cont.)

DB-5 WINDOW DEFINITION


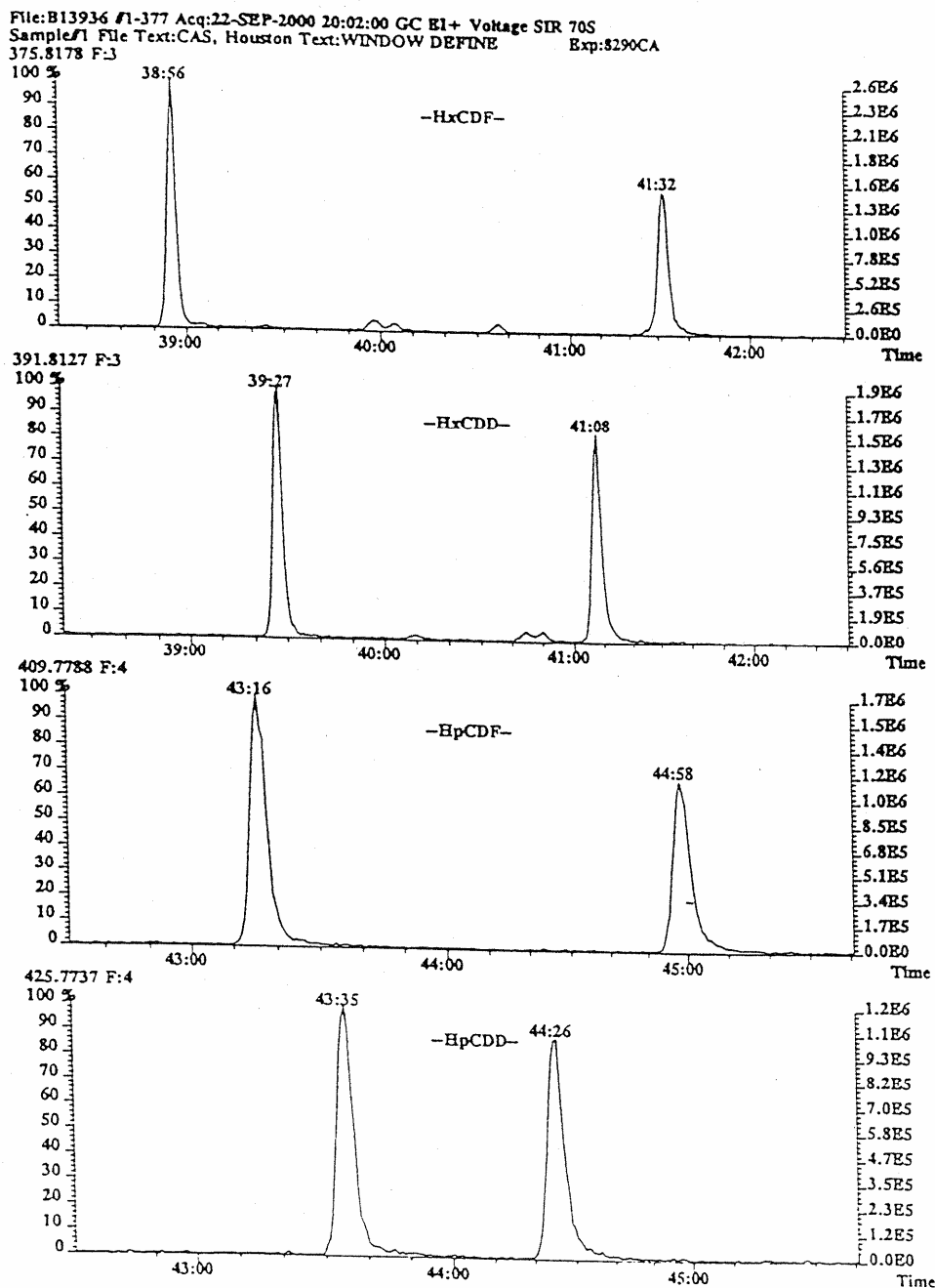


FIGURE 6
DB-225 COLUMN PERFORMANCE

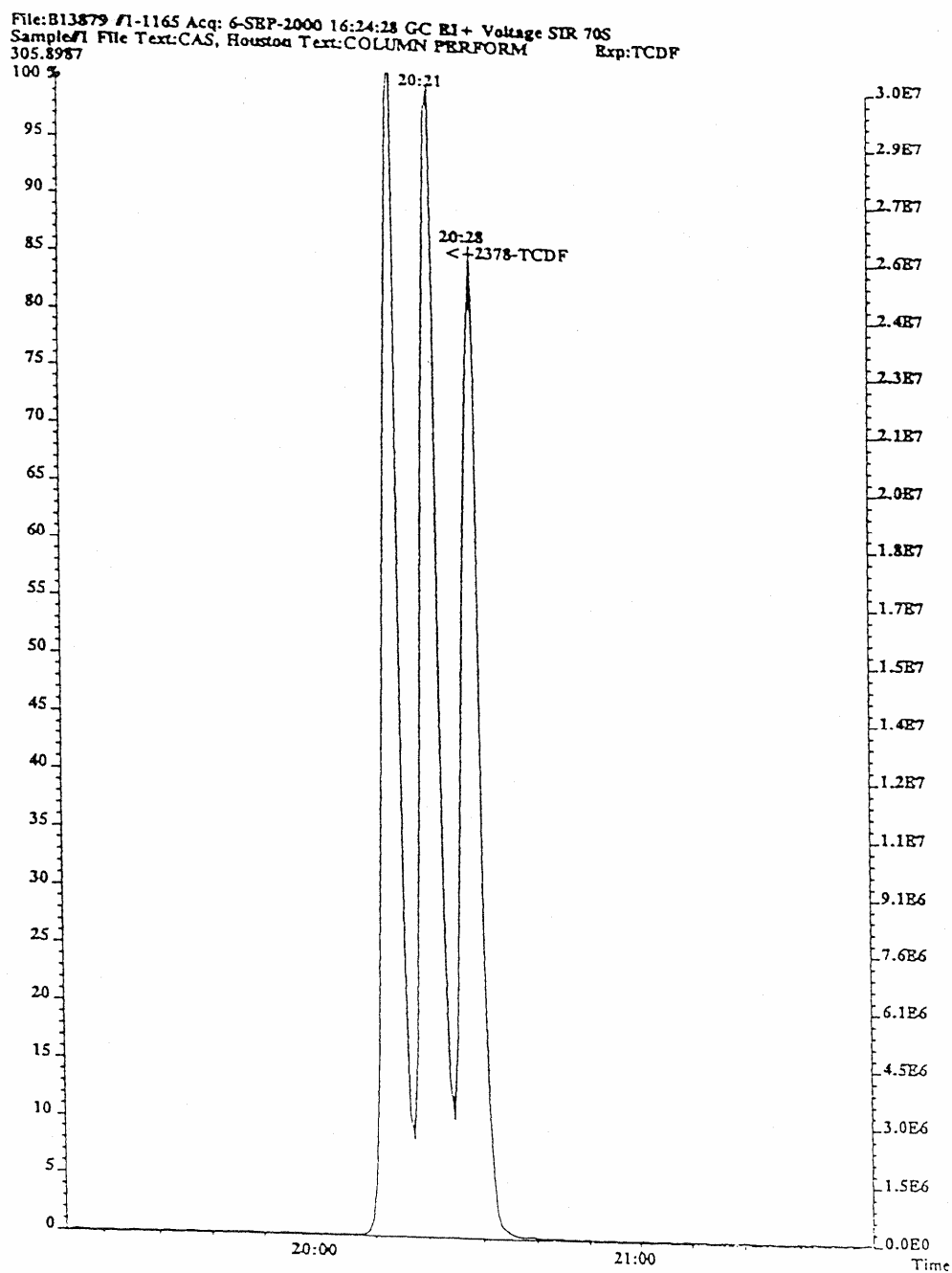


STANDARD OPERATING PROCEDURE
for
METALS DIGESTION,
SOILS, SEDIMENTS, AND SLUDGE FOR ICP ANALYSIS
FOR INDIANA PINES SITE

SOP No.: MET-3050pines

Revision: 0

September 28, 2004


Approved by:

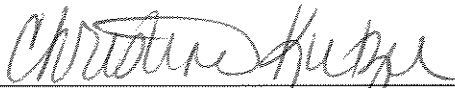
Supervisor


Date


QA Coordinator


Date


Laboratory Manager


Date

© Columbia Analytical Services, Inc., 2004
1 Mustard Street, Suite 250
Rochester, New York 14609

Annual review of this SOP has been performed
and the SOP still reflects current practice.

Initials: _____ Date: _____

Initials: _____ Date: _____

Initials: _____ Date: _____

NON-CONTROLLED COPY
Will Not Be Updated

Table of Contents**Page**

1. Scope and Applicability.....	3
2. Summary of Method.....	3
3. Definitions.....	3
4. Health and Safety Warnings.....	3
5. Cautions.....	4
6. Interferences.....	4
7. Personnel Qualifications.....	4
8. Equipment and Supplies.....	4
9. Procedure.....	5
9.1. Sample Collection.....	5
9.2. Sample Handling and Preservation.....	5
9.3. Sample Preparation	5
9.4. Sample Analysis.....	6
9.5. Troubleshooting.....	6
9.6. Data Acquisition, Calculations, and Data Reduction Requirements.....	6
10. Data and Records Management.....	6
11. Quality Control and Quality Assurance.....	7
12. References.....	7
Attachments	
Table 1 Spike Concentrations.....	8
Digestion Log Benchsheet.....	9

1 SCOPE AND APPLICABILITY

This SOP uses EPA SW-846 Method 3050B for the digestion of soils, sludges, or sediments for analysis by ICP. As stated in the EPA method, "this method is not a total digestion technique for most samples. It is a very strong acid digestion that will dissolve almost all elements that could become environmentally available." By design, elements bound in silicate structures are not normally dissolved by this procedure as they are not usually mobile in the environment." This SOP was written specifically for the Indiana Pines Site.

2 SUMMARY OF METHOD

A representative aliquot of sample is digested in nitric acid and hydrogen peroxide. Hydrochloric acid is used as a final reflux acid.

3 DEFINITIONS

- 3.1 **Laboratory Duplicates** - Two aliquots of the same sample taken in the laboratory and analyzed separately with identical procedures. Analyses of duplicates indicates precision associated with laboratory procedures, but not with sample collection, preservation, or storage procedures.
- 3.2 **Laboratory Control Sample Soil (LCSS)** - An aliquot of a soil to which known quantities of the method analytes are added. The LCSS is analyzed exactly like a sample, and its purpose is to determine whether the methodology is in control and whether the laboratory is capable of making accurate and precise measurements.
- 3.3 **Matrix Spike** - An aliquot of an environmental sample to which known quantities of the method analytes are added in the laboratory. The matrix spike is analyzed exactly like a sample, and its purpose is to determine whether the sample matrix contributes bias to the analytical results.
- 3.4 **Preparation Blank (PB)** - An aliquot of reagent water or other blank matrices that are treated exactly as a sample including exposure to all glassware, equipment, solvents, reagents, and internal standards that are used with other samples. The PB is used to determine if method analytes or other interferences are present in the laboratory environment, reagents, or apparatus.
- 3.5 **Digestion Batch** - A digestion batch is no more than 20 samples of the same matrix digested as a unit per day.

4 HEALTH AND SAFETY WARNINGS

Nitric and Hydrochloric acids are extremely corrosive. Care should be taken while working with these chemicals. Personal protective equipment including safety glasses (with side shields), gloves, and lab coat shall be worn when handling samples or reagents.

5 CAUTIONS

Antimony is easily lost by volatilization. Do not boil the digestate.

6 INTERFERENCES

Use more sample for those samples with high moisture content to meet detection limits.

7 PERSONNEL QUALIFICATIONS

At a minimum, personnel must have attained at least a 2-year degree in a science-related field and have successfully completed an Initial Demonstration of Capability and the Training Plan Form (attached). Training and Demonstration of Capability are in accordance with NELAC 2002 standard.

8 EQUIPMENT AND SUPPLIES

- 8.1 Eppendorf Pipettors
- 8.2 Funnels
- 8.3 Mortar and pestle
- 8.4 Tongue depressors
- 8.5 Filter paper
- 8.6 Hot Block Digestor with ETR-3200 Controller by Environmental Express, LTD.
- 8.7 Graduated block digestor cups
- 8.8 Block Digestor Filters.
- 8.9 CPI MOD Block Digestor

- 8.10 Reagent water - ASTM Type II deionized water.

- 8.11 Concentrated nitric acid (Baker Instra-Analyzed 69-70%): Store at room temperature in the dark in the original container or in glass. Expires per manufacturer's indications or one year from receipt if no indication is given.

- 8.12 Concentrated hydrochloric acid (Baker Instra-Analyzed 36.5-38%): Store at room temperature in the original container or in glass. Expires per manufacturer's indications or one year from receipt if no indication is given.

- 8.13 Hydrogen peroxide (30%) - H_2O_2 . Purchased commercially. Should be demonstrated to be free of impurities at levels which would interfere with sample determinations. Store at room temperature in the original container. Expires upon manufacturer's indications or 1 year from receipt if no indication is given.

- 8.14 ERA Soil Laboratory Control Sample (LCSS) - Concentrations and Performance Acceptance Limits distributed through vendor. Store at room temperature. Expires upon manufacturer's indications or 1 year from receipt if no indication is given.

- 8.15 Metals spiking solutions – Purchased commercially. See Table 1. Store at room temperature. Stocks expire upon manufacturer's indications or 1 year from receipt, whichever is sooner. Solutions prepared from stocks expire 6 months from preparation.

9 PROCEDURES

9.1 Sample Collection – Collect samples in purchased, certified clean glass or plastic.

9.2 Sample Handling and Preservation – Analyze samples within 6 months of sample collection. Store samples in a refrigerator or at room temperature. Sample receiving, handling, storage, and custody procedures are in accordance with NELAC 2002 Standard.

9.3 Sample Preparation

- 9.3.1 Set the temperature on the Block Digestor to a temperature that brings the sample temperature to 90-95°C without boiling.
- 9.3.2 The Hot Block is on a timer which can be set to turn on and off whenever necessary. To set timer press the timer button and choose the days M-F (Monday through Friday). Then choose the hour and minutes to start and stop the Block Digestor.
- 9.3.3 Label graduated hot block digestor sample cups with appropriate sample IDs for digestion.
- 9.3.4 Mix the sample thoroughly to achieve homogeneity using a tongue depressor or the mortar and pestle.
- 9.3.5 Weigh (to the nearest 0.01g) 1.00g to 1.50g of sample into labeled digestor sample cup. For sludges and sediments that have a high moisture content, use more sample. The goal is to use about 1g of dry weight sample. At this point add the appropriate spiking solutions (see Table 1) directly onto the designated spike sample prior to addition of reagents.
- 9.3.6 Unless otherwise specified by project requirements, the addition of acid should be as follows: Add 10ml of 1:1 HNO₃ and 1.5 mL of 1:1 HCl, cover with reflux cap and reflux for 15 minutes. The sample temperature should be 90-95°C. Allow the sample to cool, then add 5ml of concentrated HNO₃, cover and reflux for 30 minutes. Repeat the addition of 5ml of HNO₃ and reflux to 5 mLs. Do not allow the sample to go to dryness. CAUTION: Do not boil. Antimony is easily lost by volatilization.

- 9.3.7 Cool the sample and add 2ml of DI and 3ml of 30% H_2O_2 . Cover and heat to start the peroxide reaction. Care must be taken to ensure that losses do not occur due to excessive effervescence. Heat until effervescence subsides and cool the sample cup.
- 9.3.8 If the effervescence does not subside, add 3 mLs of hydrogen peroxide with warming to each of the samples (including blanks and LCSs) in the batch. If necessary, continue to add 30% H_2O_2 in 1ml aliquots with warming until the effervescence is minimal, or until the general sample appearance is unchanged. Do not add more than 10ml of 30% H_2O_2 .
- 9.3.9 Add 10 mL 1:1 HCL.
- 9.3.10 Cover and reflux the samples for 15 minutes without boiling. Allow to cool.
- 9.3.11 Rinse filters with 1:1 nitric acid and DI.
- 9.3.12 All samples are diluted to 100 mLs with DI. Quantitatively transfer the digestate to a graduated cylinder by pouring the sample through a prepared filter into the cylinder and rinsing the beaker and reflux cap with DI into the filter. Rinse the filter with DI. Bring to volume with DI. Pour into a labeled B-cup.
- 9.4 **Sample Analysis** – Give digested samples and a copy of the prep sheet to the ICP analyst. Analyze according to MET-6010Bpines.
- 9.5 **Troubleshooting** - All hoods in the Metals Prep Lab are wiped down once a week with DI water. The tops of all digestion hot plates are wiped down daily.
- 9.6 **Data Acquisition, Calculations and Data Reduction Requirements**
- Digestion logs are used to record all sample volumes, spike volumes, etc. The Manufacturer's lot number for the reagents used are added to the digestion log (see attached digestion log benchsheet).

10 DATA AND RECORDS MANAGEMENT

- 10.1 **Responsibilities** - It is the responsibility of the analyst to perform the analysis according to this SOP and to complete all documentation required for data review. Analysis and interpretation of the results are performed by personnel in the laboratory who have demonstrated the ability to generate acceptable results utilizing this SOP. Final review and sign-off of the data is performed by the department supervisor or designee.
- 10.2 Data will be reviewed after ICP analysis according to MET-6010Bpines.

11 QA/QC REQUIREMENTS

- 11.1 Each day, digest one laboratory control sample (LCS) per digestion batch, or per 20 samples, or per EPA SDG group, whichever is more frequent. Use the appropriate solid laboratory control sample (LCSS) for soils analysis.
- 11.2 Each day, digest one blank per digestion batch, or per 20 samples, or per EPA SDG group, whichever is more frequent. Use D.I. water and follow the digestion procedures.
- 11.3 Each day, prepare one duplicate and one spiked sample with each digestion batch, or per twenty samples, or per EPA SDG group, whichever is more frequent. At times, specific samples will be assigned as duplicates of spikes depending on client requirements.
- 11.4 Matrix spikes are prepared by adding the appropriate volume of spiking solution (See Table 1).
- 11.5 See MET-6010Bpines for applicable QC limits and corrective action.

12 REFERENCES

“Test Methods For Evaluating Solid Waste, Physical/Chemical Methods”. EPA SW846, Third Edition, December 1996.

NELAC, 2002 Standard.

Table 1 Spiking Concentrations for LCS and MS Samples

SPIKE SOLUTION A		1.00ml Spk A		to Final Vol of 100ml	
<i>Metal</i>	<i>Conc. (ug/mL)</i>		<i>Metal</i>	<i>Conc. (ug/mL)</i>	
AL	200		NI	50	
AS	4		SE	1	
BA	200		AG	5	
BE	5		TL	200	
CD	5		V	50	
CR	20		ZN	50	
CO	50		B	100	
CU	25		CA	200	
FE	100		MG	200	
PB	50		NA	2000	
MN	50		K	2000	

SPIKE SOLUTION B		1.00ml Spk B		to Final Vol of 100ml	
<i>Metal</i>	<i>Conc. (ug/mL)</i>		<i>Metal</i>	<i>Conc. (ug/mL)</i>	
SB	50		TI	50	
MO	50		-	-	

INDIVIDUAL METALS		0.10ml Spk. to Final Volume of 100ml		INDIVIDUAL METALS		0.5ml Spk. to Final Volume of 100ml	
<i>Metal</i>	<i>Conc. (ug/mL)</i>			<i>Metal</i>	<i>Conc. (ug/mL)</i>		
SE	1000			SN	1000		

Batch ID:

6010B/846 // 200.7/136 // ASP/CLP4.1

Spike Witness / Lot Approval:

Date:

Analyst:

Prep Method: SW846 3050 // CLP

Batch Temp:

Report Type: Routine // ASP // Pkg5

Digest: Initial // Redigest of:

Submission / Order #	Initial Wgt. (g)	Final Vol (ml)	Initial Color / Texture	Final Color / Clarity	Metals	Spike Vol (ml)
1						
2						
3						
4						
5						
6						
7						
8						
9						
10						
11						
12						
13						
14						
15						
16						
17						
18						
19						
20						
21						
22						
23						
24						

Color / Clarity Key:

Color: C = Colorless ; Y = Yellow ; B = Brown

BL = Black ; G = Grey ; W = White

Clarity: CDY = Cloudy ; CLR = Clear ; OP = Opaque

Texture: F = Fine ; M = Medium ; CS = Coarse ; NAQ = Non Aqueous

Spiking Standards / Reagent Lot #:

Spike A,B: Spike #4: _____

TCLP Spk: TCLP Ba: _____

Se Std: Sn Std: _____

HNO3: HCL: _____

H2O2: LCSS: _____

Comments:

01

SOP NO.: MET-6010BPINES

Revision: 1

Date: 9/29/04

Page: 1 of 37

STANDARD OPERATING PROCEDURE

for

DETERMINATION OF METALS AND TRACE ELEMENTS BY INDUCTIVELY COUPLED PLASMA ATOMIC EMISSION SPECTROMETRY (ICP) FOR INDIANA PINES SITE

SOP No.: MET-6010BPINES

Revision: 1

September 29, 2004

Approved by:



Department Supervisor




Date



Laboratory Manager



Date

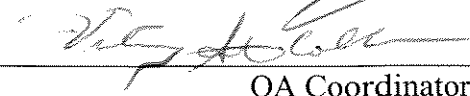

QA Coordinator



Date

© Columbia Analytical Services, Inc., 2004

1 Mustard Street, Suite 250

Rochester, New York 14609

Annual review of this SOP has been performed
and the SOP still reflects current practice.

Initials: _____ Date: _____

Initials: _____ Date: _____

Initials: _____ Date: _____

NON-CONTROLLED COPY

Will Not Be Updated

Table of Contents	<u>Page</u>
1. Scope and Applicability.....	3
2. Summary of Method.....	3
3. Definitions.....	4
4. Health and Safety Warnings.....	5
5. Interferences.....	6
6. Personnel Qualifications.....	6
7. Equipment and Supplies.....	6
8. Procedure.....	8
8.1. Calibration and Standardization.....	8
8.2. Sample Collection.....	8
8.3. Sample Handling and Preservation.....	9
8.4. Sample Preparation	9
8.5. Sample Analysis.....	10
8.6. Troubleshooting.....	11
8.7. Data Acquisition, Calculations, and Data Reduction Requirements.....	12
8.8. Computer Hardware and Software.....	13
9. Data and Records Management.....	13
10. Quality Control and Quality Assurance.....	13
11. References.....	15
Attachments	
Training Plan Form.....	16
Table 1 List of Analytes and Practical Quantitation Limits.....	17
Table 2 Recommended Wavelengths and Instrument Specifications.....	18
Table 3 Operating Conditions.....	19
Appendix I Standard Preparation Instructions.....	20

1. SCOPE AND APPLICABILITY

- 1.1. This SOP uses EPA SW-846 Method 6010B for the determination of trace elements, including metals, in solution using Inductively coupled plasma-atomic emission spectrometry (ICP-AES). The method is applicable to all of the elements listed in Table 1. All matrices, including ground water, aqueous samples, TCLP and EP extracts, industrial and organic wastes, soils, sludges, sediments, and other solid wastes, require digestion prior to analysis.
- 1.2. Detection limits, sensitivity, and the optimum and linear concentration ranges of the elements can vary with the wavelength, spectrometer, matrix and operating conditions. The Method Reporting Limits (MRL) are listed in Table 1. The reported MRL may be adjusted if required for specific project requirements, however, the capability of achieving other reported MRLs must be demonstrated. Results may be reported to the Instrument Detection Limits (IDLs) upon request.
- 1.3. This SOP was modified specifically for the Indiana Pines site project.

2. SUMMARY OF METHOD

- 2.1. Samples are digested according to one of the proper metals digestion methods listed in SW-846.
- 2.2. This method describes multielemental determinations by ICP-AES using sequential or simultaneous optical systems and axial or radial viewing of the plasma. The instrument measures characteristic emission spectra by optical spectrometry. Samples are nebulized and the resulting aerosol is transported to the plasma torch. Element-specific emission spectra are produced by a radio-frequency inductively coupled plasma. The spectra are dispersed by a grating spectrometer, and the intensities of the emission lines are monitored by photosensitive devices. Background correction is required for trace element determination. Background must be measured adjacent to analyte lines on samples during analysis. The position selected for the background-intensity measurement, on either or both sides of the analytical line, will be determined by the complexity of the spectrum adjacent to the analyte line. In one mode of analysis the position used should be as free as possible from spectral interference and should reflect the same change in background intensity as occurs at the analyte wavelength measured. Background correction is not required in cases of line broadening where a background correction measurement would actually degrade the analytical result. The possibility of additional interferences (discussed later) should also be recognized and appropriate corrections made.

3. DEFINITIONS

- 3.1. Calibration Blank - A volume of reagent water acidified with the same acid matrix as in the calibration standards. The calibration blank is a zero standard and is used to calibrate the ICP instrument.
- 3.2. Calibration Standard (CAL) - A solution prepared from the dilution of stock standard solutions. The CAL solutions are used to calibrate the instrument response with respect to analyte concentration
- 3.3. Dissolved Analyte - The concentration of analyte in an aqueous sample that will pass through a 0.45 μm membrane filter assembly prior to sample acidification.
- 3.4. Instrument Detection Limit (IDL) - The concentration equivalent to the analyte signal which is equal to three times the standard deviation of a series of 10 replicate measurements of the calibration blank signal at the same wavelength.
- 3.5. Initial/Continuing Calibration Verification Solution (ICV/CCV) - A solution of method analytes, used to evaluate the performance of the instrument system with respect to a defined set of method criteria.
- 3.6. Internal Standard - Pure analyte(s) added to a sample, extract, or standard solution in known amount(s) and used to measure the relative responses of other method analytes that are components of the same sample or solution. The internal standard must be an analyte that is not a sample component
- 3.7. Laboratory Duplicates - Two aliquots of the same sample taken in the laboratory and analyzed separately with identical procedures. Analyses of duplicates and indicates precision associated with laboratory procedures, but not with sample collection, preservation, or storage procedures.
- 3.8. Laboratory Control Sample (LCS) - An aliquot of to which known quantities of the method analytes are added in the laboratory. The LCS is analyzed exactly like a sample, and its purpose is to determine whether the methodology is in control and whether the laboratory is capable of making accurate and precise measurements.
- 3.9. Matrix Spike - An aliquot of an environmental sample to which known quantities of the method analytes are added in the laboratory. The matrix spike is analyzed exactly like a sample, and its purpose is to determine whether the sample matrix contributes bias to the analytical results.

- 3.10. Preparation Blank (PB) - An aliquot of reagent water or other blank matrices that are treated exactly as a sample including exposure to all glassware, equipment, solvents, reagents, and internal standards that are used with other samples. The PB is used to determine if method analytes or other interferences are present in the laboratory environment, reagents, or apparatus.
- 3.11. Linear Range - The concentration range over which the instrument response to an analyte is linear.
- 3.12. Method Detection Limit (MDL) - The minimum concentration of an analyte that can be identified, measured, and reported with 99% confidence that the analyte concentration is greater than zero.
- 3.13. Plasma Solution - A solution that is used to determine the optimum height above the work coil for viewing the plasma.
- 3.14. Interference Check Solution (ICS) - A solution of selected method analytes of higher concentrations which is used to evaluate the procedural routine for correcting known interelement spectral interferences with respect to a defined set of method criteria.
- 3.15. Method Reporting Limit Standard (MRL) - Standard prepared with a known concentration of elements to check accuracy at the low end of the curve.
- 3.16. HLCCV1 - A standard prepared at the bench at a high concentration to encompass the range of the samples being analyzed. This standard is used to assess accuracy at the high end of the linear range.
- 3.17. HLCCV2 - A standard prepared slightly higher than the calibration range for metals.
- 3.18. Batch - a group of no more than 20 field samples digested or analyzed together on the same day with the same reagents.

4. HEALTH AND SAFETY WARNINGS

- 4.1. Corrosives - Because all samples and standards are diluted in 2% HNO₃ and 5% HCl, there is a danger of exposure to corrosives, sufficient care must be taken in handling these solutions. Safety glasses must be worn while preparing and handling the solutions.
- 4.2. High Voltage - The power unit supplies high voltage to the RF generator which is used to form the plasma. The unit should never be opened. Exposure to high voltage can cause injury or death.

- 4.3. UV Light - The plasma when lit is a very intense light, and must not be viewed with the naked eye. Protective lenses are in place on the instrument. Glasses with special protective lenses are available.
- 4.4. When the nature of the sample is either unknown or is known to be hazardous, acidification should be done in a well ventilated area or fume hood

5. INTERFERENCES

There are several types of interferences by the ICP's: Spectral interferences can be from an overlap of spectral lines, background points or background from line emissions of high concentration elements. Physical interferences are effects associated with the sample introduction process, example high dissolved solids buildup on the nebulizer tip. Chemical interferences caused by the sample matrix itself. IEC's aid in eliminating some of these interferences. IECs are interelement correction factors that the instrument uses to compensate for spectral overlap when analyzing samples with complex spectra. Refer to Method 6010B Section 3.0 or Method 200.7 Section 4.0 for more detail and suggested procedures to correct and adjust the instrument due to interferences.

6. PERSONNEL QUALIFICATIONS

At a minimum, personnel must have attained at least a 4-year degree (or 2-yr degree plus one year experience) in a science-related field and have successfully completed an Initial Demonstration of Capability and the Training Plan Form (attached). Training and Demonstration of Capability are in accordance with NELAC 2002 standard.

7. EQUIPMENT AND SUPPLIES

- 7.1. ICP- Perkin Elmer Optima 3000XL - Inductively coupled argon plasma emission spectrometer (ICP) equipped with the following:
 - 7.1.1. Computer-controlled emission spectrometer with background correction.
 - 7.1.2. Mass flow controller for argon nebulizer gas supply.
 - 7.1.3. Peristaltic pump.
 - 7.1.4. Autosampler.
 - 7.1.5. Argon gas supply - high purity.

7.2. Volumetric flasks, class A .

7.3. Trace metals grade chemicals shall be used in all tests.

7.3.1. Hydrochloric acid (conc), HCl. Purchased commercially. Store at room temperature. Expires three years from receipt or upon manufacturer's indications, whichever is sooner.

7.3.2. Hydrochloric acid (1:1), HCl. Add 500 mL concentrated HCl to 400 mL water and dilute to 1 liter in an appropriately sized beaker. Store at room temperature. Expires one year from preparation.

7.3.3. Nitric acid (conc), HNO₃. Purchased commercially. Store at room temperature. Expires three years from receipt or upon manufacturer's indications, whichever is sooner.

7.3.4. Nitric acid (1:1), HNO₃. Add 500 mL concentrated HNO₃ to 400 mL water and dilute to 1 liter in an appropriately sized beaker. Store at room temperature. Expires one year from preparation.

7.4. Reagent Water. All references to water in the method refer to DI Type II water unless otherwise specified. Reagent water will be interference free.

7.5. All standards are prepared from NIST traceable stock standard solutions. Manufacturers expiration dates are used to determine viability of standards. Preparatory procedures for standards and QC solutions vary between instruments due to the working ranges. All preparatory information for the QA/QC samples are provided in Appendix I.

7.5.1. Mixed Calibration Standards are prepared by combining appropriate volumes of the stock solutions in volumetric flasks. Matrix match with the appropriate acid and dilute to 100ml with water. Calibration standards should be verified using a second source quality control sample (LCS, ICV, or CCV). Calibration standards should be stored at room temperature in glass volumetric flasks with a shelf-life of 7 days.

7.5.2. Initial and Continuing Calibration Verification (ICV and CCV) Standards are prepared by combining compatible analytes at concentrations equivalent to the midpoint of their respective calibration curves. The ICV and CCV standards should be prepared from a separate source independent from that used in the calibration standards. ICV / CCV standards should be stored at room temperature in glass volumetric flasks with a shelf-life of 48 hours.

- 7.5.3. MRL Standards are prepared to contain known concentrations of elements at or near the Method Reporting Limit. MRL standards should be stored in plastic containers with a shelf-life of 6 months.
- 7.5.4. Interference Check Solutions A and AB are prepared to contain known concentrations of interfering analytes that will provide an adequate test of the correction factors. ICSA / ICSAB standards should be stored in plastic containers with a shelf-life of 6 months.
- 7.5.5. Laboratory Control Sample and Matrix Spike are purchased as custom mixes stored in plastic containers with a shelf-life of 6 months at the concentrations recommended in the method. Certificates of analysis are attached in Appendix I. Each sample, up to 100 mL, is spiked with 1.0 ml of spike solution.

7.6. Blanks

- 7.6.1. Method Blanks must contain all the reagents and in the same volumes as used in the preparation of samples. The method blanks must be carried through the complete procedure and contain the same acid concentration in the final solution as the samples.
- 7.6.2. The Calibration Blank is prepared by acidifying reagent water to the same concentrations of acid found in the standards and samples.

7.7. Reagent Receiving Log

The manufacturer, lot number, standard /reagent name, concentration, date received and expiration date are recorded in a reagent log.

8. PROCEDURE

8.1. Calibration and Standardization

Calibration is accomplished daily using 3 calibration standards and a blank for each element using the internal standard technique. See Sample Analysis section for more information.

8.2. Sample Collection

Containers may be glass or plastic. Samples are cooled with ice to be shipped to the laboratory.

8.3. Sample Handling and Preservation

- 8.3.1. Solid samples require no preservation prior to analysis other than storage at 0-6°C. Samples are analyzed within 6 months of collection.
- 8.3.2. Aqueous samples are acid preserved with (1+1) nitric acid to pH <2. Samples are analyzed within 6 months of sample collection.
- 8.3.3. Samples are checked upon receipt for all the elements listed in the Sample Acceptance Policy found in NELAC 2002 Standard.
- 8.3.4. For the determination of the dissolved elements, filter the sample through a 0.45 µm pore diameter membrane filter at the time of collection or as soon thereafter as practically possible. (Glass or plastic filtering apparatus are recommended to avoid possible contamination. Only plastic apparatus should be used when the determinations of boron and silica are critical.) Use a portion of the filtered sample to rinse the filter flask, discard this portion and collect the required volume of filtrate. Acidify the filtrate with (1+1) nitric acid immediately following filtration to pH <2.

Note: When the nature of the sample is either unknown or is known to be hazardous, acidification should be done in a well ventilated area or fume hood.
- 8.3.5. Samples received by the ICP lab as digestates contain nitric and hydrochloric acid. Digestates are stored at room temperature in plastic B-cups.
- 8.3.6. Following analysis, digestates are stored until all results have been reviewed. Digestates are diluted and disposed of through the sewer system in approximately 90 days after receipt of sample.

8.4. Sample Preparation

- 8.4.1. Digest samples prior to analysis. Refer to the following Metals Methods found in SW-846:
 - 3005A Metals Digestion, Waters, Total Recoverable and Dissolved for ICP
 - 3010A Metals Digestion, Waters for ICP
 - 3020A Metals Digestion, Waters for GFAA
 - 3050B Metals Digestion, Soils, Sediments and Sludges for ICP and GFAA

8.5. Sample Analysis

- 8.5.1. Set up the instrument with proper operating parameters established as detailed below. The instrument must be allowed to become thermally stable before beginning (usually requiring at least 45 minutes of operation prior to calibration). Operating Conditions - The analyst should follow the instructions provided in Table 3.
- 8.5.2. Before using this procedure to analyze samples, there must be data available documenting initial demonstration of performance. The required data documents the selection criteria of background correction points; linear ranges, and the upper limits of those ranges; the method and instrument detection limits; and the determination and verification of interelement correction equations or other routines for correcting spectral interferences. This data must be generated using the same instrument, operating conditions and calibration routine to be used for sample analysis. These documented data must be kept on file and be available for review by the data user or auditor.
- 8.5.3. Turn on power supply for the instrument, computer, printer and light the plasma. Allow instrument to warm-up for 45-60 minutes before operation. The cooling water and the argon are on when the instrument are on.
- 8.5.4. Profile the instrument on a daily basis, and when maintenance is done to align it optically for both horizontal and vertical optimization in either mode. Aspirate a 10 ppm source of manganese (as recommended by the manufacturer). Choose the Tools menu/Spectrometer Control/Optimize X&Y. The instrument automatically adjusts the torch viewing position for maximum intensity.
- 8.5.5. Pour the 3 calibration standards, ICV/CCV standards, MRL, ICSA, and ICSAB up to 40 mL in 50 mL centrifuge tubes and add 0.80 mL of the internal standard solution. Pour all other samples, preparation blanks and laboratory control samples up to 10 mL in 15 mL centrifuge tubes and add 0.20 mL of internal standard solution. This gives an apparent concentration of 1.00 mg/L Yttrium. The Yttrium intensity is used by the instrument to ratio the analyte intensity signals for both calibration and quantitation. Cesium is used only as a stabilizer.
- 8.5.6. Internal standards can be added via pump and mixing block. This technique uses a solution of 10 mg/L Y and 10 mg/L Cs.

8.5.7. Following the calibration, analyze in the following sequence:

- ICV; ICB; MRL; ICSA; ICSAB; CCV; CCB;
- 10 environmental samples (including PBs and LCSs); CCV; CCB; repeat to the end of the run.....
- Last 10 samples; CCV; CCB, MRL, ICSA, ICSAB; HLCCV1; HLCCV2; CCV; CCB.

8.5.8. Rinse the system with the calibration blank solution before the analysis of each sample for one minute.

8.5.9. Samples which exceed the linear range of the instrument must be diluted and reanalyzed.

8.5.10. Method detection limits must be established for all wavelengths utilized for each type of matrix commonly analyzed. The matrix used for the MDL calculation must contain analytes of known concentrations within 3-5 times the anticipated detection limit. See Table 2 for approximate wavelengths. See 40 CFR Part136 Appendix B for more information.

8.6. Troubleshooting

8.6.1. All maintenance activities are recorded in a maintenance logbook kept for each instrument. Most routine maintenance and troubleshooting is performed by CAS staff. Other maintenance or repairs may, or may not require factory service, depending upon the nature of the task. Record the analytical run filename of the first acceptable run after major maintenance in the maintenance log book. Typical preventive maintenance measures include, but are not limited to, the following items:

- Cleaning the pump tubing as needed
- Empty waste container, as needed
- Cleaning the nebulizer, spray chamber, and torch, as needed
- Replace water and vacuum filters, as needed

8.7. Data Acquisition, Calculations, and Data Reduction Requirements

8.7.1. Calculations: If dilutions were performed, the appropriate factors must be applied to sample values. All results should be reported with up to three significant figures.

8.7.2. Sample Calculation (water)

$$\text{Conc. (mg/L)} = \frac{\text{Instrument Reading (mg/L)} \times \text{Final digestion volume (L)}}{\text{Initial volume (L)}}$$

8.7.3. Sample Calculation (soils)

$$\text{Conc. (mg/g)} = \frac{\text{Instrument Reading (mg/L)} \times \text{Final digestion volume (L)}}{\text{Initial mass (g)} \times \text{Percent Solids expressed as a decimal}}$$

8.7.4. **Matrix Spike Recovery** is calculated to determine accuracy for matrix and blank spikes using the following equation:

$$\text{Accuracy (\%REC)} = \frac{A - B}{C} \times 100$$

Where

A = Analyte total concentration from spiked sample

B = Analyte concentration from unspiked sample

C = Concentration of spike added

8.7.5. **Precision** is measured through the use of replicate sample analyses within the same batch and is expressed as the relative percent difference (RPD) between the replicate measurements.

$$\text{RPD} = \frac{|D1 - D2|}{(D1 + D2)/2} \times 100$$

Where

D1 = Original Result

D2 = Duplicate Result

Report each analyte concentration to the proper significant figures in mg/L or µg/L as required.

8.8. Computer Hardware and Software

Each ICP uses a Gateway GP5-233 running the ICP WinLab v.1.42.
Metals Analytical Review and Reporting System (MARRS) v.3.2.44
StarLIMS v.6.11.a

9. DATA AND RECORDS MANAGEMENT

- 9.1. **Responsibilities** - It is the responsibility of the analyst to perform the analysis according to this SOP and to complete all documentation required for data review. Analysis and interpretation of the results are performed by personnel in the laboratory who have demonstrated the ability to generate acceptable results utilizing this SOP. Final review and sign-off of the data is performed by the department supervisor or designee.
- 9.2. **Data Flow** – Samples are entered by the Project Manager into StarLIMS on a Personal Computer running on a Novell Network. On the day that the samples are received the samples appear on a daily log printed from this computer system. The Metals Prep analyst prepares a benchsheet, digests the samples and turns the samples and digest sheet over to the ICP analyst. The samples are analyzed for metals of interest using ICP software. The results are transferred to MARRS (for reporting package work) and StarLIMS for validation, reporting, and invoicing.
- 9.3. **Data Review** - Data will be reviewed by the ICP analyst and a qualified peer using a Data Review Checklist (attached) and validated by a supervisor.

10. QUALITY CONTROL AND QUALITY ASSURANCE

- 10.1. Instrument values are based on duplicate readings. Precision between the emission readings shall not exceed 20 %RSD. If RSD values exceed 20%, the sample reanalyzed and reported.
- 10.2. Preparation Blanks must be analyzed at least one PB with each batch of 20 or fewer samples of the same matrix. PB values must not exceed the MRL. Fresh aliquots of the samples must be prepared and analyzed again for affected analytes after the source of the contamination has been corrected and acceptable PB values have been obtained. If detections are greater than the MRL, the batch needs to be redigested if sample concentration is less than 5 times the concentration found in the prep blank. If the sample concentration is less than the MRL the sample does not require redigestion.
- 10.3. HLCCV1 – High standard used in curve and analyzed once during daily analysis. Should agree within 10% of the true value. If HLCCV1 is > 10% different the analysis is judged to be out of control and the source of the problem should be identified and resolved before continuing analysis.

- 10.4. HLCCV2 – standard slightly higher than calibration for some metals. Analyzed once during daily analysis. Should agree within 10% of the true value. If out of control, client data above the HLCCV1 should be re-analyzed.
- 10.5. ICV/CCV - Calibration Verification Standards must immediately follow each calibration, after every tenth sample, and at the end of the sample run. Initial Calibration Verification must verify that the instrument is within 10%. Continuing Calibration Verification standards must confirm the calibration within $\pm 10\%$ throughout the analyses. If the recovery of an analyte falls outside the required control limits, the analysis is judged to be out of control, and the source of the problem should be identified and resolved before continuing analysis. Recalibrate the instrument.
- 10.6. The results of the calibration blank (CCB) must be less than the MRL. If not, terminate the analysis, correct the problem, recalibrate, and reanalyze the samples effected.
- 10.7. Dilute and reanalyze samples that exceed the linear calibration range or use an alternate, less sensitive line for which quality control data is already established.
- 10.8. Analyze matrix spiked and duplicate samples at a frequency of one per matrix batch (max. 20 samples). Matrix spiked and duplicate samples are brought through the entire sample preparation and analytical process.
 - 10.8.1. The spiked sample or spiked duplicate sample recovery is to be within $\pm 25\%$ of the actual value or within the documented historical acceptance limits for each matrix. Sample concentrations greater than four times the spike concentration are not valid and shall not be evaluated. If the matrix spike does not meet these criteria, analyze a Post Digestion Spike.
 - 10.8.2. A control limit of $\pm 20\%$ RPD shall be used for original and duplicate samples greater than or equal to 5X the CRDL. A control limit of \pm the CRDL shall be used if either the sample or duplicate value is less than 5 times the CRDL. CRDL values are given in Table 1.
- 10.9. Laboratory Control Sample – verify sample preparation and analysis using reagent water spiked with a known amount of analytes of interest. Results should be within $\pm 20\%$. Outlying recoveries may indicate loss of analyte due to digestion procedures or laboratory contamination. If an LCS is found to be out of the specified limits, recalibrate and reanalyze. If the LCS remains out of the specified limits, redigestion of the entire batch should occur if the recovery is less than 80%. If the LCS recovery is greater than 120% redigest all positive results (greater than the MRL).

- 10.10. MRL standard- A standard at or near the MRL is analyzed at the beginning and end of each analytical run but not before the ICV. There are no limits in the 6010B method, but the CAS guideline used is +/- 50% of the true value. If the limits are not met the analysis is stopped and the instrument is recalibrated.
- 10.11. Interference Check Samples- The ICSA and ICSAB need to be run consecutively at the beginning and end of each analytical run. Results from the ICSA solution shall be monitored for false positive detections of analytes not present in the mix. The analyte recoveries for the AB solution must fall within 20% of the true value otherwise the run must be stopped, recalibrated and reanalyzed unless analytes are not detected in the associated samples or interferent elements are not present.
- 10.12. Serial Dilution Test If the analyte concentration is sufficiently high (minimally, a factor of 50 times above the IDL), an analysis of a 1:5 dilution should agree within $\pm 10\%$ of the original determination. If not, a chemical or physical interference effect should be suspected and data may be flagged accordingly.
- 10.13. Post Digestion Spike Addition: Typically if a matrix spike does not yield acceptable results, a post-digestion spike may be added to a portion of a prepared sample, or its dilution, and should be recovered to within 75% to 125% of the known value. The spike addition should produce a minimum level of 10 times and a maximum of 100 times the IDL. If the spike is not recovered within the specified limits, a matrix effect has been confirmed.
- 10.14. Instrument Performance
- InterElement Correction Factors (IEC) are analyzed annually, or as needed.
 - Linear Ranges (LR) are run biannually and must be $\pm 5\%$ of true value.
 - Instrument Detection Limits (IDL) are analyzed quarterly, or as needed.
 - Method Detection Limits (MDL) are analyzed annually.

11. REFERENCES

- *Test Methods For Evaluating Solid Waste, Physical/Chemical Methods*. USEPA SW-846, 3rd Edition, December 1996.
- *Methods For the Determination of Metals in Environmental Samples Supplement I*. USEPA/600/R-94/111, May 1994.
- 40 CFR Part 136 Appendix B
- NELAC 2002 Standard

Metals Instrument Analysis Training Plan

Procedure: _____

SOP: _____ Revision: _____ Date: _____

Trainee: _____

1. Read SOP Trainer: _____ Trainee: _____ Date: _____
2. Demonstrated understanding of the scientific basis of the analysis including:
-the chemical and physical principals behind the measurement by the instrument
Trainer: _____ Trainee: _____ Date: _____
3. Demonstrated familiarity with related SOPs
-ADM-BATCHSEQ -ADM-PCAL -ADM-SIGFIG
-ADM-DATAENTRY -ADM-DIL -ADM-SPSR
-ADM-MDL -ADM-DREV -ADM-TRANDOC
Trainer: _____ Trainee: _____ Date: _____
4. Observe performance of SOP
-standard and reagent prep and documentation – including pipet used
-instrument power up and warm-up
-instrument set-up, daily maintenance and checks
-use and loading of autosampler
-sample analysis including:
-calibration
-sample dilution
-software command of instrument
-use of QC samples and QC criteria
-common troubleshooting
-instrument logbook use
-data reduction, reporting, and review
Trainer: _____ Trainee: _____ Date: _____
5. I have read, understood and agree to perform the most recent version of the SOP:
Signature: _____ Date: _____
6. Perform SOP with supervision
- including all items in 4.
Trainer: _____ Trainee: _____ Date: _____
7. Independent performance of the SOP
-all of the item listed in 4
-IDC (4 mid-range standards performed before client samples are analyzed)
-attach IDC certificate, raw data, and summary spreadsheet.
Trainer: _____ Trainee: _____ Date: _____

TABLE 1

<i>Analyte</i>	<i>MRL Water mg/L</i>	<i>MRL Soil ug/g</i>	<i>Typical IDL ug/L</i>
Silver	0.010	1.00	0.632
Aluminum	0.100	10.0	6.57
Arsenic	0.0100	50.0	6.89
Boron	0.200	20.0	37.3
Barium	0.0200	2.00	12.2
Beryllium	0.0050	0.500	0.26
Calcium	0.500	50.0	167
Cadmium	0.0050	0.500	0.489
Cobalt	0.0500	5.00	3.03
Chromium	0.0100	1.0	1.81
Copper	0.0200	2.00	3.02
Iron	0.100	5.00	44.1
Potassium	2.00	100	857
Lithium	0.200	20.0	23.9
Magnesium	0.500	50.0	124
Manganese	0.0100	1.0	1.78
Molybdenum	0.0250	2.50	3.08
Sodium	0.500	50.0	193
Nickel	0.0400	4.00	3.92
Lead	0.00500	5.00	1.29
Antimony	0.0600	10.0	3.72
Selenium	0.00500	50.0	12.5
Silicon	1.00	100	68.7
Strontium	0.100	10.0	5.38
Tin	0.500	100	15.8
Titanium	0.0500	5.00	3.15
Thallium	0.0100	30.0	7.77
Vanadium	0.0500	5.00	2.74
Zinc	0.0200	1.0	2.47

Table 2**Recommended Wavelengths and Instrument Specifications**

Suggested wavelengths are listed below:

Analyte	Wavelength
Ag Silver	328.068
Al Aluminum	308.215
B Boron	249.773
Ba Barium	233.527
Be Beryllium	234.861
Ca Calcium	430.253
Cd Cadmium	226.502
Co Cobalt	228.616
Cr Chromium	267.716
Cu Copper	324.754
Fe Iron	238.863
Li Lithium	610.364
Mg Magnisium	279.079
Mn Manganese	257.610
Mo Molybdenum	202.030
Na Sodium	330.237
Ni Nickel	231.604
Pb Lead	220.353
Sb Antimony	206.833
Si Silicon	252.851
Sn Tin	189.933
Sr Strontium	421.552
Ti Titanium	334.941
V Vanadium	292.402
Zn Zinc	206.191
Y Yittrium	371.030

Other wavelengths may be substituted if they can provide the needed sensitivity and are corrected for spectral interference. Because of differences among various makes and models of spectrometers, specific instrument operating conditions cannot be provided. The instrument operating conditions herein are recommended based upon manufacturer's instrument manuals.

Table 3
Operating Conditions

Current Method Operating Conditions are as follows, these conditions may vary to optimize the instrument for different analyses:

Parameter	Radial Plasma	Axial Plasma
Resolution	Fixed	Fixed
Purge Gas Flow	Normal	Normal
Read Time (min/max sec.)	5/20	5/50
Replicates	2	2
Plasma (L/min)	15	15
Aux. (L/min)	0.5	0.3
Nebulizer Flow (L/min)	0.72	0.56
Power (watts)	1300	1450
Viewing Height (mm)	15	15

SOP NO.: MET-6010BPINES

Revision: 1

Date: 9/29/04

Page: 20 of 37

APPENDIX I

PREPARATION PROCEDURES FOR STANDARDS AND QC

RADIAL OPTIMA #1 CALIBRATION STANDARD #1 (Standard is prepared weekly or as necessary)

Metal	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	Letter ID	Nitric Acid Lot#	Hydrochloric Acid Lot #	Expiration Date	Pipet ID
<i>PQL Std</i> 1	AL	200	0.100	100	0.200	2% HNO3		A				
	BE	5			0.0050	5% HCl		B				
	CO	50			0.0500			C				
	MG	5000			5.00			D				
	SE	5			BELOW			E				
	V	50			0.0500			F				
	SB	60			0.0600			G				
	CD	5			0.0050			H				
	CU	25			0.0250			I				
	MN	15			0.0150			J				
	AG	10			0.0100			K				
	ZN	20			0.0200			L				
	AS	10			BELOW			M				
	CA	5000			5.00			N				
	FE	100			0.100			O				
	NI	40			0.0400			P				
	NA	5000			5.00			Q				
	BA	200			0.200			R				
	CR	10			0.0100			S				
	PB	5			BELOW			T				
	K	5000			5.00			U				
	TL	10			BELOW			V				
<i>PQL Std</i> 2	B	200	0.100		0.200			W				
	MO	25			0.0250			X				
	SN	500			0.500			Y				
	TI	50			0.0500			Z				
	PB	100	0.050		0.055			AA				
	AS	1000	0.010		0.110			BB				
	SE	1000	0.010		0.105			CC				
	TL	1000	0.010		0.110			DD				
	SR	1000	0.010		0.100			EE				
<i>Single std</i>	1/10-	100	0.050		0.055			AA				
	AS	1000	0.010		0.110			BB				
	SE	1000	0.010		0.105			CC				
	TL	1000	0.010		0.110			DD				
	SR	1000	0.010		0.100			EE				

RADIAL OPTIMA #1- CALIBRATION STANDARD #3 / HLCCV1 (Standard is prepared weekly or as necessary)
(CALIBRATION STANDARD #2 IS A 1/5 DILUTION OF THIS STANDARD)

Metal	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	Letter ID	Nitric Acid Lot#	Hydrochloric Acid Lot #	Expiration Date	Pipet ID
Cal Std 1		5000	4.00	200	100	2% HNO ₃		A				
		5000			100	5% HCl		B				
		5000			100			C				
		5000			100			D				
Cal Std 2		100	2.00		1.00			E				
		100			1.00			F				
		150			1.50			G				
		400			4.00			H				
Cal Std 3		200			2.00			I				
		2000	2.00		20.0			J				
		2000			20.0			K				
		50			0.500			L				
Cal Std 4		500			5.00			M				
		250			2.50			N				
		1000			10.0			O				
		500			5.00			P				
Single Metals		100	8.00		4.00			Q				
		50			2.00			R				
		50			2.00			S				
		50			2.00			T				
Single Metals		100			4.00			U				
		1000	2.00		10.0			V				
		1000	2.00		10.0			W				
		1000	1.00		5.00			X				
Single Metals		1000	1.00		5.00			Y				
		1000	1.00		5.00			Z				
		1000	1.00		5.00			AA				
		1000	1.00		5.00							

RADIAL OPTIMA #1 ICV/CCV STANDARD (Standard is prepared daily.)

Metal	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	Letter ID	Nitric Acid Lot#	Hydrochloric Acid Lot #	Pipet ID
Cal Std 1	CA	5000	2.00	200	50.0	2% HNO ₃		A			
	MG	5000			50.0	5% HCl		B			
	K	5000			50.0			C			
	NA	5000			50.0			D			
Cal Std 2	AG	100	1.00		0.500			E			
	CR	100			0.500			F			
	MN	150			0.750			G			
	NI	400			2.00			H			
Cal Std 3	ZN	200			1.00			I			
	AL	2000	1.00		10.0			J			
	BA	2000			10.0			K			
	BE	50			0.250			L			
Cal Std 4	CO	500			2.50			M			
	CU	250			1.25			N			
	FE	1000			5.00			O			
	V	500			2.50			P			
Single Elements	AS	100	4.00		2.00			Q			
	CD	50			1.00			R			
	PB	50			1.00			S			
	SE	50			1.00			T			
Single Elements	TL	100			2.00			U			
	SB	1000	1.00		5.00			V			
	SN	1000	1.00		5.00			W			
	B	1000	0.500		2.50			X			
Single Elements	MO	1000	0.500		2.50			Y			
	TI	1000	0.500		2.50			Z			
	SR	1000	0.500		2.50			AA			

RADIAL OPTIMA #1 - HLCCV2 (Standard is prepared weekly or as necessary.)

Metal	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	Letter ID	Nitric Acid Lot#	Hydrochloric Acid Lot #	Expiration Date	Pipet ID
Cal Std 2	AG	100	2.00	100	2.00	2% HNO ₃		A				
	CR	100			2.00	5% HCl		B				
	MN	150			3.00			C				
	NI	400			8.00			D				
	ZN	200			4.00			E				
Cal Std 3	AL	2000	2.00		40.0			F				
	BA	2000			40.0			G				
	BE	50			1.00			H				
	CO	500			10.0			I				
	CU	250			5.00			J				
Single Metals	FE	1000			20.0			K				
	V	500			10.0			L				
	MO	1000	1.00		10.0			M				
	PB	1000	0.800 1.00		10.0			N				
	TI	1000	1.00		10.0			O				
								P				
								Q				
								R				
								S				
								T				
								U				
								V				
								W				
								X				
								Y				
								Z				

RADIAL OPTIMA #1 CRI STANDARD

Element	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	ID Letter	Nitric Acid Lot #	Hydrochloric Acid Lot #	Expiration Date	Pipet ID
CRDL STD		Multi	0.500	500	Multi	5% HCL 2% HNO3		A				
AG		20			0.0200			B				
AS		20			BELOW			C				
BE		10			0.0100			D				
CD		10			0.0100			E				
CR		20			0.0200			F				
CO		100			0.100			G				
CU		50			0.0500			H				
MN		30			0.0300			I				
NI		80			0.0800			J				
PB		6			BELOW			K				
SB		120			0.120			L				
SE		10			BELOW			M				
TL		20			BELOW			N				
V		100			0.100			O				
ZN		40			0.0400			P				
AS		1000	0.050		0.120			Q				
PB		1000	0.050		0.106			R				
SE		1000	0.050		0.110			S				
TL		1000	0.050		0.120			T				

RADIAL OPTIMA #1 ICSAB STANDARD

Element	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	ID Letter	Nitric Acid Lot #	Hydrochloric Acid Lot #	Expiration Date	Pipet ID
Int. A Sol'n		Multi	100	1000	Multi			A				
AL		5000			500			B				
CA		5000			500			C				
FE		2000			200			D				
MG		5000			500			E				
Int. B Sol'n		Multi	10.0		Multi			F				
AG		100			1.00			G				
BA		50			0.500			H				
BE		50			0.500			I				
CD		100			1.00			J				
CO		50			0.500			K				
CR		50			0.500			L				
CU		50			0.500			M				
MN		50			0.500			N				
NI		100			1.00			O				
PB		100			1.00			P				
V		50			0.500			Q				
ZN		100			1.00			R				
								S				
								T				
								U				

RADIAL OPTIMA #1 MRL STANDARD

Element	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	ID Letter	Nitric Acid Lot #	Hydrochloric Acid Lot #	Exp. Date	Pipet ID
Cal #1		5000	0.20	1000	1.00	5% HCL		A				
		5000			1.00	2%HNO3		B				
		5000			1.00			C				
		5000			1.00			D				
Cal #2		100	0.10		0.0100			E				
		100			0.0100			F				
		150			0.0150			G				
		200			0.0200			H				
		400			0.0400			I				
		2000	0.10		0.200			J				
Cal #3		2000			0.200			K				
		1000			0.100			L				
		500			0.050			M				
		500			0.050			N				
		250			0.025			O				
		50			0.00500			P				
Cal #4		50	0.10		0.0050			Q				
		100			BELOW			R				
		50			BELOW			S				
PQL #2		200	1.00		0.200			T				
		25			0.0250			U				
		500			0.500			V				
		50			0.050			W				
		1000	0.060		0.060			X				
Single Sds		1000	0.100		0.100			Y				
		1000	0.300		0.310			Z				
		1000	0.500		0.510							
		1000	0.500		0.505							
		1000	0.050		0.055							

[illegible]

AXIAL OPTIMA #2- CALIBRATION STANDARD #3 / HLCCV1 (Standard is prepared weekly or as necessary)
(CALIBRATION STANDARD #2 IS A 1/5 DILUTION OF THIS STANDARD)

Metal	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	Letter ID	Nitric Acid Lot#	Hydrochloric Acid Lot #	Expiration Date	Pipet ID
Cal Std 1	CA	5000	2.00	200	50.0	2% HNO ₃		A				
	MG	5000			50.0	5% HCl		B				
	K	5000			50.0			C				
	NA	5000			50.0			D				
Cal Std 2	AG	100	2.00		1.00			E				
	CR	100			1.00			F				
	MN	150			1.50			G				
	NI	400			4.00			H				
	ZN	200			2.00			I				
Cal Std 3	AL	2000	2.00		20.0			J				
	BA	2000			20.0			K				
	BE	50			0.500			L				
	CO	500			5.00			M				
	CU	250			2.50			N				
	FE	1000			10.0			O				
	V	500			5.00			P				
Cal Std 4	AS	100	4.00		2.00			Q				
	CD	50			1.00			R				
	PB	50			1.00			S				
	SE	50			1.00			T				
	TL	100			2.00			U				
Single Metals	SB	1000	2.00		10.0			V				
	SN	1000	2.00		10.0			W				
	B	1000	1.00		5.00			X				
	MO	1000	1.00		5.00			Y				
	TI	1000	1.00		5.00			Z				

AXIAL OPTIMA #2 ICV/CCV STANDARD (Standard is prepared daily.)

	Metal	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	Letter ID	Nitric Acid Lot#	Hydrochloric Acid Lot #	Pipet ID
Cal Std 1	CA		5000	1.00	200	25.0	2%HNO3		A			
	MG		5000			25.0	5%HCl		B			
	K		5000			25.0			C			
	NA		5000			25.0			D			
Cal Std 2	AG		100	1.00		0.500			E			
	CR		100			0.500			F			
	MN		150			0.750			G			
	NI		400			2.00			H			
Cal Std 3	ZN		200			1.00			I			
	AL		2000	1.00		10.0			J			
	BA		2000			10.0			K			
	BE		50			0.250			L			
Cal Std 4	CO		500			2.50			M			
	CU		250			1.25			N			
	FE		1000			5.00			O			
	V		500			2.50			P			
Cal Std 4	AS		100	2.00		1.00			Q			
	CD		50			0.500			R			
	PB		50			0.500			S			
	SE		50			0.500			T			
Single Elements	TL		100			1.00			U			
	SB		1000	1.00		5.00			V			
	SN		1000	1.00		5.00			W			
	B		1000	0.500		2.50			X			
Single Elements	MO		1000	0.500		2.50			Y			
	TI		1000	0.500		2.50			Z			

AXIAL OPTIMA #2- HLCCV2 (Standard is prepared weekly or as necessary)

Metal	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	Letter ID	Nitric Acid Lot#	Hydrochloric Acid Lot #	Expiration Date	Pipet ID
Cal Std 2		100	2.00	100	2.00	2% HNO ₃		A				
		100			2.00	5% HCl		B				
CR		150			3.00			C				
MN		400			8.00			D				
NI		200			4.00			E				
ZN		2000	2.00		40.0			F				
AL		2000			40.0			G				
BA		50			1.00			H				
BE		500			10.0			I				
CO		250			5.00			J				
CU		1000			20.0			K				
FE		500			10.0			L				
V		100	4.00		4.00			M				
AS		50			2.00			N				
CD		50			below			O				
PB		50			2.00			P				
SE		100			4.00			Q				
TL		1000	1.00		10.0			R				
MO		1000	0.800		10.0			S				
PB								T				
								U				
								V				
								W				
								X				
								Y				
								Z				

AXIAL OPTIMA #2 CRI STANDARD

Element	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	ID Letter	Nitric Acid Lot #	Hydrochloric Acid Lot #	Expiration Date	Pipet ID
CRDL STD		Multi	0.500	500	Multi	5% HCL 2% HNO3		A				
AG		20			0.0200			B				
AS		20			0.0200			C				
BE		10			0.0100			D				
CD		10			0.0100			E				
CR		20			0.0200			F				
CO		100			0.100			G				
CU		50			0.0500			H				
MN		30			0.0300			I				
NI		80			0.0800			J				
PB		6			0.00600			K				
SB		120			0.120			L				
SE		10			0.0100			M				
TL		20			0.0200			N				
V		100			0.100			O				
ZN		40			0.0400			P				
								Q				

ICSA STANDARD

Element	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	ID Letter	Nitric Acid Lot #	Hydrochloric Acid Lot #	Expiration Date	Pipet ID
Int. A Sol'n		Multi	100	1000	Multi	5% HCL		A				
AL		5000			500	2% HNO3		B				
CA		5000			500			C				
FE		2000			200			D				
MG		5000			500			E				
								F				
								G				
								H				
								I				
								J				
								K				
								L				
								M				
								N				
								O				
								P				
								Q				
								R				
								S				
								T				
								U				
								V				

AXIAL OPTIMA #2 ICSAB STANDARD

Element	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	ID Letter	Nitric Acid Lot #	Hydrochloric Acid Lot #	Expiration Date	Pipet ID
Int. A Sol'n		Multi	100	1000	Multi			A				
AL		5000			500			B				
CA		5000			500			C				
FE		2000			200			D				
MG		5000			500			E				
Int. B Sol'n		Multi	10.0		Multi			F				
AG		20			0.200			G				
BA		50			0.500			H				
BE		50			0.500			I				
CD		100			1.00			J				
CO		50			0.500			K				
CR		50			0.500			L				
CU		50			0.500			M				
MN		50			0.500			N				
NI		100			1.00			O				
PB		5			0.0500			P				
V		50			0.500			Q				
ZN		100			1.00			R				
AS		10			0.100			S				
SB		60			0.600			T				
SE		5			0.0500			U				
TL		10			0.100			V				

OPTIMA INTERNAL STANDARD

Element	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	ID Letter	Nitric Acid Lot #	Hydrochloric Acid Lot #	Expiration Date	Pipet ID
Y		1000	50.0	500	100	2% HNO ₃		A				
Cs		1000	50.0	500	100	5% HCL		B				
								C				
								D				
								E				
								F				
								G				
								H				
								I				
								J				
								K				
								L				
								M				
								N				
								O				
								P				
								Q				
								R				
								S				

AXIAL OPTIMA #2 MRL STANDARD

Element	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	ID Letter	Nitric Acid Lot #	Hydrochloric Acid Lot #	Exp. Date	Pipet ID
Cal #1		5000	0.20	1000	1.00	5% HCL		A				
		5000			1.00	2% HNO3		B				
		5000			1.00			C				
		5000			1.00			D				
		100	0.10		0.0100			E				
Cal #2		100			0.0100			F				
		150			0.0150			G				
		200			0.0200			H				
		400			0.0400			I				
		2000	0.10		0.200			J				
Cal #3		2000			0.200			K				
		1000			0.100			L				
		500			0.050			M				
		500			0.050			N				
		250			0.025			O				
		50			0.00500			P				
		50	0.20		0.0100			Q				
Cal #4		100			0.0200			R				
		50			0.0100			S				
		200	1.00		0.200			T				
PQL #2		25			0.0250			U				
		500			0.500			V				
		50			0.050			W				
		1000	0.060		0.060			X				
Single Std								Y				
								Z				

173

MISCELLANEOUS STANDARDS

Type Or Std	Metal	CAS Lot #	Conc. (ppm)	Vol. (mls)	Final Vol. (mls)	Final Conc. (ppm)	Matrix	Analyst/ Date	Letter ID	Nitric Acid Lot#	Hydrochloric Acid Lot#	Expiration Date	Pipet ID
Cal Std	Li	Li:Si						SD 9/21/04	A	M1780055W	M1780054U	9/27/04	M17
	Li	M1780053Y	1000	0.050	500	0.100	2% HNO ₃ 5% HCl	SD 9/23/04	B	M1780055W	M1780054U	10/4/04	M17
	Si	M1780053X	1000	0.500		0.100			C				M17
	Sr	M1780051N	1000	0.050		0.100			D				
									E				
Cal Std	Li	M1780053Y	1000	2.00	200	10.0	2% HNO ₃ 5% HCl	SD 9/21/04	F	M1780055W	M1780054U	9/27/04	M18
	Si	M1780053X	1000	2.00		10.0		SD 9/23/04	G	M1780055W	M1780054U	10/4/04	M18
	Sr	M1780051N	1000	1.00		5.00			H				
									I				
									J				
Alloy					100		2% HNO ₃ 5% HCl	SD 9/21/04	K	M1780055W	M1780054U	9/27/04	M18
	Li	M1780053Y	1000	2.00		20.0			L				
	Si	M1780053X	1000	2.00		20.0			M				
	Sr	M1780051N	1000	1.00		10.0			N				
									O				
ICV/CCV	Li	M1780053W	1000	1.00	200	5.00	2% HNO ₃ 5% HCl	SD 9/21/04	P	M1780055W	M1780054U	9/27/04	M18
	Si	M1780053V	1000	1.00		5.00		SD 9/23/04	Q	M1780055W	M1780054U	9/28/04	M18
	Sr	M1780050V	1000	0.500		2.50		SD 9/27/04	R	M1780055W	M1780054U	9/28/04	M18
								SD 9/28/04	S	M1780055W	M1780054U	9/29/04	M18
									T				
									U				
									V				
									W				

STANDARD OPERATING PROCEDURE

for

**DETERMINATION OF MERCURY IN SOLID OR SEMISOLID WASTE BY COLD VAPOR
ATOMIC ABSORPTION SPECTROMETRY
FOR INDIANA PINES SITE**

SOP No.: MET-7471APines

Revision: 0


September 23, 2004

Approved by: 
Supervisor

9/23/04
Date


QA Coordinator

9/23/04
Date

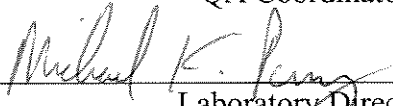
Laboratory Director

9/23/04
Date

© Columbia Analytical Services, Inc., 2004
One Mustard St, Suite 250
Rochester, NY 14609

Annual review of this SOP has been performed
and the SOP still reflects current practice.

Initials: _____ Date: _____

Initials: _____ Date: _____

Initials: _____ Date: _____

NON-CONTROLLED COPY
Will Not Be Updated

Is

Table of Contents**Page**

1. Scope and Applicability.....	3
2. Summary of Method.....	3
3. Definitions.....	3
4. Health and Safety Warnings.....	4
5. Cautions.....	4
6. Interferences.....	5
7. Personnel Qualifications.....	5
8. Equipment and Supplies.....	5
9. Procedure.....	7
9.1. Calibration and Standardization.....	7
9.2. Sample Collection.....	7
9.3. Sample Handling and Preservation.....	7
9.4. Sample Preparation	7
9.5. Sample Analysis.....	7
9.6. Troubleshooting.....	8
9.7. Data Acquisition, Calculations, and Data Reduction Requirements.....	8
9.8. Computer Hardware and Software.....	8
10. Data and Records Management.....	8
11. Quality Control and Quality Assurance.....	9
12. References.....	10

Attachments

Mercury Prep Log.....	11
Metals Training Plan Form.....	12
Metals Data Review Checklist.....	13

1. SCOPE AND APPLICABILITY

This SOP uses EPA SW-846 Method 7471A to determine the concentration of mercury in soils, sediments, bottom deposits and sludge-type materials. The range of the method is 0.2 to 10 ug/L. The range may be extended above or below the normal range by increasing or decreasing the sample size. This SOP was modified specifically for the Indiana Pines site project.

2. SUMMARY OF METHOD

A known portion of a soil sample is transferred to a hot block cup. It is digested in diluted potassium permanganate solution and oxidized for thirty minutes at 95°C. Mercury in the digested water sample is reduced with stannous chloride to elemental mercury and measured by the conventional cold vapor atomic absorption technique.

3. DEFINITIONS

- 3.1. **Calibration Blank** - A volume of reagent water acidified with the same acid matrix as in the calibration standards. The calibration blank is a zero standard and is used to auto-zero the instrument.
- 3.2. **Calibration Standard** - A solution prepared from the dilution of stock standard solutions. The CAL solutions are used to calibrate the instrument response with respect to analyte concentration.
- 3.3. **Laboratory Duplicates** - Two aliquots of the same sample taken in the laboratory and analyzed separately with identical procedures. Analyses of duplicate sample indicates precision associated with laboratory procedures, but not with sample collection, preservation, or storage procedures.
- 3.4. **Laboratory Control Sample (LCS)** - An aliquot of an ERA soil sample with a known concentration. The LCS is analyzed exactly like a sample, and its purpose is to determine whether the methodology is in control and whether the laboratory is capable of making accurate and precise measurements.
- 3.5. **Matrix Spike (MS)** - An aliquot of an environmental sample to which a known quantity of the method analyte is added in the laboratory. The MS is analyzed exactly like a sample, and its purpose is to determine whether the sample matrix contributes bias to the analytical results. The background concentrations of the analytes in the sample matrix must be determined in a separate aliquot and the measured values in the LFM corrected for background concentrations.
- 3.6. **Preparation Blank (PB)** - An aliquot of reagent water or other blank matrices that are treated exactly as a sample including exposure to all glassware, equipment, solvents, reagents, and internal standards that are used with other samples. The PB is used to determine if the method analyte or other interferences are present in the laboratory environment, reagents, or apparatus.

- 3.7. **Linear Dynamic Range (LDR)** - The concentration range over which the instrument response to an analyte is linear.
- 3.8. **Method Detection Limit (MDL)** - The minimum concentration of an analyte that can be identified, measured, and reported with 99% confidence that the analyte concentration is greater than zero.
- 3.9. **Standard Addition** - The addition of a known amount of analyte to the sample in order to determine the relative response of the detector to an analyte within the sample matrix. The relative response is then used to assess either an operative matrix effect or the sample analyte concentration.
- 3.10. **Batch** - Unit of samples prepared together on the same day, not to exceed 20 samples.

4. HEALTH AND SAFETY WARNINGS

The toxicity and carcinogenicity of each reagent used in this method has not been fully established. Each chemical should be regarded as a potential health hazard and exposure to these compounds should be minimized by good laboratory practices. Normal accepted laboratory safety practices should be followed during reagent preparation and instrument operation. Always wear safety glasses or full-face shield for eye protection when working with these reagents.

All contact with mercury should be avoided. Mercury vapor is especially toxic, causing severe respiratory tract damage. Chronic exposure to mercury through any route can produce central nervous system damage. May cause muscle tremors, personality and behavior changes, memory loss, metallic taste, loosening of the teeth, digestive disorders, skin rashes, brain damage and kidney damage. Can cause skin allergies and accumulate in the body. Repeated skin contact can cause the skin to turn gray in color. A suspected reproductive hazard; may damage the developing fetus and decrease fertility in males and females.

5. CAUTIONS

- Because of the extreme sensitivity of the analytical procedure and the presence of mercury in a laboratory environment, care must be taken to avoid extraneous contamination. Sampling devices, sample containers and plastic items should be determined to be free of mercury; the sample should not be exposed to any condition in the laboratory that may result in contamination from airborne mercury vapor.
- Samples with high organic content may required additional permanganate. Shake and add additional permanganate solution, if necessary, until the purple color persists for at least 15 minutes. Ensure that equal amounts of permanganate are added to all samples, standards and blanks

6. INTERFERENCES

- 6.1. Interferences have been reported for soils containing sulfide, chloride, copper and tellurium. Organic compounds, which have broad band UV absorbance (around 253.7 nm), are confirmed interferences. The concentration levels for interferants are difficult to define.
- 6.2. Low level mercury sample preparation, digestion, and analysis may be subject to environmental contamination if preformed in areas with high ambient backgrounds where mercury was previously employed as an analytical reagent in analyses such as chemical oxygen demand (COD).

7. PERSONNEL QUALIFICATIONS

At a minimum, personnel must have attained at least a 2-year degree in any subject and have successfully completed an Initial Demonstration of Capability after training using the Training Plan Form (found on the CAS Intranet). Training and Demonstration of Capability are in accordance with NELAC 2002 Standard.

8. EQUIPMENT AND SUPPLIES

- 8.1. Perkin Elmer FIMS - Atomic Absorption Spectrophotometer equipped with a vapor generator, quartz absorption cell and mercury hollow cathode lamp.
- 8.2. 50mL hot block cups and caps
- 8.3. 100 mL B-Cups and caps
- 8.4. Hot Block capable of maintaining a digestion temperature of 90-95°C.
- 8.5. Pipettes and graduated cylinders.
- 8.6. Mercury stock solution (1,000 mg/L) - Purchased. Store at room temperature. Dispose per manufacturer's expiration date.
- 8.7. Intermediate Stock Solution (10 mg/L) - Prepare a 1/100 dilution of the 1000mg/L Stock Solution in a volumetric flask and dilute with DI water. Acidify with 0.5 ml of concentrated HNO_3 . Store at room temperature for up to 1 week.
- 8.8. Working Solution (100 $\mu\text{g/L}$) - Prepare a 1/100 dilution of the 10mg/L Intermediate Stock Solution in a volumetric flask and dilute with DI water. Acidify with 0.5 ml of concentrated HNO_3 . Prepare fresh each day analysis is performed.
- 8.9. Calibration Standards - Prepare 0, 0.2, 0.5, 1.0, 2.0, 5.0, 10.0 $\mu\text{g/L}$ calibration curve. Transfer 0, 0.1, 0.25, 0.5, 1.0, 2.5, 5.0 mL aliquots of the 100 $\mu\text{g/L}$ working solution to a series of labeled hotblock cups. Add the appropriate amount of reagent water to bring each cup to a final volume of 5 ml. Add 5 ml of aqua regia. Loosely cap each cup. Prepare 2 blank standards to ensure sufficient volume for the analysis. The CRDL standard is prepared as the 0.2 standard.

- 8.10. ASTM Type II water
- 8.11. Concentrated Nitric Acid – Metals Grade, purchased commercially. Expires as per manufacturer's indications.
- 8.12. Concentrated Sulfuric Acid – Metals Grade, purchased commercially. Expires as per manufacturer's indications.
- 8.13. Aqua regia: Prepare immediately before use by carefully adding three volumes of concentrated HCl to one volume of concentrated HNO_3 .
- 8.14. 5% w/v Potassium Permanganate Solution - Dissolve 50 g of KMnO_4 in 1 L of reagent water. Store at room temperature for up to 6 months.
- 8.15. 12% w/v Sodium chloride-hydroxylamine chloride solution - Dissolve 120 g of NaCl and 120 g of hydroxylamine hydrochloride ($\text{NH}_2\text{OH}\cdot\text{HCl}$) in 1 L of reagent water. (Hydroxylamine sulfate $(\text{NH}_2\text{OH})_2 \text{H}_2\text{SO}_4$ may be used in place of hydroxylamine hydrochloride.) Store at room temperature for up to 6 months.
- 8.16. 1.1% Stannous chloride + 3% HCl solution - Add 11.0 g of $\text{SnCl}_2\cdot 2\text{H}_2\text{O}$ to 1 L of 3% HCl. Prepare daily.
- 8.17. The calibration blanks (ICB and CCB), prepared daily, must contain all reagents in the same concentrations and in the same volume as used in preparing the calibration solutions.
- 8.18. The preparation blank (PB) is prepared in the same manner as the calibration blank and is carried through the entire preparation scheme with each batch of samples to be analyzed.
- 8.19. With each batch of samples to be analyzed, prepare a laboratory control sample (LCS) by weighing a 0.60g portion of an ERA soil standard and place in the bottom of a 50 mL hotblock cup. The LCS must be carried through the entire sample preparation scheme.
- 8.20. Initial / Continuing Calibration Verification Standard (ICV/CCV) 3.0 $\mu\text{g/L}$ - Prepare an intermediate stock solution and working solution of 10 mg/L and 100 $\mu\text{g/L}$ using a different stock source than the calibration standards. Transfer 1.5 ml of 100 $\mu\text{g/L}$ solution (prepared daily) to a 50 ml hotblock cup. Add 3.5 ml of reagent water and 5 ml of aqua regia. Prepare 2 CCVs to ensure sufficient volume for the analysis.
- 8.21. The matrix spike sample (MS) is prepared by fortifying a 0.6g sample with 0.5 ml of 100 $\mu\text{g/L}$ CCV standard in a hotblock cup. Carry through the entire digestion and instrument procedure as a routine sample.

9. PROCEDURE

9.1. Calibration and Standardization

Calibration Standards for the initial calibration must be prepared with each daily analysis. A blank and 5 standards is required. The correlation coefficient for each calibration must ≥ 0.995 .

9.2. Sample Collection

Samples are to be collected in purchased, certified clean glass or plastic sample jars.

9.3. Sample Handling and Preservation

9.3.1. Maintain at 0-6°C from receipt until analysis.

9.3.2. Digested and analyze samples within 28 days of collection. Once digested, samples are analyzed as soon as possible.

9.3.3. Sample handling, storage, and custody procedures are in compliance with NELAC 2002 Standard.

9.4. Sample Preparation

9.4.1. Weigh 0.6g portion of a representative sample (approx. 0.2g portions from three areas of the sample) and place in the bottom of a hot block cup. Add 5 ml of reagent water and 5 ml of aqua regia. Loosely cap the sample cup.

9.4.2. Heat in the hotblock for 2 minutes at 95°C. Cool, then add 25 ml of reagent water and 15 ml of 5% potassium permanganate solution. Mix thoroughly and place in the hotblock for 30 minutes at 95°C.

9.4.3. *Note: Samples with high organic content may required additional permanganate. Shake and add additional permanganate solution, if necessary, until the purple color persists for at least 15 minutes. Ensure that equal amounts of permanganate are added to all samples, standards and blanks.*

9.4.4. Cool and add 3.0 ml of 12% sodium chloride/hydroxylamine hydrochloride solution. Add 25 ml of reagent water and the samples are now ready to be analyzed. The stannous chloride solution is added automatically by the vapor generator.

9.5. Sample Analysis

9.5.1. Analyze the standards and samples using the Perkin Elmer Flow Injection Mercury System. See Operations Manual for details.

- 9.5.2. Sample concentrations exceeding the Linear Range require sample dilution. Dilutions should be performed so that the instrument concentration will fall in the mid-range of the calibration curve.

9.6. Troubleshooting

All maintenance activities are recorded in a maintenance logbook kept for each instrument. CAS staff performs most routine maintenance and troubleshooting. Other maintenance or repairs may, or may not require factory service, depending upon the nature of the task. Typical preventive maintenance measures include, but are not limited to, the following items:

- Check gases and tubing, daily
- Check optic tubes and filter membrane for moisture before analysis

9.7. Data Acquisition, Calculations, and Data Reduction Requirements

Calculations:

From the prepared calibration curve compute sample values by comparing response with the standard curve. Calculate the mercury concentration in the sample in mg/Kg by using the formula:

$$\text{mg/Kg} = \text{Vol. (ml)/sample Wt(g)} \times 1\text{mg}/1000\text{ug} \times 1\text{L}/1000\text{ml} \times 1000\text{g}/1\text{Kg} \times C \times \text{dilution}$$

C = concentration of Hg in digestate, in ug/L

9.8. Computer Hardware and Software

- Personal Computer running Perkin Elmer AA Winlab for Window v.2.50
- Metals Analytical Review and Reporting System (MARRS) v.3.2.44
- StarLIMS v.6.11.a

10. DATA AND RECORDS MANAGEMENT

- 10.1. **Repsonsibilities** - It is the responsibility of the analyst to perform the analysis according to this SOP and to complete all documentation required for data review. Analysis and interpretation of the results are performed by personnel in the laboratory who have demonstrated the ability to generate acceptable results utilizing this SOP. Final review and sign-off of the data is performed by the department supervisor or designee.
- 10.2. **Data Flow** – Samples are entered by the Project Manager into StarLIMS on a Personal Computer running on a Novell Network. On the day that the samples are received the samples appear on a daily log printed from this computer system. The Metals Prep analyst prepares a benchsheet, digests the samples and turns the samples and digest sheet over to the ICP analyst. The samples are analyzed for metals of interest using AA software. The results are transferred to MARRS (for reporting package work) and StarLIMS for validation, reporting, and invoicing.

- 10.3. **Data Review** – Data will be reviewed by the instrument analyst and a qualified peer using a data review checklist (attached).

11. QUALITY CONTROL AND QUALITY ASSURANCE

- 11.1. Preparation Blanks must be analyzed at least once per batch of 20 or fewer samples. PB values must not exceed the MRL (method reporting limit). If the PB is out of control, fresh aliquots of the samples must be prepared and analyzed again for affected analytes after the source of the contamination has been corrected and acceptable PB values have been obtained.
- 11.2. Laboratory Control Samples assess laboratory performance against the required control limits. The control limit range is specific for each lot and is recorded on a certificate from the manufacturer. If the recovery of mercury falls outside the required control limits, the analysis is judged to be out of control, and the source of the problem should be identified and resolved before continuing analysis. Redigestion and analysis is required until acceptable LCS recovery is performed.
- 11.3. Calibration Verification Standards must immediately follow each calibration, after every tenth sample, and at the end of the sample run. Initial Calibration Verification must verify that the instrument is within $\pm 10\%$. Continuing Calibration Verification standards must confirm the calibration within $\pm 10\%$ throughout the analyses. If the recovery of mercury falls outside the required control limits, the analysis is judged to be out of control, and the source of the problem should be identified and resolved before continuing analysis. Reanalysis of any sample(s) associated with the outlying ICV or CCV standards is required. All samples must be bracketed with acceptable ICV and CCV standards.
- 11.4. Sample Matrix Accuracy and Precision are assessed based upon MS and Duplicated performance. Refer to Appendix C of the Quality Assurance Manual for frequency and QC criteria per method of analysis. If the MS is out of control and the LCS is in control, assume matrix interference and flag the associated data.
- 11.5. Method Detection Limit (MDL) - A mercury MDL must be determined annually using 7 replicates of a fortified blank solution at a concentration of 2-3 times the estimated detection limit. Practical Quantitation Limits (PQLs) are calculated from the MDL by multiplying the MDL by a factor of at least 3. The PQLs are generally used as CAS Reporting Limits. To determine the MDL, refer to 40 CFR Part 136 Appendix B.

12. REFERENCES

- *Test Methods For Evaluating Solid Waste, Physical/Chemical Methods*. USEPA SW-846, 3rd Edition, September 1994.
- *Methods For the Determination of Metals in Environmental Samples Supplement I*. USEPA/600/R-94/111, May 1994
- EPA Contract Laboratory Program, Statement of Work for Inorganic Analysis, SOW No. ILM04.0.
- Analytical Services Protocol (ASP), New York State Department of Environmental Conservation, December 1995.
- 40 CFR Part 136 Appendix B
- NELAC 2002 Standard.

Columbia Analytical Services

Analysis For: Hg EPA Method: 7470 / 7471 / 245.1 / 245.5 Analyst: _____
 In Temp: In____/ Out____ Report: Routine / ASP / Pkg.5 Date Prepped: _____
 Time In: _____ File Name: _____ Date Analyzed: _____
 Time Out: _____

Client / Submission #	Sample Number	Initial Wgt/Volume (g/ml)	Final Volume (ml)	Client / Submission #	Sample Number	Initial Wgt/Volume (g/ml)	Final Volume (ml)
1				31			
2				32			
3				33			
4				34			
5				35			
6				36			
7				37			
8				38			
9				39			
10				40			
11				41			
12				42			
13				43			
14				44			
15				45			
16				46			
17				47			
18				48			
19				49			
20					ICB	100ml DI Water	100
21					Std 0	100ml DI Water	100
22					Std 0.2*	0.20ml of 0.1ppm	100
23					Std 0.5*	0.50ml of 0.1ppm	100
24					Std 1.0*	1.00ml of 0.1ppm	100
25					Std 2.0*	2.00ml of 0.1ppm	100
26					Std 5.0*	5.00ml of 0.1ppm	100
27					Std 10.0*	10.0ml of 0.1ppm	100
28					ICV/ICV**	3.00ml of 0.1ppm	100
29					LCSW/MS**	1.00ml of 0.1ppm	100
30					CRDL*	0.20ml of 0.1ppm	100

Lot # of Reagents Used:

HNO3: _____ H2SO4: _____ HCL: _____ K2S2O8: _____ KMnO4: _____

SnCl2: _____ NaCl: _____ NH2OH-HCL: _____ LCSS CAS Lot#: _____

LCSS ERA Lot#: _____

*Source Standard: (Vendor/Lot #) _____

**ICV Standard: (Vendor/Lot #) _____

*(10ppm stock Lot#: _____)

**(10ppm stock Lot#: _____)

*(0.1ppm working std Lot#: _____)

**(0.1ppm working std Lot#: _____)

Comments/Problems: _____

Batch Name: _____

Metals Instrument Analysis Training Plan

Procedure: _____

SOP: _____ Revision: _____ Date: _____

Trainee: _____

1. Read SOP Trainer: _____ Trainee: _____ Date: _____

2. Demonstrated understanding of the scientific basis of the analysis including:
-the chemical and physical principals behind the measurement by the instrument
Trainer: _____ Trainee: _____ Date: _____

3. Demonstrated familiarity with related SOPs
-ADM-BATCHSEQ -ADM-PCAL -ADM-SIGFIG
-ADM-DATAENTRY -ADM-DIL -ADM-SPSR
-ADM-MDL -ADM-DREV -ADM-TRANDOC
Trainer: _____ Trainee: _____ Date: _____

4. Observe performance of SOP
-standard and reagent prep and documentation – including pipet used
-instrument power up and warm-up
-instrument set-up, daily maintenance and checks
-use and loading of autosampler
-sample analysis including:
-calibration
-sample dilution
-software command of instrument
-use of QC samples and QC criteria
-common troubleshooting
-instrument logbook use
-data reduction, reporting, and review
Trainer: _____ Trainee: _____ Date: _____

5. I have read, understood and agree to perform the most recent version of the SOP:
Signature: _____ Date: _____

6. Perform SOP with supervision
- including all items in 4.
Trainer: _____ Trainee: _____ Date: _____

7. Independent performance of the SOP
-all of the item listed in 4
-IDC (4 mid-range standards performed before client samples are analyzed)
-attach IDC certificate, raw data, and summary spreadsheet.
Trainer: _____ Trainee: _____ Date: _____



METALS DEPARTMENT DATA QUALITY CHECKLIST

Data File: _____

Run Date: _____

Instrument: _____

Methods Used:

ICP- 200.7 // 6010B // ASP/CLP // NIOSH
 GFAA- EPA 200 Series // SW846 // ASP/CLP
 CVAA- EPA 200 Series // SW846 // ASP/CLP

Batch ID / Metals Reviewed: _____

Yes	No	NA		Yes	No	NA
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	1. Holding Times met method requirements? ICP / GFAA- 6mths from sampling to analysis Hg- 28 days from sampling to analysis (26 days from VTSR)	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	2. ICAL met method requirements? Correlation Coefficient > or = 0.995 ICP High Check = 95-105%	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	3. ICV acceptable? ICP: 200.7 = 95-105% ; NIOSH / 6010B / ASP/CLP = 90-110% GFAA: EPA 200 Series / SW846 / ASP/CLP = 90-110% Hg: EPA 200 Series = 95-105% ; SW846 / ASP/CLP = 90-110%	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	4. CCVs acceptable? Analyzed per 10 samples? ICP: 200.7 / 6010B / ASP/CLP / NIOSH = 90-110% GFAA: EPA 200 Series = 90-110% ; SW846 / ASP/CLP = 80-120% Hg: EPA 200 Series = 90-110% ; SW846 / ASP/CLP = 80-120%	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	5. CCBs acceptable? Analyzed per 10 samples? Concentrations < RL	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	6. Method Blank results < RL?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	7. LCS recoveries within QC limits? ICP: 200.7 = 85-115% ; 6010B / ASP/CLP / NIOSH = 80-120% GFAA: EPA 200 Series = 85-115% ; SW846 / ASP/CLP = 80-120% Hg: EPA 200 Series = 85-115% ; SW846 / ASP/CLP = 80-120% LCSS (soil) Certificate of Analysis QC limits per manufacturer	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	8. All sample concentrations within LR?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	9. MS recoveries within QC limits? ICP: 200.7 = 70-130% ; 6010B / ASP/CLP = 75-125% GFAA: EPA 200 Series / SW846 / ASP/CLP = 75-125% Hg: EPA 200 Series = 70-130% ; SW846 / ASP/CLP = 75-125%	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	10. Duplicate RPD within QC limits? 20% for RPD shall be used for samples > or = 5 times the RL. RL shall be used for samples < 5 times the RL.	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	11. Is GFAA Post Digest Spike within 85-115%?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	12. Dilution factors verified and calculated correctly?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	13. Bench Sheet complete, initials, date, and time:	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	•Are standards and reagents traceable?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	•Is unused space on the sheet crossed out?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

Analyst: _____

Date: _____

Peer Review: _____

Date: _____

COMMENTS:

**Comments must be provided for any items noted above as "No"

STANDARD OPERATING PROCEDURE
for
**METALS DIGESTION,
SOILS, SEDIMENTS, AND SLUDGE FOR ICP AND GFAA ANALYSIS**

SOP No.: MET-3050B

Revision: 3

January 15, 2003

Approved by: _____
Supervisor

_____ Date

_____ QA Coordinator

_____ Date

_____ Laboratory Manager

_____ Date

© Columbia Analytical Services, Inc., 2003

1 Mustard Street, Suite 250

Rochester, New York 14609

Annual review of this SOP has been performed
and the SOP still reflects current practice.

Initials: _____ Date: _____

Initials: _____ Date: _____

Initials: _____ Date: _____

DOCUMENT CONTROL

NUMBER: _____

Initials: _____ Date: _____

1 SCOPE AND APPLICATION

Method 3050 is an acid digestion procedure used to prepare matrices such as soils, sludges, or sediments for analysis by ICP or graphite furnace atomic absorption.

2 METHOD SUMMARY

A representative aliquot of sample is digested in nitric acid and hydrogen peroxide. Hydrochloric acid is used as a final reflux acid for ICP analyses. Nitric Acid is used as the final reflux acid for most Graphite Furnace analyses.

3 DEFINITIONS

- 3.1 **Laboratory Duplicates** - Two aliquots of the same sample taken in the laboratory and analyzed separately with identical procedures. Analyses of duplicates and indicates precision associated with laboratory procedures, but not with sample collection, preservation, or storage procedures.
- 3.2 **Laboratory Control Sample Soil (LCSS)** - An aliquot of a soil to which known quantities of the method analytes are added by an outside vendor. The LCSS is analyzed exactly like a sample, and its purpose is to determine whether the methodology is in control and whether the laboratory is capable of making accurate and precise measurements.
- 3.3 **Matrix Spike** - An aliquot of an environmental sample to which known quantities of the method analytes are added in the laboratory. The matrix spike is analyzed exactly like a sample, and its purpose is to determine whether the sample matrix contributes bias to the analytical results.
- 3.4 **Preparation Blank (PB)** - An aliquot of reagent water or other blank matrices that are treated exactly as a sample including exposure to all glassware, equipment, solvents, reagents, and internal standards that are used with other samples. The PB is used to determine if method analytes or other interferences are present in the laboratory environment, reagents, or apparatus.
- 3.5 **Digestion Batch** - A digestion batch is no more than 20 samples of the same matrix digested as a unit per day.

4 INTERFERENCES

- 4.1 See appropriate analysis SOP for applicable interferences

5 SAFETY

Nitric and Hydrochloric acids are extremely corrosive. Care should be taken while working with these chemicals. Personal protective equipment including safety glasses (with side shields), gloves, and lab coat shall be worn when handling samples or reagents.

6 SAMPLE COLLECTION, PRESERVATION AND STORAGE

For non-aqueous samples, glass or plastic sample containers are acceptable. Samples are analyzed within 6 months of sample collection. Additional sample handling policies and procedures are in SMO-GEN.

7 APPARATUS AND EQUIPMENT

- 7.1 250 and 100 mL beakers
- 7.2 Ribbed watch glasses
- 7.3 Hot plates
- 7.4 Graduated cylinders
- 7.5 Eppendorf Pipettors
- 7.6 Funnels
- 7.7 Mortar and pestle
- 7.8 Tongue depressors
- 7.9 Filter paper
- 7.10 Hot Block Digester with ETR-3200 Controller by Environmental Express, LTD.
- 7.11 Graduated block digester ribbed watch glasses
- 7.12 Block Digester Filters.
- 7.13 CPI MOD Block Digester

8 PREVENTIVE MAINTENANCE

- 8.1 All hoods in the Metals Prep Lab are wiped down once a week with DI water. The tops of all digestion hot plates are wiped down daily.

9 STANDARDS, REAGENTS, AND CONSUMABLE MATERIALS

- 9.1 Reagent water - ASTM Type II deionized water. Reagent water must be interference free.
- 9.2 Concentrated nitric acid (Baker Instra-Analyzed 69-70%): Acid should be demonstrated to be free of impurities at levels which would interfere with sample determinations. Store at room temperature in the dark. Expires per manufacturer's indications or one year from receipt, whichever is sooner.
- 9.3 Concentrated hydrochloric acid (Baker Instra-Analyzed 36.5-38%): Acid should be demonstrated to be free of impurities at levels which would interfere with sample determinations. Store at room temperature. Expires per manufacturer's indications or one year from receipt, whichever is sooner.

- 9.4 Hydrogen peroxide (30%) - H_2O_2 . Purchased commercially. Should be demonstrated to be free of impurities at levels which would interfere with sample determinations. Store at room temperature. Expires upon manufacturer's indications or 1 year from receipt, whichever is sooner.
- 9.5 ERA Soil Laboratory Control Sample (LCSS) - Concentrations and Performance Acceptance Limits distributed through vendor. Store at room temperature. Expires upon manufacturer's indications or 1 year from receipt, whichever is sooner.
- 9.6 Metals spiking solutions - Purchased commercially. See Table 1. Store at room temperature. Stocks expires upon manufacturer's indications or 1 year from receipt, whichever is sooner. Solutions prepared from stocks expire 6 months from preparation.

10 RESPONSIBILITIES

- 10.1 It is the responsibility of the analyst to perform the analysis according to this SOP and to complete all documentation required for data review. Analysis and interpretation of the results are performed by personnel in the laboratory who have demonstrated the ability to generate acceptable results utilizing this SOP. Final review and sign-off of the data is performed by the department supervisor or designee.

11 PROCEDURES

11.1 HOT PLATE

- 11.1.1 Mix the sample thoroughly to achieve homogeneity using a tongue depressor or the mortar and pestle.
- 11.1.2 Weigh (to the nearest 0.01g) 1.00g to 1.50g of sample into a 250 or 100 mL beaker. For sludges and sediments that have a high moisture content, use more sample. The goal is to use about 1g of dry weight sample. At this point add the appropriate spiking solutions (see Table 1) directly onto the designated spike sample prior to addition of reagents.
- 11.1.3 Unless specified by project or state requirements, the addition of acid should be as follows: Add 10ml of 1:1 HNO_3 , cover with a ribbed watch glass and reflux for 15 minutes. The sample temperature should be 90-95°C. Allow the sample to cool, then add 5ml of concentrated HNO_3 , cover and reflux for 30 minutes. Repeat the addition of 5ml of HNO_3 and reflux to 5 mLs. Do not allow the sample to go to dryness. CAUTION: Do not boil. Antimony is easily lost by volatilization.

11.1.4 Cool the sample and add 2ml of DI and 3ml of 30% H_2O_2 . Cover and heat to start the peroxide reaction. Care must be taken to ensure that losses do not occur due to excessive effervescence. Heat until effervescence subsides and cool the beaker.

11.1.5 If the effervescence does not subside, add 3 mLs of hydrogen peroxide with warming to each of the samples (including blanks and LCSs) in the batch. If necessary, continue to add 30% H_2O_2 in 1ml aliquots with warming until the effervescence is minimal, or until the general sample appearance is unchanged. Do not add more than 10ml of 30% H_2O_2 .

11.1.6 If the sample is being prepared for analysis by ICP, add 10 mL 1:1 HCL. If the sample is being prepared for analysis by Graphite Furnace no HCl is added.

11.1.7 Cover and reflux the ICP samples for 15 minutes without boiling. Allow to cool.

11.1.8 Prepare filters by rinsing with 1:1 nitric acid and DI.

11.1.9 All samples are diluted to 100 mLs with DI. Quantitatively transfer the digestate to a graduated cylinder by pouring the sample through a prepared filter into the cylinder and rinsing the beaker and watch glass with DI into the filter. Rinse the filter with DI. Bring to volume with DI.

11.2 HOT BLOCK DIGESTOR

11.2.1 Set the temperature on the Block Digestor to a temperature that brings the sample temperature to 90-95°C without boiling.

11.2.2 The Hot Block is on a timer which can be set to turn on and off whenever necessary. To set timer press the timer button and choose the days M-F (Monday through Friday). Then choose the hour and minutes to start and stop the Block Digestor.

11.2.3 Label graduated hot block digester sample cups with appropriate sample IDs for digestion.

11.2.4 Mix the sample thoroughly to achieve homogeneity using a tongue depressor or the mortar and pestle.

11.2.5 Weigh (to the nearest 0.01g) 1.00g to 1.50g of sample into labeled digester sample cup. For sludges and sediments that have a high moisture content, use more sample. The goal is to use about 1g of dry weight sample. At this point add the appropriate spiking solutions (see Table 1) directly onto the designated spike sample prior to addition of reagents.

- 11.2.6 Unless specified by project or state requirements, the addition of acid should be as follows: Add 10ml of 1:1 HNO₃ and for ICP only add 1.5 mL of 1:1 HCl, cover with reflux cap and reflux for 15 minutes. The sample temperature should be 90-95°C. Allow the sample to cool, then add 5ml of concentrated HNO₃, cover and reflux for 30 minutes. Repeat the addition of 5ml of HNO₃ and reflux to 5 mLs. Do not allow the sample to go to dryness. CAUTION: Do not boil. Antimony is easily lost by volatilization.
- 11.2.7 Cool the sample and add 2ml of DI and 3ml of 30% H₂O₂. Cover and heat to start the peroxide reaction. Care must be taken to ensure that losses do not occur due to excessive effervescence. Heat until effervescence subsides and cool the sample cup.
- 11.2.8 If the effervescence does not subside, add 3 mLs of hydrogen peroxide with warming to each of the samples (including blanks and LCSs) in the batch. If necessary, continue to add 30% H₂O₂ in 1ml aliquots with warming until the effervescence is minimal, or until the general sample appearance is unchanged. Do not add more than 10ml of 30% H₂O₂.
- 11.2.9 If the sample is being prepared for analysis by ICP, add 10 mL 1:1 HCL. If the sample is being prepared for analysis by Graphite Furnace no HCl is added.
- 11.2.10 Cover and reflux the ICP samples for 15 minutes without boiling. Allow to cool.
- 11.2.11 Prepare filters by rinsing with 1:1 nitric acid and DI.
- 11.2.12 All samples are diluted to 100 mLs with DI. Quantitatively transfer the digestate to a graduated cylinder by pouring the sample through a prepared filter into the cylinder and rinsing the beaker and reflux cap with DI into the filter. Rinse the filter with DI. Bring to volume with DI. Pour into a labeled B-cup.

12 QA/QC REQUIREMENTS

- 12.1 Each day, digest one laboratory control sample (LCS) per digestion batch, or per 20 samples, or per EPA SDG group, whichever is more frequent. Use the appropriate solid laboratory control sample (LCSS) for soils analysis.
- 12.2 Each day, digest one blank per digestion batch, or per 20 samples, or per EPA SDG group, whichever is more frequent. Use D.I. water and follow the digestion procedures.
- 12.3 Each day, prepare one duplicate and one spiked sample with each digestion batch, or per twenty samples, or per EPA SDG group, whichever is more frequent. At times, specific samples will be assigned as duplicates of spikes depending on client requirements.

12.4 Matrix spikes are prepared by adding the appropriate volume of spiking solution (See Table 1).

12.5 See appropriate analysis SOP for applicable QC limits and corrective action.

13 DATA REDUCTION AND REPORTING

13.1 Digestion logs are used to record all sample volumes, spike volumes, etc. The Manufacturer's lot number for the reagents used are added to the digestion log (see attached digestion log benchsheet).

13.2 Reporting and method performance is discussed in the appropriate analysis SOP. Data review is discussed in ADM-DREV.

14 METHOD PERFORMANCE

Reporting limits are based upon an MDL study performed according to ADM-MDL and filed in the MDL binders in the QA office.

15 WASTE MANAGEMENT AND POLLUTION PREVENTION

15.1 Reagents are prepared upon an as-needed basis in small quantities. Minimum sample volumes are used during analysis.

15.2 Acidic waste is poured down the drain with copious amounts of water.

15.3 Samples with analyte concentrations exceeding TCLP regulatory limits are disposed of as hazardous waste. Others are dumped down the drain with plenty of water. See SMO-SPLDIS.

16 CORRECTIVE ACTION FOR OUT OF CONTROL DATA

If data is produced that is out of control, the samples are to be re-analyzed with in-control QA whenever possible. See corrective actions in Section 12 of this SOP and in the applicable Figures in Section 12 of the Quality Assurance Manual.

17 CONTINGENCIES FOR HANDLING OUT OF CONTROL OR UNACCEPTABLE DATA

If data is produced that is out of control and is not to be re-analyzed due to sample volume restrictions, holding times, or QC controls can not be met, follow the procedures in Section 15 of the Quality Assurance Manual.

18 REFERENCES

“Test Methods For Evaluating Solid Waste, Physical/Chemical Methods”. EPA SW846, Third Edition, December 1996.

19 TRAINING OUTLINE

- 19.1 Read current SOP and applicable methodologies. Demonstrate a general understanding of the methodology and chemistry. Follow policies in ADM-TRANDOC.
- 19.2 Observe Sample Preparation.
- 19.3 Participate in the methodology, documentation, and data reduction with guidance.
- 19.4 Complete a Training Plan Form for the procedure.
- 19.5** Show Initial Demonstration of Capability (IDC) by independently preparing and digesting four LCSs, or equivalent, according to the test method either concurrently or over a period of days. If recovery is within acceptable limits, complete IDC certification form, and Training Plan forms and file with QA. Continued capability shall be demonstrated annually using PE results, a single blind, or a new 4 replicate study.

20 METHOD MODIFICATIONS

None

21 INSTRUMENT-SPECIFIC ADDENDUM

Not Applicable

22 ATTACHMENTS

Table 1 Spike Concentrations
Digestion Log Benchsheets
SW846 Method 3050 Flow Chart

23 CHANGES FROM PREVIOUS REVISION

- Added Hot Block digestion procedures (11) and associated items to Apparatus and Equipment (7)
- Added sections 14, 16, 17, and 20 for NELAP compliance
- Changed the amount of time to reflux sample from 10-15 minutes to just 15 minutes after the first addition of acid (11).

Table 1 Spiking Concentrations for LCS and MS Samples

SPIKE SOLUTION A		1.00ml Spk A to Final Vol of 100ml	
<i>Metal</i>	<i>Conc. (ug/mL)</i>	<i>Metal</i>	<i>Conc. (ug/mL)</i>
AL	200	NI	50
AS	4	SE	1
BA	200	AG	5
BE	5	TL	200
CD	5	V	50
CR	20	ZN	50
CO	50	B	100
CU	25	CA	200
FE	100	MG	200
PB	50	NA	2000
MN	50	K	2000

SPIKE SOLUTION B		1.00ml Spk B to Final Vol of 100ml	
<i>Metal</i>	<i>Conc. (ug/mL)</i>	<i>Metal</i>	<i>Conc. (ug/mL)</i>
SB	50	TI	50
MO	50	-	-

INDIVIDUAL METALS		0.10ml Spk. to Final Volume of 100ml		INDIVIDUAL METALS		0.5ml Spk. to Final Volume of 100ml	
<i>Metal</i>	<i>Conc. (ug/mL)</i>			<i>Metal</i>	<i>Conc. (ug/mL)</i>		
SE	1000			SN	1000		

SPIKE #4 Furnace Spike		1.00ml Spk #4 to Final Vol of 100ml	
<i>Metal</i>	<i>Conc. (ug/mL)</i>	<i>Metal</i>	<i>Conc. (ug/mL)</i>
AS	4	SB	10
PB	2	TL	5
SE	1	CU	0.5

CAS-Rochester ICP Soil Digestion Log

6010B/846 // 200.7/136 // ASP/CLP4.1

Batch ID:

Analyst: _____

Date: _____

Spike Witness / Lot Approval: _____

Prep Method: SW846 3050 // CLP

Batch Temp:

Digest: Initial // Redigest of: _____

Report Type: Routine // ASP // Pkg5

Submission / Order #	Initial Wgt. (g)	Final Vol (ml)	Initial Color / Texture	Final Color / Clarity	Metals	Spike Vol (ml)
1						
2						
3						
4						
5						
6						
7						
8						
9						
10						
11						
12						
13						
14						
15						
16						
17						
18						
19						
20						
21						
22						
23						
24						

Spiking Standards / Reagent Lot #:

Spike A,B: _____ Spike #4: _____
 TCLP Spk: _____ TCLP Ba: _____
 Se Std: _____ Sn Std: _____
 HNO3: _____ HCL: _____
 H2O2: _____ LCSS: _____

Color / Clarity Key:

Color: C = Colorless ; Y = Yellow ; B = Brown

BL = Black ; G = Grey ; W = White

Clarity: CDY = Cloudy ; CLR = Clear ; OP = Opaque

Texture: F = Fine ; M = Medium ; CS = Coarse ; NAQ = Non Aqueous

01

Comments: _____



CAS-Rochester Furnace Soil Digest Log

Analyst: _____ Date: _____ Spike Witness: _____ Approval: _____

Mod: SW846 3050 // CLP

Digest: Initial Digestion // Redigestion of _____

Report: Routine // ASP // Pkg.5

6010B/846 // 200.7/136 // ASP/CLP

	Client/ Order #	Initial wgt(g)	Final vol(ml)	Initial Color/Clarity	Final Color/Clarity	Metals	Spike vol(ml)
1							
2							
3							
4							
5							
6							
7							
8							
9							
10							
11							
12							
13							
14							
15							
16							
17							
18							
19							
20							
21							
22							
23							
24							
25							
26							

Spiking Standards Lot # or Prep Date: _____

Spike #4 _____

Spike Lot #:

LCSS _____


HNO₃ _____H₂O₂ _____

Comments/Problems: _____

0091

SW846 Method 3050 Flow Chart

Non-CLP Soils, Sediments and Sludges

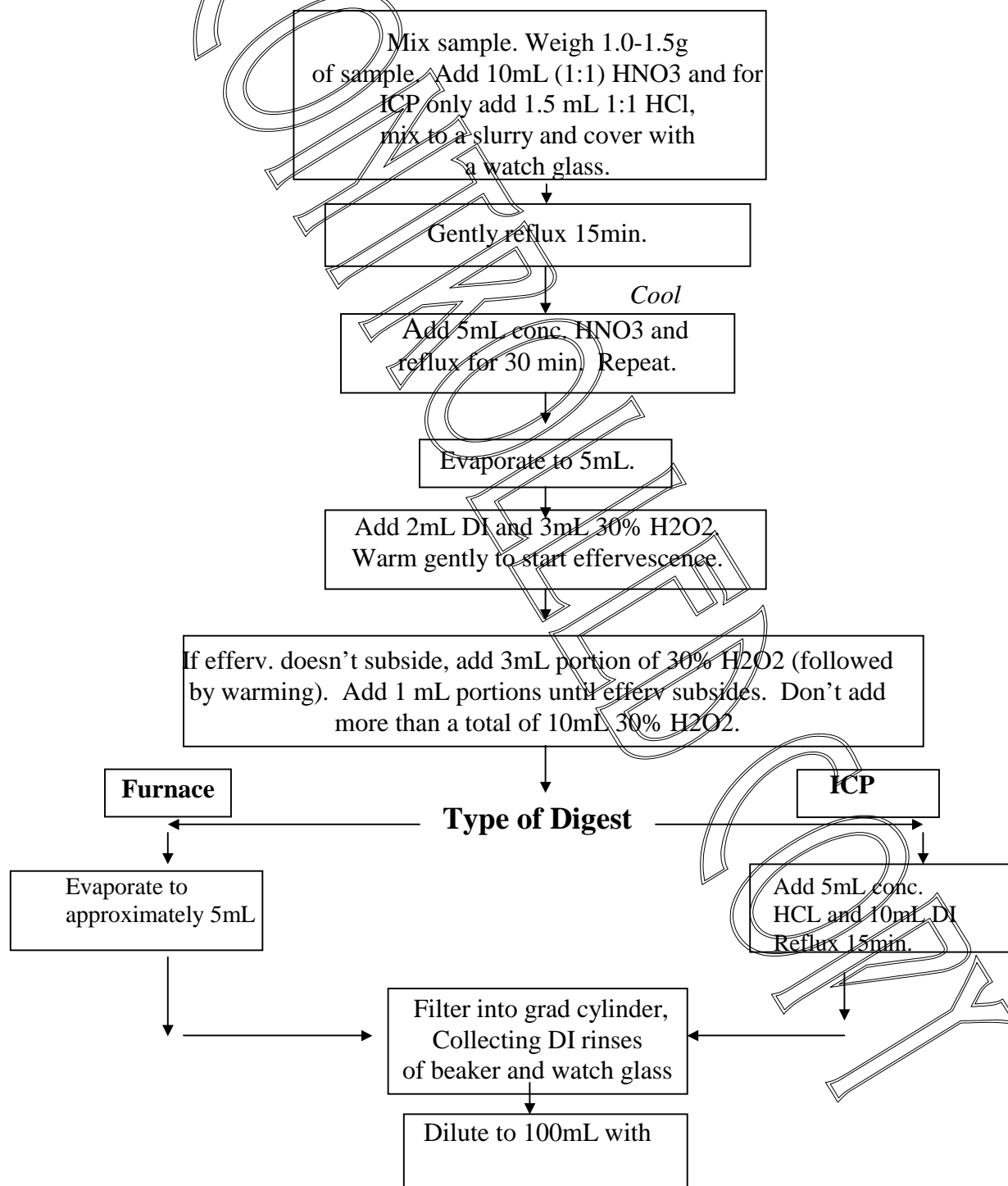

DI water.

UNCONTROLLED COPY

SOP No.: MET-GFAA

Revision: 3.0

Date: ~~11/01/01~~ 9/27/01 ^{9/27/02}

Page 1 of 11

STANDARD OPERATING PROCEDURE

DETERMINATION OF TRACE METALS BY GRAPHITE FURNACE ATOMIC ABSORPTION SPECTROMETRY (GFAA)

SOP Code: MET-GFAA

Revision: 3.0

September 27, 2001

Approved by:

Christopher Kuzel
Department Manager

11/1/01
Date

Lisa Ryan
Quality Assurance Coordinator

11/1/01
Date

Michael K. Perry
Laboratory Director

11/1/01
Date

© Columbia Analytical Services, Inc.
1 Mustard Street, Suite 250
Rochester, New York 14609

Annual review of this SOP has been performed
and the SOP still reflects current practice.

Initials: DCB Date: 12/20/02

Initials: DCB Date: 12/19/03

Initials: DCB Date: 12/14/04

DOCUMENT CONTROL

NUMBER: MET-002

Initials: SpC Date: 4/2/01

1.0 SCOPE AND APPLICATION

- 1.1 This procedure describes the procedure for the analysis of soil, sludge and water digestates by graphite furnace atomic absorption (GFAA) spectrometry. Typically, this procedure is applicable to the analytes and EPA methods listed in Table 1. Other elements may be determined when reference is made to the applicable published method. All analytical methods used are in accordance with EPA methods from SW-846 and EPA 200 Series, the EPA Contract Laboratory Program (CLP) statement of work (SOW), and the NYSDEC Analytical Services Program (ASP).
- 1.2 The Practical Quantitation Limits (PQL) are listed in Table 1. The reported PQL may be adjusted if required for specific project requirements, however, the capability of achieving other reported PQLs must be demonstrated. Results may be reported to the Instrument Detection Limits (IDLs) upon request. IDLs are updated quarterly. The Method Detection Limits (MDLs) are updated annually and are available upon request.

2.0 METHOD SUMMARY

- 2.1 Prior to analysis, samples must be digested using appropriate sample preparation methods. A representative aliquot of sample is prepared as described in the applicable digestion SOP. Refer to the following Metals Digestion SOPs:

MET-3005A	Metals Digestion, Waters, Total Recoverable and Dissolved for ICP
MET-3010A	Metals Digestion, Waters for ICP
MET-3020A	Metals Digestion, Waters for GFAA
MET-3050B	Metals Digestion, Soils, Sediments and Sludges for ICP and GFAA
MET-CLP	Metals Digestion, Waters and Soils for CLP

- 2.2 The digestate is analyzed for the element(s) of interest, using GFAA conditions (See Instrument Specifications by Metal Manual in the GFAA lab) for the element to be determined. Absorbance is measured as a function of element concentration.
- 2.3 For GFAA analyses by CLP procedures, see the applicable CLP SOW or ASP.

3.0 DEFINITIONS

- 3.1 **Analytical Sequence** - Samples are analyzed in a set referred to as an analytical sequence. The sequence begins with instrument calibration followed by analysis of sample digestates interspersed with analysis of calibration verification standards.
- 3.2 **Initial Calibration Verification (ICV)** - ICV solutions are made from a stock solution which is different from the stock used to prepare calibration standards and is used to verify the validity of the standardization.

- 3.3 Matrix Spike (MS)** - In the matrix spike analysis, predetermined quantities of standard solutions of certain analytes are added to a sample matrix prior to sample digestion and analysis. The purpose of the matrix spike is to evaluate the effects of the sample matrix on the methods used for the analyses. Percent recoveries are calculated for each of the analytes detected.
- 3.4 Duplicate Sample (DUP)** - A laboratory duplicate. The duplicate sample is a separate field sample aliquot that is processed in an identical manner as the sample proper. The relative percent difference between the samples is calculated and used to assess analytical precision.
- 3.5 Method Blank** - The method blank is an artificial sample designed to monitor introduction of artifacts into the process. The method blank is carried through the entire analytical procedure.
- 3.6 Continuing Calibration Verification Standard (CCV)** - A standard analyzed at specified intervals and used to verify the ongoing validity of the instrument calibration.
- 3.7 Instrument Blank (CCB)** - The instrument blank (also called continuing calibration blank) is a volume of blank reagent of composition identical to the digestates. The purpose of the CCB is to determine the levels of contamination associated with the instrumental analysis.

4.0 INTERFERENCES

Interferences are dealt with through the use of matrix modifiers (commonly Ni and Pd) and post digestion spikes. Detailed discussion of interferences may be found in the applicable EPA method.

Interferences from contaminated reagents must be eliminated. The purity of acids must be established by the laboratory as being high enough to eliminate the introduction of contamination above the Method Detection Limit.

5.0 SAFETY

Normal precautions as per the CAS EH&S Manual are to be followed. In addition, because acids are used in the procedure, there is a danger of exposure to corrosives. Sufficient care must be taken in handling acidic solutions. Safety glasses must be worn while preparing and handling the solutions. Gloves and a laboratory coat should be worn while handling samples, acids, and sample digestates.

6.0 SAMPLE COLLECTION, CONTAINERS, PRESERVATION, AND STORAGE

- 6.1 Either glass or plastic sample containers may be used.
- 6.2 All aqueous samples are preserved with nitric acid to a pH of <2.
- 6.3 Soil and aqueous samples for GFAA analyses by CLP procedures are stored at 0-6 °C. from time of receipt until digestion. Soil and aqueous samples for GFAA analyses by other procedures are stored at ambient temperature from time of receipt until digestion.
- 6.4 Holding time for samples for GFAA analyses is 6 months (sample collection to digestion).
- 6.5 Samples are received in the GFAA analysis laboratory as 0.5% nitric acid digestates. Sample digestates are stored in labeled plastic B-cups or Hot Block digestion vessels.

7.0 APPARATUS AND EQUIPMENT

- 7.1 Graphite furnace atomic absorption spectrophotometer (AAS). See Appendix A for element-specific instrument parameters.
- 7.2 Hollow Cathode Lamp (HCL) or Electrodeless Discharge Lamp (EDL) for each metal analyzed by this procedure.
 - 7.2.1 Electrodeless Discharge Lamp power supply.
- 7.3 100-1000uL Eppendorfs
- 7.4 2 ml Beaker cups compatible with the AAS autosampler.
- 7.5 Volumetric flasks of suitable precision and accuracy.

8.0 PREVENTIVE MAINTENANCE

All maintenance activities are recorded in a maintenance logbook kept for each instrument. Most routine maintenance and troubleshooting is performed by CAS staff. Other maintenance or repairs may, or may not require factory service, depending upon the nature of the task. Typical preventive maintenance measures include, but are not limited to, the following items:

- Cleaning the quartz windows
- Changing the graphite tubes

- Cleaning the capillary tube
- Change filter in Separator Trap annually or if necessary
- Inspection and cleaning of electrodes, shroud, and cells as needed.

9.0 STANDARDS, REAGENTS, AND CONSUMABLE MATERIALS

9.1 Concentrated Nitric Acid – Metals Grade or higher to eliminate the introduction of contamination above the method detection limit.

9.2 Matrix Modifiers

Palladium Modifier Dilute 10 mls Palladium Nitrate (1%) plus 1.0 ml Magnesium Nitrate (2%) with 100.0 ml DI. Expires within 6 months at room temperature. Used for arsenic, selenium, and thallium.

Ammonium Dihydrogen Phosphate (purchased) Dilute 1.0 g to 100 mls DI. Expires 6 months at room temperature. Used for lead analysis.

9.3 Standards

9.3.1 1000 ppm Stock Standards (AA Grade) - Commercially available certified solutions.

9.3.2 The GFAA Calibration Stock Standard is made by pipetting the volume of the 1000 ppm As, Pb, Se, and Tl stock standards shown in the table below, plus 0.50 ml of concentrated nitric acid into a 100 ml volumetric flask and diluting to volume with DI water. Prepare Calibration Stock Standard weekly.

Analyte	ml of Stock Std. (1000 ppm)	Add 0.50 ml conc. HNO ₃ and dilute to	Final Concentration (mg/L)
Arsenic	0.50	100 ml	5.0
Lead	0.45	100 ml	4.5
Selenium	0.50	100 ml	5.0
Thallium	0.50	100 ml	5.0
Lead - DW	0.30	100 ml	3.0

9.3.3 The GFAA Calibration Working Standard is prepared by diluting 1.0 ml of the GFAA Calibration Stock Standard and 0.5 ml of concentrated nitric acid to 100 ml with DI water. Final concentrations range from 0.030 to 0.050 mg/L. Prepare calibration working standard fresh each day. *LMR 8/23/02*

insert table > see attached

9.3.4 The GFAA Initial and Continuing Calibration Verification (ICV and CCV) Stock Standard is made by pipetting the volume of the 1000 ppm As, Pb, Se, and Tl

9.3.3 Continued...

"Prepare the calibration working standard fresh each day." This standard is loaded into the autosampler (typically location 38) so that the instrument may auto-dilute the working standard to prepare the calibration standards at four concentration levels intended for the initial curve (blank and 4 standards). The instrument makes the following dilutions according to the analyte being analyzed:

Element	Working Standard (ug/L)	Volume of Working Std. (uL)	Final Volume (uL)	Final Concentration (ug/L)
Arsenic	50	4	20	10
	50	8	20	20
	50	12	20	30
	50	20	20	50
Lead	30	3	30	3
	30	9	30	9
	30	15	30	15
	30	30	30	30
Antimony	50	6	30	10
	50	12	30	20
	50	18	30	30
	50	30	30	50
Selenium	50	3	30	5
	50	9	30	15
	50	18	30	30
	50	30	30	50
Thallium	50	4	20	10
	50	8	20	20
	50	12	20	30
	50	20	20	50
Copper	10	4	20	2
	10	8	20	4
	10	12	20	6
	10	20	20	10

stock standards shown in the table below, plus 0.50 ml of concentrated nitric acid into a 100 ml volumetric flask and diluting to volume with DI water. Prepare the ICV / CCV Stock Standard from a second source, weekly.

Analyte	ml of Stock Std. (1000 ppm)	Add 0.50 ml conc. HNO ₃ and dilute to	Final Concentration (mg/L)
Arsenic	0.25	100 ml	2.5
Lead	0.20	100 ml	2.0
Selenium	0.25	100 ml	2.5
Thallium	0.25	100 ml	2.5
Lead - DW	0.14	100 ml	1.4

9.3.5 The GFAA CCV Working Standard is prepared by diluting 1.0 ml of the GFAA CCV Stock Standard and 0.5 ml of concentrated nitric acid to 100 ml with DI water. Final concentrations range from 0.014 to 0.025 mg/L. Prepare CCV working standard fresh each day.

9.3.6 The Continuing Calibration Blank (CCB) is prepared by diluting 5.0 ml of concentrated HNO₃ to 1000 ml with DI water.

10.0 RESPONSIBILITIES

It is the responsibility of the analyst to perform the analysis according to the instructions in this SOP and to complete all documentation required for data review. Analysis and interpretation of the results are only to be performed by personnel in the laboratory who have demonstrated the ability to generate acceptable results utilizing this SOP. This demonstration is in accordance with the training program of the laboratory. Final review and sign-off of the data is performed by the department supervisor/manager or designee.

11.0 PROCEDURE

11.1 Instrument Operation and Data Acquisition Procedure

11.1.1 Instrument performance specifications are specified in the operations manual located in the GFAA Lab. Refer to these when setting the parameters for acquisition and select the set of parameters applicable to the element being measured. For Arsenic, Lead, Selenium, and Thallium, the gas type is 95% argon-5% hydrogen. The modifiers used are 500 ppm Pd and 500 ppm Mg(NO₃)₂.

- 11.1.4 Turn on computer and instrument. Install appropriate HCL or EDL lamp and turn on by going to the Align Lamps screen. Allow the Hollow Cathode Lamp to warm up atleast 15 minutes and the Electrodeless Discharge Lamp atleast one hour. This warm up time is used to insure maximum optimization is achieved.
- 11.1.5 Wipe down the autosampler tip before each analysis. Inspect and clean, if necessary, the electrodes, shroud, lens and cells. If the instrument has a fume extractor be sure the Separator Trap is filled with deionized water.
- 11.1.6 Select Automatic Run and choose the appropriate program for analysis.
- 11.1.7 Go to Align Lamp screen and optimize lamp alignment.
- 11.1.8 Insert sample analysis sequence into Instrument Run Log. Create Sample Information File (sample labels) from logbook. Open the Method Editor and go to the Checks page. Assign a post digestion spike to the proper samples by setting parameters where is asks "Perform Recovery Measurements".
- 11.1.9 Load autosampler wheel with standards, blank, modifier, and sample digestates being analyzed, being sure to identify samples by their autosampler position in the instrument log book.
- 11.1.10 Start automatic run.
- 11.1.11 Monitor the analytical run for calibration and sample abnormalities.
- 11.1.12 Dilutions for samples with results over the calibration range are to be made manually and added to the autosampler.

11.2 Calibration

- 11.2.1 The analysis begins with the analysis of calibration standards. Analyze an instrument blank and four calibration standards. The correlation coefficient for each calibration shall be checked to determine that the coefficient is equal to or greater than 0.995.
- 11.2.2 Following calibration, analyze an ICV standard. The resulting value must be within 95-105% of the true value for 200 Series metals and 90-110% of the true value for SW-846 and ASP / CLP4.1. If not, prepare new standards and recalibrate the system.

11.3 Sample Analysis

- 11.3.1 Following calibration, analyze samples and QC samples in an analytical sequence. Refer to the *SOP for Analytical Batches and Analytical Sequences*.
- 11.3.2 Prepare and analyze post-digestion spiked samples. For ASP / CLP4.1 analyses, all samples must be post-spiked and analyzed. For routine analyses, one sample per submission / job number is spiked and evaluated. If the spike recovery is outside the control limits of 85-115%, all the samples in the batch are post-spiked and analyzed. If spike recoveries are within the acceptable limits, analysis is continued with no further spiking. If the post-digestion spike is <40% recovery the sample is diluted by a factor of 5-10 and reanalyzed. If the post-digestion spike is between 40-85% and the sample concentration is less than half of the spike the data is reported. For CLP4.1 a "W" flag will be on the Form I when this occurs.
- 11.3.3 Method of Standard Additions (MSA) analysis is required when an outlying post-digestion spike recovery is between 40-85% or >115% AND the sample absorbance or concentration is less than 50% of the spike. Quantitation is biased to some unknown interference that has been confirmed and dilution and reanalysis does not improve performance, therefore MSA must be performed by analyzing the sample plus 3 spikes at 50, 100, and 150% of the sample concentration (single injection is only required). Plot a linear regression curve of concentration vs. absorbance and quantitate sample concentration from the curve. Refer to the CLP SOW or NYSDEC ASP for proper qualifications of sample data.

12.0 QA/QC REQUIREMENTS

- 12.1 All GFAA sample analyses shall be performed with duplicate burns and the average reported. Duplicate burns should not exceed 20%RSD to maintain precision throughout the run. If RSD exceeds 20%, reanalyze once; if continues to be >20%, dilute 1:2 and reanalyze to avoid interference.
- 12.2 The correlation coefficient for each calibration must be equal to or greater than 0.995. The software produces the calibration curve point by point and does not reject a calibration that has a correlation coefficient < 0.995. The run must be stopped by the operator if the correlation coefficient fails.
- 12.3 Analyze CCV standards and CCBs no less frequently than every ten samples in the analytical sequence.
- 12.3.1 For CCVs, the resulting value must be 90-110% of the true value for CLP analyses and 80-120% for routine analyses. If not, recalibrate the system and reanalyze samples run since the last acceptable CCV.

12.3.2 For CCBs, the resulting value must be less than the PQL. Check the CCB result for carryover. Re-analyze if the CCB result is above the PQL.

- 12.4 Each sample preparation batch must have a method blank associated with it. The method blank result should be < PQL. If not, redigest the batch of samples.
- 12.5 A laboratory control sample (LCS) is digested one per batch, or per 20 samples. The LCS recovery criteria is listed in Appendix C of the Quality Assurance Manual. If the LCS fails the acceptance criteria, redigest the batch of samples.
- 12.6 Post-digestion spike recovery acceptance limits are 85-115%. If outside these limits for the Preparation Blank stop analysis, correct problem and reanalyze. If routine sample recovery fails all samples in the corresponding submission / job # must be spiked. If post-spike recovery is > 40%, samples may be run by the method of standard additions to prevent further dilution see sections 11.3.2 and 11.3.3.
- 12.7 A duplicate sample is digested one per batch, or per 20 samples (one per 10 if 200 Series is requested). Frequency and QC criteria are listed in Appendix C of the Quality Assurance Manual. If the RPD is greater than the limit, determine if the sample is non-homogenous. Redigest if necessary, otherwise data may be flagged with a "*" for job specific QC samples.
- 12.8 A matrix spiked sample is digested one per batch, or per 20 samples. Frequency and QC criteria are listed in Appendix C of the Quality Assurance Manual. If outside acceptance limits, redigest if necessary, otherwise the data may be flagged with a "N" for job specific QC samples. If the sample concentration is >4x the spike level, no action is required.
- 12.9 Additional QC measures include annual determination of method detection limits. Refer to ADM-MDL for procedure and requirements.

13.0 DATA REDUCTION AND REPORTING *refer to ADM-DREV for data review*

13.1 Results for aqueous samples are calculated as follows and are reported in mg/L: *procedures.*

$$\text{mg/L (sample)} = C * x (\text{Digestion Dilution Factor}) x (\text{Post-Digestion Dilution Factor}) \div 1000$$

13.2 Results for soil and solid samples are calculated as follows and are reported in mg/Kg.

$$\text{mg / Kg (Sample)} = C^* x \text{Post Digestion Dilution Factor} x \frac{\text{Digestion Vol. (ml)}}{\text{Sample wt. (g)}} x \frac{1\text{mg}}{1000\text{g}} x \frac{1\text{L}}{1000\text{ml}} x \frac{1000\text{g}}{1\text{Kg}}$$

where C * is the concentration of the analyte as measured at the instrument in µg/L.

14.0 WASTE MANAGEMENT AND POLLUTION PREVENTION

- 14.1 It is the laboratory's responsibility to comply with all federal, state, and local regulations governing waste management, particularly the hazardous waste identification rules and land disposal restrictions, and to protect the air, water, and land by minimizing and controlling all releases from fume hoods and bench operations. Compliance with all sewage discharge permits and regulations is also required.
- 14.2 Excess, unused sample and testing byproducts are disposed following the procedures in the *SOP for Waste Disposal*.

15.0 REFERENCES

Test Methods for Evaluating Solid Waste, Physical/Chemical Methods, EPA SW-846, 3rd Edition; (September 1986) and Updates I (July 1992), II (September 1994), IIA (August 1993), IIB (January 1995), III (December 1996).

Methods for Chemical Analysis of Water and Wastes, EPA-600/4-79-020, (Revised March 1993).

Methods for the Determination of Metals in Environmental Samples, EPA/600/4-91/010 (June 1991) and Supplement I, EPA/600/R-94/111 (May 1994).

EPA Contract Laboratory Program, Statement of Work for Inorganic Analysis, SOW No. ILM04.0.

Analytical Services Protocol (ASP), New York State Department of Environmental Conservation, December 1995.

16.0 TRAINING OUTLINE

- 16.1 Read current SOP and applicable methodologies. Demonstrate a general understanding of the methodology and chemistry.
- 16.2 Observe Sample Preparation and Analysis.
- 16.3 Participate in the methodology, documentation, and data reduction with guidance.
- 16.4 Instrument Operation and Maintenance, if applicable.
- 16.5 Demonstrate Competency by performing the analysis independently. Analyze a known proficiency or standard four times to establish Initial Demonstration of Capability. If recovery is

within acceptable limits, complete training form and certificate and file with QA. Continuing Demonstration of Capability is required on an annual basis, refer to ADM-TRANDOC.

17.0 INSTRUMENT-SPECIFIC ADDENDUM

See Operations Manual in GFAA Lab.

18.0 ATTACHMENTS

Table 1 Summary of Parameters, Methodology, and Reporting Limits

19.0 CHANGES FROM PREVIOUS REVISION

- Added reference to the Instrument Specifications per Metal Manual to Section 2.2.
- Changed Independent to Initial in Section 3.2
- Added reference to Hot Block digestion vessels in Section 6.5
- Added EDL Lamp and EDL power supply to Section 7
- Inserted how often the Calibration and ICV / CCV Standard Stocks should be prepared in Section 9.3.2 and 9.3.4
- Sections 11.0, 12.0 of previous SOP were revised to include more detail and referenced Appendix C of Quality Assurance Manual for QC criteria and frequency requirements.
- Deleted all references to analyzing Antimony by GFAA.

Table 1

Parameters, Methodology, and Reporting Limits

Parameter	Methodology	Practical Quantitation Limit (PQL)	
		Water (mg/L)	Soil (ug/g)
Arsenic	206.2/7060A/CLP	0.0050	0.50
Lead	239.2/7421	0.0050	0.50
Lead	CLP	0.0030	0.30
Lead in Drinking Water	239.2	0.0010	
Selenium	270.2/7740/CLP	0.0050	0.50
Thallium	279.2/7841/CLP	0.010	1.0

*Contract Laboratory Program (CLP) references refer to ILM04.1 and/or 6010B-CLPM (ASP 1995).

SOP No.: MET-GFAA

[illegible]

MEMO

From the QA Department

To: Department Supervisor
From: Lisa Reyes / Vicky Collom
Date:
RE: SOP Annual Review/ Newly Revised SOP

SOP for Metals Analysis by Graphite Furnace
SOP number MET-GFAA
Revision number 3
Date: 11/1/01

☐ It has been about a year since this SOP was last reviewed and it now needs to be reviewed again. Please review the SOP with your staff and complete the applicable section below.

☒ This SOP has been revised to include the changes summarized in Section 19.0. Please review this SOP with your staff and complete the section below.

☐ I have reviewed this SOP with the following personnel and it still reflects current practice.

Signature (supervisor): _____ Date: _____

We have read, understood, and agree to perform the most recent version of this SOP or test method.

Signature	Date	Signature	Date
<i>John M. O'Leary</i>	<i>11/1/01</i>		
<i>Luigi McCarty</i>	<i>11-01-01</i>		
<i>Paul C. Zed</i>	<i>11/1/01</i>		
<i>Christopher R. Zed</i>	<i>11/1/01</i>		

☐ This SOP needs to be revised and updated to reflect current practices. This will be completed by _____.
(date)

Signature (supervisor): _____ Date: _____

*Please indicate the outcome of your review and return this memo to me within 5 working days.
Thank you for your help!!*

AV

STANDARD OPERATING PROCEDURE

**TOTAL SULFUR FOR ION CHROMATOGRAPHY
FOR INDIANA PINES SITE**

MET-ICSPINES

Revision 0

September 23, 2004

Approved By:

Christine Kuzel
Supervisor

9/23/04
Date

Lisa Reyes
QA Manager

9/23/04
Date

Michael K. Perry
Laboratory Manager

9/23/04
Date

COLUMBIA ANALYTICAL SERVICES, INC.

1 Mustard Street, Suite 250

Rochester, NY 14609

©Columbia Analytical Services, Inc. 2004

Annual review of this SOP has been performed
and the SOP still reflects current practice.

Initials: _____ Date: _____

Initials: _____ Date: _____

Initials: _____ Date: _____

NON-CONTROLLED COPY

Will Not Be Updated

Table of Contents

	<u>Page</u>
1. Scope and Applicability.....	3
2. Summary of Method.....	3
3. Definitions.....	3
4. Health and Safety Warnings.....	3
5. Cautions.....	4
6. Interferences.....	4
7. Personnel Qualifications.....	4
8. Equipment and Supplies.....	4
9. Procedure.....	5
9.1. Sample Collection.....	5
9.2. Sample Handling and Preservation.....	5
9.3. Sample Preparation	5
9.4. Sample Analysis.....	6
9.5. Troubleshooting.....	6
9.6. Data Acquisition, Calculations, and Data Reduction Requirements.....	6
10. Data and Records Management.....	6
11. Quality Control and Quality Assurance.....	7
12. References.....	7

Attachments

Digest Sheet.....	8
Training Plan Form.....	9

1. SCOPE AND APPLICABILITY

This procedure is used to determine the concentration of oxidizable sulfur in a sample using peroxide digestion and ion chromatography. This SOP describes the sample preparation step for the analysis and refers to the determinative procedure used for ion chromatography. The procedure is applicable to most sample matrices including water, wastewater, soils, and miscellaneous solids. The PQL for soils is 200 mg/Kg. This SOP was modified specifically for the Indiana Pines site project.

2. SUMMARY OF METHOD

A portion of the sample is digested using a heated peroxide solution. The resulting digestate is filtered and analyzed for sulfate using ion chromatography. The sulfate result is converted to concentration of sulfur.

3. DEFINITIONS

- 3.1. Laboratory Control Sample (LCS): A laboratory blank that has been fortified with target analyte and used to determine that the analysis is in control.
- 3.2. Matrix Spike (MS) Analysis - In the matrix spike analysis, a predetermined quantity of target analyte is added to a sample matrix prior to sample preparation and analysis. The percent recovery is calculated. The MS is used to evaluate the effects of the sample matrix on the method used for the analysis
- 3.3. Duplicate Sample (DUP) - A laboratory duplicate. The duplicate sample is a separate field sample aliquot that is processed in an identical manner as the sample proper. The relative percent difference between the samples is calculated and used to assess analytical precision.
- 3.4. Method Blank / Preparation Blank (MB) - The method blank is an artificial sample composed of analyte-free water or solid matrix and is designed to monitor the introduction of artifacts into the analytical process. The blank is carried through the entire analytical procedure.
- 3.5. Batch – Up to 20 samples of the same matrix digested together on the same day.

4. HEALTH AND SAFETY WARNINGS

The toxicity or carcinogenicity of each reagent used in this method has not been precisely determined; however, each chemical and sample should be treated as a potential health hazard. Exposure should be reduced to the lowest possible level. The laboratory maintains a compilation of Material Safety Data Sheets in binders the conference room. Always wear safety glasses or a shield for eye protection, and protective clothing, and observe proper mixing when working with these reagents.

5. CAUTIONS

Boiling samples to dryness may cause combustion

6. INTERFERENCES

Samples impervious to peroxide digestion will yield results of low bias. Samples with high organic content may require additional digestions.

7. PERSONNEL QUALIFICATIONS

At a minimum, personnel must have attained at least a 2-year degree in a science-related field and have successfully completed an Initial Demonstration of Capability and the Training Plan Form (attached). Training and Demonstration of Capability are in accordance with NELAC 2002 standard.

8. EQUIPMENT AND SUPPLIES

- 8.1. 50ml Digestion Vessel for Hot Block
- 8.2. 250 mL glass beaker and ribbed watch glasses
- 8.3. Hotplate capable of maintaining a digestion temperature of 90-95°C.
- 8.4. Hot Block Digestor- Environmental Express
- 8.5. Filter Mate 2u filter paper and plunger for Environmental Express Digestion Vessel.
- 8.6. Dionex Ion Chromatograph Series 4000i, as described in GEN-300.0 SOP.
- 8.7. 10 N Sodium Hydroxide (NaOH): Dissolve 400g sodium hydroxide in distilled water, cool and dilute to 1 liter. Store at room temperature for up to 1 year.
- 8.8. 30% peroxide; purchased solution. Store at room temperature. Expires upon manufacturer's indications or in 1 year, whichever is sooner.
- 8.9. Laboratory D.I. water
- 8.10. Granular sodium sulfite, Na_2SO_3 anhydrous FW=126.04. 254390 mg/Kg (25.4%) sulfur. To be used for the LCS and for spiking the MS.

9. PROCEDURE

9.1. Sample Collection

- 9.1.1. Samples are to be collected in purchased, precleaned, certified sample containers – (plastic, glass, etc). Samples are to be cooled upon collection and shipment to lab.
- 9.1.2. The amount of sample collected should be 3 times the analytical aliquot, at a minimum.

9.2. Sample Handling and Preservation

- 9.2.1. Maintain samples at 0-6 °C upon receipt until analysis.
- 9.2.2. No specific holding time applies.
- 9.2.3. For further sample handling, storage, and custody procedures, see SMO-GEN.

9.3. Sample Preparation

- 9.3.1. Aqueous samples: Measure a 50 ml sample aliquot into a digestion vessel. Record the volume.
- 9.3.2. Soil samples: weigh out 0.5-5g of sample into a digestion vessel. Record the weight.
- 9.3.3. Add 2 drops of 10 N NaOH to each vessel, or until the sample is basic in nature.
- 9.3.4. Add appropriate standard to matrix spike and LCS aliquots.
- 9.3.5. Add 3 mL 30% peroxide to each vessel.
- 9.3.6. Bring each vessel to 50 ml with D.I. water.
- 9.3.7. Place digestion vessel in hotblock digester OR transfer contents of digestion vessel to a beaker and place on a hotplate.
- 9.3.8. Digest each sample until digestate is clear, or three times. Bring the volume of the digestate to ~ 5 - 10 mL each time taking care to not evaporate the samples to dryness. BOILING SAMPLE TO DRYNESS MAY CAUSE COMBUSTION.
- 9.3.9. Allow samples to cool. Bring soil and water samples to a final volume of 20.0 mL in the digestion vessel. Record the final volume. If particulates are present in the sample, filter using 2u FilterMate filter for Environmental Express digestion vessels. If one sample is filtered, the entire batch is to be filtered, including the MB and LCS.

- 9.3.10. Give the batch of samples and a copy of the digest sheet to Wetchem for analysis.
Document custody transfer.

9.4. Sample Analysis

The extract is analyzed for sulfate by ion chromatography (IC) using SOP GEN-300. Refer to that SOP for specific analysis instructions.

9.5. Troubleshooting and Preventive Maintenance – Wipe down all hoods in the Metals Prep Lab once a week with DI water.

9.6. Data Acquisition, Calculations, and Data Reduction Requirements

9.6.1. The PeakNet software will multiply the solution result by any dilution made at the IC and by the final volume. Divide by the initial volume or weight.

9.6.2. The IC sulfate result will be multiplied by 0.3338 to obtain the concentration of the sulfur (S is 33.38% of SO₄ by atomic weight).

10. DATA AND RECORDS MANAGEMENT

- 10.1. **Responsibilities** – It is the responsibility of the analyst to perform the analysis according to this SOP and to complete all documentation required for data review. Analysis and interpretation of the results are performed by personnel in the laboratory who have demonstrated the ability to generate acceptable results utilizing this SOP. This demonstration is in accordance with the training program of the laboratory. Final review and sign-off of the data is performed by the department supervisor/manager or designee.
- 10.2. **Data Flow** – Samples are entered by the Project Manager into StarLIMS v.6.11.a on a Personal Computer running on a Novell Network. On the day that the samples are received the samples appear on a daily log printed from this computer system. The Metals Prep analyst prepares a benchsheet (attached), digests the samples and turns the samples and digest sheet over to the IC analyst. The samples are analyzed for sulfate using Dionex PeakNet 5 Chromatography software and the results are transferred into the StarLIMS computer system for final calculation, validation, reporting, and invoicing.
- 10.3. **Data Review** - Data will be reviewed by the IC analyst and a qualified peer using a Data Review Checklist (attached to GEN-300) and validated by a supervisor.

11. QUALITY CONTROL AND QUALITY ASSURANCE

11.1. Method Blank-

- 11.1.1. Frequency - Prepare one method blank per batch of 20 samples.
- 11.1.2. Acceptance Criteria – The result of the method blank must be less than the reporting limit. If there is method blank contamination, samples which have results less than the reporting limit may be reported.
- 11.1.3. Corrective Action - If there is method blank contamination, attempt to find the source of the contamination, correct the problem, and re-digest the batch (with the exception of the samples accepted as above).

11.2. LCS –

- 11.2.1. Frequency – one per batch of 20 or fewer samples.
- 11.2.2. Acceptance criteria – The result of the LCS must be within 80-120% of the true value.
- 11.2.3. Corrective action – If the LCS is out of control limits, find and correct the problem and re-digest the batch.

11.3. Matrix Spike –

- 11.3.1. Frequency – one per batch of 20 or fewer samples of the same matrix.
- 11.3.2. Acceptance criteria – The result of the MS should be within 70-130% of the true value.
- 11.3.3. Corrective Action – If the MS is out of control limits, and the LCS is compliant, assume matrix interference and report. If the MS is out of control and the LCS is out of control, find the problem and redigest the batch.

11.4. IC QC Requirements are outlined in Section 12 of GEN-300.

12. REFERENCES

NELAC, 2002 Standard
CAS SOP for Ion Chromatography, GEN-300.

CAS - Rochester, NY: Total Sulfur for IC Digestion Log

Report Type: Routine // 6 // ASP // Pkg5

Analyst: _____ Date: _____

Prep Method: Total Sulfur for IC by Method 300

Digest: Initial // Redigest Of: _____

Batch ID: _____
Spike Witness / Lot Approval: _____
Batch Temp: _____

Submission / Order #	pH	Initial Vol. / Wt. (ml / g)	Final Vol(ml)	Initial Color / Clarity	Final Color / Clarity	Analyte Total Sulfur/Sulfate	Spike Added Vol(ml)
1							
2							
3							
4							
5							
6							
7							
8							
9							
10							
11							
12							
13							
14							
15							
16							
17							
18							
19							
20							
21							
22							
23							
24							
25							
26							

Color / Clarity Key:
Color: C = Colorless ; Y = Yellow ; B = Brown ; G = Grey
BL = Black
Clarity: CDY = Cloudy ; CLR = Clear ; OP = Opaque
Texture: F = Fine ; M = Medium ; C = Coarse ; NA = Non Aqueous

Spiking Standards / Reagent Lot Numbers:

10 N Sodium Hydroxide (NaOH): _____

30 % Peroxide: _____

Sodium Sulfite (Na2SO4): _____

COMMENTS: _____

Metals Digestion Training Plan

Procedure: _____

SOP: _____ Revision: _____ Date: _____

Trainee: _____

1. Read SOP Trainer: _____ Trainee: _____ Date: _____

2. Demonstrated familiarity with related SOPs

-ADM-BATCHSEQ	-ADM-PCAL	-ADM-SIGFIG
-ADM-DATAENTRY	-ADM-SPSR	-ADM-MDL
-ADM-TRANDOC		

Trainer: _____ Trainee: _____ Date: _____

3. Observe performance of SOP
- standard and reagent prep and documentation – including pipet use and balance use and calibration, if applicable
- digestion unit set-up
- sample prep and reagent and spike addition
- holding times
- benchsheet/logbook use
- analytical sequence, batch QC required
- time and temperature needed to digest sample, if applicable
- preventive maintenance and other troubleshooting
- digestate filtering and dilution
- digestate labelling and storage

Trainer: _____ Trainee: _____ Date: _____

4. I have read, understood and agree to perform the most recent version of the SOP:

Signature: _____ Date: _____

5. Perform SOP with supervision
- including all items in 4.

Trainer: _____ Trainee: _____ Date: _____

6. Independent performance of the SOP
- all of the item listed in 4
- IDC (4 mid-range standards performed before client samples are analyzed)
- attach IDC certificate, raw data, and summary spreadsheet.

Trainer: _____ Trainee: _____ Date: _____


STANDARD OPERATING PROCEDURE


**DETERMINATION OF SULFUR IN SOILS USING ION CHROMATOGRAPHY
AFTER ALKALINE DIGESTION
FOR INDIANA PINES SITE**


GEN-300Pines

Revision 0

September 24, 2004

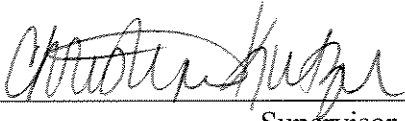Approved By:  9/24/04
Supervisor Date

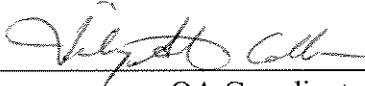 9/24/04
QA Coordinator Date

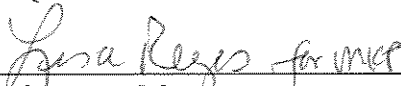 9/24/04
Laboratory Manager Date

©COLUMBIA ANALYTICAL SERVICES, INC. 2004
One Mustard St., Suite 250
Rochester, NY 14609

Annual review of this SOP has been performed
and the SOP still reflects current practice.

Initials: _____ Date: _____
Initials: _____ Date: _____
Initials: _____ Date: _____

NON-CONTROLLED COPY
Will Not Be Updated

Table of Contents	<u>Page</u>
1. Scope and Applicability.....	3
2. Summary of Method.....	3
3. Definitions.....	3
4. Health and Safety Warnings.....	4
5. Interferences.....	5
6. Personnel Qualifications.....	5
7. Equipment and Supplies.....	5
8. Procedure.....	8
8.1. Calibration and Standardization.....	8
8.2. Sample Collection.....	9
8.3. Sample Handling and Preservation.....	9
8.4. Sample Preparation	9
8.5. Sample Analysis.....	9
8.6. Troubleshooting.....	12
8.7. Data Acquisition, Calculations, and Data Reduction Requirements.....	13
8.8. Computer Hardware and Software.....	14
9. Data and Records Management.....	14
10. Quality Control and Quality Assurance.....	14
11. References.....	16
Attachments	
IC Training Plan Form.....	17
Wetchem Data Review Checklist.....	18

1 SCOPE AND APPLICABILITY

1.1 This SOP uses Method 300.0 for the analysis of sulfate by Ion Chromatography in soil samples prepared by alkaline digestion according to MET-ICS.

1.2 Range

Using the settings and calibration techniques outlined in this SOP, the upper range for sulfate is 10 ppm (solution concentration). Higher concentrations of sulfate may be determined using appropriate dilutions. Review current calibration for specific ranges.

1.3 The PQL for the current system is 0.20 mg/L

2 SUMMARY OF METHOD

Sample digested by alkaline digestion. The extract is filtered and injected into an ion chromatograph (Dionex Series 4000i). Sulfate is chromatographically separated and measured with a conductivity detector. Suppression is accomplished using an ion exchange membrane. It is assumed that all of the sulfur is converted to sulfate during the digestion. The sulfate results are converted by calculation to concentration of sulfur in the original soil sample.

3 DEFINITIONS

3.1 **Initial Calibration** - analysis of analytical standards for a series of different specified concentrations; used to define the linearity and dynamic range of the response of the system.

3.2 **Independent Calibration Verification (ICV)** - ICV solutions are made from a stock solution which is different from the stock used to prepare calibration standards and is used to verify the validity of the standardization. The ICV is analyzed immediately following the calibration standards.

3.3 **Relative Percent Difference (RPD)** – The absolute value of the difference of two values divided by the average of the same two values. Used to compare the precision of the analysis. The result is always a positive number.

3.4 **Batch** - Samples processed together as a unit, not to exceed 20 investigative samples.

3.5 **Method Detection Limit (MDL)**: a statistically derived value representing the lowest level of target analyte that may be measured by the instrument with 99% confidence that the value is greater than zero

3.6 **Method Reporting Limit (MRL)**: The minimum amount of a target analyte that can be measured and reported quantitatively. The MRL is equivalent to Practical Quantitation Level (PQL) and Estimated Quantitation Level (EQL). Typically, the MRL is calculated as five times the MDL (although this is a rule of thumb and not intended to be a strict policy of establishing the MRL for a compound).

3.7 **QA/QC Samples:** Samples added to a sample preparation batch, or an analytical batch to provide quality assurance checks on the analysis.

3.7.1 **Matrix Spike (MS)** - In the matrix spike analysis, predetermined quantities of standard solutions of certain analytes are added to a sample matrix prior to analysis. The purpose of the matrix spike is to evaluate the effects of the sample matrix on the methods used for the analyses. Percent recoveries are calculated for the analyte detected. In this method, spikes are very useful in determining proper retention times when a low concentration of an analyte is detected or expected to be adjacent to a large concentration of analyte. When a spike is used to verify retention time, calculation of recovery is not necessary.

3.7.2 **Duplicate Sample (DUP)** - A laboratory duplicate. The duplicate sample is a separate field sample aliquot that is processed in an identical manner as the sample proper. The relative percent difference between the samples is calculated and used to assess analytical precision.

3.7.3 **Continuing Calibration Verification Standard (CCV)** - A standard analyzed at specified intervals and used to verify the ongoing validity of the instrument calibration.

3.7.4 **Instrument Blank (ICB/CCB)** - The instrument blank (also called initial or continuing calibration blank) is a volume of blank reagent of composition identical to the samples (ie. not chemically preserved). The purpose of the ICB/CCB is to determine the levels of contamination associated with the instrumental analysis. The ICB is performed once, immediately after the ICV.

3.7.5 **Laboratory Control Standard (LCS)** - In the LCS or blank spike analysis, predetermined quantities of standard solutions of certain analytes are added to a blank prior to sample analysis. Percent recoveries are calculated for the analyte detected.

4 HEALTH AND SAFETY WARNINGS

- Take all appropriate safety precautions for handling reagents and samples when performing this procedure. This includes the use of personnel protective equipment, such as safety glasses, lab coat and the correct gloves.
- Handle chemicals, reagents and standards as described in the CAS safety policies, approved methods and in MSDSs where available.
- The use of pressurized gases is required for this procedure. Exercise care when moving cylinders. All gas cylinders must be secured to a wall or an immovable counter with a chain or a cylinder clamp at all times. Sources of flammable gases (e.g., pressurized hydrogen) should be clearly labeled.

- When releasing the cap on the suppressor reagent, wear a face shield and exercise caution. The container is pressurized and the reagent will emit a fine mist. Turn the cap slowly.

5 INTERFERENCES

5.1 Interferences can be caused by substances with retention times that are similar to and overlap those of the anion of interest. Large amounts of an anion can interfere with the peak resolution of an adjacent anion. Sample dilution and/or spiking can be used to solve most interference problems. The most common examples of this are:

5.1.1 Sulfite will interfere with the sulfate peak.

5.1.2 Thiosulfate can interfere if the run time of the entire chromatogram is too short.

6 PERSONNEL QUALIFICATIONS

At a minimum, personnel must have attained at least a 4-year degree (or 2-yr degree plus one year experience) in a science-related field and have successfully completed an Initial Demonstration of Capability and the Training Plan Form (attached). Training and Demonstration of Capability are in accordance with NELAC 2002 standard.

7 EQUIPMENT AND SUPPLIES

- 7.1 Analytical Balance, capable of accurately weighing to the nearest 0.0001 g.
- 7.2 Anion guard column: A protector of the separator column. If omitted from the system the retention times will be shorter. Dionex Ionpac AG4A-SC 4×50 mm (P/N 43175)
- 7.3 Anion separator column: Dionex AS14 4x250 (P/N 046124). Expires when separation between the anions of interest is no longer acceptable or upon manufacturer's indications, whichever occurs first.
- 7.4 Anion suppressor device: Dionex anion micro membrane suppressor (P/N 53946).
- 7.5 Detector-Conductivity Cell: approximately 1.25 µL internal volume.
- 7.6 Dionex PeakNet 5.1 Chromatography Workstation software or equivalent. Personal computer connected to network, capable of running the PeakNet software.
- 7.7 Calibrated MicroPipettor and tips.
- 7.8 Calibrated repipettor.

7.9 System configuration

- An automated sampler (Dionex P/N 39534)
- An analytical gradient pump (Dionex System 4000i)
- A separator column (Dionex AS14 4x250mm P/N 046124)
- A conductivity detector (Dionex P/N 40157)
- A 50 μ L sample loop
- Pump rate of 2.0 mL/min

7.10 Standards Preparation General Information

- Bring any cooled parent stocks to room temperature before use.
- All standards and reagents are to be tightly capped when not in immediate use. Protect standards and reagents from light whenever possible.

7.11 Reagent water: Distilled or deionized water, free of the anions of interest.

7.12 Stock Eluent solutions for AS14 column

7.12.1 0.5 M Sodium Carbonate Concentrate– Dissolve 26.49g Na₂CO₃ in 400 mLs DI. Bring to volume in a 500 mL volumetric flask. Expires one year. Store at room temperature.

7.12.2 0.5M Sodium Bicarbonate Concentrate – Dissolve 21.00 g of NaHCO₃ in 400 mL DI. Dilute to a final volume of 500 mLs with DI. Store at room temperature. Expires in one year.

7.13 Working Eluent Solution for AS 14 Column – 3.5 mM Sodium Carbonate / 1.0 mM Sodium Bicarbonate – Filter a sufficient volume of each of the 2 eluent reagents through 0.2 μ m syringe filters into separate dispo cups. Pipette 7.0 mL of 0.5 M Na₂CO₃ and 2.0 mL of 0.5 M NaHCO₃ into a 2 Liter volumetric flask. Dilute to volume with DI. Degas for 5 minutes with ultra high purity Helium at a rate of 1-5 bubbles per second. Store at room temperature. Expires in 1 week.

7.14 Regeneration solution (micro membrane suppressor): Sulfuric acid 0.1N. Dilute 5.6 mL of conc. sulfuric acid (H₂SO₄) to 2L with reagent grade water. Degas for 5 minutes with ultra high purity Helium at a rate of 1-5 bubbles per second. This solution is stable for one week from date of preparation. Store at room temperature in plastic.

7.15 Stock standard solutions

7.15.1 Sodium Sulfate – ACS reagent grade dried at 103-105°C for 30 mins. Store dried material in a small glass beaker. Enclose the beaker in aluminum foil to protect from light. Store the covered beaker in a desiccator. Expires in one year.

7.15.2 Sulfate (SO_4^{2-}) 1000 mg/L: Dissolve 1.479 g prepared (as above) sodium sulfate (Na_2SO_4) in reagent water and dilute to 1 L. Store at 0-6°C in amber glass for up to 1 year.

7.16 Intermediate Calibration Standards –

7.16.1 Routine Intermediate Stock – Store at 0-6°C in plastic. Expires in 6 months. Also used as LCS and MS Intermediate stock.

Analyte:	SO_4
Stock Conc (mg/L):	1000
mLs Stock:	20.0
Final Vol (mLs):	200.0
Int. Stock Conc (mg/L):	100.0

7.17 Calibration Standards - prepared in 100 mL volumetric flask as follows: record the pipette ID used in the reagent prep logbook. Make fresh weekly. Store at 0-6°C in glass or plastic.

Standard ID	mLs of Intermediate Stock	Final Volume, mLs	Working conc.
			SO_4 mg/L
Std # 9	10.0	100.0	10.0
Std # 8	8.0	100.0	8.0
Std # 7	5.0	100.0	5.0
Std # 6	2.0	100.0	2.0
Std # 5	1.0	100.0	1.0
Std # 4	0.5	100.0	0.50
Std # 3	0.2	100.0	0.20
Std # 2	0.1	100.0	0.10
Std # 1	0.0	100.0	0.00

7.18 Reference Standard Stocks:

7.18.1 Potassium Sulfate – ACS reagent grade dried at 103-105°C for 30 mins. Store dried material in a small glass beaker. Enclose the beaker in aluminum foil to protect from light. Store the covered beaker in a desiccator. Expires in one year.

7.18.2 Sulfate (SO_4) 3200 mg/L: Dissolve 5.80g dried (as above) K_2SO_4 in reagent water and dilute to 1 L. Store at 0-6 °C in amber glass for up to 1 year.

7.19 ICV/CCV Intermediate stock (12.8 mg/L) - Dilute 4.0 mL of 3200 mg/L reference stock to 1 Liter in a volumetric flask. Store at 0-6°C in plastic for up to 6 months.

7.20 ICV/CCV (6.4 mg/L) - prepare by diluting the CCV intermediate stock solution with equal parts water in small quantity (about 30 mLs DI and 30 mLs intermediate stock solution). The resulting concentrations are half of those of the intermediate solution. Prepare fresh when needed or at least once a week. Store at room temperature in plastic.

7.21 LCS (2.0 mg/L): Store at room temperature in glass or plastic for up to one week. Add 2.0 mLs of the intermediate stock solution (prepared same as the intermediate solution used for calibration standards) to DI in a 100 mL volumetric flask and bring to volume.

7.22 Matrix Spike Solution - Add 2.0 mL of the intermediate stock solution to 100 mL sample (or dilution of sample). Prepare fresh before use.

7.23 Consumable materials.

- 5 mL vials with filter caps. (Dionex P/N 038141)
- 0.2 μm syringe filters.

8 PROCEDURE

8.1 Calibration and Standardization-

8.1.1 Prepare calibration standards according to Section 7. Document preparation in standards log book. Load standards according to Autosampler Vial Loading Section. Start instrument and analyze according to sections below.

8.1.2 The initial calibration is made by linear regression. This method of quantitation uses the equation of a line ($y=mx+b$). The curve must not be forced through zero. System calibration must have correlation coefficient of 0.995 or better. Delete outlier standards. Standards must be within 10% of their true value. Method 300.0 requires a minimum of 3 standards and a blank. If the removal of outlier standards does not bring the curve into compliance, recalibrate.

8.1.3 Immediately after an acceptable calibration has been achieved, run the ICV, ICB, and an LCS. If these are compliant, continue with samples as described in the daily analytical sequence.

- 8.2 **Sample Collection** – Samples should be collected in purchased, certified clean glass or polyethylene bottles or jars.
- 8.3 **Sample Handling and Preservation** – Sulfate holding time is 28 days from collection. Samples stored at 0-6°C from receipt until analysis. Sample handling, storage, and custody procedures are in accordance with NELAC 2002 Standard.
- 8.4 **Sample Preparation** – Soil samples for Total Sulfur are digested according to MET-ICS. Further preparation of the extract is given below.
- 8.5 **Sample Analysis** –
- 8.5.1 **Prepare the Instrument** –
- 8.5.1.1 Be sure there is a current MDL and IDC for the system.
- 8.5.1.2 Check eluent and regenerant levels in containers. Fill to appropriate levels as necessary. Hand tighten caps of both jugs.
- 8.5.1.3 Remove plugs from the waste lines on the back of the instrument. Screw flow restrictor onto end of suppressor drain line (both lines are labeled).
- 8.5.1.4 Turn on Helium carrier gas (should be at approximately 17psi.) and compressed air (100 psi.) by turning the yellow handles to the up-down position and the small valves to IC#1. These are located along the column to the left of the computer.
- 8.5.1.5 Start the Dionex Gradient pump – on the bottom right half of the instrument there is a button with stop / start indicator. Press the button to light the start indicator.
- 8.5.1.6 Turn on the Conductivity Cell. In the middle of the instrument there is a CELL off/on indicator. Press the button to light the “on” indicator. Allow the system to warm up for about an hour.
- 8.5.2 **Create a schedule in the PeakNet software -**
- 8.5.2.1 While the system is warming up, determine whether an ICAL is to be run. The instruments must be calibrated if any of the following apply:
- when a new column is put in
 - when system configuration changes warrant calibration
 - every 6 months
 - when QC samples indicate the old calibration is no longer acceptable.
- 8.5.2.2 Determine which samples are to be analyzed.

8.5.2.3 Remove any standards or reagents needed from the cooler and allow to warm to room temperature before use.

8.5.2.4 Create the schedule of the day's run in the software. This may be modified later as needed, but will help with initial organization.

8.5.2.5 If a calibration is not to be run – set up the schedule to analyze samples in the following analytical sequence: CCV, CCB, LCS, 10 samples, CCV, CCB, 10 samples, CCV, CCB, LCS, etc. with a CCV/CCB set after every 10 samples and an LCS after every 20 samples and DUP/MS where appropriate (at no particular position but one set for every 10 samples). Skip the initial calibration section. Prepare the samples and load the autosampler as described below.

8.5.2.6 If a calibration is to be run – set up the schedule to analyze the calibration standards, ICV, ICB, LCS, 10 samples, CCV, CCB, 10 samples, CCV, CCB, LCS, etc. with a CCV/CCB set after every 10 samples and an LCS after every 20 samples and DUP/MS where appropriate (at no particular position but one set for every 10 samples). Continue with initial calibration section.

8.5.3 Prepare the extract for analysis–

8.5.3.1 Draw the extract up into a 10 mL pipette. Place a 0.2 μ m syringe filter on the end of the pipette and push some of the sample (only enough to make a dilution – 2 mLs is plenty) through the filter into a dispo cup.

8.5.3.2 Use the filtered extract to make an appropriate dilution.

8.5.4 Autosampler Vial Loading

8.5.4.1 Rinse all sample vials and caps to remove any debris present from the manufacture.

8.5.4.2 Once the sample or standard has been placed in the sample vial, place a vial cap in the vial and use the tool to press the cap down flush with the top of the vial.

8.5.4.3 Place the loaded vials into cassettes according to the schedule created and in compliance with the analytical sequence described below. Place the holder in the autosampler.

8.5.5 Start Instrumental Analysis

8.5.5.1 Open the run screen in the PeakNet software. Load the schedule. Select Start.

8.5.5.2 Push the “auto off-set” button on the IC unit to reset the conductivity baseline.

8.5.5.3 Press “Run” on the autosampler.

8.5.6 Evaluate sample analysis

8.5.6.1 Examine solution concentrations of target analytes in the samples. If the concentration is greater than the high calibration standard, reanalyze the sample at a dilution.

8.5.6.2 Check peak integrations.

8.5.6.2.1 Where possible, all integrations should be performed consistent with integration of the corresponding calibration standards.

8.5.6.2.2 Be sure the peaks on the chromatogram and the instrument calculated concentration make sense. Sometimes the software will attempt to integrate overrange peaks and will incorrectly assign them a concentration which would be acceptable for the dilution if it was a reasonable integration.

8.5.6.2.3 On occasion, the software integrates peaks incorrectly. The sample may be reanalyzed or the analyst may use the software to correct the integration. Any manual integration or manipulation of peaks must be consistent with the calibration standards and the QC samples.

8.5.6.3 Evaluate QC samples. All samples must be bracketed by acceptable CCVs and CCBs. See Section 10 for further discussion of QC and sample acceptance and corrective action.

8.5.7 Instrument Shut Down –

8.5.7.1 Take the daily readings. Then turn the auto offset & cell to off and the pump to stop.

8.5.7.2 Turn the gas and air off to each IC individually by turning the small valve handles perpendicular to the gas flow direction.

8.5.7.3 Vent the eluent first, leave the cap very loose, and then ASAP vent the suppressor. (Vent the suppressor by slowly opening both jugs. The

suppressor is acidic, so use care. Wear face shield and cover the jugs with a plastic bag for added protection).

8.5.7.4 Take the flow restrictor off of the suppressor drain line and plug both the eluent and suppressor drain lines.

8.5.7.5 After the last IC is shut off, turn both the gas and air yellow handles to the right.

8.6 Troubleshooting –

8.6.1 Rinsing the IC pump and valves. This should be done weekly, preferably Friday night or Saturday.

8.6.1.1 Disconnect the column from the valve. Plug the column with one of the solid plugs so that it doesn't dry out.

8.6.1.2 Attach the old column to the valve (the old column is in the IC "tool drawer" in the box on the left, behind the filters, B-cups, etc. Get the syringe then, too. It has to have the orange union fitting attached to its tip) Place the tube at the end of the column in the graduated cylinder.

8.6.1.3 Disconnect the eluent line, and plug it up, because it will continue to siphon all over you if you don't. Keep the brown-colored union fitting attached to the blue-colored tubing that leads to the pump heads.

8.6.1.4 Fill the carboy labeled "DI" about halfway with DI (rinse it once or twice first). Put the carboy back in the rack and feed the long tubing to the side of the IC. Attach the syringe to the fitting at the end of the tubing and pull the DI into the syringe to get the siphon going. When it is going, detach it from the syringe and attach it to the brown-colored union fitting attached to the blue-colored tubing that leads to the pump heads. Be sure to allow some of the water dribbling out of the DI carboy tubing to fill up any lost liquid in the brown-colored union fitting, so that you won't (hopefully) have to prime the pump.

8.6.1.5 Now you can turn on the pump. The DI should start flowing out the old column. Let it go for at least 15 minutes, after which time it can be turned off and you can go home.

8.6.1.6 As per Dionex Tech Support, this is to be done only every 6 months: While the DI is pumping through the pump & valves, lubricate the pump by opening up the pump drawer about 2 inches, exposing the pump motor housing. There is a little port in the front of the motor, with yellow grease in it. Attach the grease syringe (located in the cupboard below the IC) and squirt in 0.1 mL of grease (Dionex P/N 39440).

8.6.2 To re-configure back to operation mode:

8.6.2.1 Take off the DI carboy.

8.6.2.2 Attach the filled eluent carboy to the brown-colored union fitting after having starting the siphon, etc.

8.6.2.3 Allow the eluent to pump through the old column until you are sure that all DI has been displaced. Check with pH paper, or allow to pump >8-10 minutes.

8.6.2.4 Re-attach the valve to the guard column/analytical column.

8.6.3 Nightly: Release gas pressure in eluant/suppressor bottles and cap both waste ports. Fill in the daily log, recording Date, Column ID, Helium inlet pressure, System backpressure, Eluant pressure, Detector Background, and Reagent flow.

8.6.3.1 The incoming pressure of the Helium carrier is checked (should be approx. 17 psi.)

8.6.3.2 The system pressure is checked (usually around 1500 psi.).

8.6.3.3 The background of the detector should be around 22-24 μ s.

8.6.3.4 The flow rate of the suppressor coming from the waste line should be 3-4 mL per minute.

8.6.4 Maintenance log – Document all preventive maintenance, as well as instrument repair, in the appropriate instrument maintenance log. Most routine maintenance and troubleshooting are performed by CAS staff. Other maintenance or repairs may, or may not require factory service, depending upon the nature of the task. Any maintenance performed by outside services must also be documented – either through notes in the log or through documents provided by the service. The log entries will include the date maintenance was performed, symptoms of the problem, serial numbers of major equipment upgrades or replacements. The datafile name of the first acceptable run after maintenance is to be documented in the maintenance log.

8.7 Data Acquisition, Calculations, and Data Reduction Requirements

8.7.1 The results which are printed on the instrument report will be adjusted for any dilution made at the instrument. Further adjustment for initial weight and final volume will be made separately. The final multiplication by 0.3338 (sulfur is 33.38% of sulfate by molecular weight) will be done by StarLIMS.

8.7.2 Data will be reviewed by the analyst and a qualified peer using the Data Quality Checklist (attached) and validated by supervisor.

- 8.7.3 All sample data and QC data, including calibration verification must reference the name (date or filename) of the ICAL on the raw data report

8.8 Computer Hardware and Software

- 8.8.1 StarLIMS v.6.11.a
8.8.2 Personal Computer running Dionex PeakNet v5.1

9 DATA AND RECORDS MANAGEMENT

- 9.1 **Responsibilities** - It is the responsibility of the analyst to perform the analysis according to this SOP and to complete all documentation required for data review. Analysis and interpretation of the results are performed by personnel in the laboratory who have demonstrated the ability to generate acceptable results utilizing this SOP. Final review and sign-off of the data is performed by the department supervisor or designee.
- 9.2 **Data Flow** – Samples are entered by the Project Manager into StarLIMS on a Personal Computer running on a Novell Network. On the day that the samples are received the samples appear on a daily log printed from this computer system. The Metals Prep analyst prepares a benchsheet, digests the samples and turns the samples and digest sheet over to the IC analyst. The samples are analyzed for sulfate using PeakNet software. The results are printed and hand entered into StarLIMS. StarLIMS makes the final calculations and the results are printed for data review. When the results are approved, the StarLIMS is used for reporting, and invoicing.
- 9.3 **Data Review** - Data will be reviewed by the IC analyst and a qualified peer using a Data Review Checklist (attached) and validated by a supervisor.

10 QA/QC REQUIREMENTS

10.1 Laboratory Control Standards (LCS)

10.1.1 An LCS must be run daily and once every 20 samples.

10.1.2 The LCS must be within 10% of the true value.

10.1.3 If the LCS is outside the acceptance criteria stop the run, correct the problem and reanalyze the LCS. Exception: if the LCS recovery is high and sample results less than the reporting limit, analysis may continue and data may be reported.

10.2 Method Detection Limits (MDL)

MDLs should be performed every 6 months, when a new operator begins work or whenever there is a significant change in the background or instrument response. The result of the MDL must be less than the PQL. If it is not, correct the problem and do another MDL study or raise the PQL. See 40 CFR Part 136 Appendix B.

10.3 Initial and Continuing Calibration Verification (ICV/CCV)

- 10.3.1 An ICV is analyzed immediately after the standards. The ICV must be 90-110% of the true value or the curve may not be used.
- 10.3.2 A CCV is analyzed every 10 samples.
- 10.3.3 All CCVs must be within 10% of the true value. If the CCV is not in control, correct the problem, obtain a compliant CCV and reanalyze all samples bound by the non-compliant CCV. Recalibrate if necessary. Exception: if the CCV recovery is high and sample results are less than the reporting limit, analysis may continue and data may be reported.

10.4 Continuing Calibration Blanks (CCB)

- 10.4.1 A CCB must be analyzed every 10 samples immediately following the CCV.
- 10.4.2 All CCB's must be less than the PQL. If the CCB is above the PQL, correct the problem and obtain a compliant CCB following a compliant CCV. Reanalyze samples bound by non-compliant CCB. Recalibrate if necessary. Exception: If there is blank contamination and the sample results are less than the reporting limit, analysis may continue and data may be reported.

10.5 Matrix Spikes (MS)

- 10.5.1 A matrix spike must be analyzed once every 10 samples. Do not choose field blanks for the analysis of MS.
- 10.5.2 The matrix spike should be within the lab-generated limits of 69-120% for waters and 70-130 % for soils. If it is not, note the outlying recovery in the case narrative. If the MS is out and the LCS is in, matrix interference is assumed and the batch is acceptable. It is recommended that the MS be reanalyzed to confirm the outliers, however it is not required.

10.6 Duplicates (DUP)

- 10.6.1 A DUP must be analyzed every 20 samples. The DUP is regularly analyzed every 10 samples since the MS must be analyzed every 10 samples. Do not choose field blanks for analysis of DUP.
- 10.6.2 The acceptance criteria for a DUP is less than 20% RPD or \pm the reporting limit if the sample is less than 5 times the reporting limit.
- 10.6.3 If a DUP is outside of the acceptance criteria, reanalyze to confirm and flag with an asterisk (estimated).

11 REFERENCES

- Method 300.0, *Methods for the Determination of Inorganic Substances in Environmental Samples*, EPA/600/R-93/100 Revised August 1993.
- Method 4110 B in *Standard Methods for the Examination of Water and Wastewater*, 18th Ed., 1992.
- NELAC 2002 Standard
- 40CFR Part 136 Appendix B

Ion Chromatography Analysis Training Plan

Procedure: _____

SOP: _____ Revision: _____ Date: _____

Trainee: _____

1. Read SOP Trainer: _____ Trainee: _____ Date: _____

2. Demonstrated understanding of the scientific basis of the analysis including:
 - column separation, retention time
 - suppressor and eluent functions
 - method of detection
 - Correlation co-efficients and results based on calculated curvesTrainer: _____ Trainee: _____ Date: _____

3. Demonstrated familiarity with related SOPs

-ADM-BATCHSEQ	-ADM-PCAL	-ADM-SIGFIG
-ADM-DATAENTRY	-ADM-DIL	-ADM-SPSR
-ADM-MDL	-ADM-DREV	-ADM-TRANDOC

Trainer: _____ Trainee: _____ Date: _____

4. Observe performance of SOP
 - sample preparation (soil, water, other matrices)
 - standard and reagent prep and documentation – including pipet(s) used
 - IC start-up and warm-up procedures
 - software use, entering sample IDs in analytical sequence
 - sample dilution guidelines (1/10 or more unless F, OPO4; 1/50 or more for leachates; other nasties; use of historical data)
 - holding times
 - use and loading of vials (filter if necessary, no air in tubes)
 - use and loading of autosampler
 - sample analysis including:
 - calibration
 - software command of instrument
 - recognition of normal vs. abnormal peaks
 - linear range
 - manual integration
 - use of QC samples and QC criteria
 - IC instrument logbook use
 - data reduction, reporting, and reviewTrainer: _____ Trainee: _____ Date: _____

5. I have read, understood and agree to perform the most recent version of the SOP:

Signature: _____
Date: _____

6. Perform SOP with supervision
 - including all items in 4.Trainer: _____ Trainee: _____ Date: _____

7. Independent performance of the SOP
 - all of the item listed in 4
 - IDC (4 mid-range standards performed before client samples are analyzed)
 - attach IDC certificate, raw data, and summary spreadsheet.Trainer: _____ Trainee: _____ Date: _____



WET CHEMISTRY DATA QUALITY CHECKLIST

Analysis: _____

Date: _____

Yes	No	NA		Yes	No	NA
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	1. Holding Times met method requirements?	1. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	2. ICAL met method requirements?	2. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	3. ICV acceptable?	3. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	4. CCVs acceptable? Analyzed per 10 samples?	4. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	5. CCBs acceptable? Analyzed per 10 samples?	5. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	6. Method Blank results < RL?	6. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	7. LCS recoveries within QC limits?	7. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	8. All reported sample concentrations within LR?	8. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	9. MS recoveries within QC limits?	9. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	10. Duplicate RPD within QC limits?	10. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	11. Dilution factors verified and calculated correctly?	11. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	12. Bench Sheet complete, initials, date, and time:	12. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	•Are standards and reagents traceable?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	•Is unused space on the sheet crossed out?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	•Pipette ID referenced?	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	13. All applicable Log Books filled in?			
			14. Manual data entry to LIMS correct? Date? Time?	14. <input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

Analyst: _____

Peer Review: _____

Date: _____

Date: _____

COMMENTS:

****Comments must be provided for any items noted above as "No"**

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 1 of 14

STANDARD OPERATING PROCEDURE

for

The Determination of Method Detection Limits

SOP Code: ADM - MDL

Revision: 5

August 1, 2003

Approved by:

Lawrence J. Savolny

Quality Assurance Director

8/1/2003

Date

Sam W. ...

Chief Quality Officer

8/1/03

Date

Steve Vincent

President

8/1/2003

Date

© Columbia Analytical Services, Inc., 2003
1317 South 13th Avenue
Kelso, Washington 98626

Annual review of this SOP has been performed
and the SOP still reflects current practice.

Initials: _____ Date: _____

Initials: _____ Date: _____

Initials: _____ Date: _____

DOCUMENT CONTROL

NUMBER: _____

Initials: _____ Date: _____

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 2 of 14

Standard Operating Procedure for The Determination of Method Detection Limits

1.0 PURPOSE

This standard operating procedure (SOP) documents the procedure for the determination of method detection limits (MDLs).

2.0 APPLICABILITY

The procedure described in this SOP is designed for applicability to a wide variety of sample types ranging from reagent (blank) water or wastewater containing the analyte to solids (such as soil) containing the analyte to the analyte in a gaseous matrix. The MDL for an analytical procedure will vary as a function of sample matrix. This SOP requires a complete, specific, and well-defined analytical procedure. It is essential that all sample-processing steps of the analytical procedure are included in the determination of the MDL; that is, all the steps that a sample is processed through, from sample preparation to analytical completion, must be included in the MDL determination. The MDL obtained by this procedure is used to judge the significance of a single measurement of a future sample. This SOP for the determination of MDLs was designed for applicability to a broad variety of physical and chemical methods. To accomplish this, the procedure was made device or instrument-type independent.

3.0 DEFINITIONS

3.1 Method Detection Limit (MDL)

The MDL is the minimum concentration of a substance or analyte that can be measured and reported with 99% confidence that the analyte concentration is greater than zero and is determined from analysis of a sample in a given matrix type containing the analyte.

3.1.1 The Calculated MDL (MDL_C) is the MDL as calculated in Section 6.10 and will typically contain two or more significant figures.

3.1.2 The Reported MDL (MDL_R) is the MDL that is used for reporting purposes. MDLs for organic analytes will be reported with two significant figures.¹ MDLs for inorganic analytes will be reported with either one or two significant figures depending upon the number of significant figures in the analytes' MRL.²

3.2 Analytical Procedure

The written, step-by-step description of the operation by which samples are processed in order to obtain the concentration of an analyte in a sample.

¹ Organic analyte MDLs: see Section 6.3.2 in Reference 9.2.

² Inorganic analyte MDLs: see Section 6.4.2 in Reference 9.2

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 3 of 14

3.3 Spike Level

The spike level is the known concentration of analyte that is added to a matrix for the determination of the MDL.

3.4 Interferences

Interferences are defined as systematic errors in the measured analytical signal of an established procedure caused by the presence of known or unknown species (interferent) that hinder an accurate analysis of the target analyte(s).

3.5 Matrix³

3.5.1 When the matrix analyzed is aqueous (includes ground water, surface water, waste water, drinking water, etc.), analyte-free reagent water is to be used. When the matrix analyzed is solid (includes soil, sand, tissue, or other solid materials), analyte-free soil, sand, tissue, or a suitable material is to be used. When the matrix analyzed is gaseous (i.e., air or emissions), an analyte-free, inert gas (such as zero-grade air or ultrapure helium or nitrogen) is to be used.

3.5.2 If the analysis is performed on a matrix for which there is not available an appropriate or similar, analyte-free matrix (such as, metals analysis on soil samples), the MDL analysis will be done as prescribed by the SOP for the analysis except the sample (weight) will be omitted; that is, the analysis will be done on all the reagents but without addition of any sample.

4.0 DISCUSSION

The MDL is a property of the analytical procedure, sample matrix, and measurement system (e.g., an instrument if one is used in the analytical procedure). The MDL is a statistic. It is an estimate that includes both the systematic and random errors that are an inherent part of the analytical procedure. The MDL for a given analyte will be unique for the sample's matrix and may be different than the MDLs shown in published methods. The MDL actually achieved in a given analysis will vary depending on instrument sensitivity and matrix effects.

The relative uncertainty of an analytical measurement increases as the measured value approaches the MDL and at the MDL the uncertainty in the measured value may be 100% or greater.

³ For a list of matrices, see Section 3.3 in Reference 9.3

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 4 of 14

Regarding Method Proficiency, in Chapter One of the Third Edition of SW-846 (as updated) it states:

Procedures should be in place for demonstrating proficiency with each analytical method routinely used in the laboratory. These should include procedures for demonstrating the precision and bias of the method as performed by the laboratory and procedures for determining the method detection limit (MDL). All terminology, procedures and frequency of determinations associated with the laboratory's establishment of the MDL and the reporting limit should be well-defined and well-documented. Documented precision, bias, and MDL information should be maintained for all methods performed in the laboratory.

This SOP is based upon the procedure described in 40 CFR Part 136, Appendix B (Reference 9.1).

5.0 RESPONSIBILITIES

It is the responsibility of the laboratory manager, working with the quality assurance program manager (QA PM) and department managers and supervisors, to schedule MDL determinations as they come due. It is the responsibility of the QA PM to track the status of MDLs. Completed MDL determinations are to be reviewed by the QA PM and approved by the laboratory manager and/or the QA PM before they are implemented. The QA PM is responsible for maintaining the MDL file described in Section 8.0.

6.0 PROCEDURE

6.1 General Requirements

- 6.1.1 MDLs are to be determined for each analyte and for each matrix. This SOP describes procedures for determining MDLs for the generic matrices aqueous, solid, and gaseous. MDLs for specific matrix types may be adapted from the procedures in this SOP. See Section 3.5.1.
- 6.1.2 All sample processing steps in the analysis procedure shall be included in the determination of the MDL. MDLs shall be generated for all preparatory and cleanup procedures routinely used on samples.
- 6.1.3 An MDL study is required for PCB Aroclors 1016 and 1260 only; i.e., it is not necessary to perform an MDL study for all the PCB Aroclors, unless required by specific clients or accreditation programs.
- 6.1.4 An MDL study is not required for any analyte for which spiking solutions or quality control samples are not available; e.g., temperature.

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 5 of 14

6.2 Frequency of MDL Determination

- 6.2.1 An MDL study shall be determined initially; i.e., when the procedure is first put into production. An MDL study may also be done as part a procedure's training requirements.
- 6.2.2 An MDL study shall be performed at the frequency specified in the applicable method or as specified by an accrediting authority. For example, some state accrediting programs require annual MDL studies.
- 6.2.3 If an MDL study is not performed annually, an MDL verification check shall be performed quarterly⁴ on every instrument used to perform a particular analysis. The MDL verification check sample is spiked at approximately two times the current MDL. The MDL verification check sample shall be acceptable if it produces a response that is at least three times above the instrument's noise level. If the MDL verification check fails, additional MDL verification checks shall be performed at a higher level to set a higher MDL, or a new MDL study shall be performed.
- 6.2.4 A new MDL determination is to be performed "...each time there is a change in the test method that affects how the test is performed, or when a change in instrumentation occurs that affects the sensitivity of the analysis."⁵

6.3 Estimation of the MDL

Use one of the following guides to help estimate the MDL.

- 6.3.1 The concentration value that corresponds to an instrument signal-to-noise ratio in the range of 2.5 to 5.
- 6.3.2 The concentration equivalent of three times the standard deviation of replicate instrumental measurements of the analyte in reagent water.
- 6.3.3 That region of the calibration curve where there is a significant change in sensitivity, i.e., a break in the slope of the calibration curve.
- 6.3.4 Instrumental limitations.

⁴ The quarterly MDL verification checking procedure is based on the procedure in Reference 9.7, Section D.1.4, Clarification Box D-12.

⁵ See Reference 9.6, Section D.1.2.c) and Reference 9.7, Section D.1.4.c).

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 6 of 14

6.4 Aqueous Blank MDLs

- 6.4.1 Prepare reagent (blank) water that is as free of analyte as possible. Reagent or interference free water is defined as a water sample in which analyte and interferent concentrations are not detected at or above the MDL of each analyte of interest.
- 6.4.2 Prepare a minimum of 7 (preferably 8 to 12) analyte-spiked reagent water samples at a concentration that is 3 to 5 times the estimated MDL from Section 6.3.
- 6.4.3 Analyze the analyte-spiked reagent water samples prepared in Section 6.4.2 by processing them through the **entire** analytical procedure. Make all computations according to the directions prescribed in the analytical procedure with the final results reported in the same units as used for water samples.
Proceed to Section 6.10.

6.5 Aqueous Sample MDLs

- 6.5.1 Analyze the aqueous sample by processing it through the **entire** analytical procedure.
- 6.5.2 Calculate the analyte concentration.
 - 6.5.2.1 If the measured concentration of the analyte is in the recommended range of 3 to 5 times the estimated MDL from Section 6.3, proceed to Section 6.5.3.
 - 6.5.2.2 If the measured concentration of the analyte is less than the recommended 3 to 5 times the estimated MDL, add a known amount of analyte to bring the concentration of analyte between 3 to 5 times the estimated MDL and proceed to Section 6.5.3.
 - 6.5.2.3 If the measured concentration of the analyte is greater than 5 times the estimated MDL, either obtain another sample with a lower concentration of analyte in the same matrix, or the sample may be used as is for determining the MDL if the analyte concentration does not exceed 10 times the MDL of the analyte in reagent water. The variance of the analytical procedure changes as the analyte concentration increases from the MDL; hence the MDL determined under these circumstances may not truly reflect method variance at lower analyte concentrations. Proceed to Section 6.5.3.
- 6.5.3 Prepare and analyze a minimum of 7 (preferably 8 to 12) aliquots of the aqueous sample by processing them through the **entire** analytical procedure. Make all computations according to the directions prescribed in the analytical procedure

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 7 of 14

with the final results reported in the same units as used for water samples.
Proceed to Section 6.10.

6.6 Solid Blank MDLs

- 6.6.1 Prepare a solid material (e.g., soil, sand, tissue, Na_2SO_4 , Teflon chips, or other appropriate material) that is free of analyte.
- 6.6.2 Prepare a minimum of 7 (preferably 8 to 12) analyte-spiked solid samples at a concentration that is 3 to 5 times the estimated MDL from Section 6.3. The same weight of analyte-spiked solid is substituted for the sample weight in the analytical procedure.
- 6.6.3 Analyze the analyte-spiked solid samples prepared in Section 6.6.2 by processing them through the **entire** analytical procedure. Make all computations according to the directions prescribed in the analytical procedure with the final results reported in the same units as used for solid samples. Proceed to Section 6.10.

6.7 Solid Sample MDLs

- 6.7.1 Analyze the solid sample by processing it through the **entire** analytical procedure.
- 6.7.2 Calculate the analyte concentration.
 - 6.7.2.1 If the measured concentration of the analyte is in the recommended range of 3 to 5 times the estimated MDL from Section 6.3, proceed to Section 6.7.3.
 - 6.7.2.2 If the measured concentration of the analyte is less than the recommended 3 to 5 times the estimated MDL, add a known amount of analyte to bring the concentration of analyte between 3 to 5 times the estimated MDL and proceed to Section 6.7.3.
 - 6.7.2.3 If the measured concentration of the analyte is greater than 5 times the estimated MDL, either obtain another sample with a lower concentration of analyte in the same matrix, or the sample may be used as is for determining the MDL if the analyte concentration does not exceed 10 times the MDL of the analyte in soil. The variance of the analytical procedure changes as the analyte concentration increases from the MDL; hence the MDL determined under these circumstances may not truly reflect method variance at lower analyte concentrations.
Proceed to Section 6.7.3.

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 8 of 14

- 6.7.3 Prepare and analyze a minimum of 7 (preferably 8 to 12) aliquots of the soil sample by processing them through the **entire** analytical procedure. Make all computations according to the directions prescribed in the analytical procedure with the final results reported in the same units as used for solid samples. Proceed to Section 6.10.

6.8 Gaseous Blank MDLs

- 6.8.1 Using an appropriate sample container (e.g., Tedlar® bag or SUMMA® passivated canister) and appropriate analyte-free inert gas (such as zero-grade air or ultrapure nitrogen), prepare a minimum of 7 (preferably 8 to 12) analyte-spiked inert gas samples at a concentration that is 3 to 5 times the estimated MDL from Section 6.3.
- 6.8.2 Analyze the analyte-spiked inert gas samples prepared in Section 6.8.1 by processing them through the **entire** analytical procedure. Make all computations according to the directions in the analytical procedure with the final results reported in the same units as used for air samples. Proceed to Section 6.10.

6.9 Rejection of Replicate Sample Results

- 6.9.1 A replicate sample result may only be rejected if there is an assignable cause for not using that result. Assignable causes include, but are not limited to, replicate sample preparation error, instrument malfunction, bad injection or purge, and internal standard(s) missing or response uncharacteristically high or low. The cause for rejecting the replicate sample result must be documented in the MDL data package.
- 6.9.2 For multi-analyte analyses, if a replicate sample result is rejected for an assignable cause, results for all the analytes from that sample are to be rejected; that is, “picking and choosing” analyte results from a sample is not permitted.

6.10 Calculation of MDL_C

- 6.10.1 Determine the standard deviation, s , of the replicate sample results.

$$s = \sqrt{\frac{\sum_{n=1}^i (x_i - \bar{x})^2}{n - 1}}$$

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 9 of 14

$$\text{where } \bar{x} = \frac{\sum_{n=1}^i x_i}{n}$$

6.10.2 Multiply the standard deviation obtained in Section 6.10.1 times the appropriate one-sided 99% Student's t-statistic, which is found in the following table.

$$MDL_C = s \times \{\text{appropriate Student's t-statistic}\}$$

No. of Samples (n)	Student's t-statistic	Degrees of Freedom (n - 1)
7	3.143	6
8	2.998	7
9	2.896	8
10	2.821	9
11	2.764	10
12	2.718	11
13	2.681	12
14	2.650	13
15	2.624	14
16	2.602	15
17	2.583	16
18	2.567	17
19	2.552	18
20	2.539	19
21	2.528	20

6.11 Determination of MDL_R

The Reported MDL (MDL_R) is the calculated MDL rounded **up** to the appropriate number of significant figures. See Section 3.1.2.

6.12 Evaluation of the Quality of the MDL Study

The quality of the MDL is evaluated using the following criteria.

6.12.1 Spike Level The spike level is **too low** if the MDL_C is greater than the spike level.
The spike level is **too high** if the spike level is greater than ten times the MDL_C.

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 10 of 14

- 6.12.2 Percent Relative Standard Deviation (%RSD) The %RSD should be some value close to 20, where the %RSD is equal to the standard deviation (s) divided by the average of the spike recoveries times 100. $[\%RSD = (s \div \bar{x}) 100]$
- 6.12.3 Percent Spike Recovery The spike recovery should be approximately what is to be expected for that analyte from the analytical procedure; i.e., a 40% spike recovery for an analyte is too low if the method normally recovers 80% or more for that analyte.
- 6.12.4 MDL Quality The criteria in Section 6.12.1 must be true. At least one of the criteria in Sections 6.12.2 and 6.12.3 should also be true. If the MDL_C does not meet these criteria, then the study should be repeated, adjusting the spike level appropriately.

6.13 Instruments

If more than one instrument is used for the same analytical procedure, the replicate samples should be analyzed on each instrument to ensure there is no instrument bias. Under some specific customer contracts and for some programs (such as the Navy's Installation Restoration program), instrument-specific MDLs are required. There are two options for complying with this requirement:

1. Analyze the replicate samples on each instrument used for the analytical procedure and calculate the MDL_C for each instrument. The MDL_R will be the largest of the several MDL_C's; or
2. Analyze the replicate samples on each instrument used for the analytical procedure and calculate a single MDL_C using all the values from each instrument. A minimum of five values is needed from each instrument. For example, if two instruments are used, there would be a minimum of two times five or ten values to be used to calculate the MDL_C. Make sure to use the appropriate Student's t-statistic that corresponds to the number of values used to calculate the standard deviation. Note: This option may not be acceptable under some specific customer contracts or for some programs, such as the DOD quality systems for environmental laboratories.⁶

6.14 Review and Approval

Completed MDL determinations are to be reviewed by the supervisor of the analysis. The QA PM will review and approve the MDL determination before it is implemented.

⁶ See Reference 9.7, Section D.1.4.

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 11 of 14

6.15 Department of Defense (DoD) Requirements⁶

- 6.15.1 An MDL verification check shall always be performed immediately following an MDL study. DoD requires that the MDL check sample be spiked at approximately 2 times the current reported MDL.
- 6.15.2 If an annual MDL study is not performed, MDL verification checks shall be performed quarterly. If the quarterly MDL verification check fails, additional MDL verification checks shall be performed at a higher level to set a higher MDL, or the MDL study shall be reconducted.
- 6.15.3 For DoD, the MDL verification check sample shall be acceptable if it produces a response that lies at least 3 times above the instrument's noise level.

7.0 QUALITY ASSURANCE/QUALITY CONTROL REQUIREMENTS

7.1 Replicate Samples

No fewer than 7 replicate samples can be used; 8 to 12 replicate samples is preferred.

7.2 Analysis of the Replicate Samples

The replicate samples do not have to all be analyzed in the same analytical batch on the same day. In fact, it is preferred to spread out the replicate samples among several analytical batches analyzed on several days (to increase the contribution of the day-to-day variability). Furthermore, it is recommended that at least one MDL spike be routinely analyzed monthly and data accumulated and calculated at a later time.

7.3 MDL Quality

The MDL determination must meet the criteria in Section 6.12. If the MDLs from more than one instrument are combined as in Section 6.13, the combined MDL must meet the criteria in Section 6.12.

7.4 Matrices

MDLs shall be generated for all applicable matrices. See Section 6.1.1.

7.5 Preparatory and Clean-up Procedures

MDLs shall be generated for all preparatory and clean-up procedures routinely used on samples. See Section 6.1.2.

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 12 of 14

8.0 RECORDS

The data for the MDL determination is summarized in a table similar to the one shown in Figure 1. An Excel spread sheet similar to Figure 1 is available for this purpose. This summary and the reference to the location of the raw data are to be filed in a readily available file of MDLs. This file is to be located both in the department performing the analytical procedure and in a centralized location for MDLs from the entire laboratory. Also shown in Figure 1 are two examples illustrating how MDL data is to be summarized.

9.0 REFERENCES

- 9.1 *40 CFR Part 136, Appendix B, Definition and Procedure for the Determination of the Method Detection Limit--Revision 1.11.*
- 9.2 *SOP for Significant Figures, ADM-SIGFIG.*
- 9.3 *SOP for Sample Batches, ADM-BATCH.*
- 9.4 *Test Methods for Evaluating Solid Waste, Physical/Chemical Methods, SW-846, Third Edition, September 1986 and as amended by Updates I, II, IIA, IIB, III, and IIIA.*
- 9.5 *Standard Methods for the Examination of Water and Wastewater, APHA/AWWA/WEF, 19th Edition, 1995, Method 1030E; 20th Edition, 1998, Method 1030C.*
- 9.6 National Environmental Laboratory Accreditation Conference (NELAC), Quality Systems Standard, Appendix D, Section D.1.2.
- 9.7 Department of Defense *Quality Systems Manual for Environmental Laboratories*, Final Version 2, June 2002, Appendix D, Section D.1.4.

UNCONTROLLED COPY

SOP Code: ADM - MDL
Revision: 5
Date: August 1, 2003
Page 13 of 14

10.0 CHANGES FROM PREVIOUS REVISION

- Section 3.1.1 Section 6.10 cross reference corrected
- Section 3.3 Corrected Section number
- Section 3.4 Corrected Section number
- Section 3.5 Corrected Section number
- Section 3.5.1 Changed matrices to generic matrices aqueous, solid and gaseous
- ★ Section 6.1 Section completely revised
- ★ Section 6.2 New section – causing subsequent sections to be renumbered and section cross-references to be revised
- Sections 6.3 through 6.15 Sections renumbered and internal cross references updated
- Section 6.4 “Water” changed to “Aqueous” to be consistent with Section 3.5.1
- Section 6.5 “Water” changed to “Aqueous” to be consistent with Section 3.5.1
- Section 6.7 “Soil” changed to “Solid” to be consistent with Section 3.5.1
- ★ Section 6.12.1 MDL_R changed to MDL_C
- ★ Section 6.12.4 MDL_R changed to MDL_C
- Section 6.13 “Navy’s Installation and Restoration program” changed to “DOD quality systems for environmental laboratories” at end of paragraph 2.
- Section 7.4 Cross reference changed to Section 6.1.1
- Reference 9.7 Updated

UNCONTROLLED COPY

Revision: 5

Page 14 of 14

Figure 1

MDL Determination Summary

Instrument: SVM GC/MS No. 03

Extraction/Digestion Method: 3520CMatrix: Water / ~~Soil~~ / ~~Air~~ Units: ug/L (ppb)

Approved by: _____ Date: _____

[illegible]² Pentachlorophenol

³ Since these are organic analytes, the MDL is rounded up to two significant figures, per Section 3.1.2.

STANDARD OPERATING PROCEDURE

SAMPLE RECEIVING

SMO-GEN

Revision 2
December 8, 2004

Approved By: _____

Supervisor

_____ Date

_____ QA Coordinator

_____ Date

_____ Laboratory Manager

_____ Date

©COLUMBIA ANALYTICAL SERVICES, INC. 2004
One Mustard St., Suite 250
Rochester, NY 14609

Annual review of this SOP has been performed
and the SOP still reflects current practice.

Initials: _____ Date: _____

Initials: _____ Date: _____

Initials: _____ Date: _____

DOCUMENT CONTROL

NUMBER: _____

Initials: _____ Date: _____

1 SCOPE AND APPLICATION

This Sample Management SOP provides a key foundation to the SMO department. It explains the process of receiving samples and the steps that lead to the distribution of samples throughout the lab. By implementing an organized and thorough approach to the initial stages of the process, we can maintain an efficient and well-documented account of the samples status and condition.

2 METHOD SUMMARY

The process of receiving and distributing samples is outlined in this SOP. Upon receipt, a CR/PF is completed for each cooler. The sample information is logged into the LIMS database and then broken down and distributed within the lab.

3 DEFINITIONS

SMO	Sample Management Office
PM	Project Manager
CR/PF	Cooler Receipt & Preservation Form
QA/QC	Quality Assurance / Quality Control
LIMS	Laboratory Information Management System
COC	Chain of Custody
BREAKDOWN	The act of removing samples from coolers, labeling the samples, checking preservations, and distributing them to the correct destinations.
LOG-IN	Entering the information from the COC, CRPF, and the specifics of the job into the LIMS
MATRIX	The physical form of a sample. (water, oil, soil, solid, air)
NELAC	National Environmental Laboratory Accreditation Conference

4 INTERFERENCES

To avoid errors during log-in, detailed information shall be obtained from clients to complete chain of custody documentation.

5 SAFETY

- 5.1 The characteristics of incoming samples are often unknown. Treat all samples as potentially hazardous. SMO has first contact with samples and must be especially cautious.
- 5.2 All appropriate safety precautions for handling reagents and samples must be taken when performing this procedure. This includes the use of personnel protective equipment, such as, safety glasses, lab coat and the correct gloves.

5.3 Sodium Hydroxide (NaOH) is a strong caustic and a severe health and contact hazard. Use nitrile or latex gloves while handling pellets or preparing solutions.

5.4 Hydrochloric and Nitric Acid are used in this method. These acids are extremely corrosive and care must be taken while handling them. A face shield should be used while pouring acids. And safety glasses should be worn while working with the solutions. Lab coat and gloves should always be worn while working with these solutions.

5.5 Refer to the Safety Manual for further discussion of general safety procedures and information.

6 SAMPLE CONTAINERS, COLLECTION, PRESERVATIONS, AND STORAGE

Sample preservation and storage are discussed as part of the procedure later in this SOP.

7 APPARATUS AND EQUIPMENT

- Infrared or digital thermometer – calibrated and maintained as per ADM-DALYCK.

8 PREVENTIVE MAINTENANCE

Not Applicable

9 STANDARDS, REAGENTS, AND CONSUMABLE MATERIALS

9.1 The following reagents are purchased commercially, stored at room temperature and expire upon manufacturer's indications or 3 years from receipt if no other indication is given:

Sulfuric Acid; Instranalyzed grade
Nitric Acid; Instranalyzed grade
Hydrochloric Acid; Instranalyzed grade
Sodium Hydroxide; Lab grade

9.2 Consumable materials:
PH indicator Paper
Potassium-Iodide Starch Paper for detection of residual chlorine

10 RESPONSIBILITIES

It is the responsibility of the analyst to receive samples according to this SOP and to complete all documentation required.

11 PROCEDURE

11.1 SAMPLE ACCEPTANCE POLICY

11.1.1 The laboratory's sample acceptance policy outlines the circumstances under which samples shall be accepted or rejected. This information is made available to the sampler by an instruction sheet that accompanies each bottle set sent to the client (attached).

11.1.2 The samples received need to conform to the following acceptance criteria as per NELAC:

- Proper, full, and complete documentation (COCs), which shall include sample identification, the location, date and time of collection, collector's name, preservation type, sample type and any special remarks concerning the sample;
- Proper sample labeling to include unique identification and a labeling system for the samples with requirements concerning the durability of the labels (water resistant) and the use of indelible ink;
- Use of appropriate sample containers;
- Adherence to specified holding times;
- Adequate sample volume. Sufficient sample volume must be available to perform the necessary tests; and
- Procedures to be used when samples show signs of damage, contamination, and inadequate preservation.

The above criteria are addressed in the rest of this section. In the event of an unacceptable sample, the Project Manager is notified and they will contact the client or recommend a proper course of action. Any data from samples which do not meet the acceptance policy must be written up in the case narrative in the report.

11.2 PROCEDURES FOR SAMPLE RECEIPT

11.2.1 The CAS Cooler Receipt and Preservation Form (CR/PF) (attached) is used for the next steps to document the condition of the samples and coolers as per the Acceptance Policy.

11.2.2 Upon receipt, the coolers are examined for presence and condition of custody seals, locks, shipping bills, etc. The shipping labels are removed and placed on scrap paper and added to the receiving paper work.

11.2.3 The coolers are opened and examined for any existing hazards before subsequent processing.

CAUTION: *If samples exhibit any strong odors, or samples have been damaged, move cooler to the hood and continue processing per client.*

11.2.4 Chain of Custody (COC) forms (attached) and any other documents are located, removed and signed with date & time as received. Visually scan the COC for short holding time samples.

11.2.5 The temperature of the cooler is measured following the guidelines in this SOP. The acceptance criteria for samples is 0-6°C. If a cooler exhibits a temperature greater than 6°C, or exhibits any other anomalies (deviations from the Acceptance Policy), the anomalies are noted on the CR/PF and the CR/PF is placed on top of the COC packet.

11.2.5.1 Samples which are hand delivered immediately after collection (delivered within 4 hours of sampling) may not have had time to cool. The samples are considered acceptable by NELAC if there is evidence that the chilling process has begun. Document the presence of ice on the CR/PF with a note about the 4-hour rule.

11.2.6 Once the coolers are initially examined and observations and temperature are recorded, all of the COCs with corresponding CR/PF forms and shipping labels are submitted for review to the appropriate Project Manager. (The receiving paper work is comprised of at least 3 pages; the COC, a CR/PF, and a shipping label).

11.2.7 Rush requests and samples with short holding times are always given top priority for initial processing. CAS follows EPA guidelines for preservation and holding time as outlined in our QA/QC manual Table 7-1. A list of short holding time parameters are attached to this SOP. An additional list of holding times may be found in SMO-BPS. If short holding time samples need to be distributed immediately (before log-in and labeling), distribute the samples with the attached form for Internal Tracking. Write all short holding time samples on the white board in Wetchem.

11.2.8 In the down time between receiving and actual breakdown of samples, the coolers are stored in the SMO walk-in cooler, which also maintains a temperature of 0-6°C.

11.3 LOG-IN/LIMS

- 11.3.1 At log-in, the Project Manager enters the entire job into the LIMS system. The job is given a submission number, which consists of the lab location, year, and a unique sequential number. (R20-1508)
- 11.3.2 Each sample is given a unique order number during log-in. This number is unique to a given location or sample site. A single sample site may consist of many bottles according to the analyses requested. The individual bottles within an order number are uniquely identified with the use of a bar code placed on the bottles during sample breakdown. See SMO-ICOC.
- 11.3.3 COC's are returned to SMO after the Project Manager has approved all anomalies (See ADM-PCR) and entered the job into the LIMS database. At this point, the jobs are approved for breakdown.

11.4 SAMPLE BREAKDOWN

- 11.4.1 Sample containers are removed and organized according to chain of custody identification and analysis.
- 11.4.2 The following verifications are made as to the agreement of chain of custody information as it applies to samples and containers received:
- Sample identification, time, date
 - Number of containers received
 - Matrix
 - Correct bottles for analysis requested.
 - Correct sample volume for analysis
 - Correct preservatives for analysis – according to the labels. The actual preservation will be tested as below

Any discrepancies are reported to the Project Manager. Tests may be added or deleted so that LIMS matches the actual samples received. See the discrepancies section of this SOP (i.e. jobs can not have tests scheduled when the sample containers do not exist.)

- 11.4.3 Labels are printed from the LIMS database and placed on the sample containers.
- 11.4.4 Barcodes are generated which are unique to each container for the purpose of sample tracking. VOA bags receive one barcode and each vial within the bag must be labeled with a different number (if not already done so when preparing the bottle set)

11.4.5 Preservations are checked on the appropriate samples and recorded on the job paperwork, the preservation check log, and preservative lot numbers are recorded on the front page of the paper work (the Analytical Request).

11.4.6 To check the preservation, place a small piece of pH test paper in a dispo-cup and pour a small aliquot of the sample into a dispo-cup. Observe the color of the paper and compare to the colors on the paper dispenser to determine the pH. The preservation of samples should be as below:

HCl, HNO₃, H₂SO₄ pH < 2

NaOH pH > 12

PCB/608 pH 5-9

If the sample was not sufficiently preserved, notify the Project Manager to determine whether more preservative should be added.

11.4.7 To check the chlorine residual, place a small piece of starch paper in a dispo-cup and pour a small aliquot of the sample into a dispo-cup (this may be done on the same aliquot used to test pH). If the paper turns blue there is chlorine residual present in the sample. Note the discrepancy on the CR/PC and notify the Project Manager. The Project Manager may contact the client to determine if ascorbic acid should be added to eliminate the chlorine residual.

11.4.8 The CR/PC form is finished by the person who broke down the job.

11.4.9 All of the jobs are reviewed for completeness at the end of the day (or the following morning), and the walk-in cooler temperatures are logged into a temperature logbook.

11.5 SAMPLE DISTRIBUTION

11.5.1 After a job is labeled, the samples are distributed to the appropriate department. The samples are scanned into the appropriate storage areas as listed below.

11.5.2 CAS-Rochester currently maintains 3 walk-in coolers to refrigerate samples. Extractables share a cooler with Metals, VOAs have a cooler, and WetChem shares the cooler with SMO. Metals are maintained at room temperature and are placed on a dedicated cart in the metals department. Mercury and TCLP samples are placed in the Metals/Extractables cooler. SMO is responsible for documenting the location of the Wetchem samples in the Wetchem/SMO cooler.

11.6 SAMPLE SECURITY AND STORAGE:

11.6.1 Bar-coding is used as a means of sample tracking. Custody is maintained according to SMO-ICOC.

11.6.2 The sample coolers are secured with locks, which can be accessed by technical laboratory personnel, project managers, administrative support personnel and all senior staff members. All of the sample storage facilities are located in our building which is a secured area. Storage areas are kept clean and dry to avoid any damage or deterioration of samples while in storage.

11.6.3 Samples are held in refrigerators (if applicable) until analysis is completed and reported to the Client. Routine samples are typically held for 30 days after mailing of report and CLP samples are stored for 90 days after report has been mailed.

11.6.4 Refrigeration is maintained at a temperature of 0-6° Celsius.

11.7 DISCREPANCIES

11.7.1 Any discrepancies or concerns such as non-matching identifications, missing samples, and tests not scheduled correctly are to be verbally communicated to the Project Manager. Any action taken is recorded on the COC and/or cooler receipt form, "as per" Project Manager. The Project Manager will make any contact to the client when they deem it necessary. If a Chain of Custody is not received, the Project Manager is informed and a CAS COC is filled out.

11.7.2 In the event of broken samples, a note is entered on the Chain of Custody and/or the CR/PF accompanying the samples. Information pertaining to the sample is forwarded to the Project Manager for follow up purposes. Cleanup procedures are as follows:

- Liquids: Broken glass is handled carefully using disposable gloves and disposed of in the Glass Disposal Box. The liquid is disposed of in the SMO sink or under a hood if strong odors are evident. Any packing material is disposed of appropriately.
- Soils: The same documentation as liquids applies. Broken glass is disposed of in the Glass Disposal Box and the soil is disposed of into the garbage.

11.8 RECEIVING SAMPLE COOLERS ON WEEKENDS OR AFTER HOURS

11.8.1 The date received is the date on which the Laboratory Personnel takes possession of the samples. The client shall sign the COC as relinquished and the CAS employee shall sign in the adjacent box as received. If an employee outside of the SMO department receives the samples, the coolers or samples

are stored in a locked walk-in cooler and all paperwork is left on the lab bench in SMO for login the next working day.

11.8.2 If the samples are received by SMO staff but the samples cannot be processed through the complete log-in procedure on the date received, the receiving procedure is performed as outlined in 11.1 and 11.2. Then SMO needs to check for short holding time samples. Short holding time tests are posted in breakdown. Most of the tests have holding times of 48 hours which means samples received on Saturday need to be run before Monday. Any cooler integrity issues shall be handled on the next business day, but the samples need to be tested so that holding times are maintained. To maintain an organized system, notes are indicated on the COC and on the white board in Wetchem as to which samples have been sent for analysis. Use the attached form for Internal Tracking of Short holding time samples.

11.9 SAMPLE SHIPPING TO SUBCONTRACT LABS

11.9.1 The sample is logged in and the test code for subcontracted analysis is assigned.

11.9.2 For CAS Network Labs: Subcontracted samples are shipped to the network lab with a copy of the Internal Service Request form (ISR) and copy of the COC. The tests being subcontracted must be highlighted and the number of containers adjusted or a new COC is completed.

11.9.3 For other labs: A purchase order (PO) is filled out for work going to another contract laboratory. TAT, deliverables, etc. should be clearly specified. A new Chain of Custody form is filled out with the pertinent information so the samples can be analyzed in a fashion that meets the client's needs.

11.9.4 Samples are prepared for shipping by packing in bubble wrap, and ice. Temperature blanks, and the chain of custody are placed in shipping coolers. Custody seals are signed and dated and placed on the front of the cooler. The cooler is then sealed with packaging tape and shipped overnight through the courier system (confirm with Project Manager for am or pm delivery requirements).

11.10 THERMOMETER MEASUREMENTS

11.10.1 Unless unavailable, measure the cooler/sample temperature of an incoming cooler with the Infrared (IR) thermometer. Turn on the thermometer and point it at the temp blank or a sample (preferably clear glass or amber glass). Wipe the container with a dry paper towel. Hold the thermometer approximately 3-6 inches from the container and at least 3 inches above the counter. Temperature is rounded to the nearest whole number and recorded to the nearest whole

number on the cooler receipt form and COC. This method is preferred to the digital thermometer method.

11.10.2 If the IR thermometer is not available, a digital thermometer may be used. Place the thermometer in the temperature blank or plunge it into the packing material or place it as deep into the cooler as is practical with the lid closed. Allow to equilibrate 5 minutes. After measurement, the temperature is recorded to the nearest whole number in the appropriate space on the cooler receipt form and COC.

12 QA/QC REQUIREMENTS

12.1 Not Applicable

13 DATA REDUCTION AND REPORTING

13.1 All samples, custody documents and discrepancy forms must be clearly completed with permanent ink and filed with the project folder.

14 WASTE MANAGEMENT AND POLLUTION PREVENTION

14.1 Not applicable.

15 REFERENCES

15.1 *Test Methods for Solid and Hazardous Waste Physical and Chemical Analyses*, USEPA SW846, December 1996.

15.2 *NYSDEC Analytical Services Protocol*, October 1995.

15.3 NELAC Standard, Chapter 5, July 2002.

16 TRAINING OUTLINE

- Read this SOP.
- Follow policies in ADM-TRANDOC.
- Observe performance of this SOP. Follow Breakdown Training Plan Form.
- Perform this SOP with guidance.
- Perform this SOP independently and have a trained analyst check the trainee's work. If work is acceptable, complete Training Plan Form and file with QA.

17 INSTRUMENT-SPECIFIC ADDENDUM

Not Applicable

18 ATTACHMENTS

- 18.1 Sampling Instructions
- 18.2 Short Holding Time Parameters List
- 18.3 Cooler Receipt and Preservation Form
- 18.4 COC: Chain of Custody
- 18.5 Internal Tracking Form for Short Holding Time Samples.

19 CHANGES FROM PREVIOUS REVISION

- 3 – Added Definition of Log-In
- 5 – Modified Safety for consistency with other SOPs
- 9 – Added pH paper and KI Starch paper to supplies
- 11.2.4 Added need to check COC for shorties
- 11.2.5 Eliminated that discrepancies are highlighted in yellow on COC
- 11.2.7 Added procedure if shorties need done before log-in
- 11.3.2 Added reference to SMO-ICOC
- 11.3.3 Added reference to ADM-PCR
- 11.4.2 Moved preservative requirements out of this section and into (new) 11.4.6+.7
- 11.4.6 +.7 Added detail of how to check preservatives
- 11.5 Simplified wording – eliminated the holding cart. Added that SMO responsible for documenting the location of samples in cooler
- 11.6 Eliminated wording about paper chains and phasing in barcoding
- 11.8.2 Changed white board from SMO to WC
- 11.9.4 Changed from UPS to “courier”
- 11.10.1 Changed the use of the IR gun distance from 7-12 inches to 3-6 inches. Added need to wipe container with a dry paper towel when using IR gun.
- 15 Changed NELAC reference from 1999 to 2002.

SAMPLING INSTRUCTIONS

1. Please Use Caution! Some bottles contain Preservatives such as Acids and Bases that are **CORROSIVE!** These bottles are marked with colored stickers to inform the sampler but gloves should be worn at all times during sampling for your own protection. (preservative information is on the flip side for emergencies)
2. Please overfill VOA vials and stoppered bottles to eliminate headspace that interferes in the analysis of the samples – but do not overfill too much or the preservative will be flushed out.
3. Please use all of the containers that we provide for a sampling location. This will ensure that the proper containers and volumes are returned for the tests that were requested.
4. Vials labeled TRIP BLANK or TB are included in each cooler containing vials used for testing volatile organics. These vials are a control to determine if the cooler was exposed to contamination while in route or during sampling. Please indicate on the Chain of Custody if you wish to have these vials tested.

Returning Coolers Checklist – Please Read!

Incomplete information will result in a phone call and hence delayed processing of your samples.

1. **Labeling** all bottles or soil jars is essential. We look for a **location ID**, a **date**, a **time**, **initials**, **preservation**, and **sample type** on all bottles to verify a location against the Chain of Custody. (Remember that ice water can smudge or remove your markings. a permanent marker on a dry bottle works best)
2. The **Chain of Custody** should include: the **client information and sampler's signature** in the top left box, the **location ID's** of the sample with the **analysis and preservation** indicated in the center section, **turn-around-time and reporting** information in the bottom center section, and **sign-off** (signature/date/time) to the courier on the bottom left section.
3. **Custody stickers** are very important for CLP or ASP package work and are encouraged for all coolers. The sticker should be placed over the lid and body of the cooler and taped securely. These stickers provide an added level of security, and if they are broken when the coolers arrive at CAS, we will contact you.
4. **Packaging coolers** well is a key to avoid resampling. Beware of glass and make sure they are in the bubble bags that we provide. Also a snug fit is suggested. Extra paper or cardboard (especially with amber liters) is encouraged. **Note** that when the ice melts, the bottles can move inside of the cooler. **Any Leaking Coolers in shipment will be considered HAZARDOUS** by the courier (UPS). Please seal the cooler with packing tape and/or place samples and ice in a plastic bag in the cooler. A leaking cooler will likely result in a delay at UPS and missed holding times at the LAB.
5. **Receiving Temperature at CAS** is a key element to the validity of your results to withstand scrutiny in a court of law. Our Data is only considered valid if the samples have a temperature of 6 degrees Celsius or less. A temperature blank is include in all coolers and should be returned in ice with the other samples. Ice should be bagged **or** put the ice and samples in a large plastic bag and tie it off to reduce leakage. Submersion is the best way to cool the samples but watch out for labels falling off or smudging.
6. **Your supplies** (like coolers or icepacks) will gladly be returned if we have complete return-address information! Please document **your return-address on the cooler** or ice packs in the form of a sticker, or permanent marker so we can return your supplies.
7. **Ship samples** using overnight service or deliver within 24 hours of sampling time. If shipping for **Saturday**, Check mark the Saturday-delivery Box and you must affix several "Saturday" stickers to the cooler.

SHORT HOLDING TIMES

BREAK DOWN:

	<u>HOLDING TIME</u>	<u>WHEN</u>
pH	ASAP	Run daily
Conductivity	48hr from sample date	Run daily
Color	48hr from sample date	Run as needed
Turbidity	48hr from sample date	Run as needed
Sett. Solids	48hr from sample date	Run Same day

WET CHEM:

BOD	Wed./Friday	48hr from sample date
NO3	ASAP (IC)	48hr from sample date
NO2	ASAP (lachar)	48hr from sample date
Residual Chlorine	ASAP	
Odor	ASAP	
Dissolved Oxygen	ASAP	
Coliform test	ASAP (Bev)	
Orthophosphate	ASAP	48hr from sample date
Chrome Hex	ASAP	24hr from sample date
Surfactants	ASAP	48hr from sample date
Ferrous Iron	ASAP	24hr from sample date
TSS,TDS,TS,TVS	(ck sample date 7 days)	
Sulfite	ASAP	
Sulfides	(ck sample date 7 days)	

Cooler Receipt And Preservation Check Form

Project/Client _____ Submission Number _____

Cooler received on _____ by: _____ COURIER: CAS UPS FEDEX CD&L CLIENT

- | | | | | |
|----|--|-----------------|-------|-------|
| 1. | Were custody seals on outside of cooler? | YES | NO | |
| 2. | Were custody papers properly filled out (ink, signed, etc.)? | YES | NO | |
| 3. | Did all bottles arrive in good condition (unbroken)? | YES | NO | |
| 4. | Did any VOA vials have significant air bubbles? | YES | NO | N/A |
| 5. | Were Ice or Ice packs present? | YES | NO | |
| 6. | Where did the bottles originate? | CAS/ROC, CLIENT | | |
| 7. | Temperature of cooler(s) upon receipt: | _____ | _____ | _____ |

Is the temperature within 0° - 6° C?: Yes Yes Yes Yes Yes

If No, Explain Below No No No No No

Date/Time Temperatures Taken: _____

Thermometer ID: 161 or IR GUN Reading From: Temp Blank or Sample Bottle

If out of Temperature, Client Approval to Run Samples _____

Cooler Breakdown: Date: _____ by: _____

- | | | | |
|----|---|-----|-----|
| 1. | Were all bottle labels complete (i.e. analysis, preservation, etc.)? | YES | NO |
| 2. | Did all bottle labels and tags agree with custody papers? | YES | NO |
| 3. | Were correct containers used for the tests indicated? | YES | NO |
| 4. | Air Samples: Cassettes / Tubes Intact Canisters Pressurized Tedlar® Bags Inflated | | N/A |

Explain any discrepancies: _____

		YES	NO	Sample I.D.	Reagent	Vol. Added
pH	Reagent					
12	NaOH					
2	HNO ₃					
2	H ₂ SO ₄					
Residual Chlorine (+/-)	for TCN & Phenol					
5-9**	P/PCBs (608 only)					

YES = All samples OK

NO = Samples were preserved at lab as listed

PC OK to adjust pH

**If pH adjustment is required, use NaOH and/or H₂SO₄

VOC Vial pH Verification (Tested after Analysis) Following Samples Exhibited pH > 2			

Other Comments:



ATTACHMENT A-2

GENERAL ENGINEERING LABORATORIES, LLC STANDARD OPERATING PROCEDURES

GL-RAD-A-013, Rev. 10	The Determination of Gamma Isotopes
GL-RAD-A-045, Rev. 1	The Isotopic Determination of Plutonium, Uranium, Americium, Curium, and Thorium
GL-SR-E-001, Rev. 17	Sample Receipt, Login and Storage

VERIFY THE VALIDITY OF THIS SOP EACH DAY IN USE

**STANDARD OPERATING PROCEDURE
FOR
THE DETERMINATION
OF
GAMMA ISOTOPES**

(GL-RAD-A-013 REVISION 10)

APPLICABLE TO METHODS:
EPA 600/4-80-032 Method 901.1 (Modified)
DOE EML HASL-300 (Modified)

PROPRIETARY INFORMATION

This document contains proprietary information that is the exclusive property of General Engineering Laboratories, LLC (GEL). No contents of this document may be reproduced or otherwise used for the benefit of others except by express written permission of GEL.



TABLE OF CONTENTS

1.0	THE DETERMINATION OF GAMMA ISOTOPES.....	3
2.0	METHOD OBJECTIVE, PURPOSE, CODE AND SUMMARY	3
3.0	METHOD APPLICABILITY	3
4.0	DEFINITIONS.....	3
5.0	METHOD VARIATIONS	4
6.0	SAFETY PRECAUTIONS AND WARNINGS	4
7.0	INTERFERENCES	5
8.0	APPARATUS, MATERIALS, REAGENTS, EQUIPMENT, AND INSTRUMENTATION	5
9.0	SAMPLE HANDLING AND PRESERVATION	5
10.0	SAMPLE PREPARATION.....	6
11.0	PREPARATION OF STANDARD SOLUTIONS AND QUALITY CONTROL STANDARDS	7
12.0	INSTRUMENT CALIBRATION AND PERFORMANCE	8
13.0	ANALYSIS AND INSTRUMENT OPERATION	8
14.0	EQUIPMENT AND INSTRUMENT MAINTENANCE	8
15.0	DATA RECORDING, CALCULATION, AND REDUCTION METHODS	8
16.0	QUALITY CONTROL REQUIREMENTS	10
17.0	DATA REVIEW, APPROVAL, AND TRANSMITTAL	11
18.0	RECORDS MANAGEMENT.....	12
19.0	LABORATORY WASTE HANDLING AND WASTE DISPOSAL	12
20.0	REFERENCES.....	12

1.0 THE DETERMINATION OF GAMMA ISOTOPES**2.0 METHOD OBJECTIVE, PURPOSE, CODE AND SUMMARY**

- 2.1 This standard operating procedure provides the necessary instructions to conduct the analysis for Gamma Isotopes in water, soil, urine and miscellaneous matrices.
- 2.2 Water samples are counted in Marinelli beakers. Soil samples are sealed in aluminum cans, which are counted immediately if Ra-226 is not desired. If Ra-226 is desired, the sealed can is set aside to allow secular equilibrium between Rn-222 and Bi-214. Quantification is done by the abundance of the 609 KeV Bi-214 line.
- 2.3 This method has been modified from the source method EPA 600/4-80-032 "Prescribed Procedures for Measurement of Radioactivity in Drinking Water," August 1980, Method 901.1, and the Department of Energy (DOE) EML Procedures Manual source method for Gamma PHA in soils and sediments, HASL-300. For all matrices, similar principles of radiochemical concentration and counting are used.
- 2.4 This method has been modified on the basis of GEL's Performance Based Measurement System (PBMS).

3.0 METHOD APPLICABILITY

- 3.1 Minimum Detectable Activity (MDA): The MDA is based upon sample volume, instrument background, instrument efficiency, count time and other statistical factors, as well as specific isotopic values such as abundance and half-life.
- 3.2 Method Precision: Typical relative percent difference (RPD) is less than 20%.
- 3.3 Method Bias (Accuracy): The method accuracy requirement for gamma, measured by running a laboratory control sample (LCS) with each batch, is 25% of the true value.
- 3.4 Analysts go through a partnered training program with an already certified analyst for gamma spectroscopy. The analyst receives training on reviewing of standard analytical requirement such as RPD, method bias and technical review of gamma spectra. The analyst can then become qualified to perform the analysis by passing an unknown sample analysis and correctly identifying the isotope(s). Technical training records are maintained electronically by the Quality Systems staff.

4.0 DEFINITIONS

- 4.1 Clean Line: An energy line of an isotope with no known energy lines of other isotopes within 2 KeV. (This excludes daughters that use the same line for quantification.)
- 4.2 Interfered Line: An energy line of an isotope with one or more energy lines of one or more different isotopes within 2 KeV.
- 4.3 Single and Double Escape Interference Lines: When high energy gamma lines above 1022 KeV have a large emission rate, it is possible to see single and double escape peaks caused by escape of the 511 KeV annihilation photon(s) from the detector.
- 4.4 Summation Interference: When high gamma emission rates are seen, sample summation can occur. Prominent in geometries close to detection and in low energy range (i.e., 10,000 cps at 88 KeV, 15,000 cps at 210 KeV), a summation interference can be seen at $88+88=176$ KeV, $210+210=420$ KeV, $210+88=298$ KeV.

- 4.5 False Positive: An isotope that has failed one or more of several tests including half-life, abundance, and energy tolerance (± 2 KeV)
- 4.6 Abundance Test: The test where the software verifies the presence of 75% of the total abundance of a nuclide in the system library is present. The presence of greater than 75% of the total abundance will cause a nuclide to be identified. The abundance criteria may be reduced to less than 75% for nuclides with several lower abundant photons.
- 4.7 Energy Tolerance: The test where the software checks the energy line in the spectrum to see if it is within the energy tolerance setting. (The standard setting is 2 KeV.) If it is within this setting then the line is associated with that nuclide. The energy line can be associated with more than one nuclide.
- 4.8 Half-Life Test: The test to determine if the half-life of the isotope is long enough not to have decayed away. The half-life of the sample is the time from sample date to analysis date plus 1/2 the count time. A limit of no more than eight half-life is the standard setting.
- 4.9 Key Line: The line chosen by the builder of the library to be the prominent line of the isotope. This line is used in the MDA table for purposes of calculating activity, error and MDA. For non-identified isotopes the key line energy is used as the basis of determining the region used to calculate the activity, error, and the MDA. Usually this line is the most abundant line on a line that is relatively free from interference
- 4.10 Abundance (Photon Intensity): The value, usually expressed in percent, given to a photon of specific energy which is emitted during the decay of a radionuclide. The abundance represents the probability of emission of a specific energy photon when a radionuclide is decaying (gamma/disintegration).
- 4.11 Counting Uncertainty: The error of the reported result due to the counting statistics of the instrument used for qualification.
- 4.12 Back Scatter: The detection of a count that occurs when an event interacts with counting materials, changes direction, and scatters back to the detector.

5.0 METHOD VARIATIONS

Modifications to the procedure are limited to GEL's use of additional isotopes for the daily calibration check and the inclusion of a more stringent calibration and resolution periodicity.

6.0 SAFETY PRECAUTIONS AND WARNINGS

- 6.1 Keep hands free from moving parts of canning device and Gamma shields.
- 6.2 Personnel performing this analytical procedure are trained in and follow the safe laboratory practices outlined in the Safety, Health and Chemical Hygiene Plan, GL-LB-N-001.
- 6.3 Personnel handling radioactive materials are trained in and follow the procedures outlined in GL-RAD-S-004 for Radioactive Material Handling.
- 6.4 Personnel handling biological materials are trained in and follow the procedures outlined in GL-RAD-S-010 for Handling Biological Materials.
- 6.5 If there is any question regarding the safety of any laboratory practice, **stop immediately**, and consult qualified senior personnel such as a Group or Team Leader.

7.0 INTERFERENCES

- 7.1 Some Gamma isotopes emit gamma lines that may overlap with other isotopes. If the energies of the two isotopes are within 2 KeV, the peaks may not be resolvable and will give a positive bias to the result. This problem is minimized by careful review of the peak search.
- 7.2 Soil samples may vary in density from the standard used for calibration. This may bias the results due to self-absorption of lower energy (<100 K).

8.0 APPARATUS, MATERIALS, REAGENTS, EQUIPMENT, AND INSTRUMENTATION**8.1 Ancillary Equipment**

- 8.1.1 100 cc aluminum cans with lids for soil and miscellaneous samples
- 8.1.2 Gelman Sciences PETRI dish for soil, filters and miscellaneous samples
- 8.1.3 2 L and 500 mL Marinelli beakers for water samples
- 8.1.4 Air displacement pipettes.
- 8.1.5 Can annealing tool
- 8.1.6 Graduated cylinder

8.2 Reagents, Chemicals and Standards

- 8.2.1 NIST traceable mixed gamma standard in 100cc aluminum can
- 8.2.2 NIST traceable 2.0 liter mixed gamma standard in 2 L Marinelli beaker
- 8.2.3 NIST traceable mixed gamma standard in 0.5 L Marinelli
- 8.2.4 NIST traceable mixed gamma standard in Gelman Sciences PETRI dish
- 8.2.5 Standard soil blank
- 8.2.6 NIST traceable mixed gamma standard of 13-47mm glass fiber filter composites in Gelman Sciences PETRI dish.
- 8.2.7 NIST traceable aqueous Mixed Gamma Standard: Contains Am-241, Co-60, and Cs-137 as a minimum.
- 8.2.8 NIST traceable mixed gamma standard of 1-47mm glass fiber
- 8.2.9 NIST traceable mixed gamma standard frontloaded in BG-300 Impregnated Charcoal Sample Cartridge.
- 8.2.10 Nitric Acid, reagent grade. (16M)
- 8.2.11 Hydrofluoric acid, 48%.
- 8.2.12 Hydrochloric acid, reagent grade. (12M)
- 8.2.13 Boric acid, 5%. Dissolve 50 grams of H_3BO_3 per liter of water

8.3 Instrumentation

- 8.3.1 High purity germanium detector, with associated electronics and data reduction software
- 8.3.2 Top loader balance

9.0 SAMPLE HANDLING AND PRESERVATION

- 9.1 For soil samples, 500g of sample should be collected, preferably in a plastic container to avoid breakage.
- 9.2 For water samples, 2 liters of sample should be collected in a plastic container and preserved to pH2 with Nitric acid.

10.0 SAMPLE PREPARATION**10.1 Soil sample preparation.**

10.1.1 Prepare the sample for gamma counting in accordance with SOP GL-RAD-A-021 "Soil sample preparation for the determination of radionuclides".

10.1.2 Fill the appropriate container with sample prepared from step 10.1.1 using the following steps as a guideline:

10.1.2.1 If Ra-226 analysis is required, the sample is placed in a 100cc can for in-growth.

NOTE: It is recommended that in-growth be allowed 7 days to quantify Ra-226. Longer intervals can be used at the request of the client. However, shorter in-growth periods may decrease the accuracy of the data. If there is insufficient mass of sample to fill the 100cc can, contact the team or group leader.

10.1.2.2 All homogenized samples shall be placed in the 100cc can. Determine the net weight of the sample. If the net weight is less than 55 grams or greater than 190 grams, contact the team or group leader to determine the appropriate counting container. Record sample weight and date on sample container.

10.1.2.3 If there is insufficient sample to fill the 100cc can, place sample in the 10cc petri dish, cap and seal. Record sample weight and date on sample container.

10.1.2.4 If there is insufficient sample to fill the 10cc petri dish, perform the following digestion process:

10.1.2.4.1 Weigh out an appropriate aliquot into a labeled teflon beaker. Record this weight on the sample container.

10.1.2.4.2 Add 10 mL of concentrated nitric acid to each

10.1.2.4.3 ~~sample~~ Place samples on medium heat (~300 °F) and cover each sample with a teflon lid. Reflux all samples for 30 minutes.

10.1.2.4.4 Remove teflon lids and add 5 mL concentrated hydrochloric acid and 10 mL hydrofluoric acid to each sample. Cover samples and reflux for 120 minutes.

10.1.2.4.5 Remove teflon lids and allow samples to evaporate to dryness.

10.1.2.4.6 Add 5 mL of concentrated nitric acid and evaporate to dryness.

10.1.2.4.7 Repeat Step 10.1.2.4.6.

10.1.2.4.8 Add 5 mL of concentrated nitric acid to the dry samples. Add 1 ml of 5% boric acid. Place the samples back on the hotplate long enough so that the dried sample dissolves into solution.

10.1.2.4.9 Transfer solution to a 500 mL marinelli beaker and dilute to 500 mL. Record the original sample mass and diluted volume on sample container. Record the original sample mass on batch que sheet.

10.2 Water sample preparation

10.2.1 Mix and measure an appropriate volume into a 2 L or 500 mL Marinelli beaker and record the volume on the batch Que Sheet. If applicable, record the standard identification code, volume and expiration date on the batch Que sheet.

10.3 Urine Sample Preparation

10.3.1 Place a 24-hour urine container (or other suitable container) on a balance and tare the balance.

10.3.2 Transfer the entire volume of the sample received to the tared container and record the volume of sample received.

10.3.3 Add 8 M HNO₃ acid to the original sample container (typically 25 – 50 mL). Shake in the container and then heat in a microwave for approximately 30 seconds to remove sample residue from the sides of the sample container.

10.3.4 Add the nitric acid rinse to the 24-hour urine container and record the volume of the original sample plus acid.

10.3.5 Cap and shake the 24-hour urine container to homogenize the sample. Transfer an aliquot (typically 500 mL) of this solution to a Marinelli Beaker.

10.3.6 Record the amount of the original sample, excluding the nitric acid added, on the gamma spec que sheet.

Example: 800 mL is received and 50 mL of 8 M HNO₃ is added from the rinse of the sample container. 500 mL is transferred to the Marinelli Beaker. The recorded volume on the que sheet should be (500 mL/850 mL) x 800 mL = 470.6 mL.

10.4 Preparation of miscellaneous matrices

10.4.1 Prepare the sample in accordance with SOP GL-RAD-A-026 “Preparation of Special Matrices for the Determination of Radionuclides.”

10.4.2 Once the appropriate section of GL-RAD-A-026 has been performed, prepare the sample for gamma counting by referring to section 10.1.2 above.

11.0 PREPARATION OF STANDARD SOLUTIONS AND QUALITY CONTROL STANDARDS

11.1 A blank is performed with each batch. DI Water should be used to prepare the blank.

11.2 A duplicate should be run with each sample batch. Refer to the batch pull sheet to determine the designated batch duplicate sample.

11.3 A matrix spike sample is prepared by adding a known volume of standard directly to the designated sample. Refer to the batch pull sheet to determine the designated batch matrix spike sample.

11.4 A laboratory control sample is prepared by adding a known volume of standard directly to a Marinelli beaker with DI water.

12.0 INSTRUMENT CALIBRATION AND PERFORMANCE

12.1 Refer to "Gamma Spectroscopy System Operating Procedure" (GL-RAD-I-001) for calibration periodicity and instructions.

12.2 Refer to "Counting Room Instrument Maintenance and Performance Checks" (GL-RAD-I-010) for instructions concerning instrument maintenance.

13.0 ANALYSIS AND INSTRUMENT OPERATION

13.1 Place the sample on the detector and count the sample an appropriate amount of time in the gamma shield. See "Gamma Spectroscopy System Operating Procedure" (GL-RAD-I-001) for specific instructions on operating the gamma spectrometers.

14.0 EQUIPMENT AND INSTRUMENT MAINTENANCE

14.1 Refer to "Gamma Spectroscopy System Operating Procedure" (GL-RAD-I-001) for instructions concerning the Gamma Spectrometer.

14.2 Refer to "Counting Room Instrument Maintenance and Performance Checks" (GL-RAD-I-001) for instructions concerning instrument maintenance.

15.0 DATA RECORDING, CALCULATION, AND REDUCTION METHODS

15.1 Data Recording

Record the following information on the Gamma Que Sheet: preparation date, analyst's initials, spike isotope, spike code, spike volume, LCS isotope, LCS code, LCS volume, nominal concentration LCS, and nominal concentration MS. For each sample record the detector number, sample mass, sample date and time.

15.2 The instrument will report sample pCi/unit according to the following equation:

$$\text{Sample pCi/unit} = \frac{A * d}{2.22 * E * V * B * Ct * ABS}$$

Where:

A = net peak area (counts)

ABS = relative absorption factor

B = abundance (gammas/disintegration)

E = counting Efficiency (counts/gamma)

V = sample volume (grams or liters)

Ct = sample count time (minutes)

d = decay factor = $\frac{1}{e^{-\lambda t}}$

15.3 Counting uncertainty is calculated according to the following equation:

$$\text{pCi/unit} = Ac * 1.96 \sqrt{\left(\frac{ef - er}{E}\right)^2 + \left(\frac{pk - er}{pk}\right)^2 + \left(\frac{ab - er}{A}\right)^2 + \left(\frac{sy}{100}\right)^2 + (\text{Decay})}$$

Where:

Ac = Activity from 15.2

$$\text{Decay} = \left(\frac{T_{1/2\text{err}}}{T_{1/2}} \right)^2 * \left[\frac{\lambda_{\text{Er}}}{1 - e^{-\lambda_{\text{Er}}}} - \lambda(T_s + \text{Er}) - 1 \right]$$

- 15.4 The method MDA in pCi/g or pCi/L are calculated according to the following equations:

$$\text{MDA (pCi/unit)} = \frac{d * \left(2.71 + 4.66 \sqrt{\text{cpm}_b * \text{ct}} \right)}{2.22 * E * V * B * \text{ct}}$$

Where:

A = net peak area (counts)

ABS = relative absorption factor

B = abundance (gammas/disintegration)

E = counting Efficiency (counts/gamma)

V = sample volume (grams or liters)

ct = sample count time (minutes)

d = decay factor = $d = \frac{1}{e^{-\lambda t}}$

- 15.5 The absorption factor is calculated by the following equations:

$$I_1 = \frac{\ln((\text{SScpm} - \text{Scpm})/\text{ECcpm})}{(((\text{SScpm} - \text{Scpm})/\text{ECcpm}) - 1)}$$

$$I_0 = \frac{\ln((\text{SSTcpm} - \text{STcpm})/\text{ECcpm})}{(((\text{SSTcpm} - \text{Scpm})/\text{ECcpm}) - 1)}$$

$$\text{ABS} = \frac{I_1}{I_0}$$

Where:

SScpm = sample plus the source cpm at the region of interest

Scpm = sample cpm at the region of interest

ECcpm = source cpm on the empty can at the region of interest

ln = natural logarithm

SSTcpm = standard plus the source cpm at the region of interest

STcpm = standard cpm at the region of interest

- 15.6 The VAX operating system will report the following information with each completed sample:

15.6.1 The nuclide identification report

15.6.2 The minimum detectable activity report

15.6.3 The peak search report.

- 15.7 The following criteria are used to accept a reported gamma isotope from the NID report:

- 15.7.1 The peak FWHM should be less than 3 KeV.
- 15.7.2 The activity of a non-target isotope will not be reported unless the result is greater than the minimum detectable activity and the result is greater than the three sigma uncertainty..
- 15.7.3 The energy tolerance should be between 2 and 3 KeV.
- 15.7.4 The sensitivity setting should be between 0.1 and 3. The default setting is 3.
- 15.7.5 Start channel on peak search should be approximately 50 and end channel should be 4096.
- 15.7.6 The confidence level setting should be 5.
- 15.7.7 These settings should not be changed without approval from a group leader.
- 15.8 The following guidelines are used to accept unidentified lines on the peak search after environmental background subtraction:
 - 15.8.1 The line matches the natural fingerprint of the Uranium-238 or Thorium-232 decay chains (i.e. 63, 75, 93, 239, 295, 352, 511, 609, 1120, etc.).
 - 15.8.2 The line matches as a summation peak from two other lines in the spectrum.
 - 15.8.3 The line has a net area of less than 20.
 - 15.8.4 The line matches as a escape peak from an identified nuclide which emits photons greater than 1022 KeV.

16.0 QUALITY CONTROL REQUIREMENTS

16.1 Analyst and Method Verification

Refer to "Analyst and Analytical Methods Validation Procedures" (G-RAD-D-003) for instructions concerning the validation of analysts and analytical methods.

16.2 Method Specific Quality Control Requirements

- 16.2.1 A method blank will accompany each batch of 20 or less samples. The reported value should be less than or equal to the CRDL for all target isotopes. Matrix spikes are prepared by spiking a portion of the QC sample with Cs-137 (as a minimum).
- 16.2.2 For water samples only, a matrix spike (MS) should be run with every batch of 20 samples. The recovery of the spike should fall between 75 and 125%. The recovery is calculated as follows:

$$\%REC = \frac{\text{spike(pCi/g)} - \text{sample(pCi/g)}}{\text{spikedamount(pCi/g)}} * 100$$

or:

$$\%REC = \frac{\text{spike(pCi/L)} - \text{sample(pCi/L)}}{\text{spikedamount(pCi/L)}} * 100$$

NOTE: Performing a matrix spike on a soil sample would result in direct contamination of the sample, therefore, only water samples require an MS.

- 16.2.3 A sample duplicate should be run with every batch of 20 or less samples. The relative percent difference (RPD) between the sample and the duplicate should be $\leq 20\%$. The RPD is calculated as follows.

$$RPD = \frac{\text{high sample (pCi/g)} - \text{low sample (pCi/g)}}{\text{Average (pCi/g)}}$$

or:

$$RPD = \frac{\text{high sample (pCi/L)} - \text{low sample (pCi/L)}}{\text{Average (pCi/L)}}$$

- 16.2.4 A laboratory control spike (LCS) should be run with every batch of 20 samples or less. The recovery of the spike should fall between 75 and 125%. The LCS should contain Cs-137 as a minimum. Some clients may request a mixed gamma standard. For soils, a mixed gamma expired calibration source may be used as an LCS. For liquids and filters, spike a blank sample with Cs-137 as a minimum.

- 16.2.5 The recovery is calculated as follows:

$$LCS = \frac{\text{observed_pCi/g}}{\text{known_pCi/g}} * 100$$

or:

$$LCS = \frac{\text{observed_pCi/L}}{\text{known_pCi/L}} * 100$$

16.3 Actions required if the Quality Control Requirements Are Not Met

If any of the above criteria cannot be satisfied, the analyst should inform the group leader and initiate a non-conformance report as outlined in "Documentation of Nonconformance Reporting and Dispositioning, and Control of Nonconforming Items" (GL-QS-E-004).

17.0 DATA REVIEW, APPROVAL, AND TRANSMITTAL

- 17.1 The data is transmitted from the laboratory personnel to the reporting personnel as outlined in "Data Review and Validation Procedures" (GL-RAD-D-003):

- 17.1.1 Visually check the que sheet, spreadsheet, raw data and data report to make sure the information has been transcribed correctly.

- 17.1.2 Review the raw data to see if there are any hits not on the requested list. If there are, report to the client by adding the information into LIMS.

A true identification or a "hit" is any isotope greater than 10 pCi/L or 5 pCi/g on the identified nuclide list. The error must also be less than 40% of the result and not have interference by another isotope or have a very short half-life.

- 17.1.3 Check to see that the required detection limit (RDL) is met if required.

- 17.1.4 Check hits to see if they are true hits and not an interference or a false positive.

Identifications are classified into two categories: false positives (interference), and true identification (hit). The false positives are rejected by checking the abundance test results for the isotope and by checking last results for the half-life. The result is considered interference and rejected by checking to see if there are any clean lines in sample spectrum for the isotope. If none exist, then the identification is rejected. If the key line has a possible interference and secondary lines do not confirm the activity calculation, the identification is rejected. Isotopes that pass these criteria are accepted as true identifications. The above tests and criteria are standard and will be followed unless directed otherwise by contract, specification or instructions.

17.1.5 Complete the batch checklist.

18.0 RECORDS MANAGEMENT

- 18.1 Each analysis that is performed on the instrument is documented in the run log according to "Run Logs" (GL-LB-E-009).
- 18.2 All raw data printouts, calculation spreadsheets and batch checklists are filed with the sample data for archival and review.

19.0 LABORATORY WASTE HANDLING AND WASTE DISPOSAL

- 19.1 After analysis, return sample containers to storage as outlined in "Verifying the maintenance of sample integrity" (GL-LB-E-012).
- 19.2 Radioactive waste is disposed of as outlined in the Laboratory Waste Management Plan (GL-LB-G-001).

20.0 REFERENCES

- 20.1 USEPA. Prescribed Procedures for Measurement of Radioactivity in Drinking Water. Method 901.1, August 1980.
- 20.2 Canberra Nuclear Genie System Spectroscopy, Applications and Display User's Guide. Vol. I and II, May 1991.
- 20.3 EML procedures manual. HASL-300-Ed.25, 1982.

VERIFY THE VALIDITY OF THIS SOP EACH DAY IN USE

STANDARD OPERATING PROCEDURE
FOR
THE ISOTOPIC DETERMINATION OF PLUTONIUM,
URANIUM, AMERICIUM, CURIUM, AND THORIUM

Applicable to:
EML HASL-300 E-U-04 (Modified)
EML HASL-300

(GL-RAD-A-045-REVISION 1)

PROPRIETARY INFORMATION

This document contains proprietary information that is the exclusive property of General Engineering Laboratories, LLC (GEL). No contents of this document may be reproduced or otherwise used for the benefit of others except by express written permission of GEL.



VERIFY THE VALIDITY OF THIS SOP EACH DAY IN USE

STANDARD OPERATING PROCEDURE
FOR
SAMPLE RECEIPT, LOGIN AND STORAGE
(GL-SR-E-001 REVISION 17)

PROPRIETARY INFORMATION

This document contains proprietary information that is the exclusive property of General Engineering Laboratories, LLC (GEL). No contents of this document may be reproduced or otherwise used for the benefit of others except by express written permission of GEL.



TABLE OF CONTENTS

1.0	STANDARD OPERATING PROCEDURE FOR SAMPLE RECEIPT, LOGIN AND STORAGE	3
2.0	PURPOSE	3
3.0	DISCUSSION	3
4.0	DEFINITIONS	3
5.0	SAFETY, HEALTH AND ENVIRONMENTAL HAZARDS	4
6.0	PROCEDURES	4
7.0	RECORDS MANAGEMENT	9
8.0	REFERENCES	10
APPENDIX 1: STORAGE AND PRESERVATION		11
APPENDIX 2: SAMPLE RECEIPT REVIEW SHEET		14
APPENDIX 3		15

1.0 STANDARD OPERATING PROCEDURE FOR SAMPLE RECEIPT, LOGIN AND STORAGE

2.0 PURPOSE

To describe the routine operational procedures for the receipt, login and storage of samples received by General Engineering Laboratories, LLC (GEL).

3.0 DISCUSSION

3.1 Sample custody is a pre-planned mechanism for tracking a sample from the collection of the sample in the field through the release of the finished analytical data to the client. At the collection site, the sample containers are filled with sample and the Chain of Custody form is initiated. The sample collector fills out the form, which includes the name of the client, the requested analysis parameters, sample location, the date and time of collection, sampling technique, preservatives used, and any comments or remarks that may be useful in the analytical work or data interpretation that will follow. Proper sample receipt, login and storage assure accurate chain of custody.

3.2 Custody is defined as:

- 3.2.1 Being in your physical possession, or
- 3.2.2 Being in your view, after being in your possession, or
- 3.2.3 Being locked up after being in your possession, or
- 3.2.4 Being in a designated secure area

3.3 Upon arrival at the laboratory, sampling personnel, delivery service and carriers relinquish the samples to the sample management group. Each sample container receives a unique sample identifier that is assigned electronically by LIMS (Laboratory Information Management System). LIMS tracks the status and location of each sample container, and serves as the database for analytical results.

4.0 DEFINITIONS

4.1 ALPHALIMS: Laboratory Information Management System.

4.2 Chain of Custody (COC): A written record of sample transfer and possession.

4.3 Custody Seal: Security seals that are attached to sample containers and/or bottles that are used to detect unauthorized tampering.

4.4 Holding Time: The period of time between sample collection and preparation or analysis.

4.5 Labeled Package: A package containing radioactive material labeled with a Radioactive White-I, Radioactive Yellow-II or Radioactive Yellow-III label as specified in US Department of Transportation Regulations, 49 CFR 172.403 and 172.436-440.

4.6 Matrix: The physical appearance or make-up of a sample (groundwater, drinking water, wastewater, soil, sludge, etc.) as determined by the client or Project Manager.

4.7 Material Safety Data Sheet (MSDS): A document that may accompany samples of known chemical characteristics. (See our "Safety, Health and Chemical Hygiene Plan" for more information on MSDSs.)

4.8 Preservative: Additives that are introduced to a sample at the time of collection to help retard chemical and biological changes that may occur.

- 4.9 Turn Around Time (TAT): A numeric designation to the degree of attention a sample should receive. This designation is used to convey the client's requested data delivery dates to the laboratory.
- 4.10 Sample Delivery Group (SDG): One or more samples (typically not to exceed 20 samples) from a specific client that are reported by the laboratory at the same time.
- 4.11 Sample Receipt Review (SRR): A form used to document a sample's arrival and the condition of its arrival at the laboratory.
- 4.12 Sample: Any item that has been submitted for analysis to GEL.

5.0 SAFETY, HEALTH AND ENVIRONMENTAL HAZARDS

- 5.1 All samples must be handled with care during the login process. Wear protective gear such as gloves, aprons, safety glasses and laboratory coats when handling all samples. Some samples may be accompanied by MSDSs that contain vital information on potential hazards. The sample description and client labels may also give this information.

NOTE: Gloves and protective eyewear must be worn when handling samples. Lab coats should be worn when handling any samples but are only required to be worn when handling radioactive or hazardous samples. (Refer to the "Safety, Health and Chemical Hygiene Plan.")

- 5.2 If there is a spill of a known hazard (based on historical results, MSDSs, and/or sample description), immediately contact the Group Leader, Laboratory Waste Manager, or Radiation Safety Officer as appropriate.
- 5.3 All sample management personnel are required to read and understand GEL's "Safety, Health and Chemical Hygiene Plan," which is found on GEL's Intranet.

6.0 PROCEDURES

6.1 Sample Package Receipt

- 6.1.1 All sample packages submitted to GEL are received by sample management personnel. Samples are received from a number of carriers including GEL field staff, GEL couriers, individual clients, and public and private shipping companies.
- 6.1.2 Upon arrival, all sample packages will be inspected for integrity. Note any unusual physical damage, signs of leakage, or evidence that custody seals have been tampered with. If the package appears to be leaking or has any unusual odor, place it under the fume hood and notify the Group Leader, Laboratory Waste Manager, Radiation Safety Officer, or Project Manager as appropriate before continuing.
- 6.1.3 All sample packages will also be screened for external contact radiation exposure. This screening is performed to determine the possible presence of radionuclides that may require special handling. If a radioactive "labeled" package is received, or any package exceeds 0.5 mrem/hr on contact, the RSO group should be notified, and the package is segregated in the GEL sample receiving area where the RSO or designee will unpack the package following the procedures described in GL-RAD-S-007 for "Receiving of Radioactive Samples."

- 6.1.4 Bioassay and Low Level Mercury (LLHG) sample packages are initially received and segregated in the GEL login area. Following package screening, they are then transported to the bioassay or LLHG login area for inspection, login and storage. Bioassay and LLHG receiving staff perform the same sample inspection, login and storage procedures using the Sample Receipt Review (SRR) in Appendix 2 or 3 as applicable, except as noted in Section 6.1.7.
- 6.1.5 Packages indicating <0.5 mrem/hr should be further segregated to identify samples intended to be received under the authority of GEL's Radioactive Material License. In addition to "Labeled" radioactive material, radioactive material is any material that meets the following criteria:
 - 6.1.5.1 Any Material received that was shipped as DOT Hazard Class 7 Limited Quantity- Excepted Package.
 - 6.1.5.2 Any material shipped and received that is marked as radioactive (i.e. conventional trefoil, yellow and magenta tape, etc.), or has otherwise been declared radioactive by the consignor in the accompanying documents. This material may be intended for receipt under the authority of a radioactive materials license although it was shipped under DOT exemption for radioactive material.
- 6.1.6 All discrepancies noted during receipt and inspection shall be recorded using a Sample Receipt Review (SRR) form in Appendix 2.
 - 6.1.6.1 As required client specific Sample Receipt Review forms may be created by the Project Management Group. These checklists are created because additional sample management comments and checks are required in order to meet quality objectives established for these project samples.
- 6.1.7 Open all shipping containers (excluding bioassay & LLHG samples) under the high volume exhaust duct located in Sample Receiving. (NOTE: It is only necessary to open Bioassay & LHG samples under a fume hood when the integrity of the containers is suspected/determined to be compromised.)
 - 6.1.7.1 All samples received (excluding bioassay) must be screened for radioactivity using a Geiger-Muller pancake probe. Results for the highest reading samples are to be noted on the SRR form. The Radiation Safety Group shall be notified when readings for any individual non-radioactive sample exceed 2x area background.
- 6.1.8 The COC should accompany all samples received by the Sample Management Group. The COC documentation includes sample identification (e.g., MW-1; Lagoon 17; #1234567), sampling date and time, sample collector, and requested parameters to be tested. If this documentation is not present, the Sample Management Group upon receipt shall initiate the COC. Identify this initiation by printing "INITIATED ON RECEIPT" on the COC form.

- 6.1.9 Compare the sample labels to the Chain of Custody; compare sample descriptions, collection dates, collection times, number of containers and any other available information. Note any discrepancies of the COC and the Sample Receipt Review form (SRR), and inform the Project Manager. Sign and date (including time) the COC in the appropriate box.
- 6.1.10 Analytical procedure may require preservation of the sample to ensure that changes in the samples chemistry or biology do not occur. The two predominant preservation techniques used are changing the pH of the sample and cooling the sample to about 4°C. It is important to check and document the holding time, preservation and temperature of the samples upon arrival to the laboratory. The correct methods of sample storage, chemical preservation, and maximum holding times are shown in Appendix 1. Those samples determined to be non-conforming shall be documented and the Project Manager notified.

Verify and document pH preservation using the following procedure:

- 6.1.10.1 Open the container.
- 6.1.10.2 Pour an aliquot of the original sample into a secondary container. Immerse a pH strip into the secondary container to take the measurement.
- 6.1.10.3 Observe the pH as indicated on the pH strip, and discard pH strip and secondary container.
- NOTE:** Never reuse a pH strip or one that has been contaminated.
- 6.1.10.4 Document results of the preservation verification on the appropriate line of the Sample Receipt Review form. (See Appendix 2 for example of a Sample Receipt Review form.)

NOTE: If the pH of the sample is determined to be non-conforming, place the sample on hold and notify the Project Manager. The Project Manager will call the client for further direction. If direction is given to adjust the preservation, continue processing the sample, and preserve the sample with the appropriate preservative (Appendix 1) recording the lot # of preservative used on the SRR. After adding the appropriate preservative to the sample, wait 2 minutes and perform steps 6.1.10.1 - 6.1.10.4 again. The preserved sample should now be placed on the preservation adjustment hold shelf located in the main cooler. This ensures that the 16-hour holding time for metals samples and the 24-hour holding time for radiochemistry samples is met following preservation or adjustment. Document this on the SRR: "SAMPLE PRESERVED UPON ARRIVAL."

- 6.1.10.6 Following is a list of tests that require pH verification upon arrival:

TEST	pH
Ammonia	<2
COD	<2
Cyanide	>12

Hardness	<2
Hydrazine	<2
Metals	<2
Nitrate/Nitrite	<2
Phenols	<2
Phosphates, Total	<2
Radiochemistry (all except Tritium, C-14, Rn-222, I-129, I-131)	<2
Sulfide	>9
TKN	<2
TOC/TIC/DOC	<2

NOTE: The pH of all aqueous sample fractions, preserved and unpreserved, shall be checked during sample login for the following DOE Albuquerque (DOE-AL) installations: (Exceptions to the pH check are Rn-222, tritium, iodine, VOC, TOX, oil and grease, and urine samples.)

- Los Alamos National Laboratories
- Mound Plant
- Pantex Plant
- Sandia National Laboratories, Albuquerque
- Sandia National Laboratories, Livermore

6.1.10.7 Sample receipt temperature is verified and documented upon using the following procedure:

- 6.1.10.7.1 Open the sample cooler.
- 6.1.10.7.2 Remove the Temperature Validation Container (TVC) if provided.
- 6.1.10.7.3 Open the TVC and immerse a thermometer with a valid calibration into the TVC.
- 6.1.10.7.4 Allow the thermometer reading to equilibrate, and read the thermometer result while it is still immersed in the TVC.
- 6.1.10.7.5 Alternately the receipt temperature can be measured with an infrared temperature (IR) gun with a valid calibration by selecting the TVC or another sample within the shipment.
- 6.1.10.7.6 Record the observed reading on the Sample Receipt Review form (Appendix 2), as well as on the COC if a space is provided: i.e. "TEMP. 4° UPON ARRIVAL."

- 6.1.10.7.7 Temperature verification results of $4^{\circ} \pm 2^{\circ}\text{C}$, are considered conforming for those samples listed in Appendix 1 for 4° .
- 6.1.10.7.8 If the initial temperature verification results are determined to be non-conforming, select another sample container from the shipment for temperature verification and re-perform steps.
- 6.1.10.7.9 Record the verification temperature on the SRR as well as on the COC if a space is provided. Label the temperature as a verification temperature (i.e., $\text{VT}=7^{\circ}\text{C}$).
- 6.1.10.7.10 If another container is not available within the shipment to verify the temperature, the secondary temperature verification is not performed and duly noted.
- 6.1.10.8 Project Managers may specify via client specific SRRs that aqueous organic analysis sample containers (excluding volatile 40-mL vials) be checked for the presence of chlorine residual at the time of sample receipt. If chlorine residual is present, document as such on the SRR and inform the Project Manager. The Project Manager may specify that samples with chlorine residual require the addition of $\text{Na}_2\text{S}_2\text{O}_3$.
- 6.1.11 Deliver the completed COC and SRR to the appropriate Project manager.
- 6.1.12 The Project Manager or Project Manager Assistant “logs” the data from the COC and SRR into ALPHALIMS. Once samples are logged into the system, unique bar code labels are generated for each sample container.
- 6.1.13 The bar code labels are ready to be affixed to the appropriate containers.
- 6.1.14 Sample bar code labels are color-coded as follows:
 - 6.1.14.1 Yellow and magenta for radioactive samples.
 - 6.1.14.2 Solid white for Federal Division non-radioactive samples
 - 6.1.14.3 White and green for Industrial Division samples.
- 6.1.15 Compare the sample description on the printed GEL bar code label to the client sample bottle label before attaching labels to containers. You should not cover the client’s label or any other information provided by the client or sample collector.
- 6.1.16 If the sample is a solid submitted for volatiles analysis and a single container is provided, a designation is generated on the barcode label indicating, “Volatiles must aliquot sample first.” It is then stored in the appropriate volatile cooler until removed for volatiles testing. Once the volatiles lab takes its required aliquot the container will be marked with the analyst’s initials and the date completed. The sample container will

then be placed in the appropriate walk in cooler and released for other laboratory analyses. Note exception in Section 6.2.2.1.

6.2 Sample Storage and Security

6.2.1 Once the samples have been properly labeled, the samples are placed in the appropriate storage areas. The storage areas are located within the laboratory area of the building. Access to the laboratory is limited to those with security clearance identification badges. Entrance into the laboratory is electronically monitored. All visitors to the building must sign in at the reception area where they will receive "Visitor" identification badges. Visitors must be escorted while they are in the laboratory.

The samples are scanned into the electronic tracking system. Containers are loaded into the system by container type and size (i.e., 1000 mL – nalgene), preservative (i.e., H_2SO_4), and storage area of destination.

6.2.2 Samples are placed in numerical order in the appropriate storage locations throughout the facility.

6.2.2.1 Samples requiring analysis of volatile organics shall be segregated from other samples by placing them in either the radioactive or non-radioactive coolers, which are located in the Volatiles area and maintained at $4^\circ \pm 2^\circ\text{C}$.

NOTE: Samples requiring volatile analyses known to contain high concentrations of organic solvents or hydrocarbons should not be stored in the volatiles coolers. Place these samples in either the general use walk-in cooler.

6.2.2.2 Samples requiring radiochemical analyses only (except radon) are stored, in numerical order, in ambient storage. Radioactive and Non Radioactive sample are segregated in these storage areas.

6.2.2.3 Samples required cold preservation (other than volatile organics samples) are stored, in numerical order, in general use walk-in coolers, which are maintained at $4^\circ \pm 2^\circ\text{C}$. Radioactive and non-radioactive samples are segregated in these storage areas.

6.2.3 The Sample Management Group monitors cooler temperature twice daily, every working day. Calibrated thermometers are located in each walk-in cooler and readings are taken once in the morning and once in the afternoon, no less than two hours apart. Contact the Group Leader if temperatures fall outside of acceptance ranges. Document all non-conformances and corrective actions in the temperature logs.

7.0 RECORDS MANAGEMENT

7.1 The Sample Receipt Review form is attached to the Chain of Custody and forwarded to the Project Manager.

7.2 Cooler temperature logs are reviewed each month by Quality Systems. At the end of the year, completed logs are forwarded to Quality Systems for archiving.

8.0 REFERENCES

- 8.1 Example Standard Operating Procedures for Contract Laboratory Program (CLP), National Enforcement Investigations Center (NEIC), Contract Evidence Audit Team (CEAT-TechLaw), EPA Contract 68-01-6838, 1986.

APPENDIX 1: STORAGE AND PRESERVATION**SAMPLE STORAGE AND PRESERVATION REQUIREMENTS**

Parameter	Container¹	Preservation	Holding Time²
Inorganics			
Acidity	P,G	4 ⁰ C	14 days
Alkalinity	P,G	4 ⁰ C	14 days
Demand (BOD)	P,G	4 ⁰ C	48 hours
Bromide	P,G	None	28 days
Chemical Oxygen Demand (COD)	P,G	4 ⁰ C, H ₂ SO ₄ to pH<2	28 days
Chlorine by Bomb	P,G	None	None
Chloride	P,G	None	28 days
Color	P,G	4 ⁰ C	48 hours
Conductivity	P,G	4 ⁰ C	28 days
Corrosivity by pH	P	None	Immediate
Corrosivity to Steel	P	None	None
Cyanide amenable to chlorination	P,G	4 ⁰ C, NaOH to pH>12, 0.6g ascorbic acid ³	14 days ⁴
Cyanide, total	P,G	4 ⁰ C, NaOH to pH>12, 0.6g ascorbic acid ³	14 days ⁴
Dissolved Oxygen	G (bottle and tap)	None	Immediate
Fixed and Volatile Solids	P,G	4 ⁰ C	7 days
Flashpoint	P,G	None	None
Fluoride	P	None	28 days
Hardness	P,G	HNO ₃ to pH<2, H ₂ SO ₄ to pH<2	6 months
Heating Value	P	None	None
Hydrazine	G	HC1 to pH<2	Immediate
Percent (%) Moisture	P	4 ⁰ C	None
Ammonia Nitrogen	P,G	4 ⁰ C, H ₂ SO ₄ to pH<2	28 days
Nitrate	P,G	4 ⁰ C	48 hours
Nitrite	P,G	4 ⁰ C	48 hours
Nitrate/Nitrite	P,G	4 ⁰ C, H ₂ SO ₄ to pH<2	28 days
Total Kjeldahl and Organic Nitrogen	P,G	4 ⁰ C, H ₂ SO ₄ to pH<2	28 days
Odor	G	4⁰C, Zero headspace	Immediate
Oil and Grease	G	4 ⁰ C, HC1 or H ₂ SO ₄ to pH<2	28 days
Orthophosphate	P,G	Filter immediately, 4⁰C	48 hours
Total Phenols	G	4 ⁰ C, H ₂ SO ₄ to pH<2	28 days
pH	P,G	None	Immediate
Total Phosphorus	P,G	4 ⁰ C, H ₂ SO ₄ to pH<2	28 days
Residual Chlorine	P,G	None	Immediate
Salinity	P	None	28 days
Specific Gravity	P	4 ⁰ C	7 days
Sulfate	P,G	4 ⁰ C	28 days
Sulfide	P,G	4 ⁰ C, add ZNAce and NaOH to pH>9	7 days
Sulfite	P,G	None	Immediate
Sulfur by Bomb	G	None	None
Surfactants	P,G	4 ⁰ C	48 hours
Settleable Solid	P,G	4 ⁰ C	48 hours
Total Dissolved Solid	P,G	4 ⁰ C	7 days
Total Solid	P,G	4 ⁰ C	7 days
Total Suspended Solid	P,G	4 ⁰ C	7 days
Volatile Solid	P,G	4 ⁰ C	7 days
Total Organic Carbon	P,G	4 ⁰ C, HC1 or H ₂ SO ₄ to pH<2	28 days
Total Organic Halides	G	4 ⁰ C, H ₂ SO ₄ to pH<2	28 days
Total Petroleum Hydrocarbons	G	4 ⁰ C, H ₂ SO ₄ to pH<2	28 days
Turbidity	P,G	4 ⁰ C	48 hours
Metals (except chromium VI and mercury)	P	4 ⁰ C, HNO ₃ to pH<2	6 months
Chromium VI - Aqueous	P	4 ⁰ C	24 hours

Sample Receipt, Login and Storage			
SOP Effective 11/01/92 Revision 17 Effective June 2004		GL-SR-E-001 Rev 17 Page 12 of 15	
Chromium VI - Solids	P	4°C	7 days for extraction
Mercury - Wastewater and Drinking water	P,G	4°C,HNO ₃ to pH<2	28 days
Mercury - Others	G	4°C,HNO ₃ to pH<2	28 days
<u>Bacteriology</u>			
Coliform, fecal	P,G	4, 0.008% Na ₂ S ₂ O ₃ ³	6 hours
Standard Plate Count	P,G	4°C, 0.008% Na ₂ S ₂ O ₃	24 hours
Coliform, total - Wastewater	P,G	4°C, 0.008% Na ₂ S ₂ O ₃	6 hours
Coliform, total - Groundwater	P,G	4°C, 0.008% Na ₂ S ₂ O ₃	24 hours
<u>Organics</u>			
Base/Neutral and Acid Extractables - Water	Amber G, teflon-lined cap	4°C 0.008% sodium thiosulfate solution	7 days for extraction 40 days after extraction for analysis
Base/Neutral and Acid Extractables - Solid and Waste	G, teflon-lined cap	4°C	14 days for extraction 40 days after extraction for analysis
Base/Neutral and Acid Extractables - Concentrated Waste	G, teflon-lined cap	None	7 days for extraction 40 days after extraction for analysis
BTEX - Solid and sludge	G, teflon-lined septum	4°C	14 days
BTEX - Water	G, teflon-lined septum	4°C, 0.008% Na ₂ S ₂ O ₃ , zero headspace	14 days
TPH-GRO	G, teflon-lined cap	4°C, HCl to pH s, zero headspace	14 days
TPH-DRO	G, teflon-lined cap	4°C	14 days
Volatiles - Groundwater	G, teflon-lined cap	4°C, HCl to pH s, zero headspace	14 days
Chlorinated Herbicides - Water	Amber G, teflon-lined cap	4°C 0.008% sodium thiosulfate solution	7 days for extraction 40 days after extraction for analysis
Chlorinated Herbicides - Solid and Waste	G, teflon-lined cap	4°C	14 days for extraction 40 days after extraction
Volatiles - Drinking Water	G, teflon-lined cap	4°C, 0.008% Na ₂ S ₂ O ₃ ³ , zero headspace	14 days
Volatiles (excluding 2 chloroethylvinylether) - Wastewater	Encore Sampler	4°C, zero head-space, HCl to pH 2	14 days
Volatiles - Wastewater	G, teflon-lined cap	4°C, 0.008% Na ₂ S ₂ O ₃ ³ , zero headspace	7 days
Volatiles - Solid and Sludge -	Encore Sampler	4°C	14 days
Volatiles - Concentrated Waste	G, teflon-lined septum	None	14 days
Industrial Solvents	G, teflon-lined septum	4°C	None
Organochlorine Pesticides and PCBs	Amber G, teflon-lined cap	4°C 0.008% sodium thiosulfate solution	7 days for extraction 40 days after extraction for analysis
PCBs in Oil	G, teflon-lined cap	None	7 days for extraction 40 days after extraction for analysis
Dioxin	G, teflon-lined cap	4°C	7 days for extraction 40 days after extraction for analysis
Total Petroleum Hydrocarbon	G, teflon-lined septum	4°C	14 days
Coliform, total - Drinking water	P,G	4°C, 0.008% Na ₂ S ₂ O ₃	30 hours
<u>Radiochemistry</u>			
Carbon-14 - Water and Soil	P	4°C	6 months
Gamma Isotopes - Water	P	HNO ₃ to pH-2	6 months
Gamma Isotopes - Soil	P	None	6 months
Gross Alpha and Beta - Water	P	HNO ₃ to pH-2	6 months
GENERAL ENGINEERING LABORATORIES, LLC. 2040 Savage Road Charleston SC 29417 This document is controlled when an original Set ID number is stamped on the front cover page (1). Uncontrolled documents do not bear an original Set ID number			

Gross Alpha and Beta - Soil	P	None	6 months
Iodine-129 - Water and Soil	P	None	6 months
Iodine-131 - Water	P	None	6 months
Neptunium - Water	P	HNO ₃ to pH-2	6 months
Neptunium - Soil, Vegetation, and Air Filters	P	None	6 months
Plutonium - Water	P	HNO ₃ to pH-2	6 months
Plutonium - Soil, Vegetation, and Air Filters	P	None	6 months
Thorium - Water	P	HNO ₃ to pH-2	6 months
Thorium - Soil, Vegetation, and Air Filters	P	None	6 months
Uranium - Water	P	HNO ₃ to pH-2	6 months
Uranium - Soil, Vegetation, and Air Filters	P	None	6 months
Americium - Water	P	HNO ₃ to pH-2	6 months
Americium - Soil, Vegetation, and Air Filters	P	None	6 months
Curium - Water	P	HNO ₃ to pH-2	6 months
Curium - Soil, Vegetation, and Air Filters	P	None	6 months
Lead-210 - Water	P	HNO ₃ to pH-2	6 months
Nickel-59 - Water and Soil	P	None	6 months
Nickel-63 - Water and Soil	P	None	6 months
Phosphorus-32 -Water	P	HNO ₃ to pH-2	6 months
Phosphorus-32 -Soil	P	None	6 months
Polonium -Water	P	HNO ₃ to pH-2	6 months
Polonium -Soil	P	None	6 months
Promethium-147 -Water	P	HNO ₃ to pH-2	6 months
Promethium-147 -Soil	P	None	6 months
Radium-223 - Water	P	None	6 months
Radium-224 - Water	P	None	6 months
Radium-226 - Water	P	HNO ₃ to pH-2	6 months
Radium-228 - Water	P	HNO ₃ to pH-2	6 months
Radon-222 - Water	40ml volatile bottle	4°C, Zero headspace	7 days
Radon-222 - Soil	P	4°C	6 months
Strontium-89/90 -Water	P	HNO ₃ to pH-2	6 months
Strontium-89/90 -Soil	P	None	6 months
Technetium-99 -Water	P	HNO ₃ to pH-2	6 months
Technetium-99 -Soil	P	None	6 months
Total Alpha Radium -Water	P	HNO ₃ to pH-2	6 months
Total Alpha Radium -Soil	P	None	6 months
Total Uranium -Water	P	HNO ₃ to pH-2	6 months
Tritium - Water, Soil, Vegetation, and Air Filters	P	4°C	6 months
Iron 55 -Water	P	HNO ₃ to pH-2	6 months
Iron 55 -Soil	P	None	6 months
Total Uranium -Soil	P	None	6 months

¹ P = Polyethylene; G = Glass

² Samples should be analyzed as soon as possible after collection. The holding times listed are maximum times that samples may be held before analysis and be considered valid.

³ Used only in the presence of residual chlorine.

⁴ Maximum holding time is 24 hours when sulfide is present. All samples may be tested with lead acetate paper before pH adjustments in order to determine if sulfide is present. If present, remove by adding cadmium nitrate powder until a negative spot test is obtained. Filter sample and add NaOH to pH12.

APPENDIX 2: SAMPLE RECEIPT REVIEW SHEET



SAMPLE RECEIPT & REVIEW FORM

PM use only

Client:	SDG/ARCO/Work Order:
Date Received:	PM(A) Review (ensure non-conforming items are resolved prior to signing):
Received By:	

Sample Receipt Criteria	Conforming	NA	Non-Conforming	Comments/Qualifiers (Required for Non-Conforming Items)
1 Shipping containers received intact and sealed?				Circle Applicable: seals broken damaged container leaking container other (describe)
2 Samples requiring cold preservation within (4 +/- 2 C)? Record preservation method.				Circle Temp device serial # ice bags blue ice dry ice none other(describe)
3 Chain of custody documents included with shipment?				
4 Sample containers intact and sealed?				Circle Applicable: seals broken damaged container leaking container other (describe)
5 Samples requiring chemical preservation at proper pH?				Sample ID's, containers affected and observed pH:
6 VOA vials free of headspace (defined as < 6mm bubble)?				Sample ID's and containers affected:
7 Samples received within holding time?				Id's and tests affected:
8 Sample ID's on COC match ID's on bottles?				Sample ID's and containers affected:
9 Date & time on COC match date & time on bottles?				Sample ID's affected:
10 Number of containers received match number indicated on COC?				Sample ID's affected:
11 COC form is properly signed in relinquished/received sections?				
12 Air Bill ,Tracking #'s, & Additional Comments				

Radiological Information	Non-RAD	RAD	RADI	RSO RAD Receipt #
What is the radiological classification of the samples?				Comments:
Radioactivity Screening Results (maximum observed CPM)				*If > x2 area background is observed on a non-radioactive sample, contact the RSO to investigate.
PM (or PMA) review of Receiving Rad classification: _____ Initials _____ Date _____				

APPENDIX 3

BIOASSAY SAMPLE RECEIPT REVIEW

CLIENT _____

GEL COOLER _____

CLIENT COOLER _____

OTHER _____

SAMPLE REVIEW CRITERIA

YES NO

COMMENTS/QUALIFIERS

Were shipping containers received intact and sealed?			
Were chain of custody documents included?			
Were chain of custody documents completed correctly?			
Were all sample containers properly labeled?			
Were all sample containers received?			
Were samples received within holding time?			
Were the sample containers 500 mL or less?			
For KHCO - Did the sample ID and the customer number match the			
Chain of custody?			

Signature: _____

Date: _____

TABLE OF CONTENTS

1.0	STANDARD OPERATING PROCEDURE FOR THE DETERMINATION OF PLUTONIUM, URANIUM, AMERICIUM, CURIUM, AND THORIUM.....	3
2.0	METHOD OBJECTIVE AND APPLICABILITY	3
3.0	METHOD APPLICABILITY	3
4.0	DEFINITIONS.....	3
5.0	METHOD VARIATIONS	3
6.0	SAFETY PRECAUTIONS AND WARNINGS	3
7.0	INTERFERENCES.....	4
8.0	APPARATUS, MATERIALS, REAGENTS, EQUIPMENT, AND INSTRUMENTATION	4
9.0	SAMPLE HANDLING AND PRESERVATION	5
10.0	SAMPLE PREPARATION	6
11.0	PREPARATION OF STANDARD SOLUTIONS AND QUALITY CONTROL STANDARDS ...	8
12.0	INSTRUMENT CALIBRATION AND PERFORMANCE.....	8
13.0	ANALYSIS AND INSTRUMENT OPERATION	9
14.0	EQUIPMENT AND INSTRUMENT MAINTENANCE.....	9
15.0	DATA RECORDING, CALCULATION, AND REDUCTION METHODS	9
16.0	QUALITY CONTROL REQUIREMENTS	10
17.0	DATA REVIEW, APPROVAL, AND TRANSMITTAL	11
18.0	RECORDS MANAGEMENT	11
19.0	LABORATORY WASTE HANDLING AND WASTE DISPOSAL	11
20.0	REFERENCES.....	11
	APPENDIX 1.....	12
	APPENDIX 2.....	13
	APPENDIX 3.....	14

1.0 STANDARD OPERATING PROCEDURE FOR THE DETERMINATION OF PLUTONIUM, URANIUM, AMERICIUM, CURIUM, AND THORIUM

2.0 METHOD OBJECTIVE AND APPLICABILITY

- 2.1 This standard operating procedure provides the necessary instructions to conduct the analysis for isotopic plutonium, uranium, americium, and thorium in variety of matrices.
- 2.2 A sample is aliquoted and if necessary digested. Actinide elements are scavenged by co-precipitation with iron hydroxide. The precipitate is dissolved and separation of elements is accomplished through ion exchange resins and solid phase extraction. The elements are then prepared for the measurement of radioactive isotopes by co-precipitation with Neodymium fluoride. The Neodymium fluoride precipitate is trapped on a filter, mounted on a stainless steel disk and placed in a partially evacuated chamber for measurement of isotopic alpha emission.
- 2.3 This method has been modified on the basis of GEL's Performance Based Measurement System (PBMS).

3.0 METHOD APPLICABILITY

- 3.1 Method Detection Limit: Typical minimum detectable activity (MDA) for samples analyzed for plutonium, uranium, and thorium is 1 pCi/L or 1 pCi/g.
- 3.2 Method Precision: Typical relative percent difference (RPD) is 20%.
- 3.3 Method Bias (Accuracy): Acceptable criteria for method accuracy, measured by running with each batch a laboratory control sample, is $\pm 25\%$ of true value.
- 3.4 Analysts are trained and certified to run this analysis after the analyst has completed a batch with acceptable duplicate and laboratory control sample, as well as completed an unknown sample within $\pm 25\%$ of true value. Analyst training records are kept on hand in the Human Resources department.

4.0 DEFINITIONS

- 4.1 National Institute of Standards and Technology (NIST). For the purpose of this method, the national scientific body responsible for the standardization and acceptability of analyte solutions.
- 4.2 Type II water: Deionized (DI) water.
- 4.3 ALPHA LIMS: Laboratory Information Management System. The database system used to store and report data.

5.0 METHOD VARIATIONS

Some variation may be necessary due to special matrices encountered in the laboratory. These variations may be used with approval from a Group Leader or Team Leader. Variations to a method will be documented with the analytical raw data.

6.0 SAFETY PRECAUTIONS AND WARNINGS

- 6.1 Personnel performing this analytical procedure are trained in and follow the safe laboratory practices outlined in the Safety, Health & Chemical Hygiene Plan, GL-LB-N-001.
- 6.2 Personnel handling Radioactive Materials are trained in and follow the procedures outlined in GL-RAD-S-004 for Radioactive Material Handling.
- 6.3 Personnel handling biological materials are trained in and follow the procedures outlined in GL-RAD-S-010 for Handling Biological Materials.

- 6.4 If there is any question regarding the safety of any laboratory practice, **stop immediately**, and consult senior qualified personnel such as a Group or Team Leader.

7.0 INTERFERENCES

- 7.1 Internal tracer standards may have ingrown daughters that may interfere with the analysis. For example, Th-228 will be present in aged U-232 standard. This problem is overcome by using a self-cleaning U-232 tracer that utilizes barium sulfate to remove the Th-228.
- 7.2 Uranium and Thorium may not be run together if U-236 is used as a tracer due to the decay of U-236 to Th-232. The U-236 tracer does not currently utilize a barium sulfate clean-up to remove Th-232.
- 7.3 Short-lived radioactive progeny may ingrow on prepared filters. For example, the Ra-224 alpha peak will be present if the Th-228 parent is present. Counting samples as soon as possible after separation chemistry is completed minimizes this interference.
- 7.4 When present, Th-228 alpha energies interfere with the proper quantification of Pu-238. Steps are taken to ensure that the final plutonium counting sources are free from thorium interference.
- 7.5 Pu-236 decays to U-232 therefore Pu-236 may not be used as a tracer for this method.

8.0 APPARATUS, MATERIALS, REAGENTS, EQUIPMENT, AND INSTRUMENTATION

8.1 Ancillary Equipment

- 8.1.1 Ion exchange columns
- 8.1.2 Polypropylene centrifuge tube (50 mL)
- 8.1.3 Sample drying apparatus
- 8.1.4 Sample homogenizing apparatus
- 8.1.5 AG1X8 anion exchange resin, 100 – 200 mesh
- 8.1.6 Eichrom Technologies TRU Resin, 100 – 200 mesh
- 8.1.7 Silicon surface barrier detectors with associated electronics, vacuum chambers, and data reduction capabilities
- 8.1.8 Vacuum pump and filtration apparatus (25 mm)
- 8.1.9 Gelman 25 mm filters with 0.1 μ pore size
- 8.1.10 Gelman polypropylene 25 mm support filter
- 8.1.11 Stainless steel disks, 29 mm
- 8.1.12 Stainless steel tweezers
- 8.1.13 Hotplate

8.2 Reagents, Chemicals, and Standards

- 8.2.1 Ammonium hydroxide concentrated (14 N).
- 8.2.2 Neodymium (500 mg/L)
- 8.2.3 Carbon Colorant. Place two 47mm cellulose nitrate filters in a beaker and add 5mL concentrated H₂SO₄. Cover and heat on a hot plate until fumes of H₂SO₄ appear. Cool. Slurry the residue in DI water and dilute to 1L with DI water.
- 8.2.4 Sulfuric acid concentrated (18N)

- 8.2.5 Hydrochloric acid concentrated (12 M).
- 8.2.6 Iron Carrier (10mg/mL). Dissolve 62.7 g of $\text{Fe}(\text{NO}_3)_3 \bullet 6\text{H}_2\text{O}$ or 72.3 g $\text{Fe}(\text{NO}_3)_3 \bullet 9\text{H}_2\text{O}$ in 800 mL DI water and dilute to 1 L with DI water.
- 8.2.7 Hydrochloric acid (9 M). Dilute 750 mL of concentrated hydrochloric acid to 1 L with DI water.
- 8.2.8 Hydrochloric acid (6 M). Dilute 500 mL of concentrated hydrochloric acid to 1 L with DI water.
- 8.2.9 Hydrochloric acid (4 M). Dilute 333 mL of concentrated hydrochloric acid to 1 L with DI water.
- 8.2.10 Hydrochloric acid (2 M). Dilute 167 mL of concentrated hydrochloric acid to 1 L with DI water.
- 8.2.11 Hydrochloric acid (0.1 M). Dilute 8.3 mL of concentrated hydrochloric acid to 1 L with DI water.
- 8.2.12 Hydrogen peroxide (30%).
- 8.2.13 Hydrazine dihydrochloride (25%). Dissolve 25 g of hydrazine dihydrochloride in 75 mL of DI water.
- 8.2.14 Hydrofluoric acid concentrated (49%).
- 8.2.15 Ethyl alcohol (80%). Dilute 400 mL ethanol to 500 mL with DI water.
- 8.2.16 NIST traceable standards: Pu-242, Pu-238, Pu-239, Th-229, Th-230, Th-232, U-232, U-236, U-238, Am-243, Am-241, Cm-244
- 8.2.17 Nitric acid concentrated (16 M).
- 8.2.18 6 M Hydrochloric acid / 0.52 M Hydrofluoric acid. Dilute 500 mL concentrated hydrochloric acid and 16.8 mL concentrated hydrofluoric acid to 1 L with DI water.
- 8.2.19 9 M Hydrochloric acid / 0.05 M Ammonium iodide. Dissolve 7.24 g of ammonium iodide in 750 mL concentrated hydrochloric acid and dilute to 1 L with DI water. PREPARE DAILY.
- 8.2.20 9 M Hydrochloric acid / 0.04% Hydrogen peroxide. Add 8 drops of 30% H_2O_2 to 1 L 9M hydrochloric acid. PREPARE DAILY.
- 8.2.21 1 M Hydrochloric acid / 0.05 M Oxalic acid. Dissolve 6.3 g Oxalic acid in 83.5 mL Hydrochloric acid and dilute to 1 L with DI water.
- 8.2.22 Titanium (III) Chloride. 20% reagent

9.0 SAMPLE HANDLING AND PRESERVATION

- 9.1 Samples should be preserved to approximately pH 2 with nitric acid and collected in a plastic bottle.
- 9.2 Before beginning an analysis the analyst should check the sample pH with a pH strip. If necessary, adjust the pH with nitric acid to a pH 1 – 2. If the sample was pH adjusted let the sample sit overnight before continuing the batch.
- 9.3 If the sample has exceeded the hold time the analyst should contact the Group Leader or Team Leader before continuing with the batch.
- 9.4 Soil samples require no preservation and may be shipped in any suitable container.

10.0 SAMPLE PREPARATION

10.1 Soil Sample Preparation

- 10.1.1 If not already done, prepare the sample by performing GL-RAD-A-021 "Preparation of Soils for the Determination of Radionuclides".
- 10.1.2 It is recommended that the samples be ashed in a muffle furnace by performing GL-RAD-A-021B "Soil Sample Ashing for the Determination of Radionuclides".
- 10.1.3 For plutonium, uranium, and thorium analysis, take an appropriate aliquot and digest as specified in GL-RAD-A-015 "Digestion for Soils and Sand".
- 10.1.4 Proceed to step 10.2.5.

10.2 Aqueous Sample Preparation

- 10.2.1 Add an appropriate aliquot of sample to a labeled beaker. Add a certified dpm of appropriate tracers to each sample.

NOTE: Other sample matrices, such as vegetation, air filters, tissue etc. are run as outlined in GL-RAD-A-026 "Preparation of Special Matrices for the Determination of Radionuclides".

- 10.2.2 Add 1 mL of iron carrier.
- 10.2.3 Bring to a slight boil and add concentrated NH_4OH until turbidity persists, or $\text{pH} > 9$. Heat to boiling for 10 minutes and then allow to settle and cool.
- 10.2.4 Decant excess supernate and discard. Collect the remaining precipitate by centrifugation in a 50 mL centrifuge tube and discard the supernate.

NOTE: Exercise care in this step because finely divided material that contains the actinides may also be present in addition to the large iron hydroxide flocks.

- 10.2.5 Dissolve the precipitate from 10.2.4 or the residue from 10.1.4 in 10 mL of 9 M HCl / 0.04% H_2O_2 .

NOTE: Samples may be dissolved with 10 to 15 mL of 9M HCl and then add 1 drop of 30% H_2O_2 as an alternative to dissolving with 9M HCl/0.04% H_2O_2 .

NOTE: The load solution needs to be 10 mL due to limitations of the americium portion of this procedure. If the load solution needs to be increased then refer to GL-RAD-A-011 for the determination of americium, plutonium, and uranium.

- 10.2.6 Slurry AG 1X8 anion resin (Cl^- form 100 – 200 mesh) in a squirt bottle with DI water. Transfer the resin to a small column to obtain a settled resin bed of 2.5 mL.
- 10.2.7 Condition the column with 10 mL 9 M HCl. Catch in a drip pan for disposal.
- 10.2.8 Pass the sample solution for step 10.2.5 through the column and catch the effluent in a labeled, disposable 50 mL centrifuge tube for americium and thorium analysis.
- 10.2.9 Rinse the column with 5 mL of 9 M HCl and catch in centrifuge tube for americium and thorium analysis. Proceed to step 10.2.15 for americium and thorium analysis.
- 10.2.10 Rinse the column with 15 mL of 9 M HCl and catch in a drip pan for disposal.

- 10.2.11 Elute plutonium with 15 mL 9 M HCl / 0.05 M NH₄I catching in a clean, labeled, disposable centrifuge tube. Proceed to step 10.2.21 for plutonium micro-precipitation source preparation for alpha spectroscopy.
 - 10.2.12 Rinse the column with 15 mL of 6 M HCl / 0.52 M HF and catch in a drip pan for disposal.
 - 10.2.13 Rinse the column with 5 mL of 6 M HCl and catch in a drip pan for disposal.
 - 10.2.14 Elute uranium from the column by adding 15mL of 0.1M HCl, catching the uranium elute in a labeled, disposable 50mL centrifuge tube.
 - 10.2.14.1 Transfer sample to a clean beaker with DI water and evaporate to dryness over low heat.
 - 10.2.14.2 Dissolve sample with 4mL of 2M HCl. Transfer to a clean centrifuge tube with DI water.
 - 10.2.14.3 Proceed to step 10.2.23 for uranium microprecipitation source preparation for alpha spectroscopy.
 - 10.2.15 Precondition a 2 mL TRU column with 5 mL of 9 M HCl.
 - 10.2.16 Pass the sample solution from step 10.2.9 through the column catching in a drip pan for disposal.
 - 10.2.17 Rinse the column with 5 mL of 9 M HCl catching in a drip pan for disposal.
 - 10.2.18 Elute Americium and Curium with 20 mL of 4 M HCl catching in a clean, labeled, disposable centrifuge tube. Proceed to step 10.2.24 for americium and curium micro-precipitation source preparation for alpha spectroscopy.
- NOTE:** If Americium and Curium are not required then catch this rinse in a drip pan for disposal.
- 10.2.19 Elute Thorium with 20 mL of 1 M HCl / 0.05 M Oxalic acid catching in a clean, labeled, disposable centrifuge tube.
 - 10.2.20 Add 0.1 mL Neodymium (500 mg/L) to the solution and swirl to mix. Add 2 mL concentrated Hydrofluoric acid and swirl. Allow the solution to sit for 30 minutes then proceed to Step 10.2.28 for source preparation.
 - 10.2.21 Transfer the plutonium solution from step 10.2.11 to a clean beaker. Add 4 drops of iron carrier (10 mg/mL) and 10 mL of concentrated nitric acid. Evaporate the solution to dryness on medium to low heat. Dissolve the residue with 4 mL of 2 M HCl. Use DI water to transfer the solution to a centrifuge tube.
 - 10.2.22 Add 0.1 mL of Neodymium (500 mg/L) and swirl. Add 3 to 4 drops of 25% Hydrazine dihydrochloride and swirl to mix. Let the solution sit for 10 minutes, then add 2 mL of 49% Hydrofluoric acid. Swirl to mix. Allow to sit for 30 minutes, then proceed to step 10.2.28 for source preparation.
 - 10.2.23 To the uranium solution from step 10.2.14, add 1 mL of titanium trichloride solution and allow the sample to sit for 30 seconds. Add 0.1 mL of Neodymium (500 mg/L) and swirl to mix. Add 2 mL of concentrated hydrofluoric acid to precipitate fluorides. Allow the solution to sit for 30 minutes, then proceed to step 10.2.28 for source preparation.
 - 10.2.24 Transfer the Am/Cm elution from step 10.2.18 to a clean beaker. Add 4 drops of iron carrier and gently cook dry.

- 10.2.25 Add 10 mL of concentrated Nitric acid and 2 mL of 30% Hydrogen peroxide, cover and reflux until evolution of gas bubbles ceases. Remove the cover and gently cook dry.
- 10.2.26 Dissolve the residue in 4 mL of 2 M HCl with gentle heating. If residue does not dissolve repeat Step 10.2.25. Cool and transfer to a clean disposable centrifuge tube using 1 to 2 mL of DI water to rinse the beaker.
- 10.2.27 Add 0.1 mL Neodymium (500 mg/L) to the solution and swirl to mix. Add 2.0 mL concentrated Hydrofluoric acid and swirl. Allow the solution to sit for 30 minutes then proceed to Step 10.2.28 for source preparation.
- 10.2.28 Place a 25mm 0.1 μ m filter on the filter support screen. Wet the filter with 80% ethyl alcohol. Center the filter on the filter support screen and apply vacuum.
- 10.2.29 Place the filter funnel on the filter stem. Tighten firmly, being careful not to wrinkle the filter.
- 10.2.30 Rinse the funnel with 80% ethyl alcohol.
- 10.2.31 Add 1 mL of the carbon colorant.
- 10.2.32 Filter the fluoride precipitated solution through the filter paper. Rinse the centrifuge tube with \approx 5 mL DI water and pass through filter.
- 10.2.33 Rinse the centrifuge tube with \approx 5 mL 80% ethyl alcohol and pass through filter.
- 10.2.34 Rinse the funnel with 80% ethyl alcohol.
Caution – Directing a stream of liquid onto the filter will disturb the distribution of the precipitate on the filter and render the sample unsuitable for alpha spectrometry resolution.
- 10.2.35 Without turning off the vacuum, remove the funnel.
- 10.2.36 Turn off vacuum and remove filter. Mount filter on a labeled 29mm flat planchet. Ensure that the filter is centered and as flat as possible on the planchet.
NOTE: Care should be taken not to touch the active area of the filter with tweezers.
- 10.2.37 Place the mounted filter under a heat lamp for 5 minutes prior to alpha spectrometry measurement.
- 10.2.38 Count under vacuum on the alpha spectrometer long enough to reach requested MDA. Consult the operating manual for instruction on operating the alpha spectrometer.

11.0 PREPARATION OF STANDARD SOLUTIONS AND QUALITY CONTROL STANDARDS

Refer to “Preparation of Radioactive Standards” (GL-RAD-M-001).

12.0 INSTRUMENT CALIBRATION AND PERFORMANCE

For direction on calibration and instrument performance see “The Alpha Spectroscopy System” (GL-RAD-I-009).

13.0 ANALYSIS AND INSTRUMENT OPERATION

For analysis and instrument operation see “The Alpha Spectroscopy System” (GL-RAD-I-009).

14.0 EQUIPMENT AND INSTRUMENT MAINTENANCE

For maintenance of system see “Counting Room Instrumentation Maintenance and Performance Checks” (GL-RAD-I-010).

15.0 DATA RECORDING, CALCULATION, AND REDUCTION METHODS

15.1 The instrument will report sample pCi/unit according to the following equation:

$$\text{pCi / unit} = \frac{S_{\text{cpm}} - B_{\text{cpm}}}{2.22 * E * V * A * \text{decay} * R}$$

15.2 Counting uncertainty is propagated according to the following equation:

$$\text{pCi / unit} = A_c * 1.96 \sqrt{\left(\frac{ef_er}{E}\right)^2 + \left(\frac{pk_er}{pk}\right)^2 + \left(\frac{ab_er}{A}\right)^2 + \left(\frac{sy}{100}\right)^2 + (dk)^2}$$

15.3 The minimum detectable activity (MDA) is calculated according to the following equation:

$$\text{MDA}(\text{pCi / unit}) = \frac{2.71 + 4.65 * \sqrt{B_{\text{cpm}} * T_c}}{(2.22 * E * V * R * A * \text{decay} * T_c)}$$

Where:

$$\text{decay} = e^{\left(\frac{-\ln(2)T_d}{T_{1/2}}\right)}$$

$$R = \frac{T_{\text{cpm}} - B_{\text{cpm}}}{T_{\text{dpm}} * E}$$

$$dk = \frac{T_{1/2}^{\text{err}}}{T_{1/2}} * \left(\frac{\lambda T_r}{1 - e^{-\lambda T_r}} - \lambda (T_c + T_r) - 1 \right)$$

And where:

S_{cpm}	=	Sample counts per minute
B_{cpm}	=	Background counts per minute
E	=	Counting efficiency (decimal form)
V	=	Volume in liters, g, cfm, etc.
A	=	Isotopic abundance (decimal form)
ef_er	=	1 sigma efficiency error (decimal form)
pk_er	=	1 sigma peak error
ab_er	=	1 sigma isotopic abundance error (decimal form)
sy	=	1 sigma systematic error
pk	=	peak area
T_c	=	Sample count time in minutes
T_d	=	Time interval for radioactive decay

Tr	=	Elapsed real time in minutes
T _{cpm}	=	Tracer counts per minute
T _{dpm}	=	Tracer known disintegrations per minute
Ac	=	Sample calculated activity
T _{1/2}	=	Isotopic half life
T _{1/2err}	=	Isotopic half life error
λ	=	Isotopic decay constant
E	=	exponential function
R	=	Tracer Recovery
ln	=	natural log function

- 15.4 Record the following information on the alpha que sheet: preparation date, analyst initials, spike isotope, spike code, spike volume, LCS isotope, LCS code, LCS volume. For each sample record the detector number, sample mass, sample date, and sample time.

16.0 QUALITY CONTROL REQUIREMENTS

16.1 Analyst and Method Verification

- 16.1.1 Refer to “Analytical Methods Validation for Radiochemistry” (GL-RAD-D-002) for instructions concerning the validation of analysts and analytical methods.

16.2 Method Specific Quality Control Requirements

- 16.2.1 A method blank will accompany each batch of 20 or less samples. The reported value should be less than or equal to the CRDL for all target isotopes.
- 16.2.2 A matrix spike (MS) should be run with every batch of 20 samples. The recovery of the spike should fall between 75 and 125%. The recovery is calculated as follows:

$$\% \text{Rec} = \frac{\text{spike}(\text{pCi / unit}) - \text{sample}(\text{pCi / unit})}{\text{spiked amount}(\text{pCi / unit})} * 100$$

- 16.2.3 A sample duplicate should be run with every batch of 20 or less samples. The relative percent difference (RPD) between the sample and the duplicate should be less than or equal to 20%. The RPD is calculated as follows.

$$\text{RPD} = \frac{\text{high sample}(\text{pCi / unit}) - \text{low sample}(\text{pCi / unit})}{\text{Average}(\text{pCi / unit})} * 100$$

- 16.2.4 A laboratory control spike (LCS) should be run with every batch of 20 samples or less. The recovery of the spike should fall between 75 and 125%. The recovery is calculated as follows:

$$\text{LCS} = \frac{\text{observed_pCi / unit}}{\text{known_pCi / unit}} * 100$$

16.3 Actions required if the Quality Control Requirements are not met

- 16.3.1 If any of the above criteria cannot be satisfied, the analyst should inform the Group Leader and initiate a non-conformance report as outlined in "Documentation of Nonconformance Reporting and Dispositioning, and Control of Nonconforming Items" (GL-QS-E-004).

17.0 DATA REVIEW, APPROVAL, AND TRANSMITTAL

Refer to “Data Review, Validation and Data Package Assembly” (GL-RAD-D-003) for instructions concerning the data review process, approval, and transmittal.

18.0 RECORDS MANAGEMENT

- 18.1 All raw data, calculation spreadsheets and batch checklists are filed with the sample data and maintained as quality records in accordance with GL-QS-E-008 for Quality Records Management and Disposition.
- 18.2 Each analysis that is performed on the instrument is documented in the run log according to “Run Logs” (GL-LB-E-009).
- 18.3 All raw data printouts, calculation spreadsheets and batch checklists are filed with the sample data for archival and review.

19.0 LABORATORY WASTE HANDLING AND WASTE DISPOSAL

Radioactive samples and material shall be handled and disposed of as outlined in the Laboratory Waste Management Plan (GL-LB-G-001).

20.0 REFERENCES

- 20.1 EPA Environmental Monitoring and Support Laboratory. Las Vegas. Radiochemical Analytical Procedures for Analysis of Environmental Samples. March 1979.
- 20.2 EML Procedures Manual HASL-300, 1982.
- 20.3 DOE Methods Manual for Evaluating Environmental and Waste Management Samples, 1997 Edition, RP800, “Sequential Separation of Americium and Plutonium by Extraction Chromatography”.
- 20.4 Analytical Chemistry. Rapid Determination of Th-230 in Mill Tailings by alpha spectroscopy. UNC Geotech, Grand Junction Projects Office. Steve Donovan, Mark Hollenbach, and Mary Costello. Vol. 59, No. 21, 1987.
- 20.5 Los Alamos Health and Environmental Chemistry: Analytical Techniques. LA-10300-M Vol. 1, September 1987.

APPENDIX 1

PLUTONIUM, URANIUM, AMERICIUM, CURIUM, and THORIUM

Use a 2.5 cm³ column with 1X8 anion resin (Cl⁻ form 100-200 mesh)

COLUMN WORK

- _____ 10mL 9M HCl (Conditioning)
- _____ Load solution: 10mL 9M HCl / 0.04% H₂O₂ (Catch in C-Tube for Am/Cm)
- _____ 5mL 9M HCl (Catch in C-Tube for Am/Cm/Th then proceed to Appendix 2 for Am/Cm/Th procedure)
- _____ 15mL 9M HCl (Rinse)
- _____ **Elute** Pu: 15mL 9M HCl / 0.05M NH₄I (Catch in C-Tube then proceed to Plutonium Cook-down)
- _____ 15mL 6M HCl / 0.52M HF (Rinse)
- _____ 5mL 6M HCl (Rinse)
- _____ **Elute** U: 15mL 0.1M HCl (Catch in C-Tube)
- _____ Transfer to a clean beaker and evaporate to dryness
- _____ Dissolve with 4mL of 2M HCl and transfer to centrifuge tube with DI water.
- _____ Proceed to Appendix 3 for Uranium Precipitation

PLUTONIUM COOK-DOWN

- _____ Transfer to a clean beaker. Add 4-6 drops of Fe carrier, 10 mL of [HNO₃], and evaporate the solution to dryness on medium heat.
- _____ Dissolve with 4 mL of 2 M HCl and transfer to centrifuge tube with DI water.
- _____ Proceed to Appendix 3 for Plutonium Precipitation

APPENDIX 2

AMERICIUM / CURIUM / THORIUM CONTINUATION

Use a 2.5 cm³ column with TRU Resin, 100 – 200 mesh

AMERICIUM / CURIUM / THORIUM

- _____ 5mL 9M HCl (Conditioning)
- _____ Load Solution from Appendix 1 (Catch in drip pan)
- _____ 5mL 9M HCl (Rinse)
- _____ **Elute** Am/Cm: 20mL 4M HCl (Catch in C-Tube)
- _____ Transfer to a clean beaker. Add 4-6 drops of Fe carrier, and gently cook dry
- _____ 10mL [HNO₃] and 2mL 30% H₂O₂, reflux, then gently cook dry
- _____ Dissolve with 4mL of 2M HCl and transfer to centrifuge tube with DI water
- _____ Proceed to Appendix 3 for Am/Cm precipitation
- _____ **Elute** Th: 20mL 1M HCl / 0.05M Oxalic acid (Catch in C-Tube)
- _____ Proceed to Appendix 3 for Th precipitation

APPENDIX 3

AMERICIUM / CURIUM PRECIPITATION

- _____ 0.1mL 500mg/L Neodymium and swirl
- _____ 2mL 49% HF and swirl
- _____ Wait 30 minutes
- _____ Filter

PLUTONIUM PRECIPITATION

- _____ 0.1 mL 500mg/L Neodymium and swirl
- _____ 3 – 4 drops 25% Hydrazine dihydrochloride and swirl
- _____ Wait 10 minutes
- _____ 2mL 49% HF and swirl
- _____ Wait 30 minutes
- _____ Filter

URANIUM PRECIPITATION

- _____ 0.1 mL 500mg/L Neodymium and swirl
- _____ 1 mL Titanium chloride and swirl
- _____ Wait 30 seconds
- _____ 2mL 49% HF and swirl
- _____ Wait 30 minutes
- _____ Filter

THORIUM PRECIPITATION

- _____ 0.1mL 500mg/L Neodymium and swirl
- _____ 2mL 49% HF and swirl
- _____ Wait 30 minutes
- _____ Filter