


## Appendix C-1. PCB Sequence

```
Sequence: C:\HPCHEM\1\SEQUENCE\ISWS\PCB.SEQ
Operator: MONTE
Sequence preparation date: 03 Mar 95  03:56 PM
Data File Subdirectory: 950303H
Part of methods to run: full method
On a barcode mismatch: inject anyway
Comment:
  Preparation dates:   595 + PES STD = 05 AUG 94
                   30,204 I.S.   = 21 NOV 94
                   342 QC CHK    = 05 AUG 94
  Blowdown volume: 1 ml
  Injection: 2 ul, column: DB5, 60 m.
  Sample set: 941222F
```

Sample Log Table


| Seq. Vial | Sample | Sample Multiplier | ISTD   | Cal. | Method | Inj/ |
|-----------|--------|-------------------|--------|------|--------|------|
| ne Num.   | Name   | Amount            | Amount | Line | Name   | Vial |
| FRONT     |        |                   |        |      |        |      |
| 1         | 21     | HEXANE            | 1      |      | BAKE   | 1    |
| 2         | 21     | HEXANE            | 1      |      | MULLIN | 1    |
| 2         | 22     | 595 STD 950303    | 1      |      | MULLIN | 1    |
| 2         | 23     | 595 STD 950303    | 1      |      | MULLIN | 1    |
| 2         | 24     | 342QCCHK 950303   | 1      |      | MULLIN | 1    |
| 2         | 25     | XPFB 950227 HX    | 1      |      | MULLIN | 1    |
| 2         | 26     | XP01 950227 HX    | 1      |      | MULLIN | 1    |
| 2         | 27     | LB 950206C HX     | 1      |      | MULLIN | 1    |
| 2         | 28     | LBXT4 950206C HX  | 1      |      | MULLIN | 1    |
| 2         | 29     | LBX14 950206C HX  | 1      |      | MULLIN | 1    |
| 2         | 30     | HEXANE            | 1      |      | MULLIN | 1    |
| 2         | 31     | 595 STD 950303    | 1      |      | MULLIN | 1    |
| 2         | 32     | CMS 941222F HX    | 1      |      | MULLIN | 1    |
| 2         | 33     | LB 941222F HX     | 1      |      | MULLIN | 1    |
| 2         | 34     | BH01F 941100 HX   | 1      |      | MULLIN | 1    |
| 2         | 35     | CH01F 941100 HX   | 1      |      | MULLIN | 1    |
| 2         | 36     | IH01F 941100 HX   | 1      |      | MULLIN | 1    |
| 2         | 37     | IH01F 941117 HX   | 1      |      | MULLIN | 1    |
| 2         | 38     | IH02F 941117 HX   | 1      |      | MULLIN | 1    |
| 2         | 39     | JH01F 941100 HX   | 1      |      | MULLIN | 1    |
| 2         | 40     | MH01F 941100 HX   | 1      |      | MULLIN | 1    |
| 2         | 41     | UH01F 941100 HX   | 1      |      | MULLIN | 1    |
| 2         | 42     | VH01F 941100 HX   | 1      |      | MULLIN | 1    |
| 2         | 43     | WH01F 941100 HX   | 1      |      | MULLIN | 1    |
| 2         | 44     | HEXANE            | 1      |      | MULLIN | 1    |
| 2         | 45     | 595 STD 950303    | 1      |      | MULLIN | 1    |

## Appendix C-2. PCB 595 ng/mL Standard Chromatogram

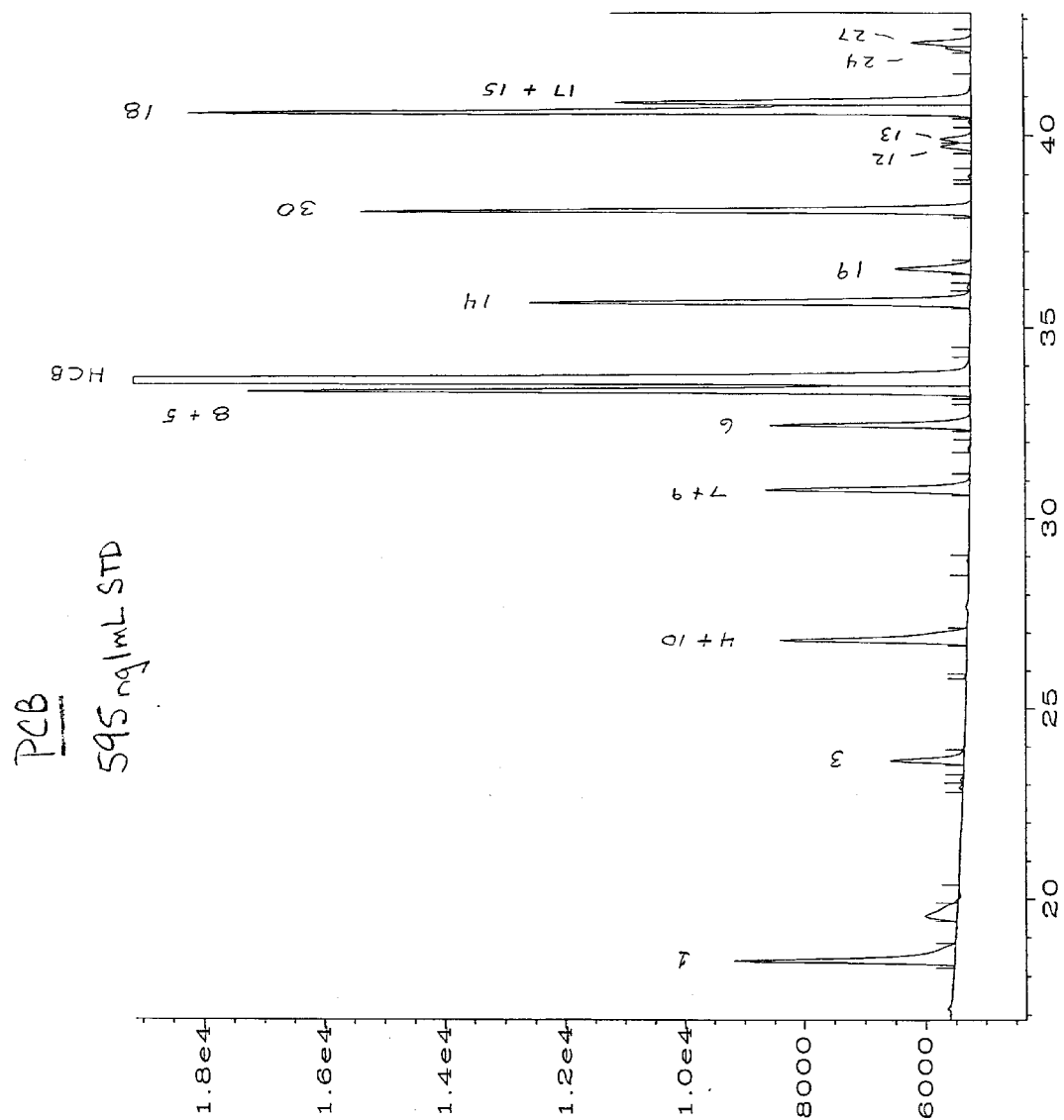


Sig. 1 in C:\HPCHEM\1\DATA\941219H\023F0201.D

# Appendix C-2. (Cont'd)

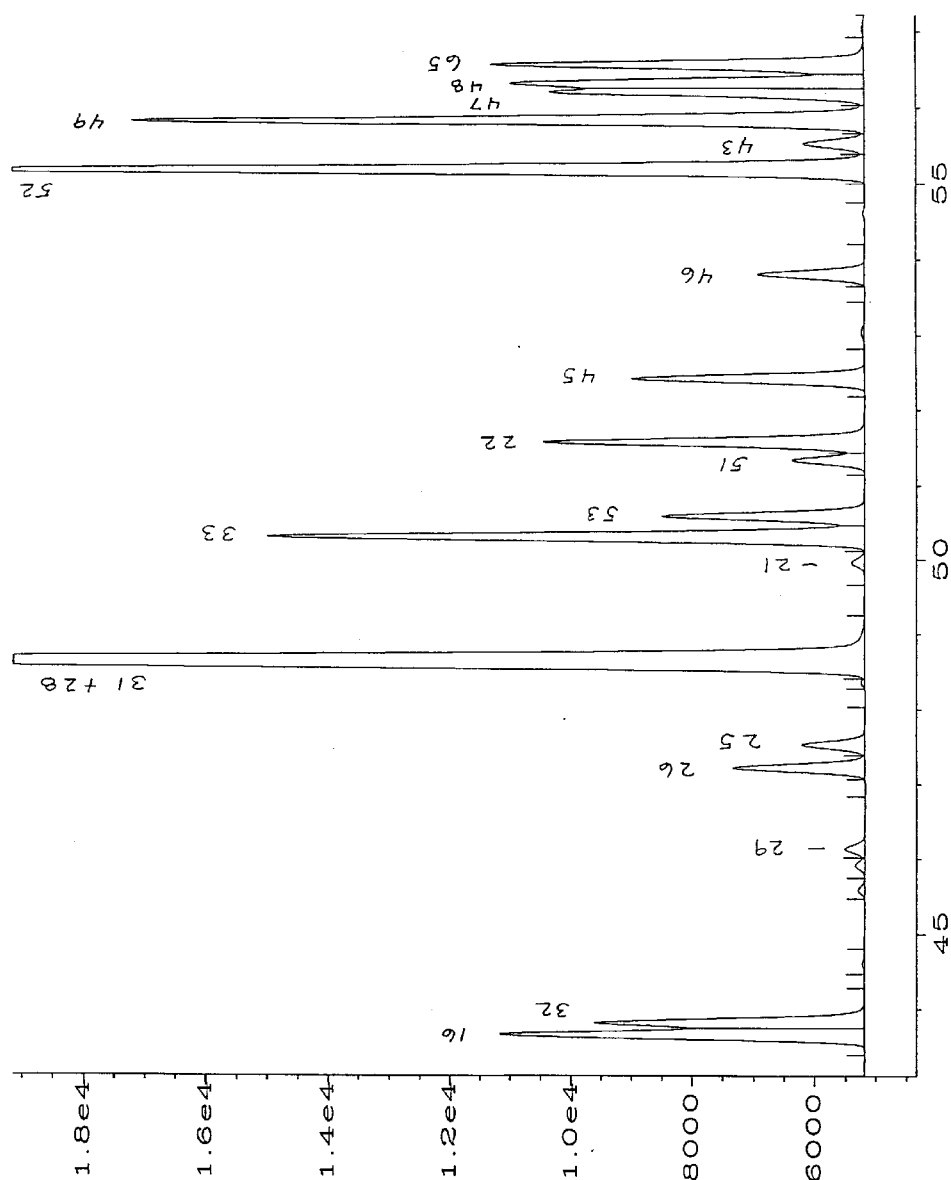


Sig. 1 in C:\HPCHEM\1\DATA\941219H\023F0201.D

## Appendix C-2. (Cont'd)

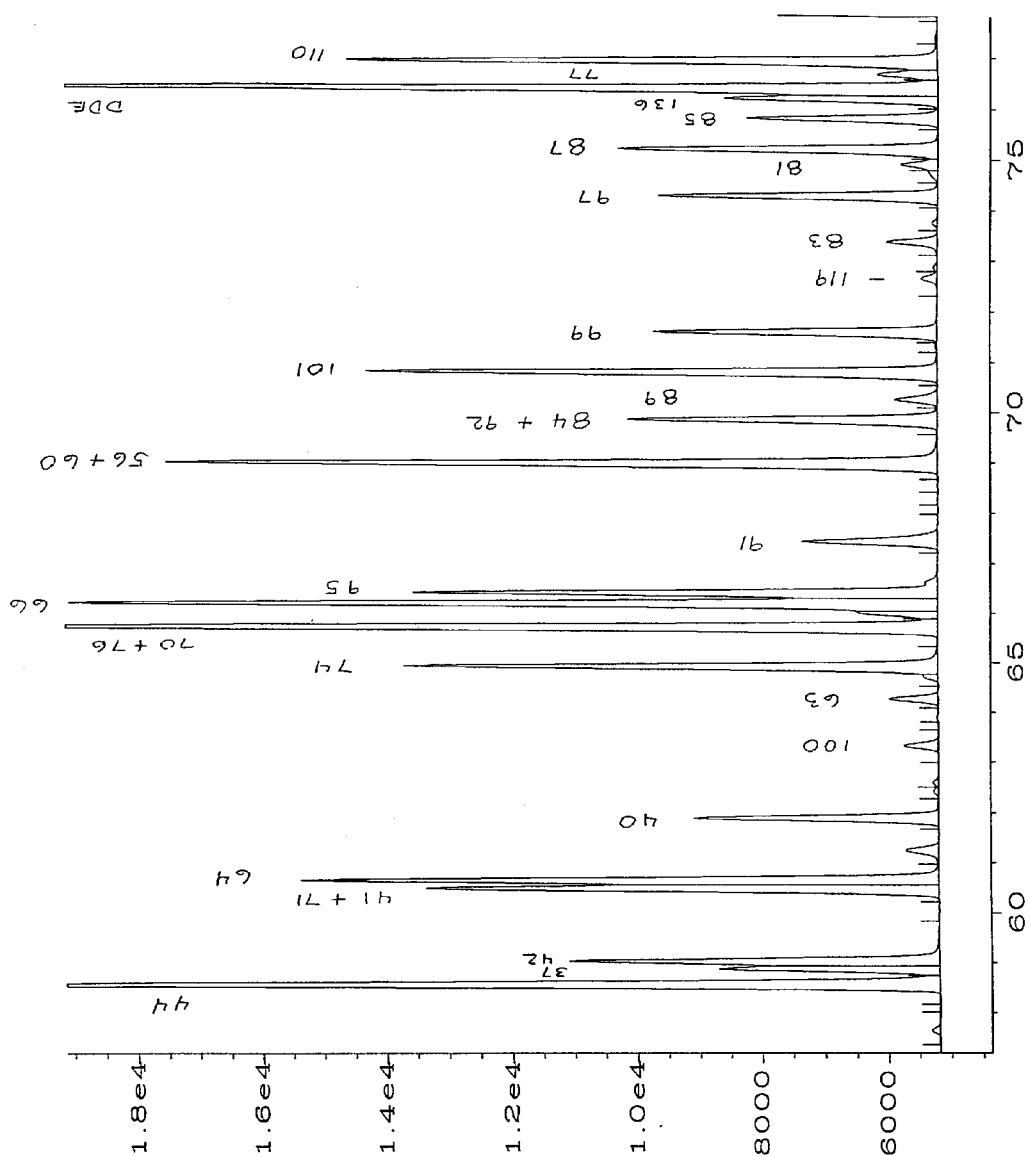


Sig. 1 in C:\HPCHEM\1\DATA\941219H\023F0201.D

## Appendix C-2. (Cont'd)

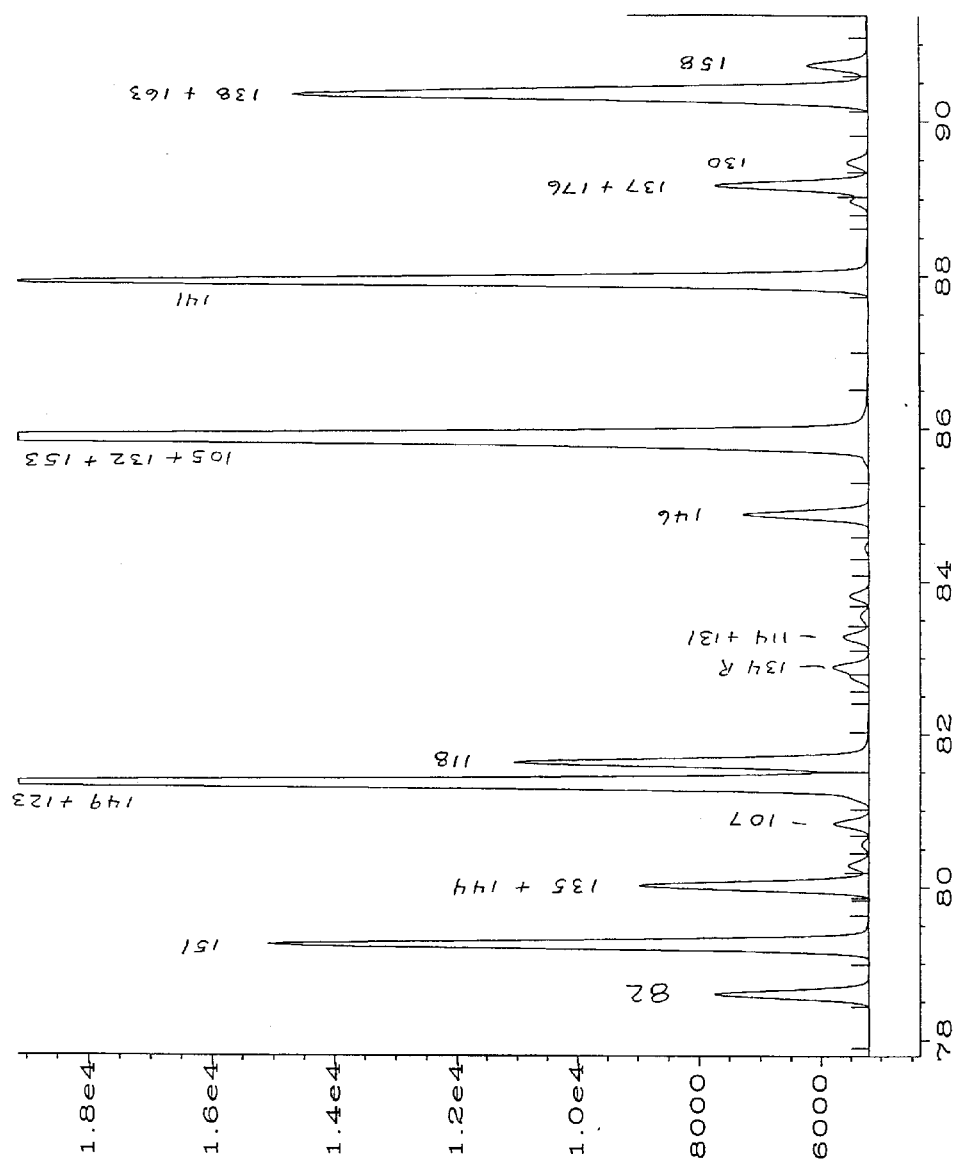


Sig. 1 in C:\HPCHEM\1\DATA\941219H\023F0201.D

## Appendix C-2. (Cont'd)

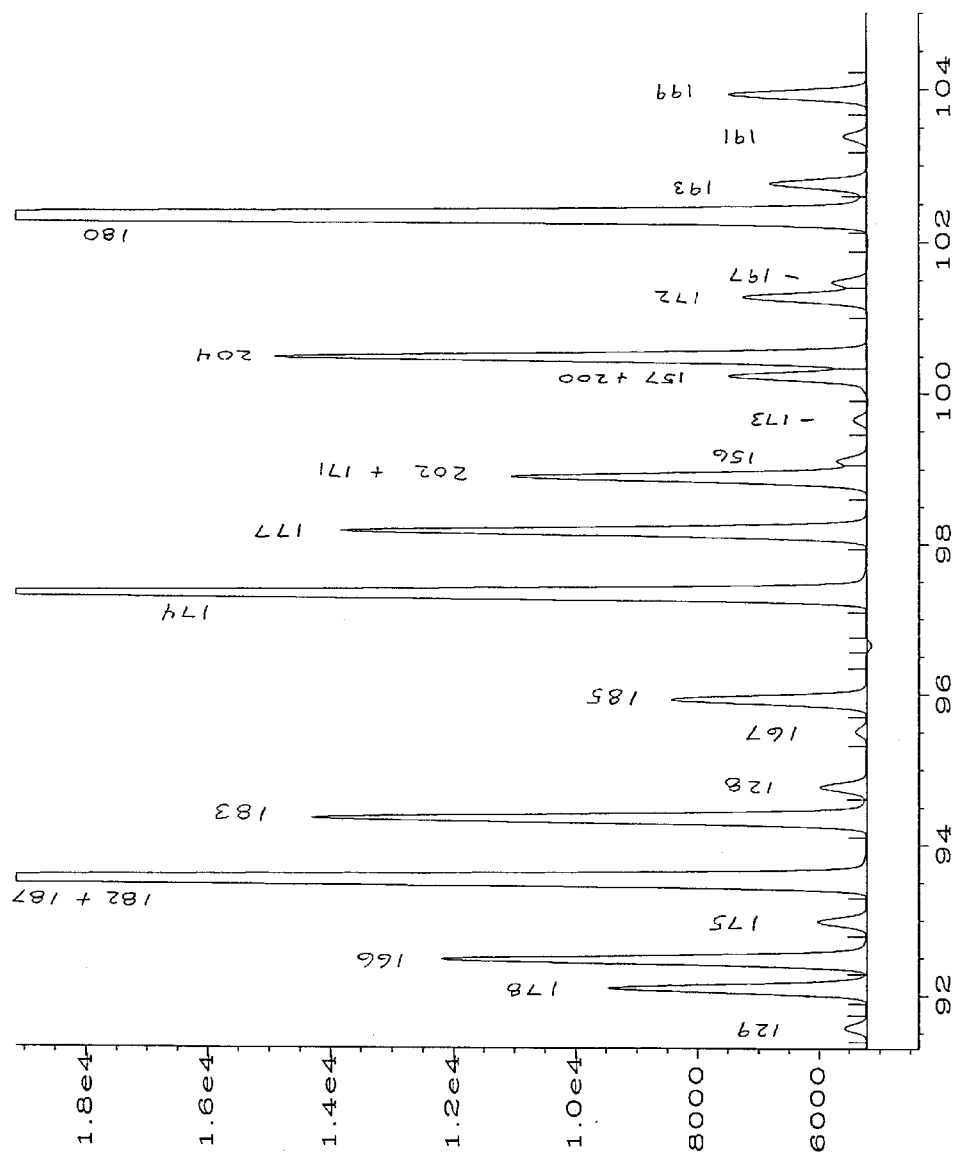

Sig. 1 in C:\HPCHEM\1\DATA\941219H\023FO201.D

## Appendix C-2. (Cont'd)

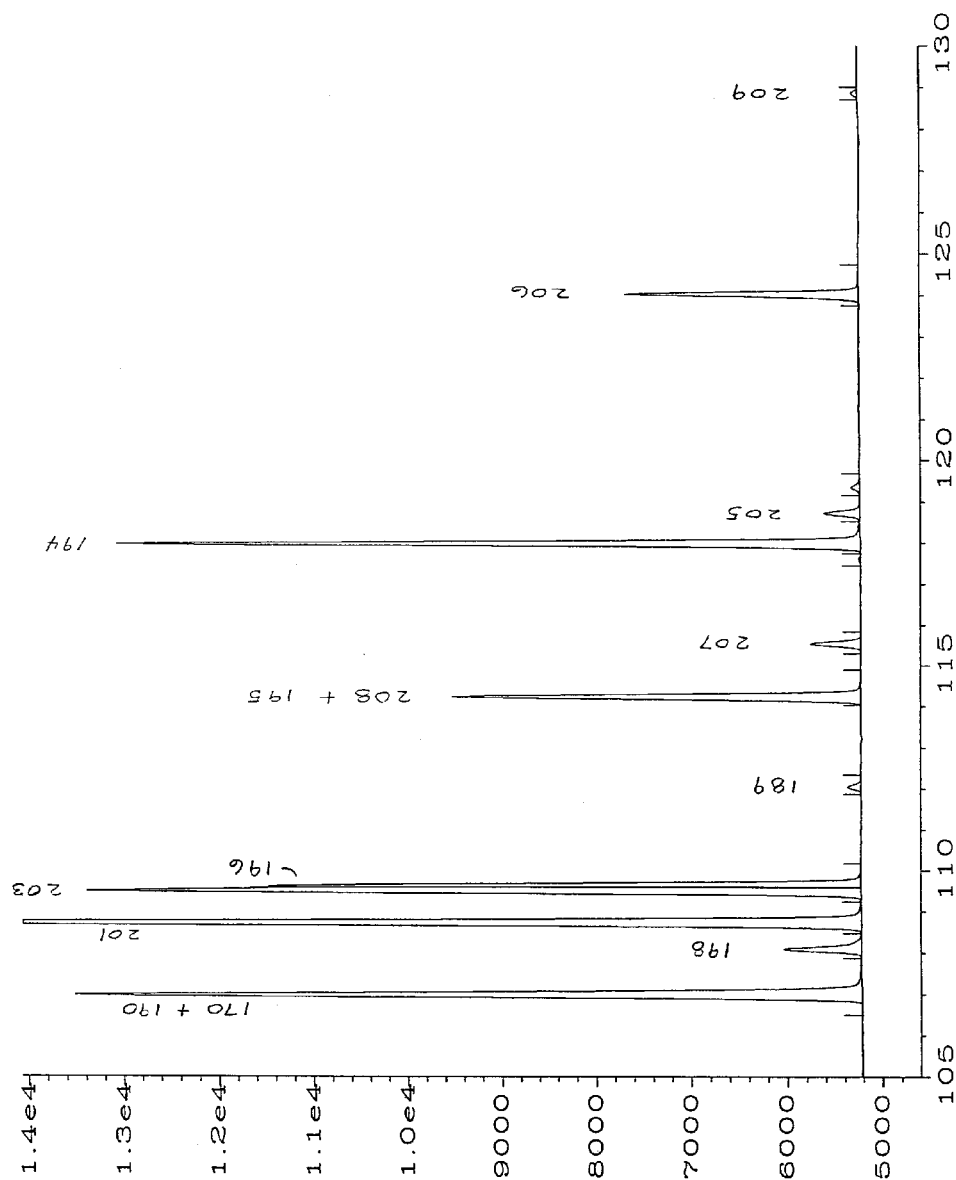

Sig. 1 in C:\HPCHEM\1\DATA\941219H\023F0201.D

## Appendix C-3. PCB Calibration Table

Method: C:\HPCHEM\1\METHODS\PCB1.MTH

| Calibration Table |        |     |        |             |          |           |              |
|-------------------|--------|-----|--------|-------------|----------|-----------|--------------|
| I                 | RT     | Lvl | ng     | Amt/Area    | Ref Istd | I#        | Name         |
| 1                 | 18.270 | 1   | 42.0   | 1.3056e-003 |          | 1 1       | (14)         |
| 2                 | 23.510 | 1   | 24.5   | 3.0529e-003 |          | 1 3       |              |
| 3                 | 26.667 | 1   | 11.9   | 3.9365e-004 |          | 1 4 + 10  |              |
| 4                 | 30.633 | 1   | 4.2    | 1.6441e-004 |          | 1 7 + 9   |              |
| 5                 | 32.320 | 1   | 6.65   | 2.8359e-004 |          | 1 6       |              |
| 6                 | 33.226 | 1   | 49.0   | 6.1031e-004 |          | 1 8 + 5   |              |
| 7                 | 33.481 | 1   | 40.0   | 7.0083e-005 |          | 1 HCB     |              |
| 8                 | 35.533 | 1   | 22.19  | 4.4475e-004 |          | 1 14      | surrogate    |
| 9                 | 36.402 | 1   | 0.98   | 8.7173e-005 |          | 1 19      |              |
| 10                | 37.938 | 1   | 9.12   | 1.3793e-004 | Ref ISTD | 1 30      | internal std |
| 11                | 39.578 | 1   | 0.595  | 1.8906e-004 |          | 1 12      |              |
| 12                | 39.772 | 1   | 0.3395 | 1.0793e-004 |          | 1 13      |              |
| 13                | 40.501 | 1   | 12.95  | 1.3403e-004 |          | 1 18      |              |
| 14                | 40.738 | 1   | 12.95  | 2.9468e-004 |          | 1 15 + 17 |              |
| 15                | 42.112 | 1   | 0.175  | 9.177e-005  |          | 1 24      |              |
| 16                | 42.259 | 1   | 0.735  | 1.0768e-004 |          | 1 27      |              |
| 17                | 43.516 | 1   | 7.0    | 1.6519e-004 |          | 1 16      |              |
| 18                | 43.678 | 1   | 6.65   | 2.2082e-004 |          | 1 32      |              |
| 19                | 45.988 | 1   | 0.1855 | 1.1757e-004 |          | 1 29      |              |
| 20                | 47.085 | 1   | 2.52   | 1.8308e-004 |          | 1 26      |              |
| 21                | 47.393 | 1   | 1.12   | 1.7135e-004 |          | 1 25      |              |
| 22                | 48.504 | 1   | 32.9   | 1.8577e-004 |          | 1 31 + 28 |              |
| 23                | 49.827 | 1   | 0.112  | 8.6874e-005 |          | 1 21      |              |
| 24                | 50.173 | 1   | 11.55  | 1.7372e-004 |          | 1 33      |              |
| 25                | 50.457 | 1   | 2.24   | 9.398e-005  |          | 1 53      |              |
| 26                | 51.210 | 1   | 0.63   | 7.6311e-005 |          | 1 51      |              |
| 27                | 51.444 | 1   | 10.15  | 2.8382e-004 |          | 1 22 (65) |              |
| 28                | 52.288 | 1   | 3.115  | 1.1769e-004 |          | 1 45      |              |
| 29                | 53.676 | 1   | 1.4    | 1.1351e-004 |          | 1 46      |              |
| 30                | 55.064 | 1   | 15.75  | 1.4022e-004 |          | 1 52      |              |
| 31                | 55.403 | 1   | 0.945  | 1.3413e-004 |          | 1 43      |              |
| 32                | 55.714 | 1   | 8.05   | 1.0043e-004 |          | 1 49      |              |
| 33                | 56.089 | 1   | 3.5    | 1.1243e-004 |          | 1 47      |              |
| 34                | 56.199 | 1   | 3.5    | 8.9794e-005 |          | 1 48      |              |
| 35                | 56.441 | 1   | 4.74   | 1.1166e-004 |          | 1 65      | surrogate    |
| 36                | 58.432 | 1   | 15.05  | 1.3076e-004 |          | 1 44      |              |
| 37                | 58.756 | 1   | 4.2    | 1.967e-004  |          | 1 37      |              |
| 38                | 58.905 | 1   | 4.9    | 1.1316e-004 |          | 1 42      |              |
| 39                | 60.351 | 1   | 8.05   | 1.4252e-004 |          | 1 41 + 71 |              |
| 40                | 60.510 | 1   | 6.3    | 9.1908e-005 |          | 1 64      |              |
| 41                | 61.758 | 1   | 3.29   | 1.204e-004  |          | 1 40      |              |
| 42                | 63.198 | 1   | 0.385  | 1.0325e-004 |          | 1 100     |              |
| 43                | 64.142 | 1   | 0.735  | 1.4191e-004 |          | 1 63      |              |
| 44                | 64.781 | 1   | 6.65   | 1.1965e-004 |          | 1 74      |              |
| 45                | 65.576 | 1   | 11.9   | 1.0377e-004 |          | 1 70 + 76 |              |
| 46                | 66.055 | 1   | 18.2   | 2.0515e-004 |          | 1 66      |              |
| 47                | 66.278 | 1   | 7.0    | 1.1703e-004 |          | 1 95      |              |
| 48                | 67.309 | 1   | 1.785  | 1.0315e-004 |          | 1 91      |              |
| 49                | 68.855 | 1   | 12.25  | 1.4375e-004 |          | 1 56 + 60 |              |
| 50                | 69.720 | 1   | 6.3    | 1.713e-004  |          | 1 92 + 84 |              |
| 51                | 70.133 | 1   | 0.35   | 6.4548e-005 |          | 1 89      |              |
| 52                | 70.667 | 1   | 6.3    | 9.7475e-005 |          | 1 101     |              |
| 53                | 71.463 | 1   | 2.59   | 8.1995e-005 |          | 1 99      |              |



## Appendix C-3. (Cont'd)

Method: C:\HPCHEM\1\METHODS\PCB1.MTH

|     |         |   |        |             |   |                  |
|-----|---------|---|--------|-------------|---|------------------|
| 54  | 72.547  | 1 | 0.098  | 5.2775e-005 | 1 | 119              |
| 55  | 73.267  | 1 | 0.525  | 8.7674e-005 | 1 | 83               |
| 56  | 74.158  | 1 | 1.96   | 6.4062e-005 | 1 | 97               |
| 57  | 74.786  | 1 | 0.56   | 1.424e-004  | 1 | 81               |
| 58  | 75.081  | 1 | 3.5    | 9.1511e-005 | 1 | 87               |
| 59  | 75.690  | 1 | 2.45   | 1.1276e-004 | 1 | 85               |
| 60  | 76.087  | 1 | 2.625  | 1.1117e-004 | 1 | 136              |
| 61  | 76.281  | 1 | 20.0   | 2.0671e-004 | 1 | DDE              |
| 62  | 76.557  | 1 | 0.805  | 1.189e-004  | 1 | 77               |
| 63  | 76.807  | 1 | 6.65   | 1.0175e-004 | 1 | 110              |
| 64  | 78.489  | 1 | 1.54   | 8.6984e-005 | 2 | 82 (166)         |
| 65  | 79.137  | 1 | 5.95   | 8.3983e-005 | 2 | 151              |
| 66  | 79.903  | 1 | 3.115  | 1.094e-004  | 2 | 135 + 144        |
| 67  | 80.720  | 1 | 0.455  | 1.1219e-004 | 2 | 107              |
| 68  | 81.247  | 1 | 9.8    | 7.5132e-005 | 2 | 123 + 149        |
| 69  | 81.512  | 1 | 4.2    | 1.0255e-004 | 2 | 118              |
| 70  | 82.777  | 1 | 0.252  | 6.6548e-005 | 2 | 134R             |
| 71  | 83.171  | 1 | 0.49   | 1.6698e-004 | 2 | 114 + 131        |
| 72  | 84.778  | 1 | 1.365  | 9.2233e-005 | 2 | 146              |
| 73  | 85.791  | 1 | 15.05  | 7.9865e-005 | 2 | 132 + 153 + 105  |
| 74  | 87.838  | 1 | 5.95   | 5.3325e-005 | 2 | 141              |
| 75  | 89.077  | 1 | 0.91   | 4.633e-005  | 2 | 137 + 176        |
| 76  | 89.377  | 1 | 0.2625 | 9.995e-005  | 2 | 130              |
| 77  | 90.273  | 1 | 9.45   | 1.0607e-004 | 2 | 163 + 138        |
| 78  | 90.627  | 1 | 0.875  | 1.0594e-004 | 2 | 158              |
| 79  | 91.476  | 1 | 0.0455 | 1.6136e-005 | 2 | 129              |
| 80  | 91.996  | 1 | 3.85   | 1.1609e-004 | 2 | 178              |
| 81  | 92.377  | 1 | 4.76   | 9.2456e-005 | 2 | 166 surrogate    |
| 82  | 92.886  | 1 | 0.7    | 1.0278e-004 | 2 | 175              |
| 83  | 93.430  | 1 | 12.6   | 6.7425e-005 | 2 | 187 + 182        |
| 84  | 94.251  | 1 | 5.95   | 8.6864e-005 | 2 | 183              |
| 85  | 94.673  | 1 | 0.35   | 5.5002e-005 | 2 | 128              |
| 86  | 95.409  | 1 | 0.1715 | 1.1549e-004 | 2 | 167              |
| 87  | 95.834  | 1 | 1.645  | 6.5655e-005 | 2 | 185              |
| 88  | 97.241  | 1 | 11.2   | 7.9148e-005 | 2 | 174              |
| 89  | 98.073  | 1 | 5.95   | 8.9665e-005 | 2 | 177              |
| 90  | 98.802  | 1 | 2.765  | 5.7744e-005 | 2 | 202 + 171        |
| 91  | 99.011  | 1 | 0.231  | 6.1043e-005 | 2 | 156              |
| 92  | 99.560  | 1 | 0.133  | 6.771e-005  | 2 | 173              |
| 93  | 100.141 | 1 | 1.365  | 7.299e-005  | 2 | 157 + 200        |
| 94  | 100.394 | 1 | 6.03   | 7.8638e-005 | 2 | 204 internal std |
| 95  | 101.182 | 1 | 1.96   | 1.1516e-004 | 2 | 172              |
| 96  | 101.376 | 1 | 0.385  | 8.7529e-005 | 2 | 197              |
| 97  | 102.248 | 1 | 21.35  | 8.7835e-005 | 2 | 180              |
| 98  | 102.667 | 1 | 1.47   | 9.8563e-005 | 2 | 193              |
| 99  | 103.291 | 1 | 0.42   | 1.1757e-004 | 2 | 191              |
| 100 | 103.846 | 1 | 1.505  | 7.8345e-005 | 2 | 199              |
| 101 | 106.907 | 1 | 5.95   | 7.8308e-005 | 2 | 170 + 190        |
| 102 | 108.001 | 1 | 0.42   | 5.5328e-005 | 2 | 198              |
| 103 | 108.642 | 1 | 14.7   | 1.2464e-004 | 2 | 201              |
| 104 | 109.427 | 1 | 7.35   | 1.0391e-004 | 2 | 203              |
| 105 | 109.532 | 1 | 7.0    | 1.6017e-004 | 2 | 196              |
| 106 | 111.955 | 1 | 0.14   | 9.8705e-005 | 2 | 189              |
| 107 | 114.145 | 1 | 2.8    | 7.1185e-005 | 2 | 208 + 195        |
| 108 | 115.463 | 1 | 0.3255 | 6.4941e-005 | 2 | 207              |
| 109 | 117.906 | 1 | 6.3    | 8.7518e-005 | 2 | 194              |

## Appendix C-3. (Cont'd)

Method: C:\HPCHEM\1\METHODS\PCB1.MTH

|     |         |   |       |             |   |     |
|-----|---------|---|-------|-------------|---|-----|
| 110 | 118.652 | 1 | 0.385 | 9.8299e-005 | 2 | 205 |
| 111 | 123.947 | 1 | 2.38  | 9.5285e-005 | 2 | 206 |
| 12  | 128.795 | 1 | 0.042 | 5.172e-005  | 2 | 209 |

### Calibration Settings

Title:

|                       |         |
|-----------------------|---------|
| Reference window:     | 0.250 % |
| Non-reference window: | 0.250 % |
| Units of amount:      | ng      |
| Multiplier:           | 1.0     |
| RF uncal peaks:       | 0.0     |
| Sample Amount:        | 0.0     |


### Sample ISTD Information

| I# | Amount |
|----|--------|
| 1  | 9.12   |
| 2  | 6.03   |

### Multilevel Information

|         |        |
|---------|--------|
| Fit:    | Linear |
| Origin: | Force  |

## Appendix C-4. PCB Sample Chromatogram

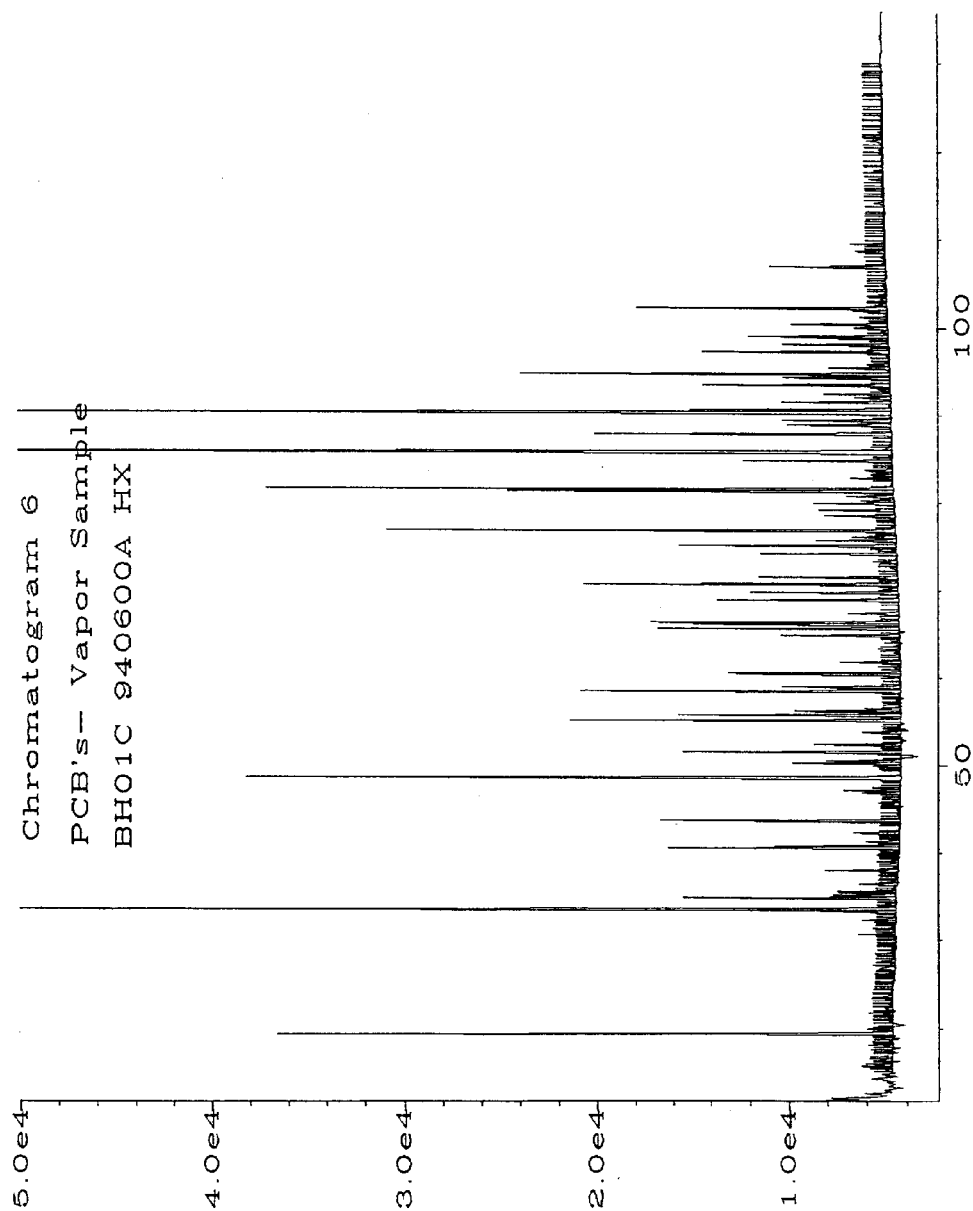

Sig. 1 in C:\HPCHEM\CHECKING\JN94CH\027F0201.D

## Appendix C-5. PCB Sample Report

### Internal Standard Report

```

Data File Name   : C:\HPCHEM\1\DATA\JN94CH\027F0201.D
Operator        : MONTE
Instrument       : ANALYZER1
Sample Name     : BH01C 940600A HX
Run Time Bar Code:
Acquired on    : 12 Nov 94 09:52 AM
Report Created on: 26 Jan 95 11:00 AM
Last Recalib on : 18 JAN 95 02:25 PM
Multiplier     : 1
Page Number    : 1
Vial Number    : 27
Injection Number: 1
Sequence Line  : 2
Instrument Method: MULLIN.MTH
Analysis Method : ISWS'PCB.MTH
Sample Amount   : 0
ISTD Amount     : 8.09
  
```

Sig. 1 in C:\HPCHEM\1\DATA\JN94CH\027F0201.D

| Ret Time | Area          | Type | Width | Ref# | ng      | Name            |
|----------|---------------|------|-------|------|---------|-----------------|
| 18.270   | * not found * |      |       | 1    |         | 1 (14)          |
| 23.510   | * not found * |      |       | 1    |         | 3               |
| 26.667   | * not found * |      |       | 1    |         | 4 + 10          |
| 30.633   | * not found * |      |       | 1    |         | 7 + 9           |
| 32.345   | 12677         | VP   | 0.142 | 1    | 8.013   | 6               |
| 33.256   | 31808         | VV   | 0.126 | 1    | 43.271  | 8 + 5           |
| 33.516   | 412806        | VP   | 0.114 | 1    | 64.487  | HCB             |
| 35.564   | 22467         | VV   | 0.114 | 1    | 22.272  | 14 surrogate    |
| 36.429   | 16004         | PP   | 0.125 | 1    | 3.110   | 19              |
| 37.976   | 26312         | PP   | 0.110 | 1-IR | 8.090   | 30 internal std |
| 39.651   | 1799          | PVA  | 0.123 | 1    | 0.758   | 12              |
| 39.812   | 1324          | VV + | 0.000 | 1    | 0.318   | 13              |
| 40.533   | 96535         | VV   | 0.122 | 1    | 28.841  | 18              |
| 40.775   | 52131         | VV   | 0.124 | 1    | 34.242  | 15 + 17         |
| 42.146   | 3918          | PV   | 0.089 | 1    | 0.801   | 24              |
| 42.291   | 17953         | VP   | 0.116 | 1    | 4.309   | 27              |
| 43.558   | 33925         | PV   | 0.105 | 1    | 12.492  | 16              |
| 43.718   | 96117         | VV   | 0.117 | 1    | 47.309  | 32              |
| 46.040   | 3301          | VP   | 0.131 | 1    | 0.865   | 29              |
| 47.130   | 21909         | VV   | 0.121 | 1    | 8.940   | 26              |
| 47.436   | 10143         | VP   | 0.116 | 1    | 3.874   | 25              |
| 48.597   | 328313        | VVA  | 0.141 | 1    | 135.949 | 31 + 28         |
| 49.908   | 12809         | VV   | 0.137 | 1    | 2.480   | 21              |
| 50.234   | 46002         | VV   | 0.128 | 1    | 17.813  | 33              |
| 50.509   | 30698         | VBA  | 0.126 | 1    | 6.431   | 53              |
| 51.268   | 10702         | PV   | 0.109 | 1    | 1.820   | 51              |
| 51.499   | 77745         | VP   | 0.107 | 1    | 49.184  | 22 (65)         |
| 52.337   | 33277         | VP   | 0.117 | 1    | 8.730   | 45              |
| 53.723   | 15803         | VP   | 0.126 | 1    | 3.998   | 46              |
| 55.119   | 129655        | PV   | 0.118 | 1    | 40.523  | 52              |
| 55.454   | 7688          | VV   | 0.118 | 1    | 2.299   | 43              |
| 55.770   | 84873         | VP   | 0.116 | 1    | 18.999  | 49              |
| 56.143   | 32633         | PV   | 0.100 | 1    | 8.178   | 47              |
| 56.256   | 41157         | VV   | 0.114 | 1    | 8.238   | 48              |
| 56.499   | 12988         | VP   | 0.113 | 1    | 3.233   | 65 surrogate    |
| 58.485   | 120051        | VV   | 0.112 | 1    | 34.990  | 44              |
| 58.807   | 26628         | VV   | 0.111 | 1    | 11.675  | 37              |
| 58.958   | 47016         | VV   | 0.116 | 1    | 11.859  | 42              |
| 60.398   | 64729         | PV   | 0.116 | 1    | 20.564  | 41 + 71         |
| 60.559   | 66442         | VV   | 0.113 | 1    | 13.612  | 64              |
| 61.805   | 22329         | BV   | 0.114 | 1    | 5.993   | 40              |

## Appendix C-5. (Cont'd)

|         |        |     |       |      |         |                  |
|---------|--------|-----|-------|------|---------|------------------|
| 63.240  | 4783   | VVA | 0.116 | 1    | 1.101   | 100              |
| 64.196  | 2246   | PV  | 0.092 | 1    | 0.711   | 63               |
| 64.831  | 47206  | VP  | 0.122 | 1    | 12.590  | 74               |
| 65.628  | 93861  | PV  | 0.118 | 1    | 21.712  | 70 + 76          |
| 66.107  | 80996  | VV  | 0.118 | 1    | 37.039  | 66               |
| 66.322  | 100355 | VV  | 0.121 | 1    | 26.180  | 95               |
| 67.353  | 23132  | VBA | 0.131 | 1    | 5.319   | 91               |
| 68.895  | 75078  | VV  | 0.127 | 1    | 24.056  | 56 + 60          |
| 69.759  | 60811  | VV  | 0.124 | 1    | 23.220  | 92 + 84          |
| 70.168  | 6948   | VV  | 0.133 | 1    | 1.00    | 89               |
| 70.710  | 126376 | VV  | 0.122 | 1    | 27.458  | 101              |
| 71.508  | 56295  | VV  | 0.121 | 1    | 10.289  | 99               |
| 72.586  | 3491   | PV  | 0.146 | 1    | 0.411   | 119              |
| 73.304  | 9134   | PV  | 0.118 | 1    | 1.785   | 83               |
| 74.194  | 53186  | VV  | 0.118 | 1    | 7.595   | 97               |
| 74.825  | 4815   | VV  | 0.111 | 1    | 1.528   | 81               |
| 75.119  | 85185  | VV  | 0.115 | 1    | 17.376  | 87               |
| 75.725  | 32140  | VV  | 0.123 | 1    | 8.078   | 85               |
| 76.121  | 17911  | VV  | 0.127 | 1    | 4.438   | 136              |
| 76.315  | 5312   | VV  | 0.117 | 1    | 2.448   | DDE              |
| 76.585  | 10803  | VV  | 0.123 | 1    | 2.863   | 77               |
| 76.846  | 205206 | VV  | 0.122 | 1    | 46.543  | 110              |
| 78.515  | 27812  | VV  | 0.121 | 2    | 4.722   | 82 (166)         |
| 79.173  | 30603  | VV  | 0.124 | 2    | 5.016   | 151              |
| 79.935  | 34593  | PV  | 0.135 | 2    | 7.386   | 135 + 144        |
| 80.750  | 19518  | VVA | 0.127 | 2    | 4.274   | 107              |
| 81.280  | 158360 | VV  | 0.125 | 2    | 23.222  | 123 + 149        |
| 81.545  | 263380 | VVA | 0.128 | 2    | 52.714  | 118              |
| 82.806  | 16415  | VV  | 0.121 | 2    | 2.132   | 134R             |
| 83.195  | 8442   | VV  | 0.123 | 2    | 2.751   | 114 + 131        |
| 84.812  | 59396  | PP  | 0.122 | 2    | 10.692  | 146              |
| 85.828  | 640338 | PVA | 0.156 | 2    | 99.814  | 132 + 153 + 105  |
| 87.851  | 122682 | PV  | 0.126 | 2    | 12.768  | 141              |
| 89.103  | 8579   | VV  | 0.112 | 2    | 0.776   | 137 + 176        |
| 89.401  | 43527  | VV  | 0.125 | 2    | 8.491   | 130              |
| 90.310  | 643704 | PV  | 0.148 | 2    | 133.262 | 163 + 138        |
| 90.659  | 82819  | VVA | 0.123 | 2    | 17.124  | 158              |
| 91.500  | 43235  | VVA | 0.121 | 2    | 1.362   | 129              |
| 92.025  | 14459  | VV  | 0.126 | 2    | 3.276   | 178              |
| 92.404  | 27226  | VV  | 0.126 | 2    | 4.913   | 166 surrogate    |
| 92.911  | 4843   | VV  | 0.130 | 2    | 0.972   | 175              |
| 93.462  | 74608  | VV  | 0.123 | 2    | 9.818   | 187 + 182        |
| 94.280  | 43770  | BVA | 0.126 | 2    | 7.421   | 183              |
| 94.696  | 150638 | VVA | 0.124 | 2    | 16.171  | 128              |
| 95.429  | 25738  | VV  | 0.128 | 2    | 5.802   | 167              |
| 95.852  | 9196   | VV  | 0.127 | 2    | 1.178   | 185              |
| 97.268  | 76068  | VV  | 0.123 | 2    | 11.751  | 174              |
| 98.100  | 42381  | PV  | 0.124 | 2    | 7.417   | 177              |
| 98.817  | 29417  | PV  | 0.127 | 2    | 3.315   | 202 + 171        |
| 99.033  | 57397  | VVA | 0.125 | 2    | 6.838   | 156              |
| 99.584  | 3408   | VV  | 0.139 | 2    | 0.450   | 173              |
| 100.160 | 3235   | VV  | 0.104 | 2    | 0.461   | 157 + 200        |
| 100.420 | 39158  | VV  | 0.123 | 2-IR | 6.010   | 204 internal std |
| 101.204 | 10858  | VV  | 0.126 | 2    | 2.440   | 172              |
| 101.391 | 1154   | VV  | 0.141 | 2    | 0.197   | 197              |
| 102.272 | 100809 | VV  | 0.120 | 2    | 17.282  | 180              |
| 102.694 | 5435   | VV  | 0.138 | 2    | 1.045   | 193              |
| 103.310 | 3085   | PV  | 0.149 | 2    | 0.708   | 191              |

## Appendix C-5. (Cont'd)

|         |          |       |   |        |           |
|---------|----------|-------|---|--------|-----------|
| 103.866 | 2834 VV  | 0.118 | 2 | 0.433  | 199       |
| 106.925 | 51283 VV | 0.135 | 2 | 7.838  | 170 + 190 |
| 107.987 | 2283 PV  | 0.141 | 2 | 0.247  | 198       |
| 108.661 | 12476 VV | 0.130 | 2 | 3.035  | 201       |
| 109.446 | 7106 BV  | 0.118 | 2 | 1.441  | 203       |
| 109.558 | 6104 VV  | 0.112 | 2 | 1.908  | 196       |
| 111.969 | 1011 VV  | 0.130 | 2 | 0.195  | 189       |
| 114.161 | 4010 PV  | 0.134 | 2 | 0.557  | 208 + 195 |
| 115.458 | 561 VV   | 0.141 | 2 | 0.0710 | 207       |
| 117.920 | 4071 PV  | 0.135 | 2 | 0.695  | 194       |
| 118.720 | 392 PV   | 0.167 | 2 | 0.0752 | 205       |
| 123.962 | 760 PV   | 0.117 | 2 | 0.141  | 206       |
| 128.801 | 153 VV   | 0.100 | 2 | 0.0155 | 209       |

| Time Reference Peak | Expected RT | Actual RT | Difference |
|---------------------|-------------|-----------|------------|
| 10                  | 37.938      | 37.976    | 0.1%       |
| 94                  | 100.394     | 100.420   | 0.0%       |

Not all calibrated peaks were found

=====

## Appendix C-6. PCB Integrator Event Report

Method: C:\HPCHEM\1\METHODS\PCB1.MTH

|                     |        | Integration Events |  |
|---------------------|--------|--------------------|--|
| Events:             | Value: | Time:              |  |
| Initial Area Reject | 200    | INITIAL            |  |
| Initial Peak Width  | 0.040  | INITIAL            |  |
| Shoulder Detection  | OFF    | INITIAL            |  |
| Initial Threshold   | -6     | INITIAL            |  |
| Integrator OFF      |        | 0.000              |  |
| Area Reject         | 100    | 1.150              |  |
| Integrator ON       |        | 15.107             |  |
| Negative Peak ON    |        | 15.107             |  |
| Baseline Now        |        | 19.190             |  |
| Baseline Now        |        | 24.907             |  |
| Baseline Now        |        | 25.960             |  |
| Baseline Hold ON    |        | 25.960             |  |
| Baseline Hold OFF   |        | 26.549             |  |
| Baseline Now        |        | 26.550             |  |
| Baseline Now        |        | 28.887             |  |
| Baseline Now        |        | 30.360             |  |
| Negative Peak OFF   |        | 30.422             |  |
| Baseline Now        |        | 30.513             |  |
| Area Sum ON         |        | 30.517             |  |
| Area Sum OFF        |        | 30.693             |  |
| Negative Peak ON    |        | 30.956             |  |
| Baseline Now        |        | 31.797             |  |
| Baseline Now        |        | 32.950             |  |
| Negative Peak OFF   |        | 32.967             |  |
| Area Sum ON         |        | 32.968             |  |
| Area Sum OFF        |        | 33.129             |  |
| Negative Peak ON    |        | 33.855             |  |
| Baseline Now        |        | 36.273             |  |
| Baseline Now        |        | 37.730             |  |
| Negative Peak OFF   |        | 39.559             |  |
| Area Sum ON         |        | 39.726             |  |
| Area Sum OFF        |        | 39.894             |  |
| Negative Peak ON    |        | 40.085             |  |
| Baseline Now        |        | 43.413             |  |
| Baseline Now        |        | 48.157             |  |
| Negative Peak OFF   |        | 48.160             |  |
| Area Sum ON         |        | 48.837             |  |
| Area Sum OFF        |        | 49.329             |  |
| Negative Peak ON    |        | 49.432             |  |
| Baseline Now        |        | 50.727             |  |
| Baseline Now        |        | 61.550             |  |
| Negative Peak OFF   |        | 63.083             |  |
| Area Sum ON         |        | 63.329             |  |
| Area Sum OFF        |        | 63.554             |  |
| Negative Peak ON    |        | 63.862             |  |
| Baseline Now        |        | 67.777             |  |
| Negative Peak OFF   |        | 74.509             |  |
| Area Sum ON         |        | 74.550             |  |
| Area Sum OFF        |        | 74.756             |  |
| Negative Peak ON    |        | 77.176             |  |
| Negative Peak OFF   |        | 79.801             |  |
| Area Sum ON         |        | 80.950             |  |

## Appendix C-6. (Cont'd)

Method: C:\HPCHEM\1\METHODS\PCB1.MTH

|                   |         |
|-------------------|---------|
| Area Sum ON       | 81.893  |
| Area Sum OFF      | 82.499  |
| Baseline Now      | 82.500  |
| Area Sum ON       | 82.522  |
| Area Sum OFF      | 82.706  |
| Negative Peak ON  | 83.975  |
| Negative Peak OFF | 85.388  |
| Area Sum ON       | 86.060  |
| Area Sum OFF      | 86.408  |
| Negative Peak ON  | 88.425  |
| Negative Peak OFF | 90.025  |
| Area Sum ON       | 90.860  |
| Area Sum OFF      | 91.300  |
| Area Sum ON       | 91.602  |
| Area Sum OFF      | 91.857  |
| Negative Peak ON  | 92.204  |
| Baseline Now      | 93.977  |
| Negative Peak OFF | 93.981  |
| Area Sum ON       | 94.384  |
| Area Sum OFF      | 94.496  |
| Area Sum ON       | 94.876  |
| Area Sum OFF      | 95.234  |
| Negative Peak ON  | 96.107  |
| Negative Peak OFF | 98.456  |
| Area Sum ON       | 99.149  |
| Area Sum OFF      | 99.440  |
| Area Sum ON       | 100.828 |
| Area Sum OFF      | 101.074 |
| Negative Peak ON  | 101.700 |
| Baseline Now      | 103.053 |
| Baseline Now      | 109.260 |
| Baseline Now      | 120.780 |
| Integrator OFF    | 130.000 |

### Calibration Settings

Title:

|                       |         |
|-----------------------|---------|
| Reference window:     | 0.250 % |
| Non-reference window: | 0.250 % |
| Units of amount:      | ng      |
| Multiplier:           | 1.0     |
| RF uncal peaks:       | 0.0     |
| Sample Amount:        | 0.0     |

### Sample ISTD Information


| I# | Amount |
|----|--------|
| 1  | 9.12   |
| 2  | 6.03   |

### Multilevel Information

Fit: Linear



## Appendix C-7. Pesticide 20 ng/mL Standard Chromatogram

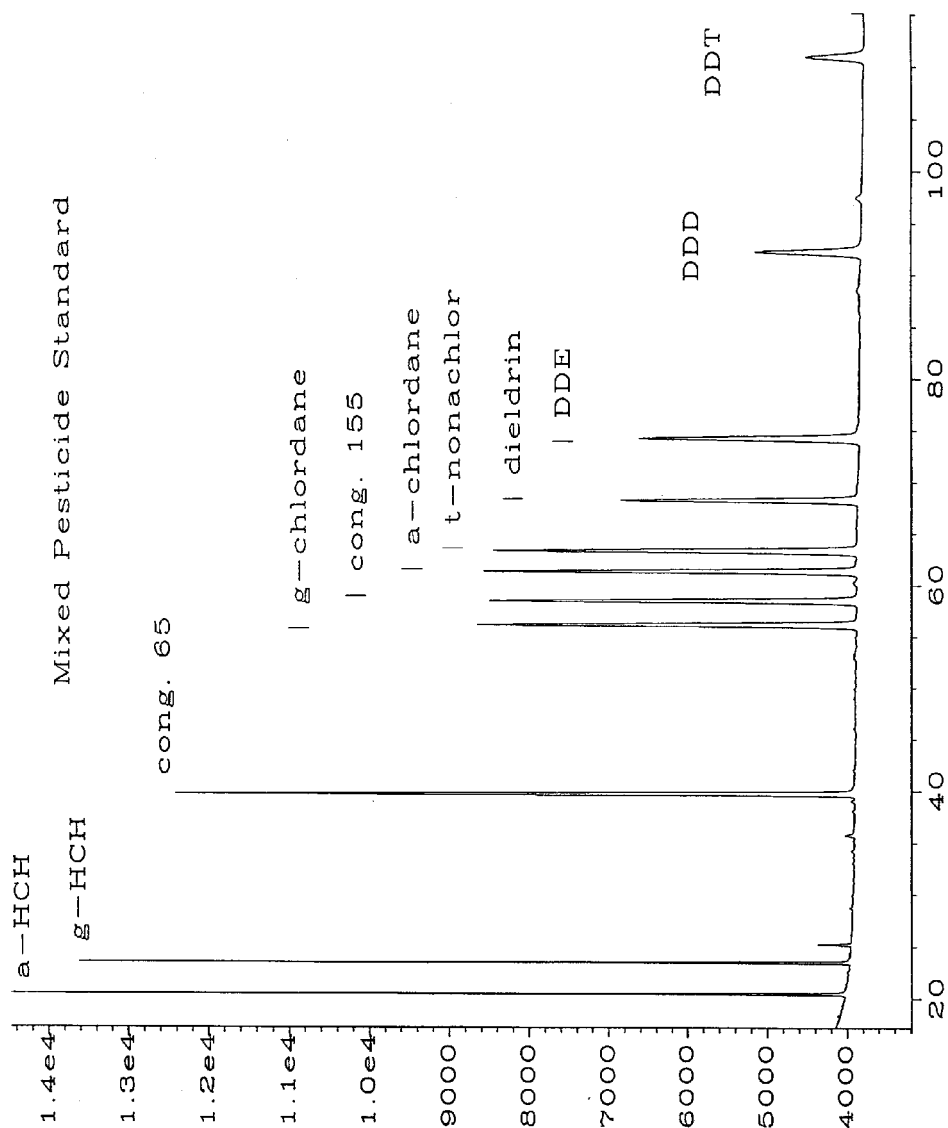


Sig. 1 in C:\HPCHEM\CHECKING\9503034F\043F0201.D

## Appendix C-8. Pesticide Calibration Table

Method: C:\HPCHEM\1\METHODS\PEST1.MTH

| Calibration Table |         |     |       |             |     |      |    |               |  |
|-------------------|---------|-----|-------|-------------|-----|------|----|---------------|--|
| Pk#               | RT      | Lvl | ng    | Amt/Area    | Ref | Istd | I# | Name          |  |
| 1                 | 20.477  | 1   | 20.0  | 3.4963e-004 |     |      | 1  | A-HCH         |  |
| 2                 | 23.517  | 1   | 20.0  | 3.9724e-004 |     |      | 1  | B-HCH         |  |
| 3                 | 39.765  | 1   | 4.74  | 4.2976e-005 |     |      | 1  | CONG 65(I.S)  |  |
| 4                 | 56.179  | 1   | 20.0  | 2.2182e-004 |     |      | 1  | G-CHLORDANE   |  |
| 5                 | 58.529  | 1   | 14.99 | 1.5123e-004 |     |      | 1  | CONG 155(I.S) |  |
| 6                 | 61.385  | 1   | 20.0  | 2.0852e-004 |     |      | 1  | A-CHLORDANE   |  |
| 7                 | 63.406  | 1   | 20.0  | 2.0738e-004 |     |      | 1  | T-NONACHLOR   |  |
| 8                 | 68.317  | 1   | 20.0  | 2.016e-004  |     |      | 1  | DIELDRIN      |  |
| 9                 | 74.299  | 1   | 20.0  | 2.6304e-004 | Ref | ISTD | 1  | DDE(I.S)      |  |
| 10                | 92.310  | 1   | 20.0  | 4.7227e-004 |     |      | 1  | DDD           |  |
| 11                | 111.059 | 1   | 20.0  | 8.9441e-004 |     |      | 1  | DDT           |  |

## Appendix C-9. Pesticide Sample Chromatogram

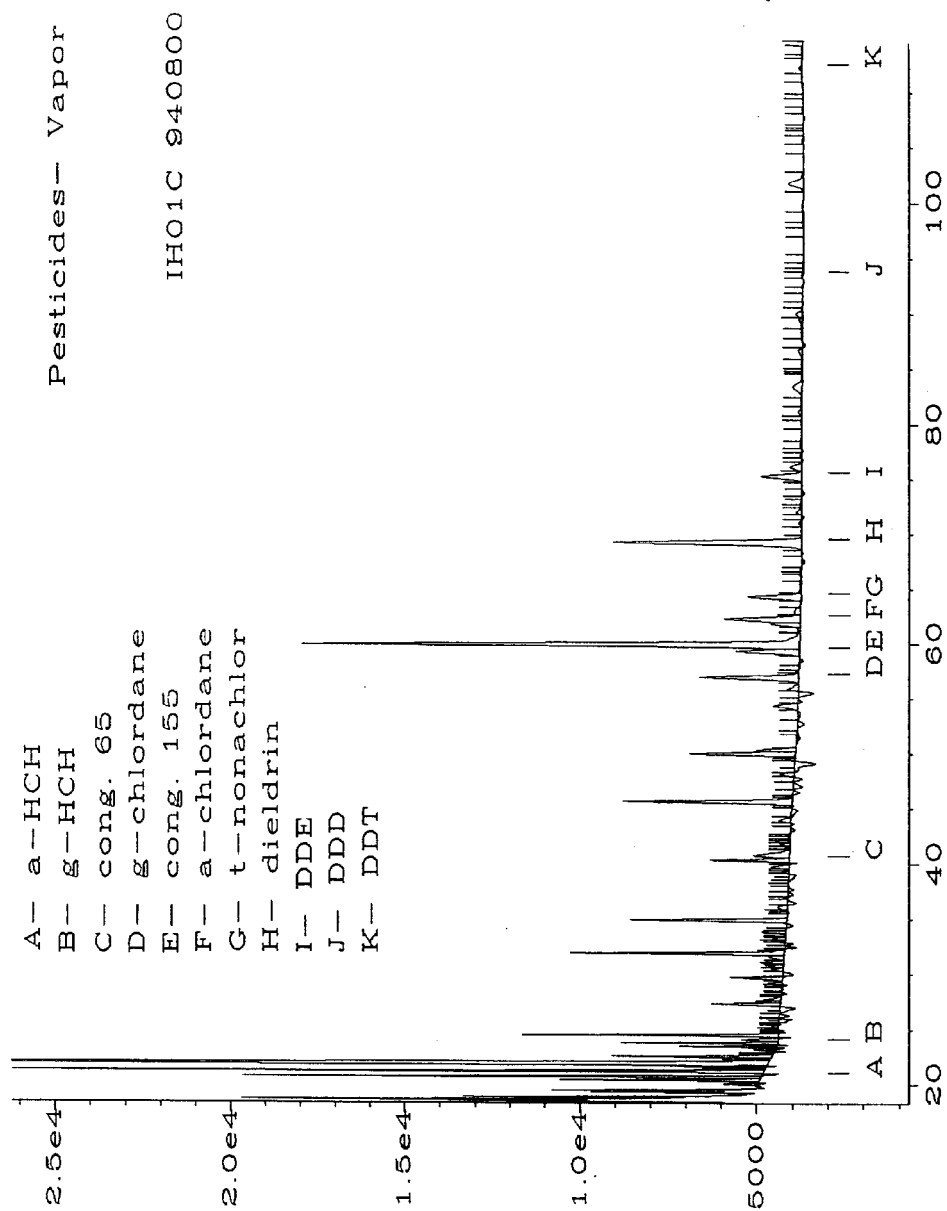

Sig. 1 in C:\HPCHEM\CHECKING\941025C4\049FO201.D

## Appendix C-10. Pesticide Sample Report

```

=====
                        Internal Standard Report
=====
Data File Name   : D:\HPCHEM\2\DATA\941025C4\049F0201.D
Operator        : MONTE
Instrument       : ANALYZER2
Sample Name     : IH01C 940800 40
Run Time Bar Code:
Acquired on     : 06 Nov 94 02:21 PM
Report Created on: 14 Feb 95 03:44 PM
Last Recalib on : 14 FEB 95 03:19 PM
Multiplier     : 1
Page Number     : 1
Vial Number     : 49
Injection Number: 1
Sequence Line   : 2
Instrument Method: ISWS'PES.MTH
Analysis Method : PEST2.MTH
Sample Amount   : 0
ISTD Amount     : 20

Sig. 1 in D:\HPCHEM\2\DATA\941025C4\049F0201.D
Ret Time   Area   Type Width Ref#   ng      Name
-----|-----|-----|-----|-----|-----|-----|
20.676     77562 UP    0.080 1    51.647 A-HCH
23.773     27066 PU    0.094 1    19.116 G-HCH
40.294     25633 PU    0.181 1    13.009 CONG 65
56.960     57957 UV    0.320 1    28.780 G-CHLORDANE
59.336     35618 PU    0.311 1    15.320 CONG 155
62.249     50086 UV    0.347 1    23.779 A-CHLORDANE
64.293     32679 PP    0.340 1    15.738 T-NONACHLOR
69.273     127351 PU    0.371 1    92.475 DIELDRIN
75.279     29322 PU    0.373 1-IR 20.000 DDE
93.453      936 PP    0.278 1     1.133 DDD
112.349     2646 PP    0.374 1     5.127 DDT

```

| Time | Reference Peak | Expected RT | Actual RT | Difference |
|------|----------------|-------------|-----------|------------|
| 9    |                | 75.220      | 75.279    | 0.1%       |

## Appendix C-11. Pesticide Integrator Event Report

Method: C:\HPCHEM\1\METHODS\PEST2.MTH

### Integration Events

| Events:             | Value: | Time:   |
|---------------------|--------|---------|
| Initial Area Reject | 100    | INITIAL |
| Initial Peak Width  | 0.040  | INITIAL |
| Shoulder Detection  | OFF    | INITIAL |
| Initial Threshold   | -6     | INITIAL |
| Integrator OFF      |        | 0.000   |
| Integrator ON       |        | 19.845  |
| Negative Peak ON    |        | 20.000  |
| Baseline Now        |        | 23.647  |
| Baseline Now        |        | 25.503  |
| Baseline Now        |        | 28.757  |
| Baseline Now        |        | 35.303  |
| Baseline Now        |        | 43.350  |
| Baseline Now        |        | 52.060  |
| Baseline Now        |        | 68.057  |
| Integrator OFF      |        | 115.000 |

### Calibration Settings

Title:

|                     |         |
|---------------------|---------|
| Reference window:   | 0.250 % |
| 1-reference window: | 0.250 % |
| Units of amount:    | ng      |
| Multiplier:         | 1.0     |
| RF uncal peaks:     | 0.0     |
| Sample Amount:      | 0.0     |

### Sample ISTD Information

| I# | Amount |
|----|--------|
| 1  | 4.74   |

### Multilevel Information

Fit: Linear  
Origin: Force